# DURONIC

Creative Consumer Electronics

## Bathroom Scales - BS901

PLEASE READ ALL THE INSTRUCTIONS BEFORE USE

[duronicon.com](http://duronicon.com)

[duronicon](https://www.youtube.com/channel/UCduronicon)

[@duroniconmedia](https://www.instagram.com/duronicon)

[Duronicon](https://www.facebook.com/duronicon)

[@duroniconmedia](https://twitter.com/duronicon)

## Contents

<b>English Manual</b>	<b>5</b>
<b>Manuel d'utilisation en Français</b>	<b>11</b>
<b>Gebrauchsanleitung auf Deutsch</b>	<b>17</b>
<b>Manual de Instrucciones en Español</b>	<b>23</b>
<b>Manuale D'istruzioni In Italiano</b>	<b>29</b>

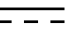
**English Manual**

**Table of Contents**

Safety Information	8
Your Scale and Its Environment	9
Efficient Use of Your Scale	10
Overview	10
Device Components	10
LCD Display	11
Initial Start-Up	11
General Instructions	11
Insert the Batteries	12
Install the App and Pair-up	12
Select Measurement Unit	14
Start Measuring	14
Daily Measurement	14
Data Management	15
Data Transmission	15
Troubleshooting	16
Error Prompt	16
When Measuring	16
When Data Transmitting	17
Specifications	18
Maintenance	19
FCC Regulations	19
Appendix	19
Table of Body Fat Level	20
Table of Body Water Level	21
Health Tips – About Body Fat	21
EMC Guidance	24
Warranty	26

## Safety and Usage Information

The warning signs and symbols are essential to ensure a correct and safe use of this product and protect you and others from injury. Please find the meanings of the warning signs and symbols, which you may encounter in the label and user manual, as follows:

	This operation guide must be read		Manufacturer
	Bluetooth		Complies with MDD 93/42 ECC Requirements
	Type of applied parts		Manufacturer date
<b>SN</b>	Serial number		ENVIRONMENT PROTECTION - Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice
	Direct current		
	Indoor use only		

### CAUTION

Our Glass Body Fat Analyzer BS901 offers you a seamless way to manage your health. Please be aware that this device is designed for adults' self-measuring and self-monitoring body fat level. Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening. This device should not be used by anyone who is acutely or chronically ill, suffering from a disease or taking medications that affect your water levels. The accuracy of readings for these patients

has not been verified. Specific medical advice should be obtained from a physician. BS901 is equipped with automatic data transmission function. It may emit electromagnetic energy so as to perform its intended function. Nearby portable and mobile RF communications equipment can affect the performance of BS901. Kindly note that the use of accessories, transducers or cables other than those specified, with the exception of transducers and cables sold by the manufacturer as replacement parts for internal components, may result in increased EMISSIONS or decreased IMMUNITY of BS901. Be aware that misuse of electrical equipments can cause electric shock, burns, fire and other hazards.

### INDICATIONS FOR USE

- Our Glass Body Fat Analyzer measures weight and uses bio-electrical impedance analysis (BIA) technology to estimate body fat, total body water percentage, bone mass, and muscle mass in generally healthy adults 18 years of age or older.
- It is intended for domestic use only.

### CONTRAINDICATIONS

1. This device is contraindicated for any female subject who may be suspected of, or is pregnant. Besides providing inaccurate readings, the affects of this device on the fetus are unknown.
2. This device is contraindicated for any person who is connected to a wearable or implantable electronic device or instrument such as a pacemaker or defibrillator.

### YOUR SCALE AND ITS ENVIRONMENT

To ensure your safety as well as the service life of your scale, please avoid using the scale under the following circumstances:

- With any other medical devices.
- On slippery floors such as a tile floor
- Jumping onto the platform immediately after bath or with wet hands
- Near a cell phone or microwave oven


Avoid storage in the following locations:

- Near to water
- Where the device may be exposed to extreme temperatures, humidity, moisture, direct sunlight, dust, or salt air
- An unstable surface
- In the reach of children

Efficient Use of Your Scale

- Place the scale on a flat, hard surface. Soft surfaces such as carpet will affect the performance of the scale.
- Step onto the scale with bare feet. Stand still and keep full contact with the electrodes until the measurement is complete.
- Start measurement at least two hours after getting up or eating.
- Avoid measurement immediately after strenuous exercise, a sauna or bath, drinking and eating.
- Always start the measurement at approximately the same time and on the same scale, on the same flat, hard surface.


Device Components

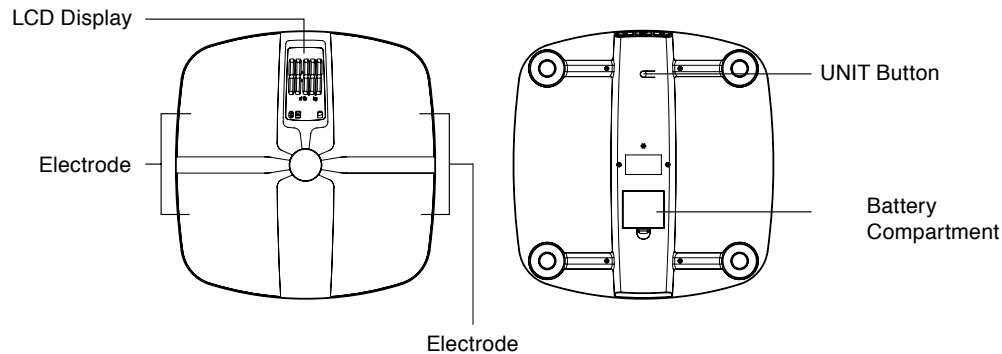List

1. Glass Body Fat Analyzer BS901
2. Four AAA-size Batteries (1.5V each)
3. User Manual

LCD Display


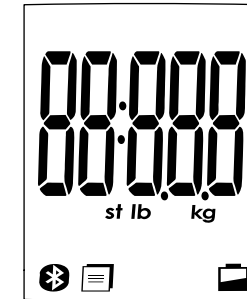
st lb	Stone		Successful Bluetooth Connection
lb	Pound		Data transmitting/pending to transmit to wireless device
kg	Kilogram		Low Battery

General Instructions

- Our Glass Body Fat Analyzer BS901 applies BIA (Bio-impedance Analysis) technology. A small amount of weak current flows through the human body so as to detect the bio-impedance and estimate body fat. The electrical current is small and may not be felt.
- This BIA technology is cheap, safe, non-invasive, toxic-free and harmless. It also simple to use and provides a lot of information.
- The current mentioned above is less than 1mA. However, please be aware that anyone with an wearable or implantable medical electronic instrument, such as a pacemaker, must avoid using this device.
- The intended use of this device is for adult and indoor use only.

Insert the Batteries

- Open the battery compartment in the back of the scale.
- Insert the batteries (4 x 1.5V AAA) into the battery compartment according to the polarity indications marked inside the compartment.
- \* The digits “8888” will be shown on the LCD.
- Close the battery door and wait until the digits “0.0kg” are shown on the LCD.

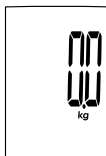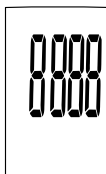CAUTION

When the symbol  appears, the device will turn off in four seconds. You should then replace the batteries. Please replace all four batteries at the same time. Do NOT mix the old batteries with the new ones. Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage. Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries.


Install the App and Pair-Up

With the advanced Bluetooth 4.0 technology, mobile or portable devices which are equipped with Bluetooth function in line with BLE Technical Specifications as well as BLP Protocol established by global organization Bluetooth SIG, are capable to receive your personal health data. Just simply install the specially-designed app and pair up your scale with your mobile or portable equipments.


The free App “TransHealth” is now available on the App Store. Simply search for the app, then download and install the app in your iPhone. Turn on Bluetooth and the App. Make sure both are ON when pair-up is proceeding. Press and hold “UNIT” button in the back of the scale to start pair-up.


Symbol  and symbol  will be shown on the LCD alternatively, indicating pair-up is proceeding.

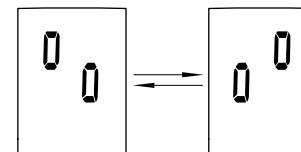If the connection is successful, the symbol [ ] will be shown on the LCD.

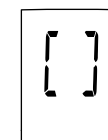If unsuccessful, the symbol “E1” will be shown on the LCD.

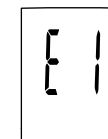
**Bluetooth Module No. : nRF8001**

Frequency Range	2.402 - 2.480 GHz	Supply Voltage	3.3V
Output Power Range	0 dB Max.	Transmitting Distance	20 Metres


### Selecting the measurement Unit

With batteries correctly installed, press “UNIT” button in the back of the scale to select measurement unit. The default measurement unit is “kilogram”. You may press the “UNIT” button to choose between kilogram, stone and pound.
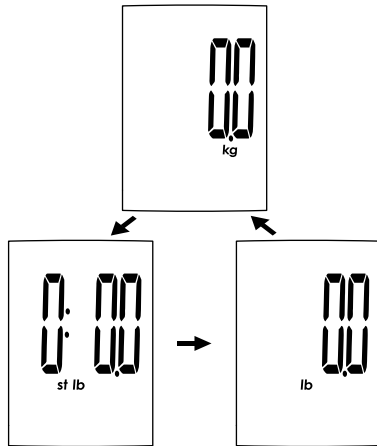
### Daily Measurement

With original SENSE ON patent technology, the BS901 will automatically switch on as you step on the platform barefooted.

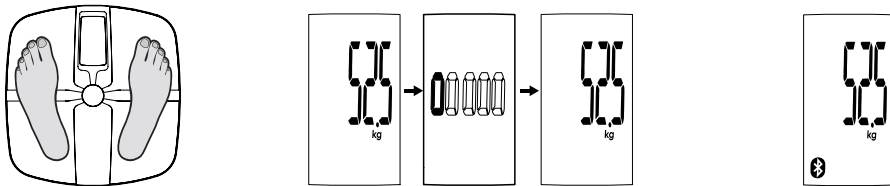
Stand still and keep full contact with the electrodes until the LCD stops displaying moving “o”.

When your scale is successfully paired with your Bluetooth device and the Bluetooth is ON, BS901 will process data transmission automatically. (Please refer to Data Transmission for more details.)

### Data Transmission

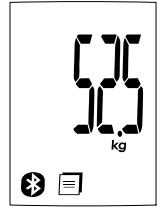When the connection is successful, the measurement data will be automatically transmitted to your device via Bluetooth.

The symbol  will disappear after successful data transmission and you may check your personal health data stored in your Bluetooth device.

If the data transmission fails, the symbol  will remain. The pending measurement data will be temporarily kept in the scale and transmitted to your Bluetooth device when next measurement is complete.

	Successful Bluetooth Connection
	Data transmitting to wireless wellness system: - If SUCCESSFUL, the symbol disappears; - If UNSUCCESSFUL, the symbol remains.

### CAUTION

Interference may occur in the vicinity of equipment marked with the following symbol .

The scale may interfere with electrical equipment in its vicinity.


To enable the data transmission function, this product should be paired to a Bluetooth device at 2.4 GHz.

How to mitigate possible interference?

1. The range between the scale and the Bluetooth device should be reasonably close, from 1 to 10 meters. Please ensure there are no obstacles between the scale and the Bluetooth device so as to obtain quality connection.
2. To avoid interference, other electronic devices (particularly those with Bluetooth transmission / Transmitter) should be kept at least 1 meter away from the scale.

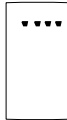Error Prompt

Error	Description	Solution
	Overload. The device will power off in four seconds.	Stop using this scale for measurement.
	Low battery. The device will power off in four seconds.	Replace all four batteries at the same time. Please purchase the authorized batteries for replacement.
	Failed to pair up your scale with your bluetooth device	Please check the items below: - Bluetooth is ON. - App Collector is ON. - Both devices are within the transmission distance of Bluetooth.

When Measuring

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; OR - Too low; OR - Huge difference between two recent measurement.	Incorrect posture	Please step on the platform barefooted and stand still.
	The device is located on a soft surface such as a carpet or rug.	Please place the device on a flat, hard surface.
	Cold body that may results in bad blood circulation.	Warm up your hands and feet to resume blood circulation and then measure again.

Problem	Root Cause	Solution
Abnormal measuring results: - Too high; OR - Too low; OR - Huge difference between two recent measurement.	Cold electrodes	Place the device in a warm room for a while and then measure again
	Either your hands or your feet are too dry	Wipe your feet with a damp cloth, keeping them slightly damp when starting measurement.
No display on LCD when the device powers on.	Batteries not yet installed	Install the batteries. (Please refer to Insert the Batteries)
	Worn batteries	Replace all four batteries in the same time. Please purchase the authorised batteries for replacement.
CANNOT proceed to analyze body fat, total body water, muscle mass and bone mass.	Step onto the platform wearing socks or shoes.	Please use barefooted and keep full contact with the electrodes during the measurement.
The device powers off automatically.	Low battery	Replace all four batteries in the same time. Please purchase the authorised batteries for replacement.

When Data Transmitting

Problem	Root Cause	Solution
Data transmission failed.	Bluetooth is OFF	Turn ON the Bluetooth via "Setting >> General >> Bluetooth".
	The App is OFF	Press the icon to turn ON your app.
	Out of range of Bluetooth transmission	Place your Bluetooth device closer to the scale.

## Specifications

Product Name	Glass Body Fat Analyzer (BS901)
Dimension	Scale: 320x330x27.5mm
	Panel: 320x330x6mm
Net Weight	Approximately 1.75kg ( Excluding the dry cells)
Display	Digital LCD with White Backlight V.A.=94mm(L) x 52mm(W)
Measurement Unit	Kilogram / Stone/ Pound
Measurement Range	5kg to 180kg / 0st: 11lb to 28: 5lb / 11lb to 396lb
Division	0.1kg / 0.2lb
Accuracy	0-50kg: ±0.3kg;      50-100kg: ±0.4kg; 100-150kg: ±0.5kg;    150-180kg: ±0.7kg
Working Environment	Temperature: 0°C to 40°C Humidity: ≤90% RH
Storage Environment	Temperature: -20°C to 60°C      Humidity: 10%RH to 93% RH
Power Source	6V (Four AAA-size Batteries)
Auto-ON	SENSE ON technology
Auto-OFF	8 to 12 seconds after clearing Otherwise, 13-17 seconds after last operation
Accessories	1. Four AAA-size Batteries 2. User Manual
Mode of Operation	Continuous Operation
Protection Against Ingress of Water	IPX0
Software Version	1.0

### About the Accuracy of This Product

This product passes strict inspection before delivery and therefore its accuracy is guaranteed by the manufacturer. Please refer to the above table for the descriptions on accuracy. This product is specially designed for body fat analysis as well as weight measurement. It should NOT be used by anyone during the process of transaction for verification of goods' weight.

## Maintenance

Please follow the following steps for correct maintenance of your scale:  
 DO use a dry soft cloth to wipe the dust.  
 DO use a dry soft cloth to dry the device.  
 DO NOT wash the device with water or immerse it in water.  
 DO NOT use propellant, abrasive or other chemicals.  
 DO NOT disassemble this device. If you have any problems please contact our customer support. (Please refer to Warranty for contact information)

## FCC Regulations

### FCC User Guide Information

#### Radio Frequency Interface Requirements - FCC

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna;
- Increase the separation between the equipment and receiver;
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- Consult the dealer or an experienced radio / TV technician for help.

Radio Transmitters (Part 15)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC RF Exposure Guidelines

Safety Information

Reducing RF Exposure - Use Properly

Only operate the device in accordance with the instructions supplied. This device complies with FCC radiation exposure limits set forth for an uncontrolled environment. In order to avoid the possibility of exceeding the FCC radio frequency exposure limits, human proximity to the antenna shall not be less than 20 cm (8 inches) during normal operation.

Table of Body Fat Level (Unit: %)

Age	Man / Sportsman				Woman / Sportswoman			
	Low	Normal	High	V. High	Low	Normal	High	V. High
20-29	<13	13.1-20	20.1-23	>23	<19	19.1-28	28.1-31	>31
30-39	<14	14.1-21	21.1-24	>24	<20	20.1-29	29.1-32	>32
40-49	<16	16.1-23	23.1-26	>26	<22	22.1-30	30.1-33	>33
50-59	<17	17.1-24	24.1-27	>27	<23	23.1-31	31.1-34	>34
60+	<18	18.1-25	25.1-28	>28	<24	24.1-32	32.1-35	>35

Table of Body Water Level (Unit: %)

Gender	Body Fat Percentage Range	Optimal Total Body Water Percentage Range
Men	4-14	70-63
	15-21	63-57
	22-24	57-55
	≥25	55-37
Women	4-20	70-58
	21-29	58-52
	30-32	52-49
	≥33	49-37

Source: Derived fr. Wang & Deurenberg: "Hydration of fat-free body mass". American Journal Clin Nutr 1999, 69833-841

Health Tips - About Body Fat

Fat is essential for the human body. It can not only store energy and protect viscera, but also regulate body temperature and maintain the normal physiological function of human body. However, too much body fat is harmful to human body. It can be accompanied by Fatty Liver, diabetes, coronary heart disease, etc.

Therefore self-measuring and self-monitoring body fat level are beneficial to your health. Since we cannot judge body fat level simply by our weight, this body fat analyzer BS901, with BIA (Bio-impedance Analysis) technology applied, is an accurate device that offers a quick and comfortable way to obtain your body fat level.

R&TTE Declaration of Conformity (DoC)

Hereby, Guangdong Transtek Medical Electronics Co., Ltd, declares that this BS901 is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. The declaration of conformity may be consulted at <http://www.transtek.cn/DoC.pdf>.

EMC Guidance

Table 1 – Guidance and MANUFACTURER’S declaration – ELECTROMAGNETIC EMISSIONS – for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacturer’s declaration – electromagnetic emissions		
The BS901 is intended for use in the electromagnetic environment specified below. The customer or the user of the BS901 should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment – guidance
RF Emission CISPR 11	Group 1	The BS901 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF Emission CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not Applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not Applicable	

Table 2 – Guidance and MANUFACTURER’S declaration – electromagnetic IMMUNITY – for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacture’s declaration – electromagnetic immunity			
The BS901 is intended for use in the electromagnetic environment specified below. The customer or the user of the BS901 should assure that it is used in such an environment .			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floor are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines	Not applicable	
Surge IEC 61000-4-5	±1 kV line(s) to line(s)	Not applicable	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT (>95% dip in UT) for 0.5 cycle	Not applicable	
	40% UT (60% dip in UT) for 5 cycles		
	70% UT (30% dip in UT) for 25 cycles		
	<5% UT (>95% dip in UT) for 5 sec		
Power frequency (50Hz) magnetic field IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE U T is the a.c. mains voltage prior to application of the test level.			

Table 4 Guidance and manufacture’s declaration – electromagnetic immunity –for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacture’s declaration – electromagnetic immunity			
The BS901 is intended for use in the electromagnetic environment specified below. The customer or the user of the BS901 should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the BS901, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> $d = 1.167 \sqrt{P}$ $d = 1.167 \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = 2.333 \sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
Radiated RF IEC 61000-4-3	3 V/m 80MHz to 2.5 GHz	3 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM

and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the BS901 is used exceeds the applicable RF compliance level above, the BS901 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the BS901.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Table 6 – Recommended separation distances between portable and mobile RF communications equipment and the ME EQUIPMENT or ME SYSTEM – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Recommended separation distances between portable and mobile RF communications equipment at the BS901			
The BS901 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the BS901 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the BS901 as recommended below, according to the maximum output power of the communications equipment			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1.167 \sqrt{P}$	80 MHz to 800 MHz $d = 1.167 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.333 \sqrt{P}$
0.01	N/A	0.117	0.233
0.1	N/A	0.369	0.738
1	N/A	1.167	2.333
10	N/A	3.690	7.377
100	N/A	11.67	23.33

For transmitters rated at a maximum output power not listed above, the recommended separation distance *d* in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where *P* is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

## Warranty

1 Year Guarantee from Shine-Mart Ltd, Trademark owners of Duronic Products.

NOTE: THESE STATEMENTS DO NOT AFFECT YOUR STATUTORY RIGHTS AS A CONSUMER

This product has been manufactured under the strictest or quality control procedures, and using the highest quality of materials, to ensure excellent performance and reliability. It will give very good and long lasting service, provided it is properly used and maintained.

The product is guaranteed for 1 years from the date of original purchase. If any defect arises due to a faulty materials or poor workmanship, the faulty product must be returned to the original place of purchase. Refund or replacement is at the discretion of that company. Duronic Products are offered with a 1 year guarantee under the following conditions:

1. The product must be returned to the retailer with original proof of purchase.
2. The product must be installed and used in accordance with the instructions contained in this operating manual.
3. It must be used only for domestic purposes.
4. It does not cover wear and tear, malicious damage, misuse, neglect, inexpert repairs or consumable parts.
5. Shine-Mart Ltd has no responsibilities for accidental or consequential loss or damage.
6. Shine-Mart Ltd are not liable to carry out any type of servicing work, under the guarantee
7. Valid in the UK only

This product is fitted with BS:1363 plug and fuse  
Imported by Shine-Mart Ltd, RM3 8SB

Information on waste Disposal for Consumers of Electrical & Electronic Equipment.

This mark on a product and/or accompanying documents indicates that when it is to be disposed of, it must be treated as Waste Electrical & Electronic Equipment, (WEEE).

Any WEEE marked products must not be mixed with general household waste, but kept separate for the treatment, recovery and recycling of the materials used.

For proper treatment, recovery and recycling; please take all WEEE marked waste to your Local Authority Civic waste site, where it will be accepted free of charge.

If customers dispose of Waste Electrical & Electronic Equipment correctly, they will be helping to save valuable resources and preventing any potential negative effects upon human health and the environment, of any hazardous materials that the waste may contain.

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1 Jahr Gewährleistung von Shine-Mart Ltd, der Markeninhaber von Duronic Produkte.
5. Shine-Mart Ltd hat keine Verantwortung für indirekte zufällige Verluste oder Folgeschaden.
6. Shine-Mart Ltd sind nicht verantwortlich Wartungsarbeiten unter die Gewährleistung zu machen
Von Shine-Mart Ltd, RM3 8SB importiert