
miroco

LIGHT THERAPY LAMP

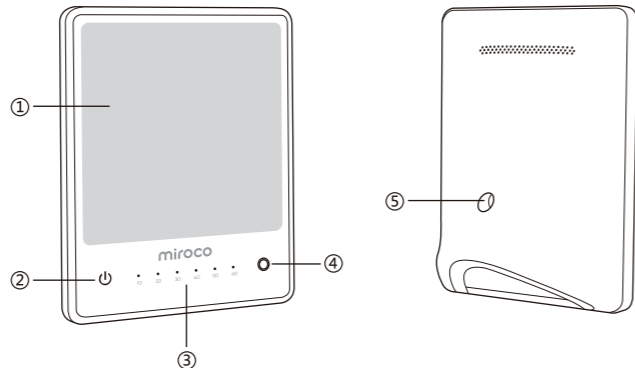EN

Product Diagram

1. Lamp
2. Power Button
3. Timer Buttons
4. Light Brightness Button
5. Power Port

DE

Produktabbildung

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EN

Features

Your new Miroco light therapy lamp is a remarkable product with many features and benefits:

- Delivers LED light to your home, office, dorm room or apartment.
- Designed with a stable stand for easy placement on a convenient spot.
- Cues your body's own energy enhancers to help you regain focus and concentration.
- UV free LED- the built-in LED light panel delivers the full brightness with no pollution, no radiation and no stroboscopic effects.
- 10,000 lux at 10cm/4" from lamp to lift up your spirit.

Important

1. Some people may suffer from light sensitivity; discontinue use immediately if any discomfort occurs.
2. Consult your health care provider before use of this product if you are under medical care or have suffered from an eye condition before.
3. Never cover the lamp or place anything on top of it during operation to avoid overheating or fire hazards.

Specifications

Model	MI-CL003 Pro
Voltage	AC 100-240V 50 / 60Hz DC 12V 1A
Power	11W Max
Storage Temperature	-20-50°C / -4-122°F

USING YOUR LIGHT THERAPY LAMP

Power On / Off

Tap the Power Button to turn on / off the lamp.

Note: If power is not disconnected, the lamp will remember the light brightness of last use when turned on again.

Adjust the Light Brightness

When the lamp is on, tap the Light Brightness Button to set the light brightness between low / semi medium / medium / high levels.

Set the Timer

1. Tap the Timer Button to set the timer from 10 minutes to 60 minutes. The timer button indicators will light up accordingly.
2. Tap and hold the set Timer Button to turn off the timer.

Warning

1. Please keep the lamp away from water or any other liquid to avoid electrical shock.
2. For indoor use only.
3. Close supervision is needed when used by or near children.
4. Keep away from heat and avoid direct sunlight.
5. Don't try to repair or dismantle the lamp by yourself.
6. Don't drop or knock the lamp on hard surfaces.

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FCC Compliance

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WEEE Compliance

Correct Disposal of This Product (Waste Electrical & Electronic Equipment)
(Applicable in countries with separate collection systems)

This marking on the product, accessories or literature indicates that the product and its electronic accessories should not be disposed of with other household waste at the end of their working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate these items from other types of waste and recycle them responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take these items for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product and its electronic accessories should not be mixed with other commercial wastes for disposal.

<https://miroco.net>

 support@mirocohome.com  <https://miroco.net>

MANUFACTURER

Manufacturer: Shenzhen NearbyExpress Technology Development Co.,Ltd.
Address: Floor 7, Building E, Galaxy World Phase II, Shenzhen, China

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