
LIGHT THERAPY LAMP
User Guide

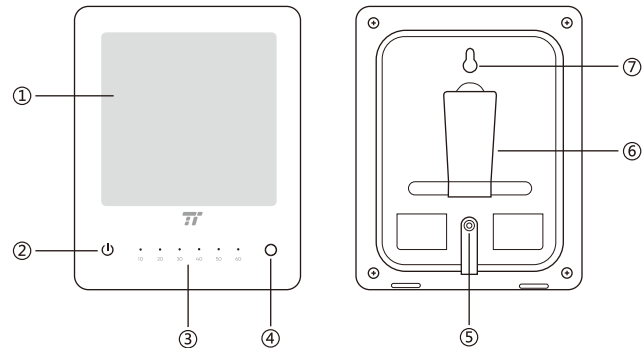EN

Product Diagram

1. Lamp
2. Power Button
3. Timer Buttons
4. Light Brightness Button
5. Power Port
6. Supporting Foot
7. Hanging Inset

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EN

Important



1. This product contains magnets. Keep a proper distance if you are wearing a pacemaker or any other implants.
2. Some people may suffer from light sensitivity; discontinue use immediately if any discomfort occurs.
3. Consult your health care provider before use of this product if you are under medical care or have suffered from an eye condition before.
4. Never cover the lamp or place anything on top of it during operation to avoid overheating or fire hazards.

Specifications

Model	TT-CL011
Voltage	AC 100 / 240V 50 / 60Hz DC 12V 1.0A
Power	13W Max
Storage Temperature	-20-50°C / -4-122°F

How to Use

1. Put the lamp on a flat and stable surface using the support foot or hang it up through the hanging inset.

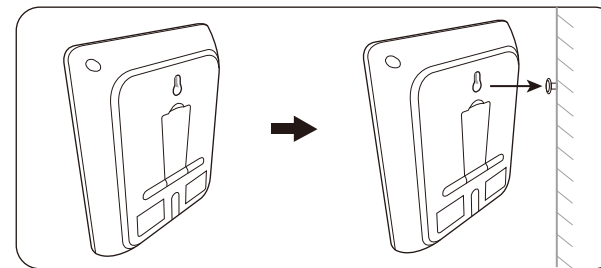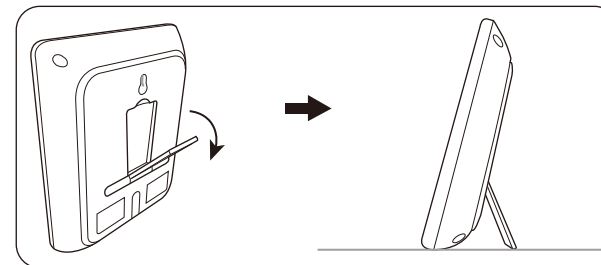2. Connect to a power source with the included power adapter. The unit will enter standby mode, and the Power Button indicator will light up.

Power On / Off

Tap the Power Button to turn on / off the lamp.

Note: If power is not disconnected, the lamp will remember the light brightness of last use when turned on again.

Adjust the Light Brightness

When the lamp is on, tap the Light Brightness Button to set the light brightness between low / medium / high levels.

Set the Timer

1. Tap the Timer Button to set the timer from 10 minutes to 60 minutes. The timer button indicators will light up accordingly.
2. Tap and hold the set Timer Button to turn off the timer.

Warning

1. Please keep the lamp away from water or any other liquid to avoid electrical shock.
2. For indoor use only.
3. Close supervision is needed when used by or near children.
4. Keep away from heat and avoid direct sunlight.
5. Don't try to repair or dismantle the lamp by yourself.
6. Don't drop or knock the lamp on hard surfaces.

TT-CL011

TT TaoTronics
www.taotronics.com

NORTH AMERICA

E-mail : support@taotronics.com(US)
support.ca@taotronics.com(CA)
Tel : 1-888-456-8468 (Monday-Friday: 9:00 – 17:00 PST)
US Distributor: SUNVALLEYTEK INTERNATIONAL INC.
Address:46724 Lakeview Blvd, Fremont, CA 94538

EUROPE

E-mail:
support.uk@taotronics.com(UK)
support.de@taotronics.com(DE)
support.fr@taotronics.com(FR)
support.es@taotronics.com(ES)
support.it@taotronics.com(IT)
EU Importer: ZBT International Trading GmbH
Address: Halstenbeker Weg 98C, 25462 Rellingen, Deutschland

ASIA PACIFIC

E-mail: support.jp@taotronics.com(JP)
Tel: 03-5542-0907 (10-17時/土日祝除く)
JP Importer: 株式会社SUNVALLEY JAPAN
〒104-0032東京都中央区八丁堀3丁目18-6PMO京橋東9F

MANUFACTURER

Manufacturer: Shenzhen NearbyExpress Technology Development Co.,Ltd.
Address: 333 Bulong Road, Shenzhen, China, 518129

Learn more about the EU Declaration of Conformity:
<https://www.taotronics.com/downloads-TT-CL011-CE-Cert.html>

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