尺寸: 100mm*120mm 折页 128g 铜版纸



Model:TT-CL018
TaoTronics Light Therapy Lamp **User Guide**

Thank you for choosing the **TaoTronics Light Therapy Lamp**. Please read this User Guide thoroughly before use and retain it for future reference.

IMPORTANT SAFETY INSTRUCTIONS



rnina

- To ensure correct and safe usage and functioning of the product, always follow these instructions carefully and keep for future reference.
- Please keep the lamp away from water or any other liquid to avoid electrical shock.
- For indoor use only.
- This appliance is not intended for use by persons (including children) with skin damage, reduced sensory or mental capabilities, or lack of experience and knowledge, unless authorized by a professional physician.
- Keep away from heat and avoid direct sunlight.
- Do not perform any repair or maintain while using this product for treatment.

- Do not try to repair or dismantle the lamp by yourself.
- Do not drop or knock the lamp on hard surfaces.
- Portable and mobile RF communication equipment may interfere with medical electrical equipment.
- Do not place the device in a location where it is difficult to operate the disconnecting device (disconnecting device: Adapter).
- Use only the supplied power adapter. Using an adapter other than the one provided by the manufacturer may result in failure to work, or even a safety incident.
- Keep the equipment out of reach of children or infants, avoid strangulation caused by the adapter cable.
- THIS IS AN ELECTRIC LAMP NOT A TOY! TO AVOID RISK OF FIRE, BURNS, PERSONAL INJURY AND ELECTRIC SHOCK IT SHOULD NOT BE PLAYED WITH OR PLACED WHERE SMALL CHILDREN CAN REACH IT

/

CAUTION

- Some people may suffer from light sensitivity, discontinue use immediately if any discomfort occurs.
- Never cover the lamp or place anything on top of it during operation to avoid overheating or fire hazards.

Consult your health care provider before use of this product if you are under medical care or have suffered from an eve condition before.

- Follow local ordinances and recycling instructions regarding disposal or recycling of the device and device components
- Do not look directly look at the light therapy lamp. It only serves as background light.
- If you have any product problems or report an unexpected situation during use, please unplug the power immediately and call the after-sales service hotline.

READ AND SAVE THESE INSTRUCTIONS

 $\overline{}$

Product Features


Your new TaoTronics Light Therapy Lamp is a remarkable product with many features and benefits:

- Delivers LED light to your home, office, dorm room or apartment.
- Foldable design for larger illumination area and easier control of power off.
- Cues your body's own energy enhancers to help you regain focus and concentration.
- UV free LED- the built-in LED light panel delivers the full brightness with no pollution, no radiation and no stroboscopic effects.
- 10,000 lux at 10cm/4" from lamp to lift up your spirit.

Package Contents

- 1 x TaoTronics Light Therapy Lamp (Model: TT-CL018)
- 1 x Power Adapter
- 1 x User Guide
- 1 x Thank You Card

Product Diagram

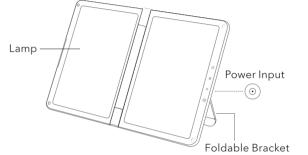


- () Power ON/OFF
- Brightness +
- · Brightness -
- 30 Minute Timer

Use Your Light Therapy Lamp

Get Started

1. Unfold the lamp.

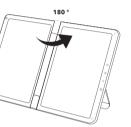

Place the lamp on a flat and stable surface using the support foot. Adjust the lamp head within proper range as per your need.



02 03 04



www.TaoTronics.com

E-mail:

support@taotronics.com


MANUFACTURER

Manufacturer: Shenzhen NearbyExpress Technology Development Co.,Ltd. Address: Floor 7, Building E, Galaxy World Phase II, Shenzhen, China



TT-CL018_V1.1_EN

 Connect the light therapy lamp to a power source with the included power adapter. The lamp is in Standby mode now with the Power indicator on.

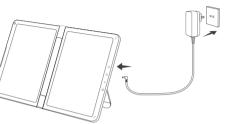



Power On / Off

- Tap () to turn on / off the lamp.
- Optionally, you can close the lamp to turn it off.



Light Brightness +

 When the lamp is on, tap and hold ☆ to increase the brightness (20%~100%).



Light Brightness -

• When the lamp is on, tap and hold 🔅 to decrease the brightness (100%~20%).



30-Minute Timer

- When the lamp is on, tap to turn on the 30-minute timer. The lamp will turn off after 30 minutes

Memory Function

 If power is not disconnected, the lamp will always remember the last set brightness once you turn it on again.

Troubleshooting

 $\overline{}$

Problem	Possible Cause	Solution
The lamp does not light up	No power supply	Correctly plug in the power adapter
	The lamp is not turned on	Refer to the using your Light Therapy Lamp section
The touch control is not responsive	Hands are wet	Tap the panel with dry hands
	Voltage fluctua- tion	Disconnect the power adapter to restart the lamp
LED failure	LED component failed	Please contact the TaoTronics Support Team
Bad smell	Power adapter failed	Please contact the TaoTronics Support Team
	Short circuit	Please contact the TaoTronics Support Team

Specifications

Model	TT-CL018
Power	14W Max
Input	DC 12V == 2A
Color Temperature	6500K±500K
Lumens	>650lm

Warranty

TaoTronics products are covered by a 12 month limited warranty from the date of its original purchase. If any problems occur, please contact TaoTronics Customer Care Centre via the contact information listed in this User Guide.

FCC Compliance

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1)This device may not cause harmful interference, and (2)This device must accept any interference received, including interference that may cause undesired operation.

08 09 10