




**MEDIC  
THERAPEUTICS**

# **HEATING PAD**

---

**User Manual**  
FOR PERSONAL USE ONLY

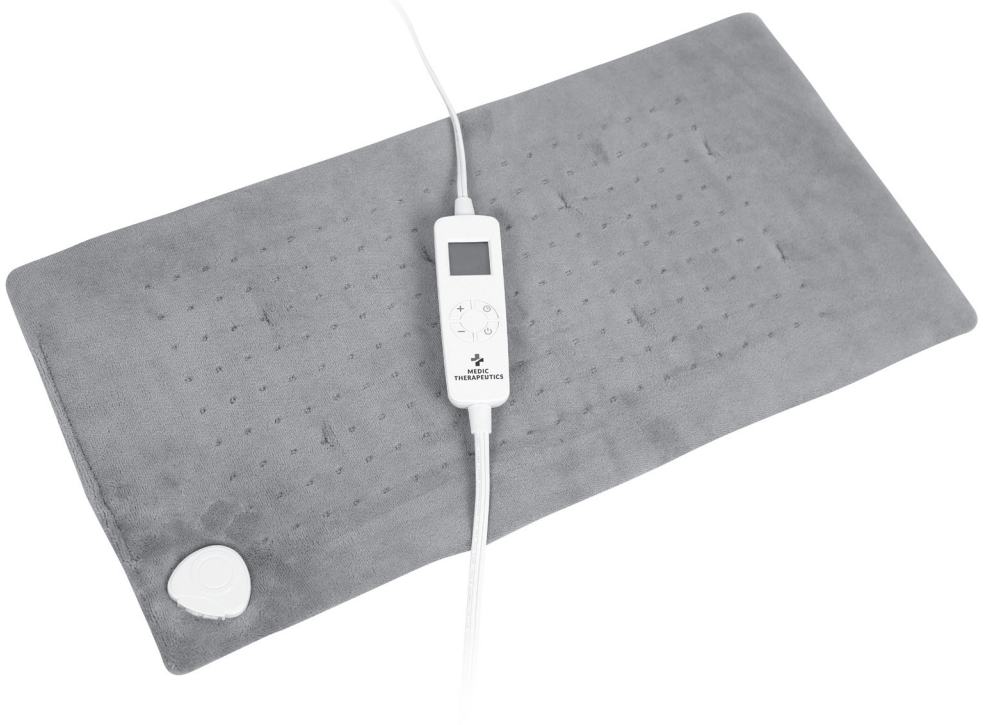

**Read the entire manual carefully  
before using this device.**

# IMPORTANT SAFETY INSTRUCTIONS

## DO NOT DESTROY

When using this or any other electrical appliance, always follow these basic safety precautions:



### DANGER

Burns will result from improper use. To reduce the risk of burns, electric shock, fire or accident of any kind, this product must be used in accordance with the following instructions.

### IMPORTANT SAFETY INSTRUCTIONS, DO NOT DESTROY

1. Do not use while sleeping.
2. Burns can occur at any setting, check the skin under the pad with frequency to avoid burns and blisters.
3. Do not use on an infant.
4. Do not use on a disabled, sleeping, or unconscious individual.
5. Do not use on an individual with poor circulation, paralysis, or diabetes.
6. Do not sit, lean or otherwise crush the pad and avoid sharp folds.
7. Always place pad on top of, not under the body.
8. Never place the pad between the body and a chair, sofa, bed or pillow.
9. Do not use on areas where the skin is sensitive.
10. Do not use in an oxygenated environment.
11. Do not use with metal pins or other metallic accessories.
12. Do not use with liniment or ointment.
13. Do not use in conjunction with any other heat-generating product or treatment.

14. Never pull on the supply cord and do not use the cord as a handle.
15. Always unplug when not in use.
16. Do not dry clean.
17. When storing, wrap the unplugged cord loosely around the pad.

**DO NOT ATTEMPT TO MODIFY THIS PLUG OR DEFEAT THIS SAFETY FEATURE IN ANY WAY.**



Please note, to reduce risk of shock, the heating pad has a polarized plug (one prong is wider than the other) and is intended to be used with a polarized outlet only. If the plug does not fit, reverse and try again. If it still does not fit, do not force, contact a qualified electrician.



**PARTS**

---

1. Heating Pad
2. Pad Connector
3. Control Connector
4. Control
5. Power Cord/Plug

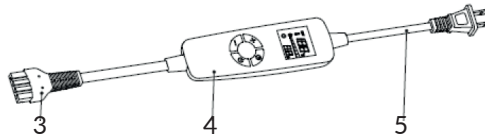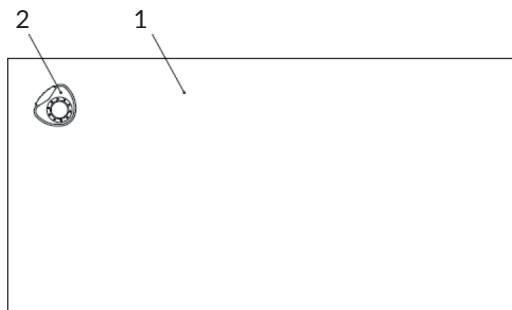


## HOW TO USE YOUR HEATING PAD

Before each use, always inspect the product to ensure there is no damage to the power cord, the Control, or the Pad. Fig. 1

The Pad Connector and Control Connector are designed to connect only one way. Align the inner prongs and then push together until Connectors are

1. Secure (see Fig. 1).
2. Plug the product into a 120-volt AC electrical outlet.

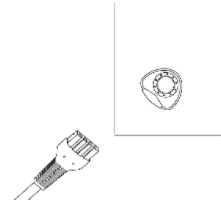


**NOTE:** If there is no power, make sure the AC electrical outlet is functional and the connection in step 1 is secure.  
Discontinue use if failure continues.

3. Use the Control to set the heat at the desired level (see Fig 2.).

**NOTE:** Heat levels are “105” (low), 115 through 145, and “155” (high)


Press the On/Off “⏻” Button on the Control. The display will illuminate, the all icon will display about 1.5 seconds, and then:

Fig. 2

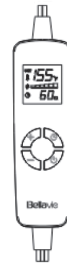

4. “105°F” will display, indicating the Control is set at the lowest heat setting, “60m” will display, indicating the Control is set at the 60 min auto off (see Fig. 3).


Fig. 3

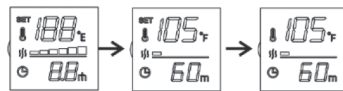


**WARNING: Burns can occur regardless of the control setting. Check skin under Pad frequently.**

5. To change the heat setting, press the “+” Button for a higher heat level or the “-” Button for a lower heat level. The desired heat level and the heat mode icon “SET” will appear for about 3 seconds, and then display the real-time temperature of the pad. (see Fig. 4).

Fig. 4

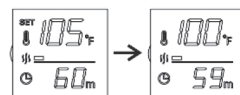

6. Display of Settings,  display as 6 setting,  display as 3 setting.

**NOTE:** If the Control shuts off, make sure the Pad and Control connectors are securely connected and re-start the product by pressing “⏻”.

7. Place the Pad on top of the body. If the Pad does not evenly heat, make sure there are no folds or bunched up areas.
8. Turn the product off by pressing the Button “⏻” and unplug it when it is not in use.


### **AUTO OFF FUNCTION:**

This product will turn off after a set duration of time. Settings are 0.5 hour (30 minutes), 1.0 hour, and 1.5 hours (90 minutes),.

**WARNING: Burns can occur regardless of the control setting; check skin under Pad frequently. Unplug the product when it is not in use.**

1. After turning the Pad on and selecting the desired heat setting (steps 4 and 5 above), press the “⏻” Button or the “⌚” Button until the desired amount of time is displayed. The display will illuminate the desired amount of time (for example “60m”), Time numerical display flickers for 3 seconds, and then the desired time setting will display steadily (see Fig. 5).

Fig. 5

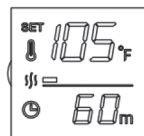

**NOTE:** When the Heating Pad is turned off, the Auto Off feature always defaults to “60m”(60min).


2. Repeat step 1 to change the Auto Off setting.

### **MOIST HEAT OPTION:**

---

1. Gently dab/pat one side of the Pad with a wet sponge/cloth or use a spray bottle to apply a fine mist.
2. Make sure the Control Connector or Pad Connector is not sprayed or dampened.

**WARNING: Do not use this product if Control or Pad Connectors are wet. Make sure the product connections are completely dry before using.**

3. Place the moistened side of the Pad on top of the body.
4. Turn the product off by pressing the On/Off Button “

## CARE AND MAINTENANCE

---

1. Disconnect the power cord from the electrical outlet.
2. Make sure the Pad Connector is disconnected from the Control Connector by pressing the tabs on the sides of the Control Connector while pulling the connectors apart.
3. Machine or hand wash the Pad using a mild detergent, cool water, and a gentle/delicate cycle.
4. Line dry or lay the Pad flat to dry. Do not tumble dry.
5. Make sure the product is completely dry before the next use.
6. To store the product, loosely roll the Heating Pad toward the Storage Straps, loop the Straps over the rolled ends, and then store it in a cool, dry, dark location after each use. Avoid sharply folding the Heating Pad.

## ELECTRICAL SPECIFICATIONS

---

Voltage	Frequency	Wattage	Amperage
120 VAC	60 Hz	85 W	0.71

## FREQUENTLY ASKED QUESTIONS

---

1. **When should a heating pad be used?**  
Whenever a hot form of temporary relief is desired for personal comfort.
2. **What temperature does this product reach?**  
Maximum temperature does not exceed the Underwriters Laboratories temperature limit of 176 °F (80 °C) when set to “6”. Temperature varies by heat setting; you may need to use a lower heat setting depending on the specific application and your personal preference.

### 3. How much time does it take to warm up to each temperature selection?

Temperature Selection	Indicator	Heat Up Elapsed Time
Low	100°F	3 min
2	110°F	4 min
3	120°F	6 min
4	130°F	8 min
5	140°F	10 min
High	150°F	15 min

### LIMITED LIFETIME WARRANTY

---

Your Medic Therapeutics Heating Pad is backed by a limited lifetime manufacturer's warranty. Medic Therapeutics will repair or replace your device at any time should it fail due to a defect in material or workmanship, subject to the certain limitations.

This limited warranty does not cover any damage that results from unauthorized or improper use, service, or repair. Further, it does not cover damage caused by accident, impact, negligence, or normal wear and tear. Should you discover your Medic Therapeutics Heating Pad is not functioning properly, please send your device to our repair center for evaluation. If your product cannot be repaired or serviced we will reserve the right to change it for a similar or newer model.

Please note that a flat rate of \$35.00 will be charged to cover a service evaluation fee and return shipping of your device. All warranty claims must be accompanied by a copy of your proof of purchase from an authorized retailer. Please send your device, proof of purchase, and a check or money order in the amount of \$35.00 made out to Medic Therapeutics to:

Medic Therapeutics Service Center  
3069 Taft Street  
Hollywood, FL 33009

Contact:  
warranty@medictherapeutics.com



**MEDIC  
THERAPEUTICS**