



HOW TO USE


Home-care Beauty Device

POINT LIFT

EP-30

Names and Functions of Components

Main unit

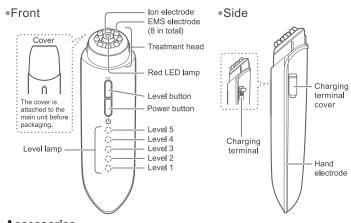

Accessories

- Cover
- (Attached to the main unit before packaging)
- AC adapter
- USB charging cable (micro USB)*
- User's Manual
- *Do not use with any other product.

The design, specifications of the device and its accessories are subject to change for a purpose of improvement without notice.

Functions

EMS

The 8 EMS electrodes stimulate the muscles electronically.

Double actions

lontophoresis and moisture pulses are simultaneously applied to the skin, supporting moisture penetration to deep inside the stratum corneum.

Auto Power OFF

Reminder buzzer

The device will emit a "beep!" sound approx. every 10 seconds.

Before Use


Precautions

Read the points below carefully before use and use the device only as directed.

- Make sure to perform a pretest before using the device for the first time.
- Remove any jewelry or accessories (necklaces, piercing, etc.) before starting treatment.
- Make sure to remove makeup before using the device.
- The following usage time and frequency are a guide.

Use the device for around 5 minutes per treatment, 2 to 3 times per week. When using again, be sure to leave an interval of 5 or more hours between treatments. **The above information is a guide only. Please observe your skin's condition and use the device above information.

Charge the unit

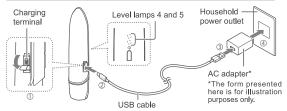

- Attach the cover to the main unit and stand the main unit in a position where the cover comes to the bottom.
- While the battery is charging, Level lamps 4 and 5 will blink alternately. When charging is complete, Level lamps 4 and 5 will light up for approx. 1 minute and then switch off.
- At an indoor temperature of approximately 25°C, it will take about 3 hours to reach full charge, which will allow the device to operate for about 1 hour Actual charging time and operating time will differ depending upon the remaining power and the charging environment.

 *If charging has not begun after about 10 seconds, check whether the USB
- charging cable or the AC adapter is not connected properly.



If there is any moisture on or around the main unit and charging terminal, be sure to wipe it dry. Charging while moist could cause damage to the battery. When charging the battery for the first time after purchase, use the device only after the battery is fully charged.

When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge displayed. In such a case, charge the device again before using it.


How to Use

Make sure to remove makeup and wash your face before using the device.




Turn on the device

If you press down on the Power button for approx. 3 seconds or longer the red LED lamp and Level lamp 1 will light up.

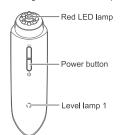


*For safety reasons, the device is set to Level 1 when turned on

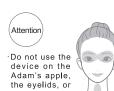



Each time you press the Level button, the device will emit a "peep!" sound and the EMS strength will switch (Levels 1 - 5). You can use the level lamp to check the level.

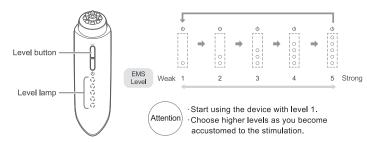



Gently apply the treatment head to your skin and move the head over the area being treated.

Apply gel or cream that you usually use on your skin in the areas over which the device is to be operated.For the neck and décolleté, use an amount of gel/cream approx. the size of a cherry (large), and for the face overall use an amount approx. the size of two cherries.

•The device will emit a "beep!" sound approx. every 10 seconds. Using the sound of the buzzer as a guide, move on to the next treatment point.





See "Moving the Unit" (reverse side) for how to move the unit.



Wipe off any gel or cream on your hands before using the main unit.

Press down firmly so that all of the EMS electrodes are in contact with your skin as you move the device over the skin.



Care for the skin when complete


After treatment has been completed, blend the gel/cream into your skin and finish treatment by applying the basic skin-care products that you usually use.

Moving the unit

Using the sound of the buzzer (approx. every 10 seconds) as a guide, move on to the next treatment point as shown in the diagram below. Press down firmly so that all of the EMS electrodes are in contact with your skin as you move the device over the skin.



Apply gel or cream to your neck and décolleté.

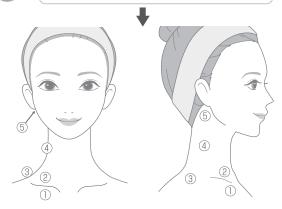


Move the device from the right side to the left.

Press the device down as if stamping on each point as you move the device over the skin.



Apply gel or cream to your entire face.

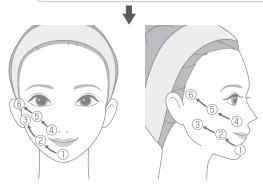

Move the device from the right side to the left.


Move the device to each point as if lifting up the skin.



- ·Apply an ample amount of gel or cream without blending it into your skin.
- ·Wipe off any gel or cream on your hands before using the main unit.





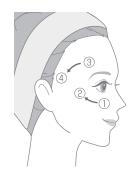

Move the device from the right side to the left.

- ·During use, if you feel that the device is not sliding smoothly over the skin, reapply the gel or cream.
- ·If the EMS stimulation feels too strong, lower the position of the device slightly towards the cheek.



Blend any remaining gel or cream into the skin on your face overall and your neck.



Smooth the gel/cream down your neck around 3 times.

After Treatment

How to clean

Clean the device after each use.



✓!\ WARNING

· When cleaning the device, make sure to turn it off and disconnect the USB cable from the main unit.

Failure to comply may cause an electric shock, short circuit, or fire.

 Do not wash in water. This device is not waterproof. May result in a failure of the device.





Wipe off any soiling using a well-wrung soft cloth or similar material.

- ·Wipe between the EMS electrodes using a cotton bud.
- ·For stubborn soiling, wipe gently with a cloth soaked in a diluted neutral detergent.



- . Do not leave the unit with gel or cream on it.
- Refrain from using an alkali cleaning agent, alcohol, thinner, benzene, nail polish remover, or other such substances.



Check out the instruction video

You can watch a video showing how to use the product on your smartphone, etc.



