

English User's Manual come with this product.
Please download the other language's User's Manual version from this link (<https://www.ya-man.co.jp/en/support/documents/>) or scan with QR code with your camera phone.
Please carefully read User's Manual to ensure safe and proper use of the device.



TO ACCESS SCAN WITH
YOUR CAMERA PHONE

User's Manual

Home Use Beauty Device

Bloom 5

YJFS16

Table of Contents

Introduction	1
Safety Precautions	1
Use in Foreign Countries	6
Names and Functions of Components.....	7
Before Use.....	9
How to Use	11

After Treatment	15
Additional Information	17
• Checklist for Possible Problems	17
• Frequently Asked Questions	19
• Product Specifications	20
After-sales Service	21



Introduction

Thank you for purchasing this YA-MAN product.

- Please read **Bloom 5 User's Manual** carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

Safety Precautions

- Read these Safety Precautions carefully before using the device in order to use it correctly.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- Do not use the device while it is in a defective state, as doing so could cause physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or ATTENTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

● Classifications of Precautions



DANGER

Improper use of this device poses a severe and likely risk of death or serious injury.



WARNING

Improper use of this device may pose a risk of death or serious injury.



ATTENTION

Improper use of this device may cause an injury or property damage*.

* 'Property damage' refers to damage affecting dwellings, household belongings, or domestic animals and pets.

● Symbols



This symbol shows you what you must not do.



This symbol shows you what you must do.



This symbol shows you what you must be careful with.



DANGER (Use and storage)

● **The following individuals should not use the device:**

Use of the device could cause electronic medical devices to malfunction.

- ◆ Do not use the device in conjunction with any medical electronic equipment, including:
pacemakers or other implanted electronic medical devices; artificial heart-lung machines or other electronic life-support medical devices; electrocardiographs or other worn electronic medical devices.



Prohibition



WARNING (Use and storage)

● **The following individuals should not use the device:**

Failure to comply may cause accidents or other problems.

- If you have existing or suspected epilepsy, heart disease, or related disorders
- Infants or those who are unable to express their intentions on their own
- If you are unable to sense temperature or stimulation
- If you have an autoimmune disease

● **Do not allow physically challenged persons to use the device unattended.**

Consult a physician prior to use under the following conditions:

- If you have an acute disease
- If you have an infectious disease
- If you have a malignant tumor
- If you have a febrile disease
- If you have high blood pressure or a blood disease
- If you have a fever
- If you have sensory impairment due to severe peripheral circulatory dysfunction (e.g. diabetes)
- If you have a nerve damage
- If you have a skin disease or atopic dermatitis
- If you are a hospital outpatient and are taking medication
- If you have a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or if you have a concern about your skin when using this device

● **Do not use in the following cases or on the following areas.**

Failure to comply may cause accidents or other problems.

- If your skin is damaged, experiencing a rash, inflammation, or other abnormalities
- If you are under the influence of alcohol
- Head
- Anywhere inside the mouth, around the pubic region, or any mucous membrane
- Eyeballs
- Upper eyelids
- Lower eyelids, outer corner of the eye (except power level eye treatment)
- Throat
- Bulging moles
- Viral warts
- Sites where you have broken a bone
- Areas of internal bleeding, swelling or scabbing after cosmetic medical procedures
- * Avoid using the device for 2 weeks after cosmetic medical procedures. After this period, consult with a physician and only use the device if the physician confirms there is no problem.
(e.g. Hyaluronic acid injections, silicone injections, Botox, thread lifts, double eyelid surgery, HIFU, etc.)
- Areas where metal was implanted by surgery (reconstructive surgery, plastic surgery, cosmetic surgery, etc.)

● **Do not use the device for a long time over areas where the skin is thin or areas above bones.**

Failure to comply may cause skin problems or physical ailments.



Prohibition



WARNING (Use and storage)

- **Do not overuse the device. For usage time and frequency, follow the instructions in the User's Manual.**

Failure to comply may cause physical injuries.

- **Before use, check to make sure the electrode and treatment head are not cracked, deformed or broken. Do not use the device if it is damaged or malfunctioning.**

Failure to comply may cause injuries or skin problems.

- **Do not use, store or charge the device in the following locations:**

- In bathrooms or near bathtubs, showers, washbasins or other containers filled with water
- A dusty location
- Outdoors or a location where there is exposure to direct sunlight
- A slope or a location where there are vibrations or shocks
- A location where there are strong radio waves or magnetism
- Near a flammable substance (such as alcohol, thinner, benzene, nail polish remover or spray).
- A location where there is exposure to oil smoke or steam
- A location that may become hot, such as near a heating device (e.g. a heater)
- A location with drastic temperature changes (e.g. near an air conditioner or heater)

Failure to comply may cause a fire, electric shocks, accidents, malfunction of the device, falls or injuries.

- **Do not use the device if the AC adapter is damaged.**

Inquiry the Global After Service

(<https://www.ya-man.co.jp/en/support/repair/>) or your local distributor and use the dedicated AC Adapter.

Do not connect anything other than the rated power supply (DC9V 2.0A).

- **Do not pull out the plug from a power outlet or pull out the plug from the device during use.**

Failure to comply may cause an electric shock or a fire due to a short circuit. Be sure to turn off the power before pulling out the plug.



Prohibition

- **Do not leave the device unattended while it is turned on.**

Failure to comply may cause a fire or other accidents.

- **Always use the dedicated AC adapter.**

Failure to comply may cause an electric shock, a short circuit or a fire.

- **Do not use the device for any purpose or in areas other than those stated in the User's Manual.**

Failure to comply may cause a malfunction of the device or injuries.

- **Children should obtain the consent of their guardians before using the device and should only use it under the supervision of a guardian or a specialist.**

Failure to comply may cause accidents or injuries.

- **Keep the device and AC adapter away from children and do not allow them to use it alone. Children shall not play with the appliance.**

Failure to comply may cause accidents or injuries.

- **Do not dispose of the device in a fire or place it near fire or flames. Do not use, charge or store in a high-temperature or high-pressure area.**

The device features a built-in rechargeable battery. Failure to comply may lead to the generation of heat, a fire or explosion.



Prohibition



Non-intended uses
are prohibited



User restrictions



Combustion
warning



WARNING (Use and storage)

- **If the device is dropped or damaged, turn off the power immediately, pull out the plug, and request a repair/inspection.**

Failure to comply may cause a fire or electric shock.



Must Comply

- **Be sure to use the device in an operating voltage range between 100~240 volts. Make sure the plug is fully inserted into the power outlet.**

Failure to comply may cause a fire, electric shock, malfunction of the device or other accidents.

* The free-repair warranty is only valid in Japan.

- **Do not disassemble, attempt to repair or modify this device in any way.**

Failure to comply may cause a fire, an electric shock or other injuries. Contact your local distributor for advice on repairing the device.



Disassembly
is prohibited

- **If you feel an earthquake, see lightning or hear thunder, turn off the power and unplug the device.**

Failure to comply may cause a fire or electric shock.



Removing
the plug

- **When disconnecting the plug from a power outlet, be sure to pull it out by grasping its head instead of the cord (do not touch the metal parts). Do not insert or disconnect the plug with wet or damp hands.**

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.

- **After using the device, when the device is not in use and when cleaning it, be sure to turn off the power and remove the plug from a power outlet.**

Make sure to disconnect the AC adapter from a power outlet after charging is completed. Failure to comply may cause a fire or other accidents.

- **Do not exceed the rated capacity of a power outlet or wiring accessory (do not connect too many plugs to one outlet).**

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device.



Overloading of
an electric circuit
is prohibited

- **Do not use if the AC adapter is damaged or the wall socket is too loose.**

Failure to comply may cause an electric shock, short circuit, fire, or malfunction of the device. Contact your local distributor for advice on repairing the device.



Use of a
damaged cable
or plug
is prohibited

- **Do not damage or modify the AC adapter. Do not apply unnecessary force, bundle or twist.**

Failure to comply may cause the cable to break, resulting in a fire or an electric shock. Straighten out any twists in the cord before storing.



Careless or
reckless use
is prohibited

ATTENTION (Use and storage)

- **Do not press the device strongly or use it intensively on a single area.**

Pressing the device stronger than necessary or using it intensively on a single area may damage the skin.

- **Do not use while driving.**

Failure to comply may cause accidents or injuries.

- **Do not allow the device to come into contact with another beauty device or electronic appliance while the power is on. Do not use with other electronic products (electric blankets/bed sheets, heated tables, etc.) at the same time.**

Failure to comply may lead to accidents due to malfunctions.

- **Do not look directly into the LED light. When using the LED light to treat the upper eyelids, lower eyelids and outer corners of the eye, be sure to close your eyes and keep the electrode head at least 5 cm away from the skin.**

Failure to comply may cause damage to your eyes or skin.

- **Use only Level 1 upon the first use and until you are accustomed to use.**

Even if you do not feel any stimulation, make only small adjustments to the intensity level. Use the device at an intensity level that suits you.

Make only gradual adjustments to the intensity level. Failure to comply may cause accidents or injuries.



Prohibition



Must Comply

- **Check carefully that your skin is not showing any abnormalities during the treatment. Stop using the device immediately if you start to feel pain, fatigue or if you see any skin reaction during use.**

Stop using the device immediately if the treatment does not agree with your skin. If symptoms persist the following day and beyond, consult a physician.

- **If you are currently undergoing dental treatment or have teeth fillings, choose a lower intensity level and apply the device to a different area.**

You may feel strong irritation or discomfort due to EMS on your cheeks or chin.

If you feel strong discomfort, consult a physician.

- **Avoid using the device on the same area for longer than 3 minutes.**

When using on the same area again, be sure to leave an interval of 24 or more hours between treatments.

Excessive use may damage the skin. The above recommendations are guidelines only; please observe your skin condition and use the device accordingly.

- **If a power failure occurs, turn off the power immediately and unplug the device from the power outlet.**

Failure to comply may cause an accident.



Must Comply



Removing the plug

Use in Foreign Countries



WARNING

- **Be sure to use the device in an operating voltage range between 100~240 volts.**

Failure to comply may cause a fire, electric shock or malfunction of the device.



Must Comply

- **Be sure to check the plug configuration and the voltage of the power outlet before use.**

Failure to comply may lead to the generation of heat, fire or malfunction of the device.




Combustion
warning

- The plug configuration and voltage may vary depending on the region and the type of accommodation. Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a consumer electronics retailer.


Names and Functions of Components

Main Unit

● Front view

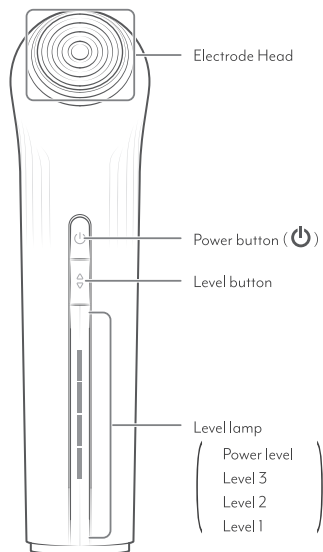

● Rear

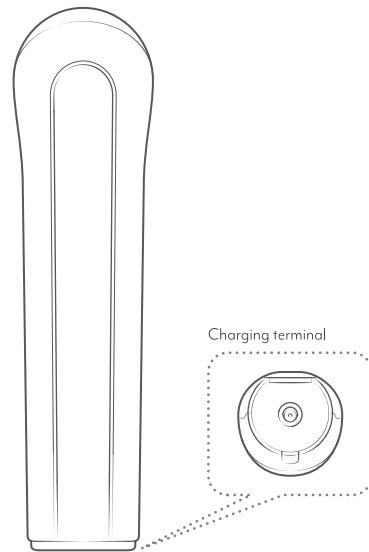

Accessories

● AC Adapter

● User's Manual

Functions

Safety functions

- Temperature sensor
Output will be reduced automatically for safety reasons when the temperature on the skin or the electrode head rises.
If the temperature becomes too high, the power will be turned off for safety. (When this happens, the electrode head flashes and beeps five times)
If this happens, let the device rest for about 30 minutes.
- Current sensor
Measures the current output value and controls the output.
- Motion sensor
The motion sensor detects movement of the device and pauses output for safety reasons when the electrode head has been used on the same area continuously.
When movement of the device slows, the motion sensor may reduce the level of warmth.

Timer / Standby

Estimated usage time is about 6 minutes for levels 1 to 3 and about 10 minutes for the power level. 3 minutes after the end of half-face treatment (for levels 1 to 3) or after a set time from the end of each treatment (for the power level), a buzzer will sound and RF and EMS output will be paused (standby function).

The Power button () and Level lamp will flash. When you are ready to perform treatment, press the Power button () to restart.

Auto Power-off

Immediately after the end of all modes, or if the device is not operated within approx. 5 minutes of going into standby mode, the device will be turned off automatically after 2 short beeps (Auto Power-off).

Charging Guide Display

The color of the power button (white, orange or red) indicates remaining battery levels.

If the rechargeable battery's remaining level is low when the power is turned on, the button will start flashing and a buzzer will sound.

The remaining battery level display will differ according to the condition of the battery and the operating environment.


Personal memory function

The device remembers the level being used and whether EMS was turned on/off when the last treatment ended, and returns to the same settings the next time it is switched on.

Before Use

Charge the device

This product can be used either cordless or by connecting it to an AC adapter.

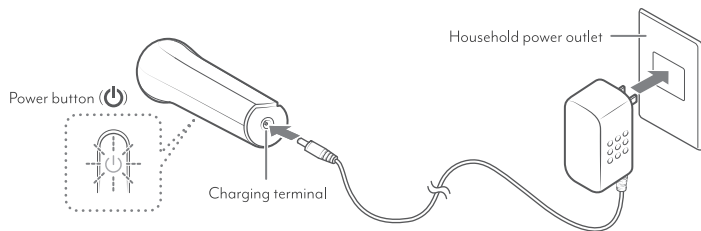

You can check the remaining battery level through the Power button (⏻).

Remaining battery level indicator	
White light	The battery is sufficiently charged.
Orange light	Low battery.
Red flashing	Charge the battery. Or connect the device to an AC adapter.

- While the battery is charging, the Power button (⏻) will blink on and off.
- When the battery has fully charged, the Power button (⏻) will be on for approx. 5 seconds and then will switch off. Disconnect the AC adapter.
- At an indoor temperature of approx. 25°C, it will take about three hours to reach a full charge, at which point the device can operate for about 30 minutes. Actual charging time and operating time will differ depending upon the remaining battery level and the charging environment.

* If charging does not begin after approx. 10 seconds, check if the AC adapter connection is loose.

* Design and specifications of the main unit and accessories are subject to change for improvements without prior notice.

ATTENTION

- This device is not waterproof. If water droplets should fall on the main unit or the charging terminal, wipe them off immediately and thoroughly. Storing or charging the device while it is moist could cause a malfunction of the device.
- When you have not used the device for a long period of time, check to see that it is functioning normally before re-starting use. Failure to comply may lead to an accident or malfunction of the device.
- Fully charge the battery before using the device for the first time after purchasing it.
- When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge displayed. In such a case, recharge the device before using it.
- Use only the included AC adapter.

Pre-testing on Skin

When using the device for the first time, always test for compatibility.



ATTENTION

- Be sure to perform a pre-test in order to prevent skin problems from occurring.
- If a rash, reddening, itchiness, stinging sensation or other abnormality occurs at the pre-testing site, immediately stop the pre-test.
- Check carefully that your skin is not showing any abnormalities during the treatment. Stop using the device immediately if the treatment does not agree with your skin.
- Remove all metal accessories (watches, necklaces, rings, earrings, etc.) when using the device. Failure to comply may cause accidents or injuries.
- The recommended ambient temperature is between 10 °C and 40 °C.

1

Wash the inner side of your arm with soap and water, then wipe dry with a towel.




2

Apply your regular skin lotion or gel (about one pearl-sized drop) to the center of the electrode head, turn the power on, then use the device for about 3 minutes while moving it across the skin (page 13).

Be sure to test at Level 1.

How to hold the device and use it on your skin

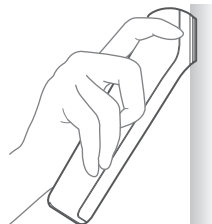

3

Rub the skin-care products remaining on your skin into the skin and check for any abnormalities.

Continue to use the product only if you experience no reddening or burning sensation on your skin within 24 hours.

Hold the main unit as shown below, moving it slowly so that the entire electrode head remains in contact with your skin.

- Be sure to move the device with the entire electrode head in contact with your skin.

How to Use

Basic method of use

ATTENTION


Make sure to remove makeup before using the device on your face.

If you use the device over makeup, the foundation may cause dark staining of the skin.

1

Turn on the device

Press and hold the Power button (⏻), the Power button, Level button, Level lamp and electrode head light up, and RF, EMS and red LED are turned on.

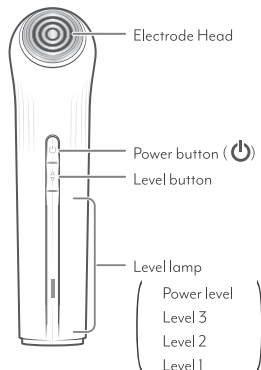

2

Select the level

Each time the Level button is pressed, levels can be adjusted. (Level 1 → Level 2 → Level 3 → Power level → Level 1 → Level 2...).

You can check the treatment level from the Level lamp.



EMS on/off switch

Press the Level button and hold down to switch EMS on/off
(Level button lit: EMS on; unlit: EMS off).

* When using the power level, EMS cannot be turned off.


ATTENTION

Use only Level 1 upon the first use and until you are accustomed to use.

3

Apply your regular skin lotion or gel to the center of the electrode head, then apply the electrode head to the skin and move it around.

The whole face is divided into 6 areas (level 1-3) or 8 areas (power level) (pages 13-14).
Apply one pearl-sized drop of your regular skin lotion or gel for each area.

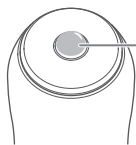

Your regular skin lotion or gel
(about one pearl-sized drop)



ATTENTION

- Do not apply your regular skin lotion or gel with your hand. Failure to comply may lead the skin-care product to stick on the device or cause it to slip and fall.
 - Please reapply skin-care products or water during treatment to ensure there is always some liquid between the electrode head and the skin. If there is an insufficient amount of skin-care products or water, the stimulation will weaken, causing skin problems.
 - If the device does not move smoothly, apply more of your regular skin lotion or gel.
 - Use the device in a position in which the entire surface of the electrode head comes in contact with your skin.
 - Be sure to move the device over the skin as you use it.
 - If it feels too hot, move it a little more quickly or lower the level.
 - Please carefully clean the parts of the device that come into contact with the skin every time you use it and keep the device clean at all times.
- Failure to comply may cause skin problems.

4

Care for the skin after the treatment

After treatment is complete, rub the remaining skin-care product into the skin and use your usual basic skin-care products to condition the skin.

After each use, refer to 'After Treatment' and remove the remaining sebum, serum and other particles stuck to the device. Failure to comply may cause a malfunction of the device or skin problems.

ATTENTION

See 'How to move the device' (pages 13-14) for information on how to move the device correctly.

- If you feel an unpleasant tingling on your skin, reduce the intensity level or move to another area.
You could also moisten the skin with water, lotion or gel before using the device. Failure to comply may cause skin damage.
- Throat, eyeball, upper eyelids, lower eyelids/outer corner of the eye (except power level eye treatment).
- Do not leave skin-care products on the electrode head while the power is turned on.




Areas you should
not treat

How to move the device in levels 1 to 3

Move the device in the directions of the arrows in order from ① to ③ below. Repeat the movement three times for each arrow. Perform the treatment for the right half face as indicated below.

①

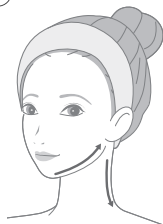


Apply your regular skin lotion or gel product



Move towards the neck as if lifting the facelines upwards

②

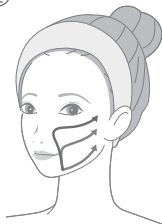

Apply your regular skin lotion or gel product



Move smoothly in the direction of the arrow

③



Apply your regular skin lotion or gel product



In a zigzagging direction from the area between your eyebrow toward your temple

Left half face

Approx.
3 minutes

Standby

Right half face

Approx.
3 minutes


- Recommended speed for moving the device: 4 cm/second
- Recommended duration of use: about 6 minutes per use (over the entire face)
- **The device will make one long beep after approx. 3 minutes and the treatment will pause. Please use this as a guide. To restart, press the Power button (⏻).**
- If you want to stop the treatment during use, press the Power button (⏻) for at least one second to turn off the power.

How to move the device in power level

Move the device in the directions of the arrows in order from ① to ④ below. Perform the treatment for the right half* In power level, the output varies depending on the area treated.

*In power level, the output varies depending on the area treated.

① Face lines treatment

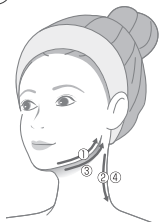


Apply your regular skin lotion or gel

↓
When the standby function is activated
Move towards the neck as if lifting the facelines upwards
After about 60 seconds



Level lamp

② Laugh lines treatment

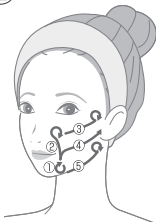

Apply your regular skin lotion or gel

↓
When the standby function is activated
Move sideways as if lifting up the skin from the mouth.
After about 110 seconds



③ Eye region treatment



Apply your regular skin lotion or gel

↓
When the standby function is activated
After moving as if pulling up the skin from the lower eyelid toward the outer corner of the eye, pull up from the area between your eyebrow toward the temple.
After about 70 seconds



④ Forehead treatment




Apply your regular skin lotion or gel

↓
When the standby function is activated
From above the eyebrows toward the hairline
After about 60 seconds



When the entire face has been treated, the device will shut off entirely.

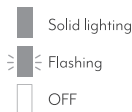

Left half face

Approx. 5 minutes

Right half face

Approx. 5 minutes

- Recommended speed for moving the device: 4 cm/second
- Recommended duration of use: about 10 minutes per use (over the entire face)
- **After a set period of time, the standby function will be activated and paused. To resume, press the power button.**
- If you want to stop the treatment during use, press the Power button () for at least one second to turn off the power.

After Treatment

How To Clean

Please carefully clean the parts of the device that come into contact with the skin every time you use it and keep the device clean at all times.

Failure to comply may cause skin problems.

1 When cleaning the device, be sure to turn it off and disconnect the AC adapter from the main unit.

Failure to comply may cause an electric shock, a short circuit or a fire.

2 Wipe off any dirt using a well-wrung soft cloth or similar material.


- For stubborn stains, wipe gently with a cloth soaked in a diluted neutral detergent.
- Remove any skin-care products remaining between the electrode rings of the electrode head with a cotton bud or cotton pad.



WARNING

- Do not wash the device in water.

This device is not waterproof. Failure to comply may cause a malfunction of the device.

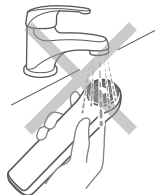

ATTENTION

- Do not leave the device with skin-care products still on it.
- Do not wipe the device with solvents such as alcohol, thinner, benzene, nail polish remover, kerosene, chemicals or polishing powder.
Failure to comply may cause damage to the surface of the device.



Storage

- Store in a clean, well-ventilated place out of the reach of children.
- Remove the AC adapter from the main unit and the power outlet when not in use for a long period of time.
- Do not leave the device in hot or humid places, such as near flames, in bright sunlight or in a bathroom.
- Cleaning and user maintenance shall not be made by children without supervision.

Maintenance and Inspections

- When you have not used the device for a long period of time, check to see that it is functioning normally before re-starting use.
Failure to comply may lead to an accident or malfunction of the device.
- If the device is malfunctioning, go to 'Checklist for Possible Problems' (page 17) and attempt to identify the problem before arranging to have the device repaired.
If the problem persists, inquiry the Global After Service (<https://www.ya-man.co.jp/en/support/repair/>) or your local distributor.
- If the device malfunctions, appears damaged or needs service, please contact the Global After Service (<https://www.ya-man.co.jp/en/support/repair/>) or your local distributor.

Additional Information

Checklist for Possible Problems

If you perceive any abnormality while using the device, verify the following points before arranging to have it repaired.

Problems	Possible Causes and Solutions
The device does not turn on.	Press and hold the Power button (⏻) for about one second.
	The device has run out of power. Charge the battery (page 9).
The device shuts off.	The device is in standby or Auto-off status (page 8). Press the Power button (⏻) to turn the power on (Standby function). If the device is not operated within five minutes of going into Standby Function, it will switch off automatically (Auto Power-off). If the device is left at low temperature for a long time, it may turn itself off. After leaving it at the recommended ambient temperature of between 10 °C and 40 °C for a while, check the charging status and press again the power button (⏻) and keep it pressed for about 1 second.
	The device will automatically switch off approx. 6 or 10 minutes after the treatment started.
The device only operates for a few minutes even after charging.	The battery level is low. Charge the battery (page 9).
	The battery has expired. Consult the Global After Service (https://www.ya-man.co.jp/en/support/repair/) or your local distributor.
	At low temperature, the power may switch off immediately. Use the device while charging it.
The device does not start charging.	Insert the included AC adapter into the main unit (page 9). Check to make sure the cable is not loose.
	Connect the included AC adapter to a wall outlet (page 9).
The radio frequency (RF) feels too hot	The device has been pressed against the same area of skin for too long. Be sure to keep the roller moving over the skin during operation.
	Move the device a little more quickly.
	Lower the level. Be sure to use Level 1 on the first use.

Problems	Possible Causes and Solutions
The radio frequency (RF) does not feel warm.	Your skin is dry. Moisturize your skin by applying your regular skin lotion or gel.
	Raise the level.
	Once started, after the end of half-face treatment (for levels 1 to 3) or after a set time from the end of each treatment (for the power level), RF and EMS functions will be paused (timer/standby function). Press the Power button (⏻) to restart.
EMS causes uncomfortable tingling or pain.	Lower the level. Be sure to use Level 1 on the first use.
	Use in a position in which the entire surface of the electrode head comes in contact with your skin.
The EMS feels weak or I cannot feel it at all.	Raise the level.
	Your skin is dry. Moisturize your skin by applying your regular skin lotion or gel.
	Depending on muscle mass, there may be some skin areas where it is difficult to feel the EMS.
	EMS output is turned off. Press and hold the Power/Level button (page 11). (In power level, EMS is always ON)
	3 minutes after the end of half-face treatment (for levels 1 to 3) or after a set time from the end of each treatment (for the power level), RF and EMS functions will be paused (timer/standby function). Press the Power button (⏻) to restart.
The Power button light turns red (or orange).	It is displaying the remaining battery level (page 8).
The level does not begin at Level 1.	This is the Personal Memory Function (page 11) operating.

Frequently Asked Questions

Q What is the recommended duration and frequency of use?

A

As a general rule, use it once a day (for 6 minutes, or 10 minutes on power level), 5 days a week.

* Please observe your skin condition and use the device accordingly.

* Check carefully that your skin is not showing any abnormalities during the treatment. Stop using the device immediately if the treatment does not agree with your skin.

Q Can I use the device continuously?

A

We recommend using the device once a day (for 6 minutes, or 10 minutes on power level).

Q How much of the skin-care product should I use each time?

A

Apply about six to eight pearl-sized drops of your regular skin lotion or gel for the entire face.

If you feel that the amount of skin-care product is insufficient, add more.

Q Can I use the device in the bathroom?

A

This device is not waterproof. Do not use the device in the bathroom.

Q My skin feels hot after use.

A

This is due to the heat of the radio frequency (RF).

Q The electrode head does not warm up.

A

The radio frequency (RF) does not warm the electrode head. Rather, it uses the power of electricity to warm your skin.

Q What happens when I increase levels?

A

As the levels increase, your skin is warmed more quickly and strongly. You will also feel the EMS more strongly.

* Use only Level 1 when starting until you are accustomed to use.

Q Can I wash the device with water?

A

This device is not waterproof. Do not wash it with water. Failure to comply may cause a malfunction of the device.

Product Specifications

Rated power supply	DC9V 2.0A
Power Consumption	Approx. 18 W (when charging)
Battery	Lithium-ion battery
Charging time	Approx. 3 hours (May vary depending on the charging environment.)
Operating Time	Approx. 30 minutes (when fully charged) (May vary with conditions of use and the EMS level.)
Outer Dimensions	Approx. W47 x D47 x H181 mm
Weight	Approx. 180 g
Materials	Main unit: ABS resin and polycarbonate Electrode head: PET resin and stainless steel
Country of origin	Japan
Accessories	AC adapter User's Manual
Recommended environmental temperature	10° C~40° C

After-sales Service

Rechargeable Battery

- This appliance contains batteries that are non-replaceable.
- The rechargeable battery is a consumable item. (not covered by the product warranty)
- If the operation time after charging has become extremely short, it could be the end of the battery's useful life. The battery can be charged approx. 300 times.
- Consult the Global After Service (<https://www.ya-man.co.jp/en/support/repair/>) or your local distributor for battery replacement.

Disposing of the device

- Follow your local government rules for proper disposal of this product.
- This product uses a lithium-ion storage battery.



Correct Disposal of this product

- This marking indicates that this product should not be disposed with other household wastes throughout the UK and EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.



How to remove the battery for disposal

The lithium-ion battery used in this product is a valuable recyclable resource.


- The battery must be removed from the appliance before it is disposed.
- Be sure to turn off the power and disconnect the adapter before removing the batteries.
- Use a screwdriver to remove the battery. Be careful not to injure yourself.
- When removing the battery, be careful not to damage it on a hard object.* Failure to comply may cause the generation of heat or a fire.
- Bring the removed battery to the appropriate collection point for disposal.
- Follow your local government rules for proper disposal or recycling of removed battery.

Operation Steps


Refer to figures 1 to 4 below to remove the battery.

- 1 Remove the outer cover using a flat-tipped object like a flat-blade screwdriver.


* It is fixed with strongly adhesive tape

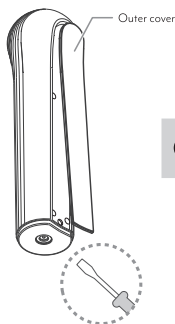

- 2 Remove the four screws on the lower case, separate the upper case from the lower case and disconnect the cable from the connector.

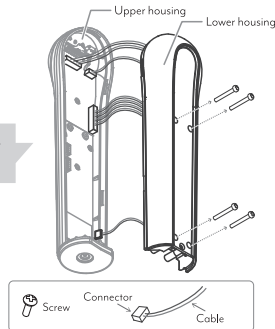


- 3 Unscrew the 3 screws on the lower housing and remove the battery compartment.

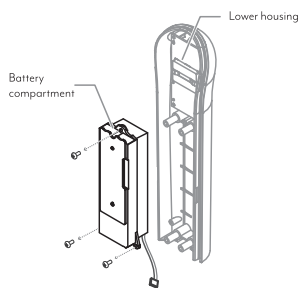

- 4 Remove the battery using a flat object like a ruler.

* It is fixed with strongly adhesive tape

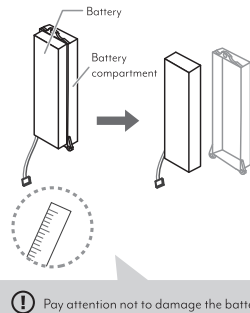

This figure is intended to be used for disposing of the product, not repairing it.
CAUTION: Disassembling the device by yourself may cause a failure.

Warranty

If you have purchased from www.global.ya-man.com,

Warranty: Please register your device at <https://global.ya-man.com/apps/warranty/registrants/new>

Contact: inquiries to contactus@ya-man.com

If you have purchased from local distributor,

Warranty: Please contact your local distributor

User's manual (WEB)

User's manual is also available at following address:

<https://www.ya-man.co.jp/en/support/documents/>