glo[®]

Face

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Thank you for trusting Glo, the Glo Face beauty and wellness device will help you take care of your skin at home to make it look luminous and firm, thanks to the combination of three technologies: facial massage, microcurrent toning and light therapy.

Be sure to READ ALL INSTRUCTIONS BEFORE USE and use this device only for its intended purpose.

INTENDED USE Glo Face is designed for face, neck and décolleté stimulation and for cosmetic and over-the-counter use.

 **WARNING:** MODIFICATIONS TO THE DEVICE ARE NOT ALLOWED, IT WILL VOID THE WARRANTY.

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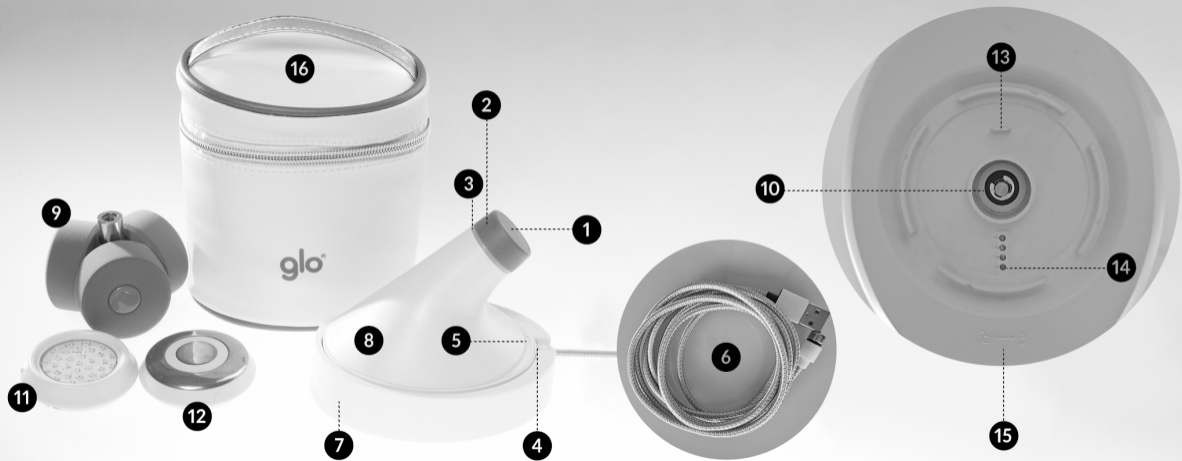
- ❶ PUSH BUTTON: It has two functions. 1 On/Off: Turns the device on and off. 2 LED therapy selection: press to change the LED light (red, green or blue).
- ❷ INTENSITY SELECTOR: Rotates to increase or decrease the intensity of microcurrents and facial massage.
- ❸ INDICATOR LIGHT: Indicates if the device is on and the selected intensity.
- ❹ CHARGING PORT: Connect the USB cable to charge the battery (1.5 hours face with 5V/2A charger). After charging, remove the USB cable.
- ❺ BATTERY INDICATOR LIGHTS: Blink during charging. Charging is complete when all 4 battery LED indicators are fully lit and have stopped flashing. The time of use will depend on each attachment, as not all of them consume the same energy. **Important: Fully charge for the first use. Only works without cord.**
- ❻ USB CHARGING CABLE
- ❼ STAND
- ❽ ACOUSTIC WARNER: helps the application routine. 1min for facial massage. 30 sec for phototherapy.
- ❾ MASSAGE: Massage head
- ❿ MAGNETIC AXIS: For easy positioning of the MASSAGE head.
- ⓫ GLOW: PBM (Photobiomodulation) head.
- ⓬ LIFT: Microcurrent head
- ⓭ PLACEMENT GUIDE: Helps to place the LIFT microcurrent head and the GLOW light therapy head in their correct position.
- ⓮ ELECTRICAL CONNECTOR: The LIFT and GLOW heads must fit with their electrical connection on the base to function.
- ⓯ LOCK SYMBOLS: help to place the head and fix it.
- ⓰ BEAUTY BAG: elegant and convenient for storage and travel.

BEFORE FIRST USE


- You can take a photo before and after the treatment with the same light and position to check your evolution.
- Remove hair from the face. This is especially important when using the MASSAGE head.
- The battery must be fully charged. Check the battery status and recharge your Glo Face if necessary.

PREPARATION

- Clean and dry your face and neck, removing all residue.
- Attach the treatment head of your choice

LIFT: Microcurrent Technology

Transmits microcurrents on the skin to tone and firm the face. Must be used with conductive gel. Apply conductive gel on your skin, and spread it evenly over the area you want to treat. It will beep every 10 seconds. Massage the skin until it is completely absorbed. **IMPORTANT NOTICE: DEVICE AND HEAD MUST BE COMPLETELY DRY BEFORE STARTING THE TREATMENT.** Avoid thyroid area.

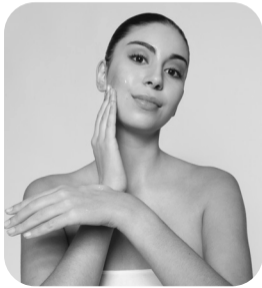- I) Press the power button for 2 seconds to activate the device. The lights below the push button indicate the intensity level of the treatment, which can be adjusted by turning the selector to adapt it to your sensitivity.
- II) Lightly press the two metal rings of the microcurrent head on the areas to be treated, slide it following the indications in the image. Check that the two metal rings are touching your skin. Slowly slide the device with upward movements until you hear a beep (10 sec), and press lightly for optimal results.
- III) When you finish your treatment, press the power button for 2 seconds to turn off the device. If you wish, you can apply your favorite facial care products afterwards.
- Remove any conductive gel residual.

GLOW: MULTI-LEDS TECHNOLOGY

Formed with 24 Multi-LEDs and 7 near infrared (NIR) - Allows you to choose between red, blue or green therapy.

- BLUE (415 nm- epidermis): reduces skin imperfections.
- GREEN (568nm-dermis): brightens dull complexion and evens skin tone.
- RED (633nm-hypodermis): reduces the signs of aging for a youthful, luminous-looking skin by stimulating collagen and elastin production.
- NIR (near infrared-830 nm): works by reducing inflammation, soothing irritation and strengthening the skin. Infrared light is not visible to the human eye, which is why infrared LEDs appear to be off, when in fact they are always on.

Important note: Phototherapy is not harmful to the eyes, but it can be uncomfortable, so we recommend closing the eyes each time you move the device from one area of the face to another. The light should not be applied directly to the eyes, even when eyes are closed. During the treatment a subtle heat emitted by LEDs might be felt.

- It can be used on clean, dry skin or apply a cosmetic beforehand.
- I) Press the power button for 2 seconds to activate the device and the light therapy head. The lights underneath the button indicate that it is on and the LEDs on the head will light up (the infrared will appear to be off because they are not perceivable by the human eye). You can change the treatment LED light by pressing the power button once, from red to green or blue.
- II) Place the head on your face, and hold it until you hear a “beep” (30 sec). Move to another area.
- III) When you finish your treatment, press the power button for 2 seconds to turn off the device.
- If desired, you can apply your favorite facial care products.




Note: This head can be applied to other areas of the body, such as the décolleté, hands, back or scalp.

MASSAGE: Relaxing Facial Massage

Soft sponges made of anti-bacterial material for a relaxing facial massage. **Important:** Make sure your hair is well gathered and removed from the face to prevent it from getting caught in the MASSAGE head. We recommend making a ponytail or using a headband.

- I) Press the power button for 2 seconds to activate the device. The lights below the push button indicate the intensity level of the massage. You can change the intensity by turning the selector from min to max. The head will start to rotate.
- II) Place the head on your face, and hold it firmly. Move it over the skin loosely until you hear a “beep” (1 min). Move to another area.
- III) When you finish your treatment, press the central button for 2 seconds to turn off the device.
- If you wish, you can apply your favorite facial care products.

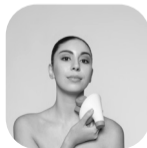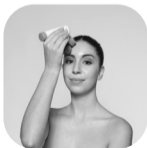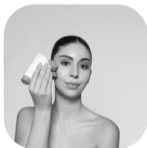END OF TREATMENT

- Turn off the device, clean it and store it in its base or in its case.
- You can charge it and leave it ready for the next use.

HOW TO CLEAN Glo FACE

- Always clean Glo Face after each use.
- Glo Face must be turned off, not connected to the USB cable. Remove the head in use in order to clean the base thoroughly.
- **IMPORTANT:** Never use cleaning products containing alcohol, acetone, petroleum.
- Clean the metal rings and the gap between them with a cloth dampened with lukewarm water. Dry it with a towel or non-abrasive cloth.
- Clean the phototherapy head with a cloth dampened with warm water. Dry it with a towel or non-abrasive cloth.
- **Massage accessory cleaning:**
 - The sponges are made of an anti-bacterial material to guarantee maximum hygiene, but after several uses it is advisable to wash them.
 - The sponges can be removed for cleaning, and should not be washed without being removed first **because the metal axis could be damaged.**
 - **WARNING!** Do not get the metal parts of the connecting axis wet.
 - I) Remove the sponges carefully to avoid damaging them.
 - II) Wash the sponges with neutral soap and lukewarm water under running tap water. Squeeze them gently several times until the water runs out clean. Allow the sponges to air dry before placing them back on the shaft.


MAINTENANCE

- The massage accessory is recommended to be replaced every 3 months. This time is indicative and depends on the use given, if it is deformed or damaged.
- When you are not going to use the device for a while, you can store it in its case. When you start using it again, charge it fully before the first use.
- All accessories are available at www.glo910.com

SPECIFICATIONS

MATERIALS: ABS, PC. **USAGE:** Up to 90 uses per charge*. **CHARGE TIME:** 1.5 h 5V/2A
MAX NOISE LEVEL: <67 dB **WEIGHT:** 170 G **BATTERY:** Li-Ion 1800mAh 3.7V
*Based on 1-minute treatments.

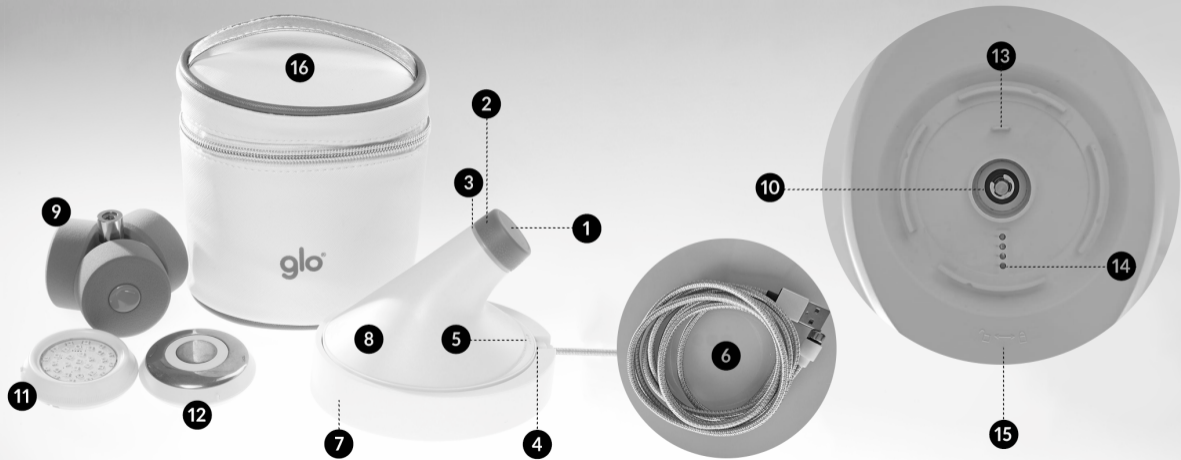Non-contractual photos: the product may differ from the photos in this manual.
These instructions are also available on our website, www.glo910.com, along with videos and more details on how to get the most out of your device. Follow us on Instagram @glo910spain to stay up to date.

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