JOLT BOLT

TECHNOLOGIE BREVETÉE SOLFT FLOW™

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CONTENTS

ADVANTAGES OF THE MESSAGE GUN JOLT BOLT page 33

INSTRUCTIONS AND OWNER'S MANUAL page 35

Instructions page 35

Charging page 36

Massage heads..... page 37

PRECAUTIONS FOR USE page 39

Warnings..... page 39

health warnings page 39

Environmental protection page 40

TECHNICAL SPECIFICATIONS page 41

FAQ & WARRANTY page 41

FAQ page 42

Warranty page 43

JOLT

**Thank you for your order,
we truly appreciate your trust in us.**

In order to use our recovery boots properly, please read the instructions, directions for use and safety information before use. More information is available on our website.

www.getjolt.com

JOLT

ADVANTAGES OF THE MESSAGE GUN JOLT BOLT


The JOLT™ BOLT Muscle Massage Gun is a portable, rechargeable vibration massage device. JOLT™ helps to increase blood flow and effectively alleviates the discomfort caused by the accumulation of lactic acid after sports and thus reduce discomfort for people suffering from muscle stiffness / aches.

JOLT™ can be used in the following cases:

- For warming up before a physical effort
- For the recovery after a sporting effort
- To relieve aches and joint stiffness
- To stimulate blood circulation
- To relieve stress and tension
- Massage of knots and trigger points
- Suitable for beginners, amateurs and professionals
- To relax and to loosen the tissues
- To reduce lactic acid build up in the muscles
- To increase joint mobility
- To fluidize edema in the swollen area
- To relieve aches and pains
- To relieve tendinitis

PRESENTATION OF DIFFERENT ATTACHMENTS

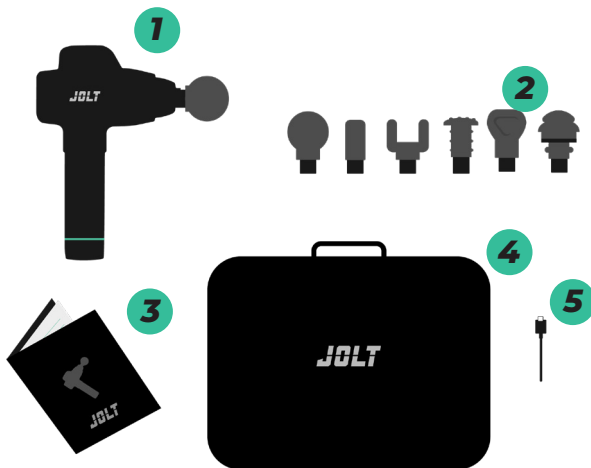
1. 1 JOLT BOLT massage gun
2. 4 massage heads
3. 1 Micro USB charger cable
4. 1 manual
5. 1 carry case

JOLT

INSTRUCTIONS AND OWNER'S MANUAL

INSTRUCTIONS

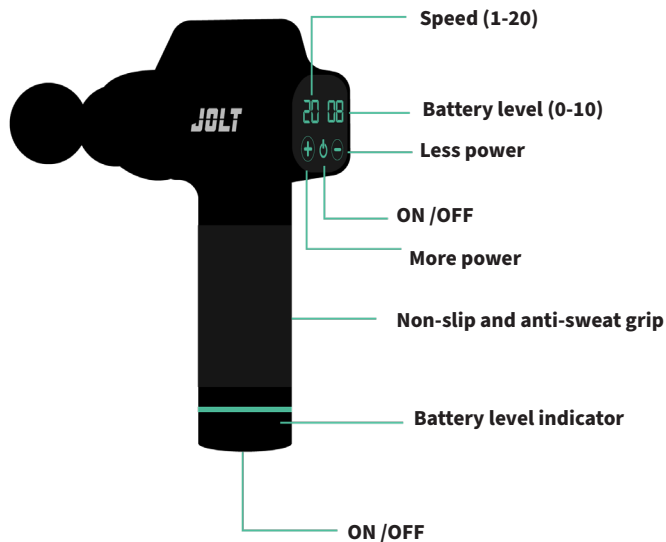1. Turn on the device using the ON/OFF button on the bottom of the battery (under the handle),
2. Press the power button on the digital display to start the gun.
3. Use the « + » and « - » buttons to adjust the speed from 1 to 20.
4. The LED on the handle will turn red when the battery level is less than 20%. If the battery level is less than 10%, the display will show «E1», the unit will not operate and you should recharge it.
5. The display goes from 01 to 10 successively while charging. The battery is fully charged when it reaches 10, and the indicator light at the bottom of the handle will turn green.
6. The unit features a stopwatch function and will automatically shut off after 15 minutes of continuous operation (for safety reasons to avoid overheating)
7. When the massage head is pressed too hard, the unit will start its automatic protection function and stop. «E4» will display and it can be used normally after powering on.

CHARGING

The charging port is on the bottom of the massage gun (Micro USB cable) Fully charge the battery for the first use and turn it off while it is charging.

1. Connect the power cable to the gun (Micro USB). The indicator will flash to show that the device is charging. Do not move the massage gun while it is charging.
2. When the indicator turns **green**, it means the battery is full.

MASSAGE HEADS

With the massage gun turned off, slowly insert the massage head into the barrel until it stops. As shown in the illustration, the position of the red line on the massage head should be aligned with the position of the red line on the gun.

To remove the massage head, you just need to pull on the head to remove it

JOLT™ BOLT can be used with 4 massage heads depending on the parts of your body you want to massage.

Round head

Its spherical design with soft coating allows to work large muscle groups such as the back of the thigh and the calf, the triceps, ... It is mainly suitable for new users who are sensitive to vibrations.

Flat head

It allows you to relax any muscle of the body for a classic massage.

Cylindrical head

The conical shape of this head allows an impacting and toning massage with a significant pressure depth. It is suitable for the meridians, soles of the feet, palms of the hands, Achilles heel, rhomboid muscles of the back, ...

U-shaped head

Its shape fits perfectly the spine and cervical/neck. It can also fit the shoulder and other bilateral muscles (do not touch directly the spine and cervical to avoid injury).

Arrow head

This head is a variation of the flat head, a soft arched coating has been added to its top in order to adapt to the roughness of your body. This head is suitable for soft tissues (tissues that support body structures and organs as muscles), internal and external muscles of the abdomen, chest and legs.

Key head

Its flat and long shape enables you to massage deep and long muscle groups. It has a similar action to the cylindrical head but covers a larger muscle area. It can be used on the hips as well.

PRECAUTIONS FOR USE

WARNING

- Do not place the device in water, keep it away from liquids or heat sources, and prevent water from entering the product's slots. They must be free of dust and debris . Store it in a dry place.
- Do not drop or use the product incorrectly.
- Do not remove the screws or attempt to disassemble it.
- Do not change or modify the product in any way.
- Do not use the device continuously for more than 1 hour and let it stand for 30 minutes before using it again.
- Only charge the device with the charging cable. Charge the device every 6 months if you are not using it for a long time
- Unplug the charger from the device after charging or before use.
- Do not use or charge the device unattended.
- Use a soft cloth for daily maintenance. Do not use engine oil, organic solvents or others.

HEALTH WARNING

Do not use JOLT™ or any other massage device without the approval of your doctor or physiotherapist or osteopath.

The use of JOLT™ is strongly discouraged in the following cases:

- Pregnancy
- Diabetes with complications such as neuropathy
- Retinal damage
- Heart disease, with pacemaker or defibrillator
- Recent surgery (less than 90 days)
- Epilepsy
- Migraines
- Joint prostheses
- Microtrauma
- Stenosis of the blood
- Tumor area
- Sensitive areas: head, face, cervical vertebrae, spine
- Positions close to the bone with less muscle coverage: humerus, back of the foot, back of the hand
- Other health concerns you may have.

To reduce the risk of personal injury, this product should be used as follows:

1. Use is exclusively for adults, use only on dry parts of the body on clean skin.
2. Use the product only on soft tissues of the body (without pain or discomfort), do not use on hard parts like skull or body bones.
3. Only use massage heads which are recommended for the muscle groups involved.
4. Whether controlled by deep or surface pressure, too much massage can cause hematoma; check the massage area regularly and stop using it as soon

ENVIRONMENTAL PROTECTION

Do not dispose of this product in a conventional trash can at the end of its life. Dispose of it in a recycling garbage can or return it to the JOLT™ Service Center. We will recycle it in an environmentally friendly way.

JOLT

TECHNICAL SPECIFICATIONS

- Model : JOLT Bolt MG-002
- Product type: Muscle massage gun
- 20 adjustable speeds
- Digital display
- 6 interchangeable massage tips (included)
- Brushless motor 55m
- Average battery life of 5-7h depending on the selected power level
- Charging time : 3h
- Rechargeable Lithium battery 2500mAh
- Charging with Micro USB charger (included)
- 3 speeds from 1200rpm to 3200rpm (revolutions per minute) / 20hz à 58hz
- Quiet Technology lower than 50 decibels
- Amplitude : 14mm
- Comfortable non-slip and anti-sweat grip
- On/Off LED indicator • Speed LED indicator
- Charging LED indicator • 950g net weight
- Massage gun size: 24cm x 17cm x 5cm
- Storage box included

FAQ & WARRANTY

FAQ

Any questions or particular concerns? Please visit our Frequently Asked Questions (FAQ) page for answers to all your questions:

MORE QUESTIONS / ANSWERS

WARRANTY

The warranty period for this product is 24 months from the date of purchase.

Wear and tear due to appearance and mechanical operation is not covered by the warranty.

The warranty does not cover:

1. Accidental damages, misuse or transport
2. Incorrect maintenance
3. Operation without following the instructions in the manual
4. Not able to provide a formal purchase certificate

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Tous droits réservés / All rights reserved