JOLT MINI

TECHNOLOGIE BREVETÉE SOLFT FLOW™

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JOLT

CONTENTS

ADVANTAGES OF THE MASSAGE GUN JOLT MINI	page 33
INSTRUCTIONS AND OWNER'S MANUAL	page 35
Instructions	page 35
Charging	page 36
Massage heads.....	page 37
PRECAUTIONS FOR USE	page 39
Warnings.....	page 39
health warnings	page 39
Environmental protection	page 40
TECHNICAL SPECIFICATIONS	page 41
FAQ & WARRANTY	page 42
FAQ	page 42
Warranty	page 43

**Thank you for your order,
we truly appreciate your trust in us.**

In order to use our recovery boots properly, please read the instructions, directions for use and safety information before use. More information is available on our website.

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
The JOLT™ Mini Muscle Massage Gun is a portable, rechargeable vibration massage device. JOLT™ helps to increase blood flow and effectively alleviates the discomfort caused by the accumulation of lactic acid after sports and thus reduce discomfort for people suffering from muscle stiffness / aches.

JOLT™ can be used in the following cases:

- For warming up before a physical effort
- For the recovery after a sporting effort
- To Relieve aches and joint stiffness
- To Stimulate blood circulation
- To Relieve stress and tension
- Massage of knots and trigger points
- Suitable for beginners, amateurs and professionals
- To Relax and to loosen the tissues
- To Reduce lactic acid build up in the muscles
- To Increase joint mobility
- To Fluidize edema in the swollen area
- To Relieve aches and pains
- To Relieve tendinitis

PRESENTATION OF DIFFERENT ATTACHMENTS

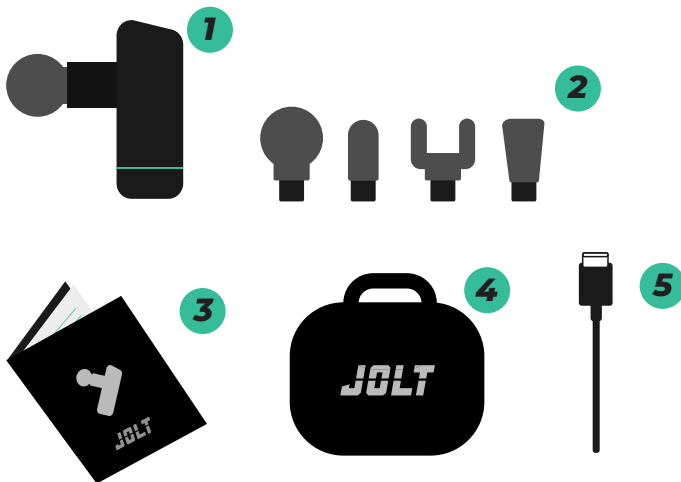
1. 1 JOLT mini massage gun
2. 4 massage heads
3. 1 USB-C charger cable
4. 1 manual
5. 1 carry case

JOLT

INSTRUCTIONS AND OWNER'S MANUAL

INSTRUCTIONS

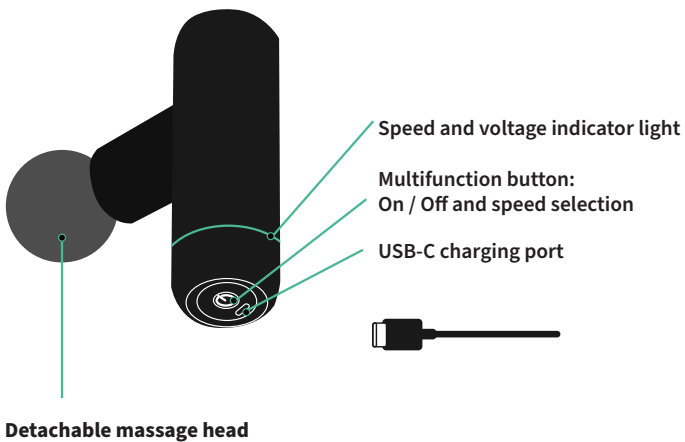1. Press and hold the multifunction button for 1 second to turn the device on, the indicator light around the handle shows that the device is on.
2. Then press the multifunction button briefly, the massage gun will start working, the indicator will indicate the selected speed (there are 3 speeds). Each push on the button will increase the vibration force by 1 speed.

3 massage speeds:

1. **Green** 1%-30% Relaxation
2. **Blue** 31%-60% Deep tissue massage
3. **Red** 61%-100% Professional mode

3. Press and hold the multifunction button for 1 second to turn off the massage gun.

CHARGING

The charging port is on the bottom of the massage gun (USB-C cable) Fully charge the battery for the first use and turn it off while it is charging.

1. Connect the power cable to the gun (USB-C cable). The indicator will flash to show that the device is charging. Do not move the massage gun while it is charging.
2. When the indicator turns **green**, it means the battery is full.

MASSAGE HEADS

With the massage gun turned off, slowly insert the massage head into the barrel until it stops. As shown in the illustration, the position of the red line on the massage head should be aligned with the position of the red line on the gun. To remove the massage head, you just need to pull on the head to remove it.

JOLT™ Mini can be used with 4 massage heads depending on the parts of your body you want to massage.

Round head

Its spherical design with soft coating allows to work large muscle groups such as the back of the thigh and the calf, the triceps, ... It is mainly suitable for new users who are sensitive to vibrations.

Flat head

It allows you to relax any muscle of the body for a classic massage.

Cylindrical head

The conical shape of this head allows an impacting and toning massage with a significant pressure depth. It is suitable for the meridians, soles of the feet, palms of the hands, Achilles heel, rhomboid muscles of the back, ...

U-shaped head

Its shape fits perfectly the spine and cervical/neck. It can also fit the shoulder and other bilateral muscles (do not touch directly the spine and cervical to avoid injury).

PRECAUTIONS FOR USE

WARNING

- Do not place the device in water, keep it away from liquids or heat sources, and prevent water from entering the product's slots. They must be free of dust and debris . Store it in a dry place.
- Do not drop or use the product incorrectly.
- Do not remove the screws or attempt to disassemble it.
- Do not change or modify the product in any way.
- Do not use the device continuously for more than 1 hour and let it stand for 30 minutes before using it again.
- Only charge the device with the charging cable. Charge the device every 6 months if you are not using it for a long time
- Unplug the charger from the device after charging or before use.
- Do not use or charge the device unattended.
- Use a soft cloth for daily maintenance. Do not use engine oil, organic solvents or others.

HEALTH WARNING

Do not use JOLT™ or any other massage device without the approval of your doctor or physiotherapist or osteopath.

The use of JOLT™ is strongly discouraged in the following cases:

- Pregnancy
- Diabetes with complications such as neuropathy
- Retinal damage
- Heart disease, with pacemaker or defibrillator
- Recent surgery (less than 90 days)
- Epilepsy
- Migraines
- Joint prostheses
- Microtrauma
- Stenosis of the blood
- Tumor area
- Sensitive areas: head, face, cervical vertebrae, spine
- Positions close to the bone with less muscle coverage: humerus, back of the foot, back of the hand
- Other health concerns you may have.

To reduce the risk of personal injury, this product should be used as follows:

1. Use is exclusively for adults, use only on dry parts of the body on clean skin.
2. Use the product only on soft tissues of the body (without pain or discomfort), do not use on hard parts like skull or body bones.
3. Only use massage heads which are recommended for the muscle groups involved.
4. Whether controlled by deep or surface pressure, too much massage can cause hematoma; check the massage area regularly and stop using it as soon

ENVIRONMENTAL PROTECTION

Do not dispose of this product in a conventional trash can at the end of its life. Dispose of it in a recycling garbage can or return it to the JOLT™ Service Center. We will recycle it in an environmentally friendly way.

JOLT

TECHNICAL SPECIFICATIONS

- Model : JOLT™ Mini
- 3 adjustable speeds
- 4 interchangeable attachments (included)
- 85W high performance brushless motor
- Average autonomy of 3 hours per charge depending on the chosen power of use
- Charging time : 3hours
- Rechargeable Lithium battery of 2500mAh
- USB-C cable charging (input 5V 2A - voltage 11,1V)
- 3 speeds from 1300rpm to 3300 rpm (Revolutions Per Minute) / 20hz to 60hz
- Quiet technology below 45 decibels
- Amplitude : 10mm
- Comfortable anti-slip grip
- ON/OFF LED indicator
- Speed LED indicator
- Charging LED indicator
- Cordless and lightweight (net weight 420g)
- Massage gun size: 14.9cm x 9cm x 4.75cm

FAQ & WARRANTY

FAQ

CAN'T TURN ON YOUR JOLT MINI?

1. Too low temperature:

If you keep the massage gun below 0°C for a long time, it will reduce the battery activity and prevent normal operation.

Solution: Place the product for two hours in a place where the temperature is above 0°C.

2. Low battery :

When the battery is low or the user does not charge it for a long period of time, the battery enters a «protection» mode.

Solution : Turn off the massage gun and plug in the charger. Wait until it is fully charged before operating it normally again.

MORE QUESTIONS / ANSWERS

WARRANTY

The warranty period for this product is 24 months from the date of purchase.

Wear and tear due to appearance and mechanical operation is not covered by the warranty.

The warranty does not cover:

1. Accidental damages, misuse or transport
2. Incorrect maintenance
3. Operation without following the instructions in the manual
4. Not able to provide a formal purchase certificate

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Tous droits réservés / All rights reserved