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**JOLT BOOTS**

## TABLE OF CONTENTS

<b>Chapter 1 : Introduction to the different accessories</b>	<b>page 44</b>
<b>Chapter 2 : Installation and operating instructions</b>	<b>page 46</b>
<b>Chapter 3 : Technical specifications</b>	<b>page 53</b>
<b>Chapter 4 : Special precautions</b>	<b>page 54</b>
<b>Chapter 5: Maintenance and Cleaning</b>	<b>page 57</b>
<b>Chapter 6 : Storage and transport</b>	<b>page 57</b>

**Thank you for your order  
and your trust.**

*In order to use our massage recovery boots properly,  
please read the instructions, directions for use  
and safety guidelines before use. You can get more  
information by visiting our website.*

*[www.getjolt.com](http://www.getjolt.com)*

**JOLT**BOOTS

## JOLT™BOOTS COMPRESSION BOOTS ADVANTAGES

JOLT™ Recovery Boots are effective after hard workouts, whether you're a professional or an amateur. You'll recover much faster while reducing the risk of inflammation and muscle soreness. JOLT™ Boots operate on the massage principle to relax your muscles and remove metabolic waste products that build up in your muscles during training. Our recovery boots provide a passive, targeted and progressive uniform massage. You'll be able to complete another session the following day with your muscles regenerated. JOLT™ Boots will help you in the medium and long term to improve your performance and your mobility.

## OPERATING PRINCIPLE

JOLT™ Boots is an air pressure system activated via the control unit to supply 6 «air chambers» per leg. The modes of operation, pressure system, duration and other parameters are set through the interface of the control unit. The air chambers are controlled by on-board electronics (sequential compression, rest, inflation and deflation times) to apply periodic pressure to the blood and human peripheral tissues.

You can use JOLT™ Boots at home on a power outlet or on the move thanks to its built-in rechargeable battery (5200mAh = up to 4 hours of battery life).

For optimal recovery, we recommend that you use your JOLT™ Boots within the first hour after your workout if possible. If that is not possible, you should do a JOLT™ Boots session as soon as possible after your workout.

# CHAPTER 1 : INTRODUCTION TO THE DIFFERENT ACCESSORIES

Our compression boots include a central unit (control box), 2 connection hoses, 1 single to double adapter hose, 2 leg sleeves and a charging cable with power adapter. Take the products and accessories out of the packaging and make sure they are complete by referring to the «list» below. Also check that the products and accessories are not damaged.

Model: JOLT BOOTS


### Packing list

### Quantité


1. Leg sleeves, connecting piece incl.	2
2. Control unit	1
3. Mains adapter	1
4. Carrying case	1
5. Single to dual adapter	1
6. Connecting hose	2
7. User Manual	1

### Control unit display details

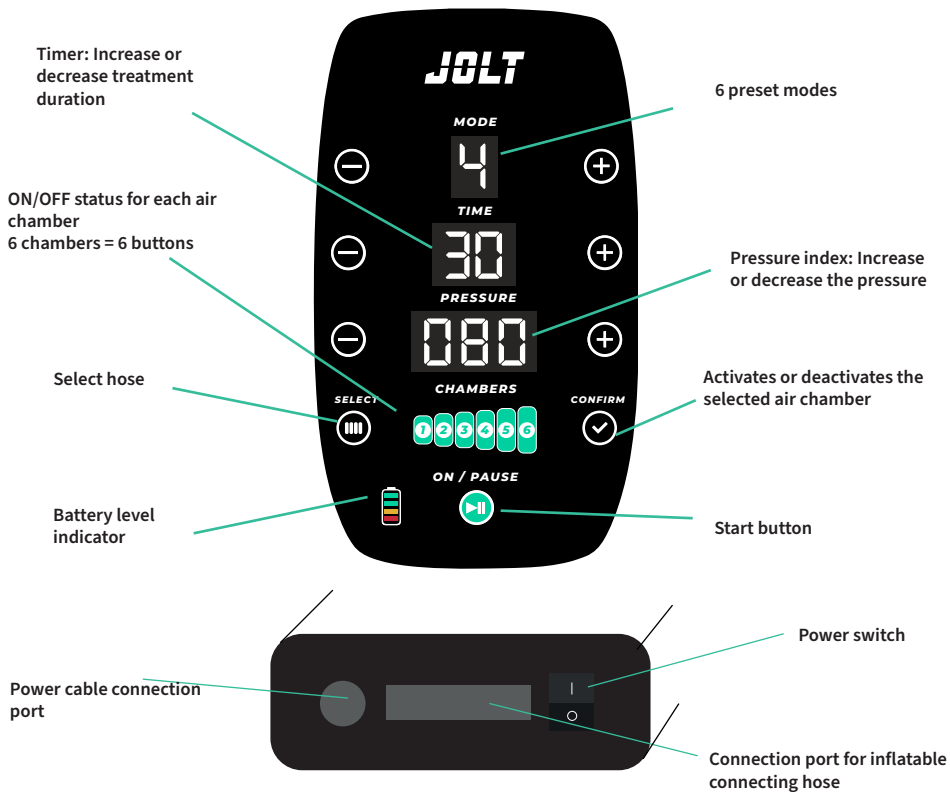
## CHAPTER 2 : INSTALLATION AND OPERATING INSTRUCTIONS


### 2.1 Connections before use

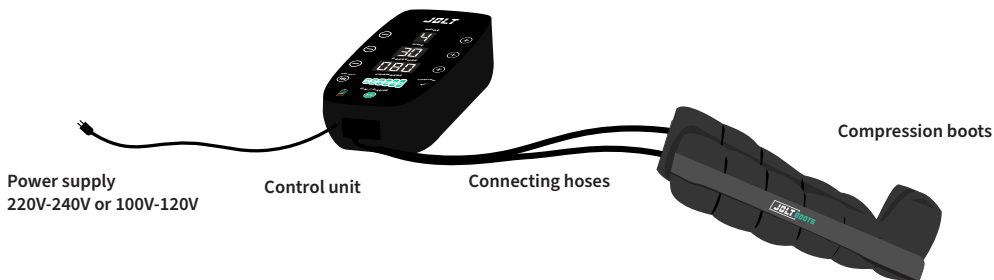**Before use :** Check the battery level by turning on the power switch on the front of the control unit (the symbol «I» turns the device on). If the control unit does not turn on, it means the battery is discharged, so plug in the charger before using.

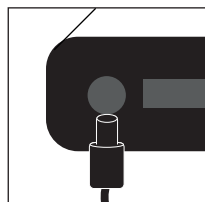
A) Hose connection for single sleeve use: connect the one end of the «connection hose» to the «connection port» of the control unit and the other end to the sleeve.

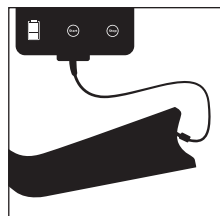**B)** Hoses connection for simultaneous use of both sleeves: First connect one end of the «single to double adapter» to the «connection port» of the control unit. Then connect the other end of the «Single to Dual Adapter» to each of the 2 « connecting hoses». Finally, connect the other end of the « connecting hoses » to the 2 sleeves.

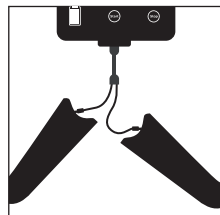**Note :** There is a connection direction when you connect the hoses and connectors. The arrow shows you the correct connection direction. Make sure that the « single to dual adapter» and the «connecting hoses» are correctly connected.

## 2.2 Use of JOLT™ compression boots

Before you start with your session, it is important to stay comfortable and relaxed. Adopt the right position ( sitting or lying ). It is recommended to slightly lift the limbs during the session so that they are at the same level as the heart, more convenient for the lymphatic and blood circulation in these areas.

### A) Put on the recovery boots

You should pull your boots up as high as you can. Only your toes should stick out at the end. This way you ensure that your entire legs benefit from the compression technology.

### B) Turn on the device

Turn on the power switch, the screen will switch on.


**Note :** If the screen does not turn on after starting, check the battery status. Please charge the battery before use or plug in the charger if the battery is discharged.

## C) Select the mode of operation

The device has 4 modes. The default mode is A. You can only switch to another mode when the device is in pause.


### Mode 1

Mode 1 operates like a flow, exerting pressure in a specific area and is therefore ideal for muscle relaxation. The muscles are gradually compressed from the feet to the chest. The use of this mode before sport should be short so as not to deprive your legs of the muscle tone they need for the exercise.

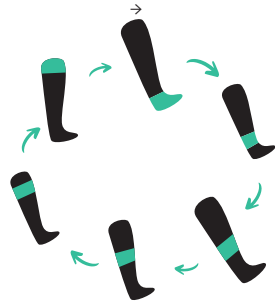
### Mode 2

Mode 2 helps you to get rid of lactic acid after an intensive workout by applying continuous and progressive pressure on your upper legs for fast recovery. Recovery is most efficient at high pressure. We suggest using your JOLT™ Boots as soon as possible after you complete your workout.
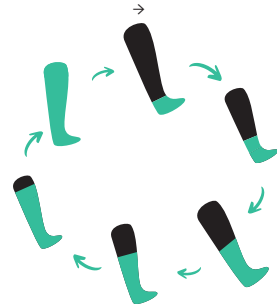
### Mode 3


It is a combination of the programs 1 and 2. Several chambers remain inflated at the same time, thus avoiding reflux. This mode also promotes lactate removal. Suitable for a faster but intensive recovery session, as your whole leg is not permanently under pressure (unlike mode 2).

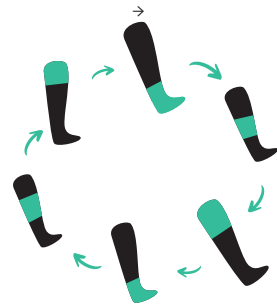
## Mode 4

This mode is for short, intense sessions. Because all chambers are inflated at the same time, the number of possible compression phases in a given time is maximized.

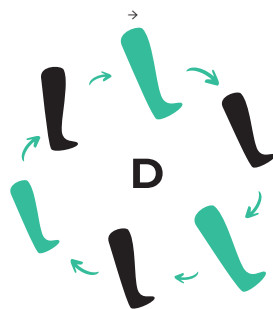
## Mode 5

Mode 5 combines mode 1 and 2. This mode starts with mode 1, once mode 1 is completed, it switches to mode 2. The cycle then repeats. The objective is to couple the use of relaxation with a deep recovery.
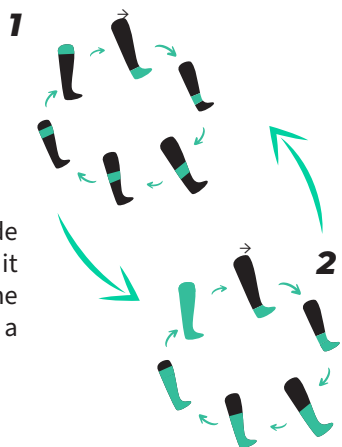
## Mode 6

Mode 6 combines mode 3 and 4. This mode starts with mode 3, once mode 3 is completed, it switches to mode 4. The cycle then repeats. The objective is to couple a fast recovery with an intense pressure.
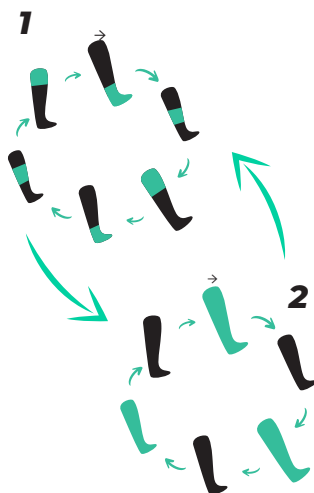
## D) Select the duration

Press the corresponding button on the control unit to set the duration of use, adjustable from 10 to 90 minutes, the default duration of use is 30 minutes. The optimal duration is between 30 and 50 minutes. Following an easy workout, you can start with 30 minutes. If you have had an intense or even very intense session, you can increase the duration of use to 50 minutes. Exceeding 1 hour of use will not necessarily bring significant additional benefits.

## E) Set the pressure

Press the pressure buttons «+» and «-» to set the pressure, adjustable from 80 mmHg (10.7 kPa) to 260 mmHg (34.7 kPa). By default, the device starts with a pressure of 80 mmHg (10.7 kPa).

## F) Chambers deactivation (optional)

You can choose to deactivate one or more chambers in case of injury/pain on the concerned area.

For example, when you press the chamber selection button, the light on the control box for chamber 1 will flash. It will not inflate. Do the same to reactivate the chamber or chambers concerned.

**Note :** All chambers are activated by default. Each chamber can be activated or deactivated individually.

## G) Starting the session

After setting the above parameters, press «Start (ON)» to start using the device, a countdown timer will start. If you wish to stop/pause the device during the session, press «STOP (PAUSE)».

## H) End of cycle

At the end of the session, the device automatically shuts off the air outlet and emits an audible alert. The display screen will reset to its original state. The boots can be removed then.




## I) Switching off the device

Turn off the power switch on the front of the unit by setting it to the «O» OFF position, then unplug the power cable to completely power off the unit.

## J) Storage and instructions after use

To store your equipment after use, you can remove ( not obligatory) the connection hoses by pressing the clips on both sides of the «connection hoses» and the « single to dual adapter».

Do not remove the individual inflation hoses (6 per legging) so often. This is only necessary if you travel and absolutely want to save space. Otherwise, you should avoid it because frequent « removing « may damage the individual inflation tubes. To remove, it is recommended that you hold the plastic cap and then pull it.

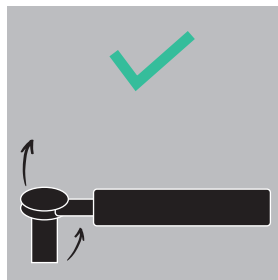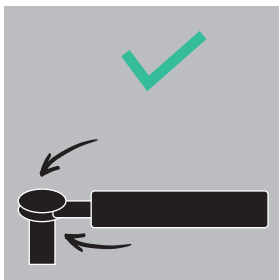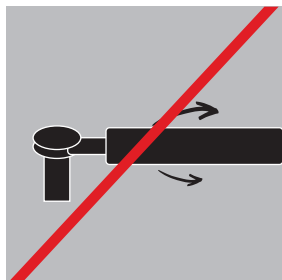**Do not pull on the hose. Disconnect by holding and removing the plastic sealing plug.**

### ⚠ Warning :

- It is important to relax during use. If you experience any discomfort, turn off the device immediately or readjust it.
- Because of different pressure tolerances and the different leg girths of every person, it is impossible to determine exactly the pressure value. However, we recommend initially setting the pressure at 80 mmHg and then gradually adjusting it according to your personal comfort.
- A pressotherapy session is generally 80-120 mmHg for 20 to 30 minutes, one or two sessions a day for 6 to 10 days.

## 2.3 Charging the battery

Plug the charger into an electrical outlet, the other end into the connection port on the control unit: the indicator light on the adapter is red: the device is charging. The indicator light turns green when the battery is fully charged.

### Warning :

Please unplug the power cable from the unit and turn it off when the battery is fully charged.

The device will turn on while recharging if left plugged in.

Should the fully charged battery remain plugged into an electrical outlet, a bidirectional power will occur which may damage both the battery and the circuit board ( affecting the device's lifespan and increasing the risk of overheating caused by overload). So it is important to unplug the power cable when the battery is fully charged.

The charger will keep heating up if it is connected to a power outlet for a long time, and therefore the various components and materials deteriorate more quickly. For longer charger life, always keep the charger unplugged from the wall outlet or turn off the power switch of the outlet.

## CHAPTER 3 : TECHNICAL SPECIFICATIONS

- **Pressure:** The pressure is adjustable from 80 mmHg (10.7 kPa) to 260 mmHg (34.7 kPa). The default pressure at startup is 80 mmHg (10.67 kPa).
- **Treatment time (session):** Each session duration can be programmed between 10 and 90 minutes. The default time is 30 minutes.
- **Sound intensity in operation:** Ambient sound  $\leq 20$  dB
- **Power consumption:** 49W
- **Normal power supply requirements:** AC 100-240 V, 50-60Hz DC 12.6 V/4 A output OR Lithium battery 12.6 V, 5200 mAh.
- **Emergency stop button:** The treatment process can be stopped at any time by pressing the emergency stop button.
- **Compression modes and number of chambers:** 6 chambers and 4 compression modes.
- **Chamber deactivation:** individually via the control unit.
- **Rechargeable battery:** Built-in lithium battery.
- **Normal conditions of use:**  
Operating temperature: between 5 °C and 40 °C  
Degree of humidity:  $\leq 85$  %  
Air pressure: between 80 kPa and 106 kPa.
- **Net weight of product :** 1,50 kg
- **Control unit size (length x width x height) :**  
210 mm x 130 mm x 86 mm.

## CHAPTER 4 : SPECIAL PRECAUTIONS

### 4.1 Warnings

- Do not put the product in water and do not use it in a wet or humid environment, which could cause an electrical short circuit or fire.
- In the event of any abnormal situation or incident during use, immediately switch off the device.
- Plug the device only into a properly grounded outlet.
- Do not kink or twist the power cable to avoid fire or electric shock.
- In case of worn, damaged or loose parts, please contact our company immediately (hello@getjolt.com)
- Do not use the unit in tilted, vibrating or impacted environments.
- Do not pull the recovery boots off while the device is running.
- Only remove the sleeves when the unit is switched off.
- Place the control unit on a flat surface.
- Do not move with the recovery boots during use.
- Please consult your doctor before using this device. This product is not intended to diagnose, treat or cure any disease.

### 4.2 Precautions for use

- Please check that the device, the connecting hoses, the single to dual adapter and the sleeves are in good condition.
- Ensure that the connecting cables are properly connected and secured before use.
- Before using the device, ensure that it is functioning properly.
- Using damaged or non-compliant electrical cables ( following a modification) may cause fire and electric shock.
- Using damaged or non-compliant electrical cables ( following a modification) may cause fire and electric shock.
- Turn off the power and unplug the power cable before moving the device.
- During transport, the device should be subjected to as little vibration as possible.

- Be careful not to exert too much pressure on the device during transport so as not to damage its overall appearance.
- The device should be placed on a stable surface. Allow at least 5 cm clearance around it to allow air circulation and efficient heat dissipation.
- Do not use the device during thunderstorms.
- To avoid injury, do not bring the inflatable connecting hose into contact with the nose, mouth, ear or any other part of the body.
- Keep the Boots away from heat sources such as stove, cigarette butts, etc. to reduce the risk of leakage and to prevent fire.
- Before use, remove jewelry, watches, and other metal objects that could cause physical injury or damage the sleeves.
- Keep pointed or sharp instruments (needles, scissors, etc.) away from the sleeve. They may puncture the sleeve and cause failure of the device.
- Do not inflate the sleeve with other tools so as not to damage it and to avoid any physical injury.
- Sharp metallic objects such as paper clips, staplers, needles or food must not be inserted into the device.
- In the event of an emergency or physical discomfort (pain, paralysis, dyspnea, dizziness, etc.), it is highly recommended to immediately press the stop button or to remove the inflatable connecting hose.
- For the elderly and people with poor arterial elasticity, the compression level can be gradually increased during the session before reaching the required pressure. Do not start with too high a pressure, but increase the intensity gradually. It is important to use the sleeve in a lying position with the body as relaxed as possible.
- It is recommended to wear underwear and socks when using the device to prevent cross-infection from the number of people using the same sleeve.
- Do not disassemble, attempt to repair or otherwise modify the design of the device. Should you need repair, please contact our customer service department at [hello@getjolt.com](mailto:hello@getjolt.com)
- Dispose of the product and its different accessories in accordance with local environmental protection regulations to prevent environmental pollution.

- Do not use the product in environments with any of the following conditions: intense sunlight, humidity, near a water source; vibrations, inclined surfaces; in a very strong magnetic field; electromagnetic waves, on-site pulse voltage; storage location of chemicals or corrosive gases.
- Do not use in the presence of flammable anesthetic mixtures with air, oxygen or nitric oxide.

**Please pay special attention to the following points to properly store the product.**

- a) Do not store the products in high temperature and humidity conditions, in direct sunlight, in dusty environments or near corrosive gases.
- b) Do not place the products near sources of water or moisture.
- c) Before you store the product, all the accessories, including the power cable should be cleaned to prevent failure the next time you use the product.

## CHAPTER 5: MAINTENANCE AND CLEANING

### 5.1. Maintenance

A) Regularly check the product and its accessories for cracks, damage, distortion or poor connection. Contact our customer service at [hello@getjolt.com](mailto:hello@getjolt.com) should this be the case.

B) Store the product and its accessories in a cool, dry place to protect them from moisture, dust, and direct sunlight. When you do not use the device for a long period of time, it is recommended to store it in its original packaging.

### 5.2 Cleaning

You can remove any dust with a soft, dry cloth or paper towel. Do not use oil, benzene, alcohol, gasoline or other chemicals to clean the device.

**Note :** Before cleaning or maintaining the pressotherapy boots, check that the power supply has been switched off.

## CHAPTER 6 : STORAGE AND TRANSPORT

### 6.1 Storage

Keep the boots and the accessories in a well-ventilated area away from direct sunlight or flammable, explosive and corrosive gases.

### 6.2 Transport

The device can be transported by any means, but must be protected from overturning, falling, strong impact, vibration, intense sunlight, rainwater or snow during transport. Do not transport the product with flammable, explosive or corrosive substances, as specified in the order contract.

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Tous droits réservés / All rights reserved