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KAPITEL 7 : REPARATUR

TABLE OF CONTENTS

Chapter 1 : Introduction to the different accessories	page 48
Chapter 2 : Installation and operating instructions	page 50
Chapter 3 : Technical specifications	page 57
Chapter 4 : Special precautions	page 58
Chapter 5: Maintenance and Cleaning	page 61
Chapter 6 : Storage and transport	page 61
Chapter 7 : Repair	page 62

**Thank you for your order
and your trust.**

*To properly use our pressotherapy arm sleeves,
please first read the instructions, directions for use
and safety guidelines. More information is available
on our website.*

www.getjolt.com

JOLT™ WINGS COMPRESSION WINGS ADVANTAGES

JOLT™ Wings Recovery Arm Sleeves are very efficient after intense workouts, for both professionals and amateurs. You can recover faster and also reduce muscle soreness and the risk of inflammation.

JOLT™ Wings is based on the massage principle to relax the muscles and get rid of the metabolic waste products accumulated in the muscles. With our pressotherapy arm cuffs, you can uniformly massage the muscles in a targeted, progressive, passive way. You can train with regenerated muscles the next day.

JOLT™ Wings help you improve your performance and increase your mobility in the mid and long term.

OPERATING PRINCIPLE

JOLT™ Wings is an air pressure system, activated via the control unit to supply 6 «air chambers» per arm. The operating modes, pressure system, time and other parameters are adjustable via the interface of the control unit. The chambers are controlled by on-board electronics (sequential compression, rest time, inflation and deflation) to exert periodic pressure on human blood and peripheral tissue. You can use JOLT™ Wings at home on a power outlet or on the go with its built-in battery (5000mAh = 4h).

CHAPTER 1 : INTRODUCTION TO THE DIFFERENT ACCESSORIES

Our compression sleeves consist of a central unit (control box), 2 connecting hoses, 2 arm sleeves and a charging cable with power adapter. Remove the product(s) and all accessories from the packaging and check that it is complete by referring to the «list» below. Also check that the products and accessories are not damaged.

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
Packing list

	Quantity
1. Arm sleeve with connecting pipe	1
2. Carrying case	1
3. Instruction manual	1


Detail of control box display

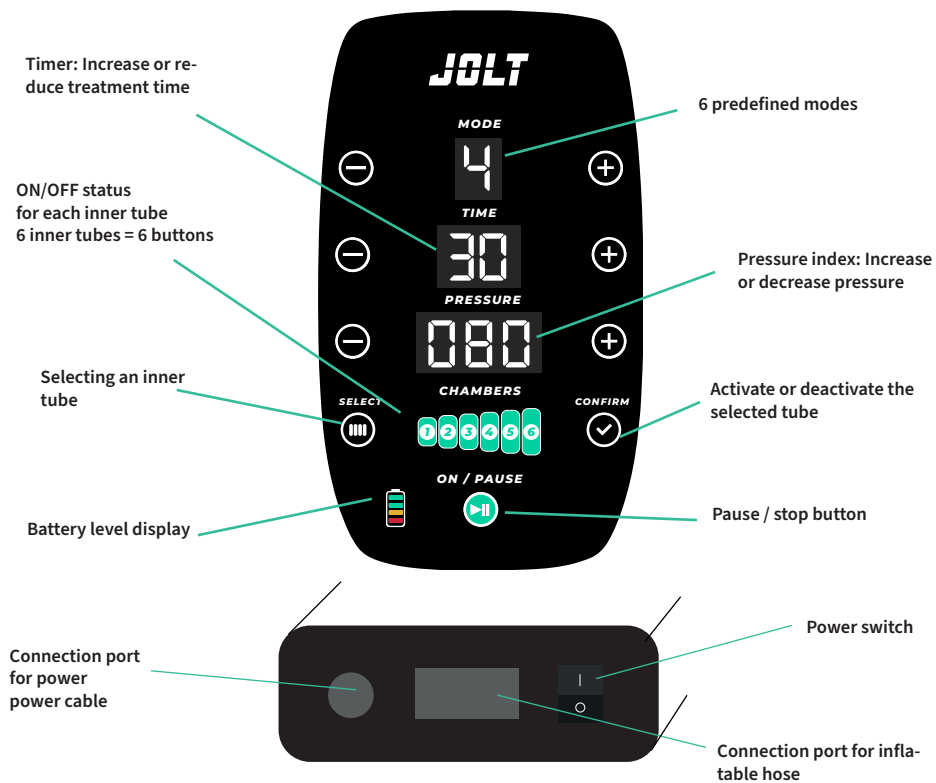
CHAPTER 2 : INSTALLATION AND OPERATING INSTRUCTIONS

2.1 Connections before use

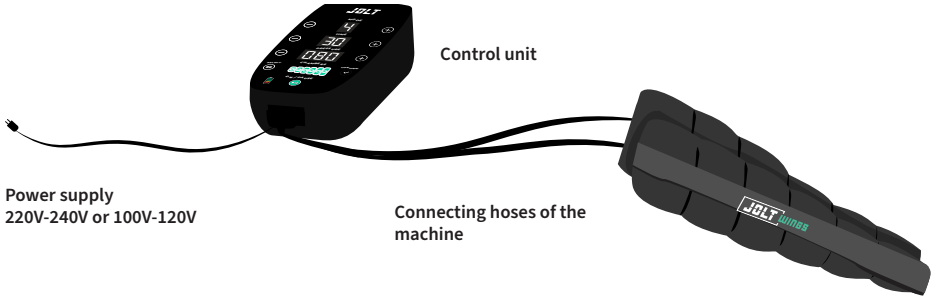
Before use: Check the battery level by pressing the power button on the control box (the power button is located on the front of the control box; the «I» symbol switches the device on). If the control box does not light up, the battery is discharged. Plug in the charger before using.

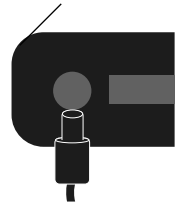A) Connecting hoses for single-leg use: connect one end of the «connecting hose» to the «connection port» on the control box, and the other end to the sleeve.

If you purchase 2 sleeves, a connector will also be supplied so that you can use your 2 sleeves simultaneously (same connection method as JOLT™ Boots).

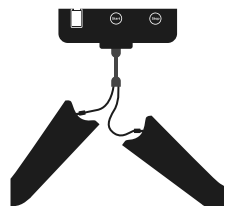B) Connecting hoses for simultaneous use of both leggings: first connect one end of the «single-to-double adapter» to the «connection port» on the control box. Then connect the other end of the «single-to-double adapter» to each of the 2 «connection hoses». Finally, connect the other end of the «connection hoses» to the 2 sleeves.

Note: When connecting pipes and connectors, the arrow indicates the direction of connection. Check that the «single-to-double adapter» and the «connecting hoses» are properly connected.

2.2 Using JOLT™ compression sleeves

Before starting your session, try to remain comfortable and relaxed. Adopt the right position (sitting or lying down). During the session, it's a good idea to lift your limbs slightly so that they're at the same level as your heart, which is more conducive to lymphatic and blood circulation in these areas.

A) Put on the recovery sleeves

Your sleeves should be pulled as high as possible. Only the ends of your hands should stick out at the end. In this way, you ensure that your whole arm benefits from the compression technology.

B) Switch on the appliance


Note: if the display does not light up after start-up, check the battery level. Charge it before use, or plug in the mains charger if it is discharged.

C) Select operating mode

There are 4 modes. By default, the device switches to mode A. You can only switch to another mode when the unit is paused.


Mode 1: Warm-up and relaxing use

Mode 1 works like a flow, exerting punctual pressure and is therefore perfectly suited to relaxing muscles. Muscles are progressively compressed from the feet to the chest. Pre-sport use should be short, so as not to deprive your legs of the muscle tone they need for the effort that follows.

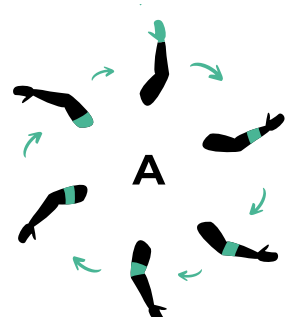
Mode 2: Deep recovery

Thanks to the continuous, progressive pressure exerted upwards on your legs, mode 2 helps eliminate lactic acid after intensive training, so you can regenerate ultra-fast. Recovery is particularly effective at high pressure. We recommend using your JOLT™ Wings as soon as possible after the end of your workout.

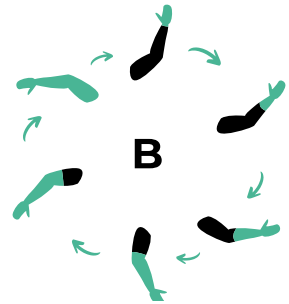Mode 3: Fast session recovery

This mode combines the flow of programs 1 and 2. Several chambers remain inflated at the same time, preventing reflux. This mode also promotes lactate elimination. Suitable for a rapid but intensive recovery session, as your whole leg is not permanently under pressure (unlike mode 2). This mode is also effective for lymphatic drainage.

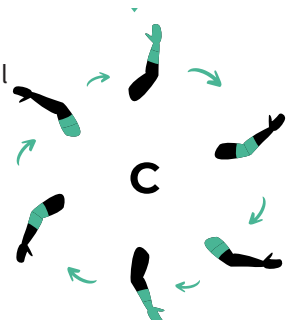Mode 4: Short & intense


This mode is designed for short, intense sessions. As all chambers are inflated at the same time, the number of possible compression phases in a given time is maximized.

Mode 5

Mode 5 combines mode 1 and 2. This mode starts with mode 1, and once mode 1 is completed, it switches to mode 2. The cycle then repeats itself. The aim is to combine relaxation with deep recovery.

Mode 6

Mode 6 combines modes 3 and 4. This mode starts with mode 3, and once mode 3 is completed, it moves on to mode 4. The cycle then repeats itself. The aim is to combine rapid recovery with intense pressure.

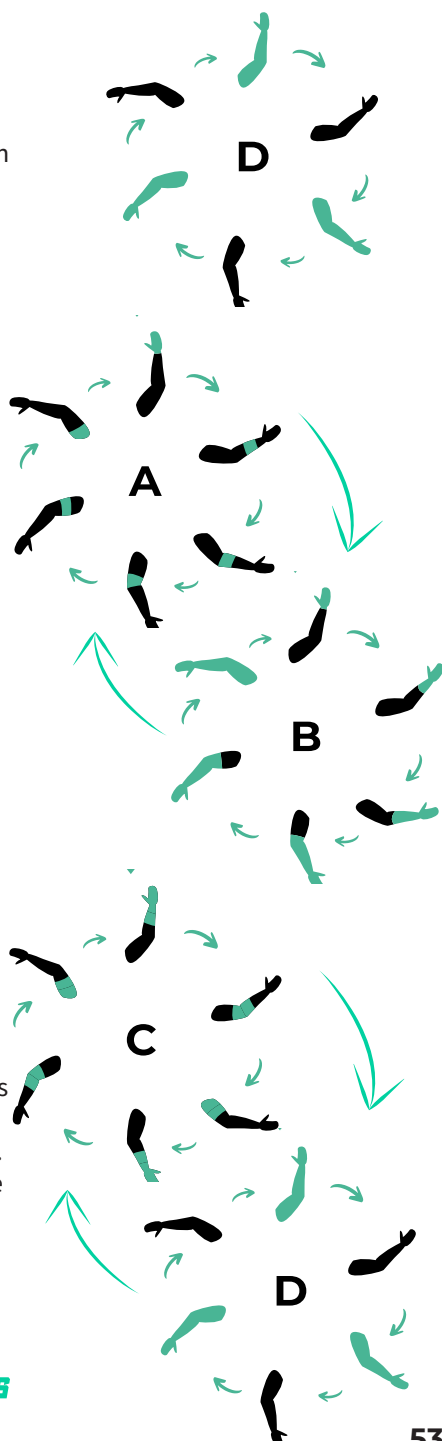D) Select the duration

Press the corresponding button on the control unit to set the duration of use. The default duration of use is 30 minutes, but is adjustable from 10 to 90 minutes. The optimal duration is between 30 and 50 minutes. After light exercise, you can start with 30 minutes. If your session was intense or very intense, you can increase the time to 50 minutes. A session longer than 1 hour will not necessarily provide significant additional benefits.

E) Set the pressure

Press the pressure buttons «+» and «-» to set the pressure. It is adjustable from 80 mmHg (11 kPa) to 220 mmHg (28 kPa). The default pressure at startup is 80 mmHg.

F) Chambers deactivation (optional)

There are six buttons (1 per chamber): 1, 2, 3, 4, 5, 6. You can choose to disable one or more chambers in case of injury/pain in the affected area.

For example, by pressing 1 (hand), the indicator light for chamber 1 on the control unit turns off, so this chamber will not inflate.

Note: All chambers are activated by default. You can disable each air chamber individually.

Example: If you disable chamber 3, it will be disabled for both cuffs simultaneously.

G) Starting the session

After completing the above settings, press the Start button (triangle logo «play») to get started with the device, a countdown timer will start. If you need to pause the device during your session, press the Pause button (2 vertical bars). If you want to stop your session, press the Stop button (square logo).

H) End of cycle

At the end of the session, the device automatically shuts off the air outlet and beeps. The screen will revert to its original state. Wait 2 minutes for the cuffs to deflate, then you can «unzip» and remove them.

I) Switching off the device




Press the power switch on the front of the unit to OFF and then unplug the power cable to shut it down.

J) Storage and instructions after use

If you want to store your equipment after use, remove the connecting hoses by pressing the clip on the connector connected to the control unit.

It is not recommended to remove the individual inflation hoses (6 along each cuff) as frequent «removal» may damage them.

If you still want to remove them, it is recommended to disconnect them by holding and removing the plastic sealing cap, not by pulling.

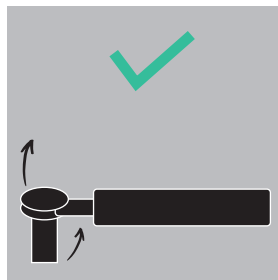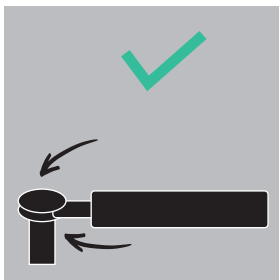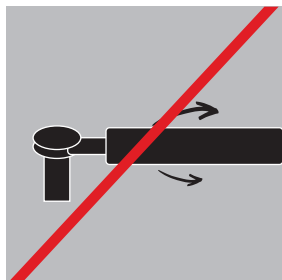Do not pull on the hose. Disconnect by holding and removing the plastic sealing plug.

⚠ Warning :

- It is important to relax during use. If you experience any discomfort, turn off the device immediately or readjust it.
- Because of different pressure tolerances and the different leg girths of every person, it is impossible to determine exactly the pressure value. However, we recommend initially setting the pressure at 80 mmHg and then gradually adjusting it according to your personal comfort.
- A pressotherapy session is generally 80-120 mmHg for 20 to 30 minutes, one or two sessions a day for 6 to 10 days.

2.3 Battery Recharge

Plug the charger into an electrical outlet, the other end into the connection port on the control unit: the indicator light on the adapter is red: the device is charging. The indicator light turns green when the battery is fully charged.

Warning :

Please unplug the power cable from the unit and turn it off when the battery is fully charged.

The device will turn on while recharging if left plugged in.

When a fully charged battery is left plugged into an outlet, a bidirectional power supply will occur, which may damage the battery and the circuit board (The life of the device will be shortened and the risk of overheating from overcharging will increase). It is therefore necessary to remove the power cable after the battery is fully charged.

When plugged into an outlet for a long time, the charger will continue to heat up, causing the various components and materials to age faster. To extend the life of the charger, always make sure that it is unplugged from the outlet or that the power switch on the outlet is turned off.

CHAPTER 3 : TECHNICAL SPECIFICATIONS

- **Pressure:** The pressure is adjustable from 80 mmHg (10.7 kPa) to 260 mmHg (34.7 kPa). The default pressure at startup is 80 mmHg (10.67 kPa).
- **Treatment time (session):** Each session duration can be programmed between 10 and 90 minutes. The default time is 30 minutes.
- **Sound intensity in operation:** Ambient sound \leq 20 dB
- **Power consumption:** 42,9W
- **Normal power supply requirements:** AC 100-240 V, 50-60Hz DC 12.6 V/4 A output OR Lithium battery 12.6 V, 5000 mAh.
- **Stop button:** The device has an emergency stop button that allows you to stop the treatment process at any time.
- **Compression modes and number of chambers:** 4 compression modes and 6 air chambers.
- **Chamber deactivation:** individually via the control unit.
- **Rechargeable battery:** Integrated lithium battery.
- **Normal conditions of use:**
Operating temperature: between 5 °C and 40 °C
Degree of humidity: \leq 85 %
Air pressure: between 80 kPa and 106 kPa.
- **Net weight of product :** 1,50 kg
- **Control unit size (length x width x height) :**
210 mm x 130 mm x 86 mm.

CHAPTER 4 : SPECIAL PRECAUTIONS

4.1 Warnings

- Do not put the product in water and do not use it in a wet or humid environment, which could cause an electrical short circuit or fire.
- In the event of any abnormal situation or incident during use, immediately switch off the device.
- Plug the device only into a properly grounded outlet.
- Do not kink or twist the power cable to avoid fire or electric shock.
- In case of worn, damaged or loose parts, please contact our company immediately (hello@getjolt.com)
- Do not use the unit in tilted, vibrating or impacted environments.
- Do not pull the recovery Wings off while the device is running.
- Only remove the sleeves when the unit is switched off.
- Place the control unit on a flat surface.
- Do not move with the recovery Wings during use.
- Please consult your doctor before using this device. This product is not intended to diagnose, treat or cure any disease.

4.2 Precautions for use

- Please check that the device, the connecting hoses, the single to dual adapter and the sleeves are in good condition.
- Ensure that the connecting cables are properly connected and secured before use.
- Before using the device, ensure that it is functioning properly.
- Using damaged or non-compliant electrical cables (following a modification) may cause fire and electric shock.
- Using damaged or non-compliant electrical cables (following a modification) may cause fire and electric shock.
- Turn off the power and unplug the power cable before moving the device.
- During transport, the device should be subjected to as little vibration as possible.

- Be careful not to exert too much pressure on the device during transport so as not to damage its overall appearance.
- The device should be placed on a stable surface. Allow at least 5 cm clearance around it to allow air circulation and efficient heat dissipation.
- Do not use the device during thunderstorms.
- To avoid injury, do not bring the inflatable connecting hose into contact with the nose, mouth, ear or any other part of the body.
- Keep the Wings away from heat sources such as stove, cigarette butts, etc. to reduce the risk of leakage and to prevent fire.
- Before use, remove jewelry, watches, and other metal objects that could cause physical injury or damage the sleeves.
- Keep pointed or sharp instruments (needles, scissors, etc.) away from the sleeve. They may puncture the sleeve and cause failure of the device.
- Do not inflate the sleeve with other tools so as not to damage it and to avoid any physical injury.
- Sharp metallic objects such as paper clips, staplers, needles or food must not be inserted into the device.
- In the event of an emergency or physical discomfort (pain, paralysis, dyspnea, dizziness, etc.), it is highly recommended to immediately press the stop button or to remove the inflatable connecting hose.
- For the elderly and people with poor arterial elasticity, the compression level can be gradually increased during the session before reaching the required pressure. Do not start with too high a pressure, but increase the intensity gradually. It is important to use the sleeve in a lying position with the body as relaxed as possible.
- It is recommended to wear underwear and socks when using the device to prevent cross-infection from the number of people using the same sleeve.
- Do not disassemble, attempt to repair or otherwise modify the design of the device. Should you need repair, please contact our customer service department at hello@getjolt.com
- Dispose of the product and its different accessories in accordance with local environmental protection regulations to prevent environmental pollution.

- Do not use the product in environments with any of the following conditions: intense sunlight, humidity, near a water source; vibrations, inclined surfaces; in a very strong magnetic field; electromagnetic waves, on-site pulse voltage; storage location of chemicals or corrosive gases.
- Do not use in the presence of flammable anesthetic mixtures with air, oxygen or nitric oxide.

Please pay special attention to the following points to properly store the product.

- a) Do not store the products in high temperature and humidity conditions, in direct sunlight, in dusty environments or near corrosive gases.
- b) Do not place the products near sources of water or moisture.
- c) Before you store the product, all the accessories, including the power cable should be cleaned to prevent failure the next time you use the product.

CHAPTER 5: MAINTENANCE AND CLEANING

5.1. Maintenance

A) Regularly check the product and its accessories for cracks, damage, distortion or poor connection. Contact our customer service at hello@getjolt.com should this be the case.

B) Store the product and its accessories in a cool, dry place to protect them from moisture, dust, and direct sunlight. When you do not use the device for a long period of time, it is recommended to store it in its original packaging.

5.2 Cleaning

You can remove any dust with a soft, dry cloth or paper towel. Do not use oil, benzene, alcohol, gasoline or other chemicals to clean the device.

Note : Before cleaning or maintaining the pressotherapy Wings, check that the power supply has been switched off.

CHAPTER 6 : STORAGE AND TRANSPORT

6.1 Storage

Keep the Wings and the accessories in a well-ventilated area away from direct sunlight or flammable, explosive and corrosive gases.

6.2 Transport

You can use any means of transport to transport the device. However, it must be protected from falling down, strong impacts, vibrations, strong sunlight, rain water or snow during transport. Avoid mixtures with flammable, explosive and corrosive substances.

CHAPTER 7 : REPAIR

Do not attempt to repair the unit by yourself in case of a failure. Please contact us as soon as possible at hello@getjolt.com

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Tous droits réservés / All rights reserved