



Cordless Water Flosser

Instructions manual

Need help? Contact us!




support@shopremi.com




1800 234 1874

Product Components

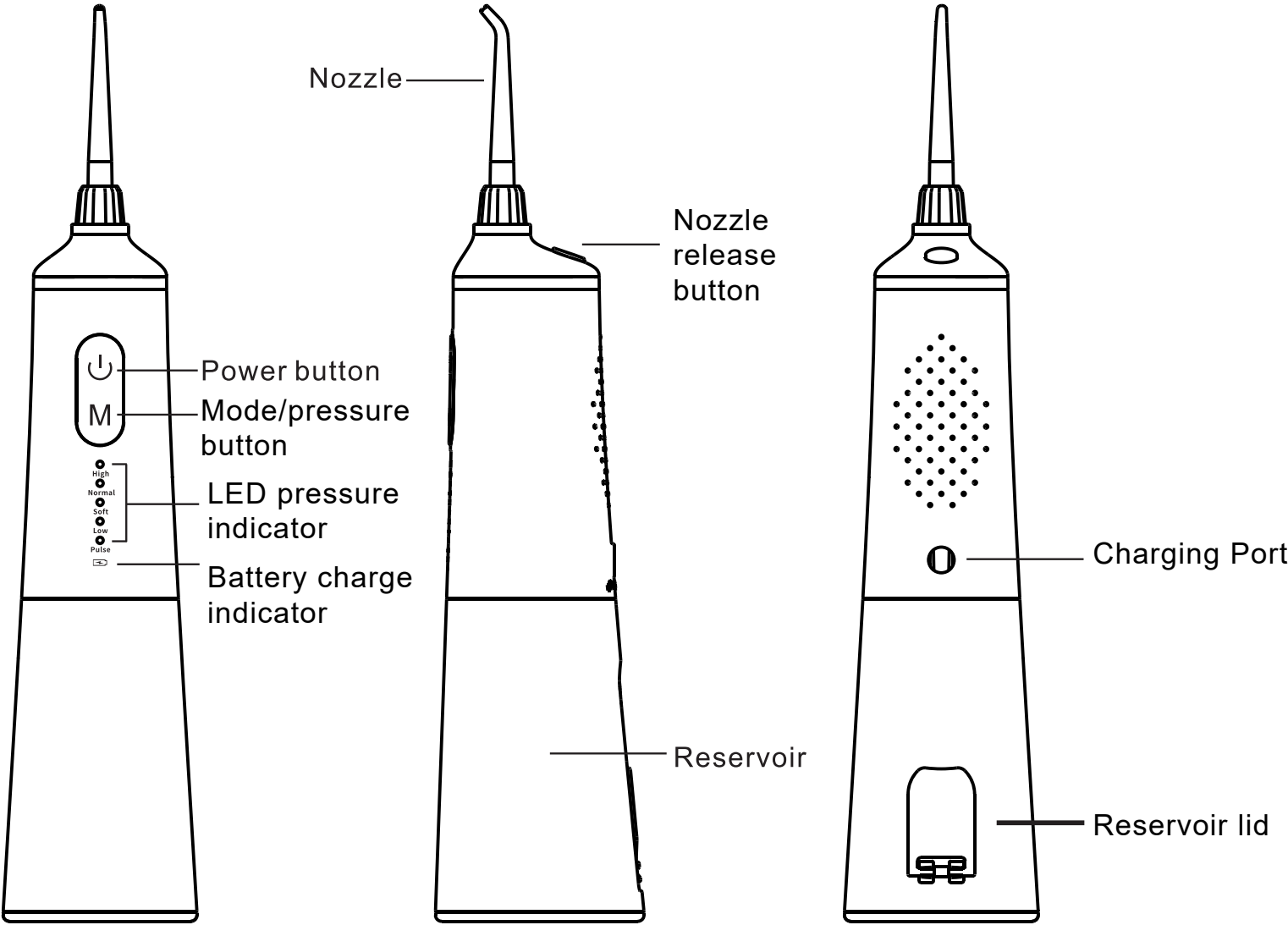

Accessories

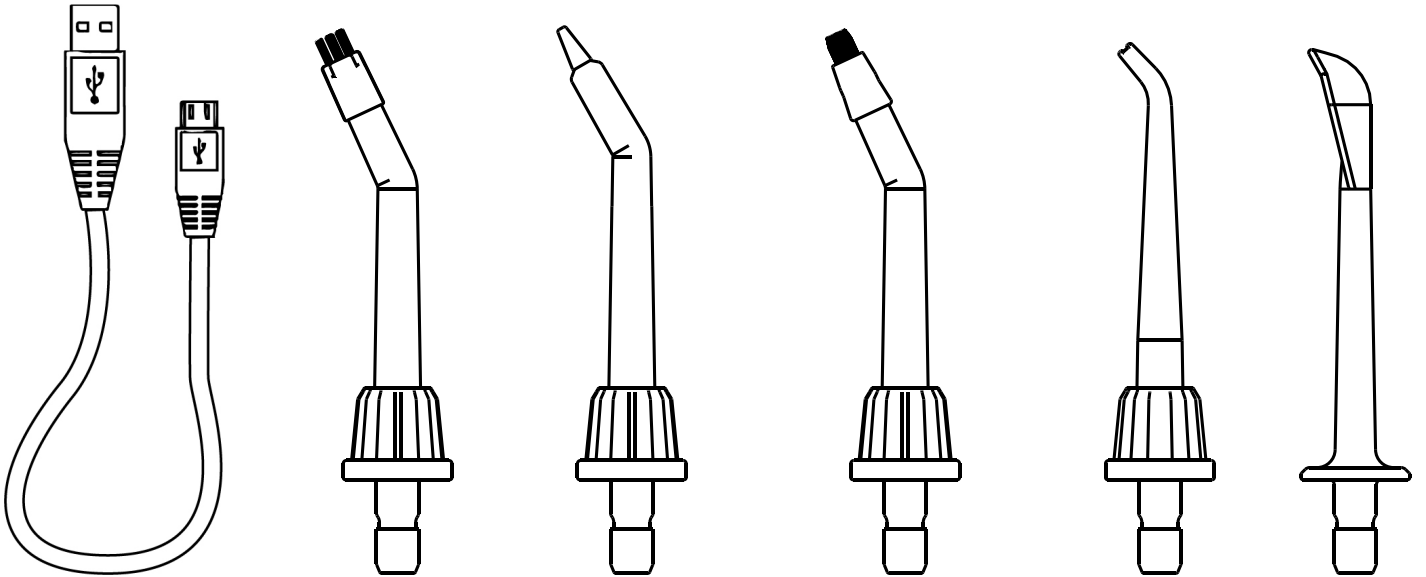

USB Charging Cable

Plaque Removal Nozzle

GumCare Nozzle

Floss n' Braces Nozzle

MultiClean Nozzle

Tongue Refresh Nozzle


Product Parameters

Product model	Remi Cordless Water Flosser / H300
Product color	Silver/Black
Working modes	High/Normal/Soft/Low/Pulse
Pulse	1200-1400 times/minute
Reservoir capacity	240ml/300ml
Charging time	About 4 hours
Adapter	DC 5V 1A
Waterproof level	IPX7
Standard	GB4706.59-2008 GB4706.1-2005


Instructions

The Remi cordless water flosser needs to be fully charged before use. Please charge until battery charge indicator changes from red or blue indicator to green light for a full charge.

1. Nozzle installation: Firmly hold base and push nozzle downward until the nozzle clicks securely into place. To remove nozzle, firmly hold base then press nozzle release button and pull nozzle away from base.

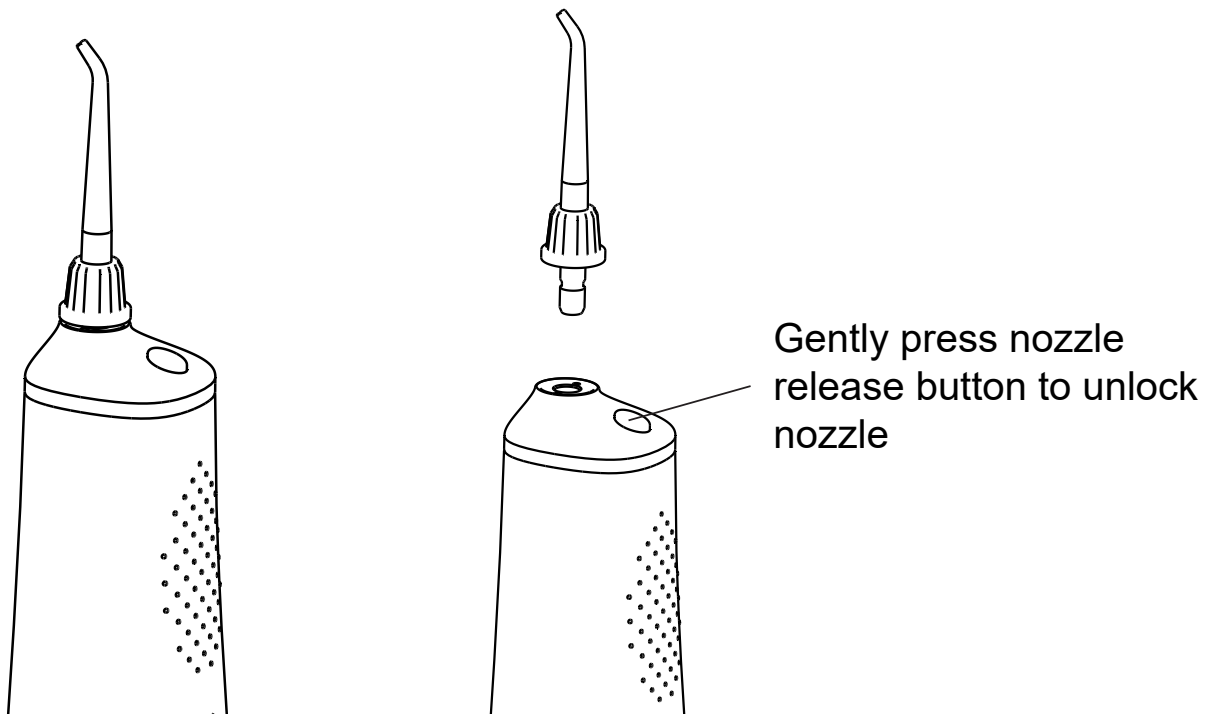

2. Open the reservoir cover and hold the flosser horizontally to fill the water tank. Cover the tank after filling the appropriate amount of water.
(Water should not exceed 104°F).

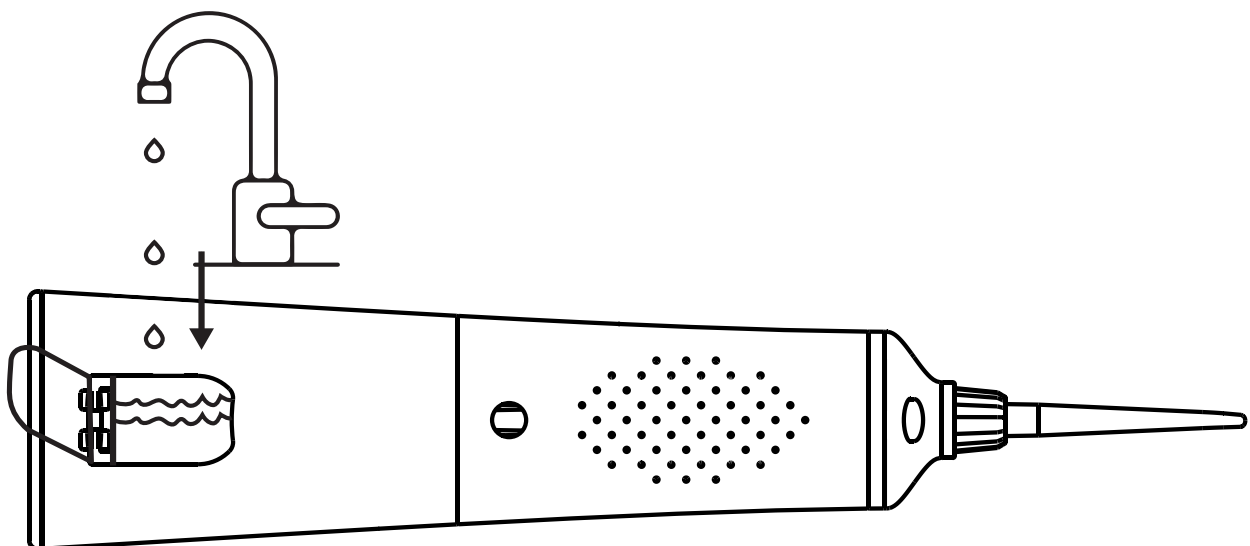

3. Press the mode button to change the level settings mode after the Remi water flosser is turned on. The mode level options are:

High Mode: Highest water pressure, perfect for deep cleaning and people who floss regularly and have excellent gum health.

Normal Mode: Strong water pressure, perfect for efficient cleaning and those who routinely use a water flosser to maintain good oral health.

Soft Mode: Gentle pressure for first time users and those with somewhat sensitive teeth and gums.

Low Mode: Very gentle pressure for first time users and those with very sensitive teeth and gums.

Pulse Mode: Mid to high range pulsating frequency, ideal for pulsating gum massage, great for those who routinely use a water flosser.

- 4. Direct the water jet flow at a 90-degree angle to teeth and/or gumline. Direct the water nozzle slowly along the teeth and gums and gumline until finished.
- 5. Direct the water flow between teeth and around any braces, bridgework, or dental work to clean thoroughly.
- 6. The Remi water flosser has a two-minute timed shutdown function. If you need to continue to use the water flosser for a more thorough cleaning please press the power button again.

Warning: It is common to experience light bleeding or see pink in your sink when you begin to first use a water flosser. This usually stops within 1-2 weeks. If it has been longer than two weeks please stop use and consult with a dental professional.

Battery Indicator

Start up power Indicator	
90% - 100%	Green light
20% - 89%	Blue light
<20%	Red light
Charging power Indicator	
100%	Green light
20%-99%	Flash blue
<20%	Flash red

Cleaning and Maintenance

Clean the product with clean water or neutral detergent instead of corrosive or abrasive detergent. Do not use hot water over 120°F to avoid damage to the product.

- **Cleaning Exterior:** Dry cloth or wet cloth can be used to clean the body, do not immerse the water flosser in water to clean.
- **Cleaning of water tank:** Dry or wet cloth can be used to clean the tank. If you do not plan to use for an extended time, please remove all water and moisture from the inside.
- **Nozzle cleaning:** can be wiped with dry cloth or wet cloth, or remove with clean water rinse and then wipe with cloth, do not bend, pull and twist the nozzle pipe.
- **Cleaning of charging plug:** Please wipe the metal part of plug regularly with cloth to avoid bad contact caused by dust accumulation.

Troubleshooting

<u>Troubleshooting</u>	<u>Cause of failure</u>	<u>Solution</u>
Low water pressure	Nozzle blockage or damage / low battery voltage	Replace nozzle or charge fully
Charging indicator abnormal or not on	Circuit fault	Return for repair
	Power plug not connected correctly	Connect the power supply correctly
Short use time after full charge	Insufficient charging time	Fully Charge
	Circuit fault	Return for repair

Important Notice

In order to reduce the risk of failure, electric shock, injury, fire, death and equipment damage, please read carefully all the instructions and warning tips in this manual before using this product.

Warning

1. Do not operate this product any way except its intended use.
2. Do not disassemble or change the parts of this product in order to avoid damage or danger.
3. Do not put hot water above 104°F into the water tank and clean the mouth. to avoid scalding or burning while using the device.
4. If the water flosser falls into the water when charging, please unplug the power plug before pulling out the water flosser.
5. If the power plug or power cord is damaged, please replace it with a full functioning cord before using the water flosser again.

Prohibited use:


1. Please do not use this product for other purposes than oral cleaning.
2. Do not press the nozzle hard on the gums or teeth to avoid injury.
3. Please do not immerse the product in water for a long time to avoid failure.
4. Please do not use wet hands to plug and unplug power plugs and charging plugs.
5. Do not damage, destroy or modify the USB power cord. Do not stretch bend or twist the USB power cord. Do not press heavy objects on the USB power cord or stick the USB power cord between objects.

Important

1. Please do not use this product for other purposes than oral cleaning.
2. Do not press the nozzle hard on the gums or teeth to avoid injury.
3. Please do not immerse the product in water for a long time to avoid failure.
4. Please do not use wet hands to plug and unplug power plugs and charging plugs.
5. Do not damage, destroy or modify the USB power cord. Do not stretch, bend or twist power cord between objects.

Attention

1. Do not use this product without nozzles installed.
2. Please do not put pure mouthwash into this product. If you need to use mouthwash, dilute it proportionally before use.
3. Persons (including children) who lack relevant experience and knowledge, such as those who do not have full civil capacity must use, clean and maintain the product under the supervision and guidance of their guardians.
4. Please put the product out of the reach of children and infants. The water flosser should not be used as a toy.
5. For patients with periodontal disease or those who have had oral surgery in the past 2 months, please consult a doctor before using it.
6. When plugging in and out of the power supply, be sure to hold the plug instead of cord.
7. Do not leave the water flosser plugged in for an extended period of time. If you do not plan to use the water flosser in the near future, unplug the device.

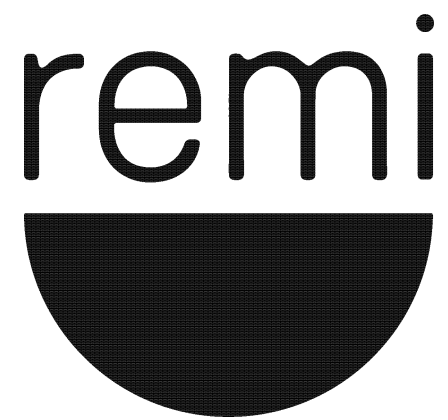

support@shopremi.com



1800 234 1874