# Orthomechanik 2.0

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# Orthomechanik 2.0

## User Manual

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

# CONTENT

<b>Structure and settings of the OM 2.0</b>	<b>3</b>
<b>First steps with your OM 2.0</b>	<b>4</b>
<b>Massage heads and their use</b>	<b>5-6</b>
<b>Our Orthomechanik-App</b>	<b>7</b>
<b>What do you get with your OM 2.0?</b>	<b>8</b>
<b>Technical specifications</b>	<b>8</b>
<b>Cleaning your OM 2.0 &amp; accessories</b>	<b>9</b>
<b>Warnings</b>	<b>10</b>
<b>Warnings for use</b>	
<b>Warnings - Handling the device itself</b>	<b>11</b>
<b>Your Warranty</b>	<b>12</b>
<b>Disposal</b>	<b>13</b>


**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Structure and settings of the OM 2.0
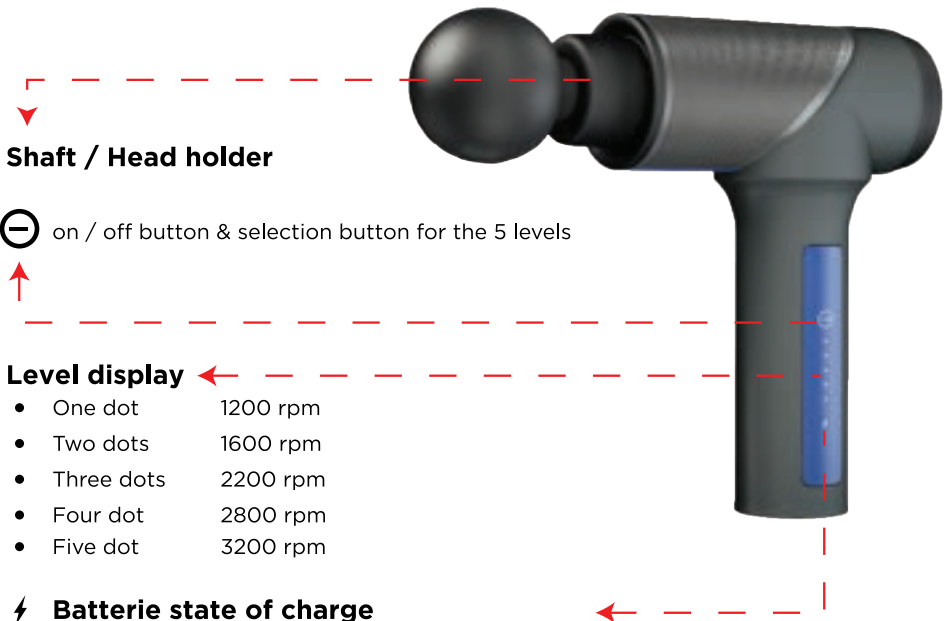
### Shaft / Head holder

⊖ on / off button & selection button for the 5 levels

### Level display

- One dot      1200 rpm
- Two dots    1600 rpm
- Three dots   2200 rpm
- Four dot     2800 rpm
- Five dot     3200 rpm

### ⚡ Batterie state of charge

- Green      100 to 70 % (2-6 hours left)
- Yellow     70 - 30% (1-2 hours left)
- Red        30 - 0% (0-1 hour left)

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## First steps with your OM 2.0

1. Fully charge your OM 2.0 before using it for the first time (approx. 5h).  
**ATTENTION:** DO NOT switch on the OM 2.0 during the charging process.
2. While you are waiting, download our OrthoApp (see QR code on flyer).
3. Select a massage head and insert it firmly into the Gun.
4. Long press and hold the (⊖) on / off button until the light comes on and the Gun starts pulsating.
5. Move the OM 2.0 slowly and not longer than 5 minutes over the respective muscle group.
6. By short pressing the (⊖) on / off button again, you switch up one level at a time (level 1-2-3-4-5-idle-1-...).
7. If you want to end the massage, long press the (⊖) on / off button again. The vibration stops and the light goes off.

**CAUTION: Make sure that no oils, creams or lotions come into direct contact with the massage gun! This may cause the paint to peel off. Any oils or creams should ONLY be used with the metal attachments.**

**CAUTION: Read our warnings on page 10 before use!**

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

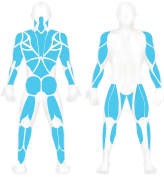
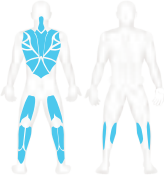
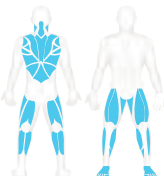
## Massage heads and their use

We recommend to use your OM 2.0 as following:

- use high frequencies for warm up before sports
- lower frequencies for relaxation

If it feels comfortable, you can use any head on any muscle. Nevertheless we recommend using our heads in the following combination:

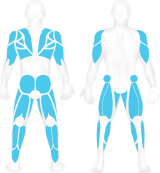
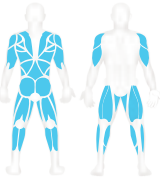
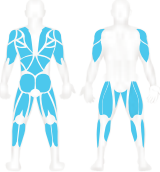
	Head	Use	Especially good for	Body part
<b>Ball Head</b>		For large muscle groups  A first mobilization for painful muscles  With high frequency for the warm-up.	upper arms, upper thigh, back, neck, shoulder muscles, forearms, gluteal muscles	
<b>Fork Head</b>		For areas where the center is to be spared  E.g. muscles next to the spine or in the area of tendons or muscle insertions	muscles around the spine, transition calf, achilles tendon, shoulder blade	
<b>Cone Head</b>		For precise & deep treatment of smaller muscle areas and trigger points	back muscles, thighs, shoulders, calves, feet	

**NEED HELP?**  
 Visit our website for videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

	Head	Use	Especially good for	Body part
<b>Flat Head Metal</b>		<p>Suitable for use on all muscle groups</p> <p>For general use with oil</p> <p>Easy to clean</p> <p>For application by several persons</p>	<p>glutes, upper arms, shoulders, calves, thighs, back</p>	
<b>Round Head Metal</b>		<p>Suitable for use on all muscle groups</p> <p>For general use with oil</p> <p>Easy to clean</p> <p>For application by several persons</p>	<p>thighs, glutes, upper arms, shoulders, calves</p>	
<b>Flat Head</b>		<p>For general use with large muscle groups</p>	<p>glutes, upper arms, shoulders, calves, thighs</p>	
<b>Soft Head</b>		<p>To approach carefully painful muscles</p> <p>In the area of palpable bone points or flat</p>	<p>neck, cervical spine, anterior tibial muscles, forearms, shoulder blades</p>	


**NEED HELP?**  
Visit our website for videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Our Orthomechanik-App

Our specialists will guide you in short videos explaining the most efficient treatment of each muscle group. In addition, we cover the most common ailments in each area and show you suitable massage techniques.

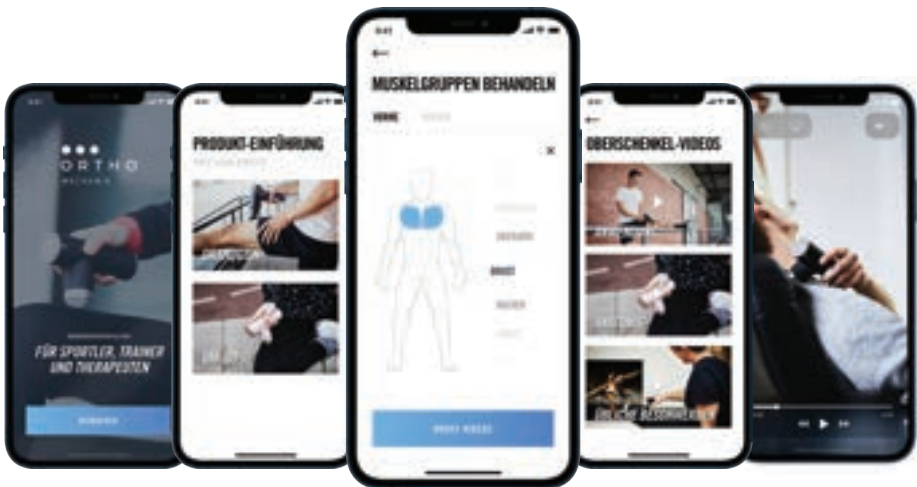And here you can  
Download the app

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**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## What do you get with your OM 2.0?

	1 x	Orthomechanik 2.0 massage gun
	1 x	Charger
	7 x	Massage heads
	1 x	Case
	1 x	Instruction manual
	1 x	Instructions for the app

## Technical specifications

I am Model	OrthoGun 2.0 mini
My battery power	3400 mAh / 42 Wh
What my motor can do	15 - 160 W
My frequencies	Level 1: 1200 U/Min Level 2: 1600 U/ Min Level 3: 2200 U/Min Level 4: 2800 U/Min Level 5: 3200 U/Min
My torque	0 ~ 7,5 kg/cm
I automatically switch off	After 10 min

**NEED HELP?**  
 Visit our website for  
 videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## **Cleaning your OM 2.0 & accessories**

### **OM 2.0 - Massage gun**

- Wipe with a dry to slightly damp cloth
- Then immediately dry thoroughly with a dry cloth

### **Attachments**

- Wipe with a damp cloth
- Disinfectants can be used
- Then dry thoroughly with a dry cloth

### **ATTENTION**

- Do not allow liquids to enter the unit or the shaft
- Do not use corrosive or sharp Detergents

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Warnings

### Warnings for use

- Use only on dry and clean body surfaces
- Best to use over clothing
- Use only on soft tissues of the body
- Do not use on head, bones or sensitive body parts (eyes, face, etc.)
- Keep fingers, hair and other body parts away from the shaft
- Pressure marks and slight bruising may occur even at low levels and with little pressure. Observe the area and talk to your doctor / physiotherapist if necessary.
- Oils, creams and liquids may only be used with the metal attachments.
- Make sure that oils, creams and liquids never come into contact with the gun or enter the shaft.

**If you have health complaints or / and belong to the following groups of people, you should talk to your attending physician or physiotherapist before treatment with the OM 2.0:**

- Pregnant women
- People with diabetes or neuropathies
- People with chronic diseases
- Wearers of pacemakers
- Wearers of artificial joints or metal plates
- If you have circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis

**Children under 8 years of age and persons with limited physical or mental abilities should only use the device under supervision.**

**CAUTION:** The massage should not feel uncomfortable at any time! If this is the case, stop the treatment and talk to your doctor!

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Warnings

### Handling the device itself

- Only use the device with the supplied power supply unit.
- Keep the OM 2.0 away from heat and liquids.
- Keep the ventilation openings free of dust and dirt.
- Do not remove any screws, do not try to open the device.
- Do not use the device for more than 1 hour at a time, let it cool down for 30 minutes before using the OM 2.0 again.
- Never use the OM 2.0 while it is charging. Disconnect the Gun from the charger before turning it on.
- Only charge the unit with the supplied charger.

**IMPORTANT:** DO NOT attempt to repair the OM 2.0 yourself if it is defective, as this will void your warranty. Instead, send us an email at [info@orthomechanik.de](mailto:info@orthomechanik.de)

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Your warranty

For the OM 2.0, OrthoMechanik GmbH provides a 30-day money-back guarantee from the date of purchase as well as a 2-year warranty on parts and workmanship.

These warranties are non-transferable and apply only to the purchaser.

Orthomechanik GmbH is entitled to request pictures, videos or other materials showing the defect.

All warranty claims must be reported to [info@orthomechanik.de](mailto:info@orthomechanik.de) immediately after the occurrence of the defect. The e-mail must contain the full name, as well as the order or invoice number.

Replacement parts and replacement products are subject to the remainder of the original warranty period or a thirty (30) day warranty period, whichever is longer.

Warranty claims CANNOT be made if the defect was caused by improper use.

NOTE: Use of unapproved third party products (e.g. charger) will immediately void any warranty claim. Third-party parts may cause the product to malfunction or not to function at all and may cause damage to the product and connected devices (e.g. charger).

NOTE: This limited warranty does not affect your statutory rights under the laws governing the sale of consumer goods.

**NEED HELP?**  
Visit our website for  
videos, blog articles and Q&A  
[www.orthomechanik.de](http://www.orthomechanik.de)

**CONTACT US**  
[info@orthomechanik.de](mailto:info@orthomechanik.de)

## Information for private households on the disposal of electronic equipment.

**Due to the 2018 regulations on the disposal of electronic equipment (WEEE), you are subject to the following regulations:**

The crossed-out waste garbage can indicates that you are required by law to dispose of this equipment separately from unsorted municipal waste.

Disposal via the residual waste bin can or the yellow bin can is prohibited.

If the products contain batteries or rechargeable batteries that are not permanently installed, they must be removed before disposal and disposed of separately as batteries.

### Your options for returning old equipment

Old equipment can be handed in at the facilities set up and made available by the public waste disposal authorities. This ensures proper disposal of the old equipment.

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We would like to point out to all end users of waste electrical and electronic equipment that you are responsible for deleting personal data on the waste equipment to be disposed of.

We are registered as a producer of electrical and/ or electronic equipment under the registration number DE53771343 at the Stiftung Elektroaltgeräten, Nordostpark 72, 90411 Nürnberg.

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