Orthomechanik 3.0

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Orthomechanik 3.0

User Manual

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

CONTENT

Structure & controls of the OM 3.0	3
First steps with your OM 3.0	4
Massage heads and their use	5-6
Our Orthomechanik-App	7
What do you get with your OM 3.0?	8
Technical Specifications	8
Cleaning your OM 3.0 & accessorie	9
Warnings	10
Warnings for use	
Warnings - Handling the device itself	11
Your Warranty	12
Disposal	13

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Structure & controls of the OM 3.0

Shaft / Head holder

⊖ On / off button & button to upshift

⊖ Button to downshift

Level display

- One dot 1200 rpm
- Two dots 1600 rpm
- Three dots 2200 rpm
- Four dot 2800 rpm
- Five dot 3200 rpm

⚡ Batterie state of charge

- Green 100 to 70 % (2-6 hours left)
- Yellow 70 - 30% (1-2 hours left)
- Red 30 - 0% (0-1 hour left)

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

First steps with your OM 3.0

1. Fully charge your OM 3.0 before using it for the first time (approx. 5h).
ATTENTION: DO NOT switch on the OM 3.0 during the charging process.
2. While you are waiting, download our OrthoApp (see QR code on flyer).
3. Select a massage head and insert it firmly into the Gun.
4. Long press and hold the (⏻) on / off button until the light comes on and the Gun starts pulsating.
5. Move the OM 3.0 slowly and not longer than 5 minutes over the respective muscle group.
6. By short pressing the top button again, you upshift one level and if you short press the lower button, you downshift one level.
7. If you want to end the massage, long press the (⏻) on / off button again. The vibration stops and the light goes off.

CAUTION: Make sure that no oils, creams or lotions come into direct contact with the massage gun! This may cause the paint to peel off. Any oils or creams should ONLY be used with the metal attachments.

CAUTION: Read our warnings on page 10 before use!

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

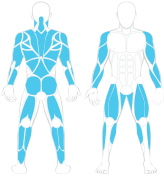
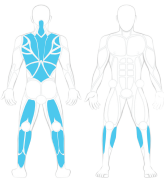
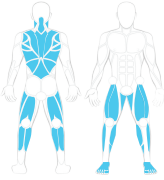
Massage heads and their use

We recommend to use your OM 3,0 as following:

- use high frequencies for warm up before sports
- lower frequencies for relaxation

If it feels comfortable, you can use any head on any muscle. Nevertheless we recommend using our heads in the following combination:

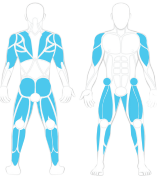
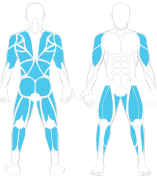
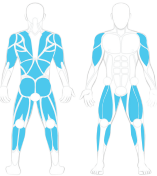
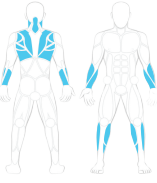
	Head	Use	Especially good for	Body part
Ball Head		For large muscle groups A first mobilization for painful muscles With high frequency for the warm-up.	upper arms, upper thigh, back, neck, shoulder muscles, forearms, gluteal muscles	
Fork Head		For areas where the center is to be spared E.g. muscles next to the spine or in the area of tendons or muscle insertions	muscles around the spine, transition calf, achilles tendon, shoulder blade	
Cone Head		For precise & deep treatment of smaller muscle areas and trigger points	back muscles, thighs, shoulders, calves, feet	

DO YOU NEED HELP?
 For Videos, Blog posts and Q&A
 visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

	Head	Use	Especially good for	Body part
Flat Head Metal		<p>Suitable for use on all muscle groups</p> <p>For general use with oil</p> <p>Easy to clean</p> <p>For application by several persons</p>	glutes, upper arms, shoulders, calves, thighs, back	
Round Head Metal		<p>Suitable for use on all muscle groups</p> <p>For general use with oil</p> <p>Easy to clean</p> <p>For application by several persons</p>	thighs, glutes, upper arms, shoulders, calves	
Flat Head		<p>For general use with large muscle groups</p>	glutes, upper arms, shoulders, calves, thighs	
Soft Head		<p>To approach carefully painful muscles</p> <p>In the area of palpable bone points or flat</p>	neck, cervical spine, anterior tibial muscles, forearms, shoulder blades	


DO YOU NEED HELP?
 For Videos, Blog posts and Q&A
 visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Our Orthomechanik-App

Our Orthomechanics - App offers you help in treating your problem muscle group as well as general tips on how to use our Guns. Our specialists are guiding you in short videos and explaining the special treatments of individual muscle groups.

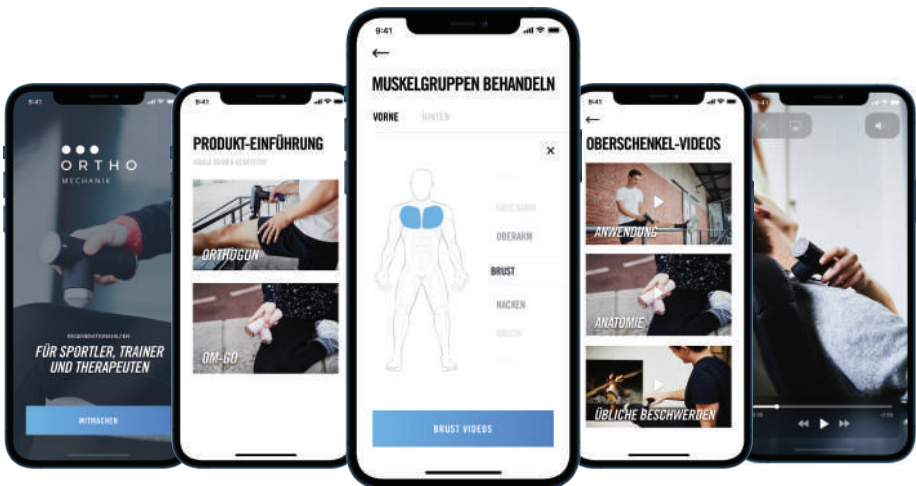And here you can
Download the app

[ORTHOMECHANIK.DE/APP](https://www.orthomechanik.de/app)

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

What do you get with your OM 3.0?

	1 x	Orthomechanik 3.0 massage gun
	1 x	Charger
	7 x	Massage heads
	1 x	Case
	1 x	Instruction manual
	1 x	Instructions for the app

Technical specifications

I am Model	Orthomechanik 3.0
My battery power	3350 mAh
What my motor can do	15 - 160 W
My frequencies	Level 1: 1200 U/Min Level 2: 1600 U/ Min Level 3: 2200 U/Min Level 4: 2800 U/Min Level 5: 3200 U/Min
My torque	0 ~ 7,5 kg/cm
I automatically switch off	After 10 min

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Cleaning your OM 3.0 & accessories

OM 3.0 - Massage gun

- Wipe with a dry to slightly damp cloth
- Then immediately dry thoroughly with a dry cloth

Attachments

- Wipe with a damp cloth
- Disinfectants can be used
- Then dry thoroughly with a dry cloth

ATTENTION

- Do not allow liquids to enter the unit or the shaft
- Do not use corrosive or sharp Detergents

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Warnings

Warnings for use

- Use only on dry and clean body surfaces
- Best to use over clothing
- Use only on soft tissues of the body
- Do not use on head, bones or sensitive body parts (eyes, face, etc.)
- Keep fingers, hair and other body parts away from the shaft
- Pressure marks and slight bruising may occur even at low levels and with little pressure. Observe the area and talk to your doctor / physiotherapist if necessary.
- Oils, creams and liquids may only be used with the metal attachments.
- Make sure that oils, creams and liquids never come into contact with the gun or enter the shaft.

If you have health complaints or / and belong to the following groups of people, you should talk to your attending physician or physiotherapist before treatment with the OM 3.0:

- Pregnant women
- People with diabetes or neuropathies
- People with chronic diseases
- Wearers of pacemakers
- Wearers of artificial joints or metal plates
- If you have circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis

Children under 8 years of age and persons with limited physical or mental abilities should only use the device under supervision.

CAUTION: The massage should not feel uncomfortable at any time! If this is the case, stop the treatment and talk to your doctor!

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Warnings

Handling the device itself

- Only use the device with the supplied power supply unit.
- Keep the OM 3.0 away from heat and liquids.
- Keep the ventilation openings free of dust and dirt.
- Do not remove any screws, do not try to open the device.
- Do not use the device for more than 1 hour at a time, let it cool down for 30 minutes before using the OM 3.0 again.
- Never use the OM 3.0 while it is charging. Disconnect the Gun from the charger before turning it on.
- Only charge the unit with the supplied charger.

IMPORTANT: DO NOT attempt to repair the OM 3.0 yourself if it is defective, as this will void your warranty. Instead, send us an email at info@orthomechanik.de

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Your warranty

For the OM 3.0, OrthoMechanik GmbH provides a 30-day money-back guarantee from the date of purchase as well as a 2-year warranty on parts and workmanship.

These warranties are non-transferable and apply only to the purchaser.

Orthomechanik GmbH is entitled to request pictures, videos or other materials showing the defect.

All warranty claims must be reported to info@orthomechanik.de immediately after the occurrence of the defect. The e-mail must contain the full name, as well as the order or invoice number.

Replacement parts and replacement products are subject to the remainder of the original warranty period or a thirty (30) day warranty period, whichever is longer.

Warranty claims CANNOT be made if the defect was caused by improper use.

NOTE: Use of unapproved third party products (e.g. charger) will immediately void any warranty claim. Third-party parts may cause the product to malfunction or not to function at all and may cause damage to the product and connected devices (e.g. charger).

NOTE: This limited warranty does not affect your statutory rights under the laws governing the sale of consumer goods.

DO YOU NEED HELP?
For Videos, Blog posts and Q&A
visit our Homepage
www.orthomechanik.de

CONTACT US
info@orthomechanik.de

Information for private households on the disposal of electronic equipment.

Due to the 2018 regulations on the disposal of electronic equipment (WEEE), you are subject to the following regulations:

The crossed-out wheeled garbage can indicates that you are required by law to send these devices for separate collection from unsorted municipal waste.

Disposal via the residual waste garbage can or the yellow garbage can is prohibited.

If the products contain batteries or rechargeable batteries that are not permanently installed, these must be removed before disposal and disposed of separately as batteries.

Your options for returning old equipment

Old devices can be returned to the facilities set up and made available by public waste disposal authorities. This ensures proper disposal of the old devices.

The following link provides an online directory of collection and return points:
<https://www.ear-system.de/ear-verzeichnis/sammel-und-ruecknahmestellen>.

We would like to point out to all end users of waste electrical and electronic equipment that you are responsible for deleting personal data on the waste equipment to be disposed of.

We are registered as a producer of electrical and/ or electronic equipment under the registration number DE53771343 at the Stiftung Elektroaltgeräten, Nordostpark 72, 90411 Nürnberg.

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