
# Infrared Heat Lamp
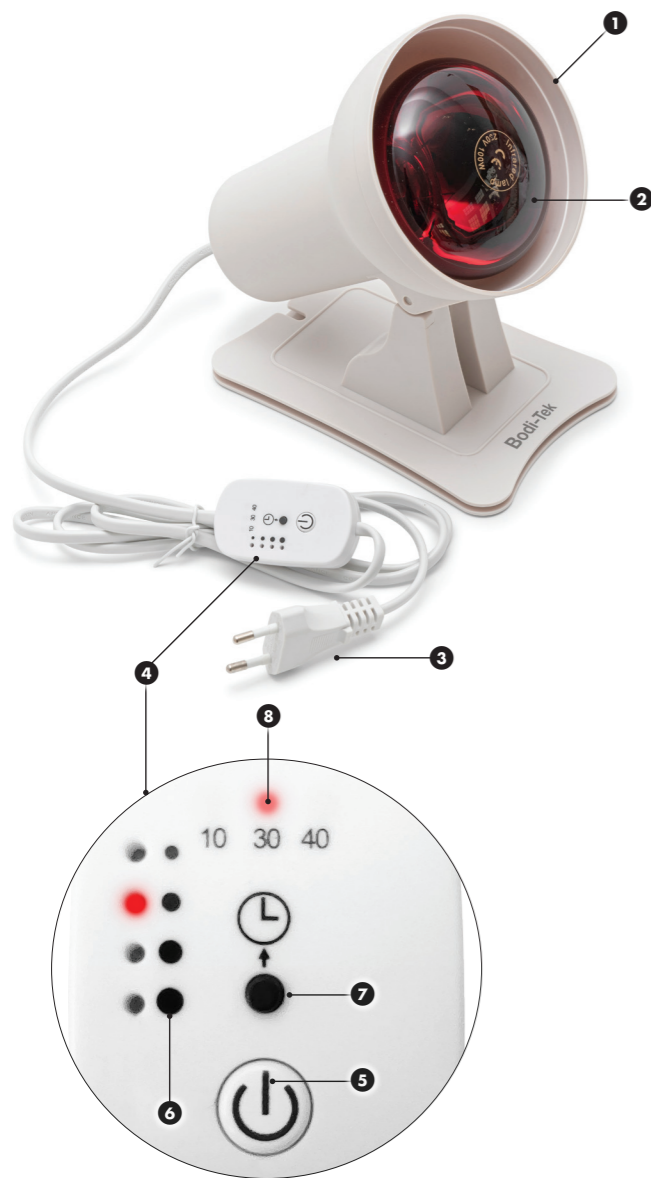
## ⚠ Warnings

Read all instructions thoroughly and retain for future reference.

- Keep out of reach from children and pets.
- This appliance can be used by adults aged 18 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do not use in the bath or shower.
- Do NOT use this appliance if it is not working properly, or if it has been dropped or damaged.
- If the power lead is damaged, it must be replaced by the manufacturer or its service agent or similarly qualified person in order to avoid a hazard. Contact Customer Care to arrange a replacement.
- Do not leave the device unattended. Turn off before you leave the room.
- This device can get very hot so avoid contact with the skin during use.
- Do not allow skin to overheat and keep a safe distance between the product and the skin. Stop use if it feels uncomfortable.
- Do not cover the appliance with clothes or any other item, as this is a fire hazard.
- This device must only be used for its intended purpose as described in this user guide.
- Do not operate the appliance without the bulb fitted or if it is damaged.
- Do not use on damaged or irritated skin.
- Do not use on the eyes or surrounding areas.

Consult your doctor before use if you have any of the following:

- If you suffer from a skin condition.
- If you are insensitive to heat due to a disease (e.g. advanced diabetes) or the use of painkillers.
- If you suffer from conditions associated with nerve demyelination such as carpal tunnel syndrome, multiple sclerosis and ulnar neuropathy.
- If you have an implanted device.
- If you are pregnant.

⚠ "Hot Surface" do not touch until cool.

⚠ CAUTION: IR emitted from this device. Do NOT stare directly at the light source.

## Know Your Unit


- |                        |                         |
|------------------------|-------------------------|
| 1 Lamp unit            | 5 ON/OFF and Heat Level |
| 2 Infrared bulb        | 6 Heat levels           |
| 3 Power cable and plug | 7 Timer                 |
| 4 Control unit         | 8 Timer indicator       |

## Introduction

Infrared light penetrates deeply into the skin delivering warming soothing heat to the underlying muscles and tissue.

## How to use


- Place the lamp on a dry and smooth heat resistant surface which is at least 30cm away from other surfaces.

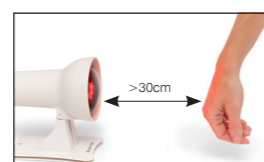1. Plug the infrared lamp into the mains supply.
2. Place the lamp in front of the area to be treated and no closer than 30cm from the skin.

3. Press the ON/OFF Heat Level button once to turn the lamp on. The lowest heat and timer setting (10 minutes) will automatically be selected and indicated by a small red light.

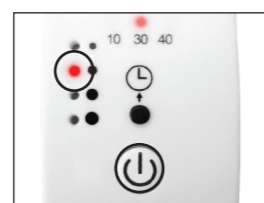
4. Press the ON/OFF and Heat Level again to increase the heat level – there are four heat levels. A small red light will indicate the heat level selected.

5. Press the Timer button to change the timer setting from 10 minutes to 30 minutes. Press again to increase from 30 minutes to 40 minutes. The lamp will automatically switch off at the end of the selected time.

- Do not use on one area of skin for more than five minutes at a time.
- If you notice any burning, swelling or discomfort, discontinue use immediately.

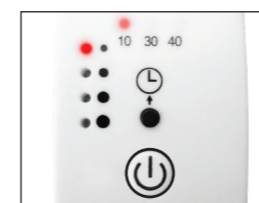6. To turn off at any point, press the ON/OFF Heat Level button repeatedly until the lamp turns off (depending on your chosen heat level, this will be a maximum of four presses).

- The infrared bulb will become very hot, take care not to touch the bulb.
- With correct use some warmth and redness of the skin is normal. This should subside quickly after use.

## After Use

Unplug from the mains supply and allow to cool for 20 minutes before handling or storing

## Cleaning & Maintenance

Do not immerse the unit in water or use any harsh cleaning products. Wipe clean with a clean, damp cloth.

## Important Information

- Do NOT stare directly at the light source.
- Replace with E27 150W Infra-red bulb

## Disposal

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. This product is classed as Electrical or Electronic Equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

## Specifications

Model: BT-LAMH  
Input: 230V ~ 50 Hz 150 W

## Guarantee

To register your two year guarantee online, visit: [www.bodi-tek.com/warranty](http://www.bodi-tek.com/warranty)

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©The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, England GL53 7ET
©The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, Inghilterra GL53 7ET
**RP prodotti nel Regno Unito:** The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England, GL53 7ET
**RP prodotti in UE:** The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands