

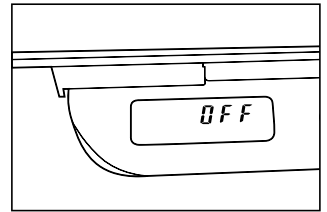
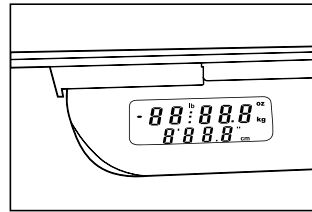
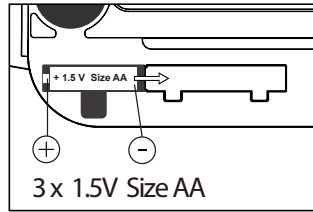



Installation

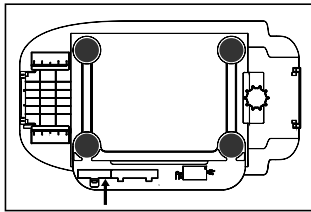

Turn on

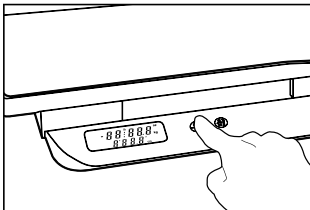


OR

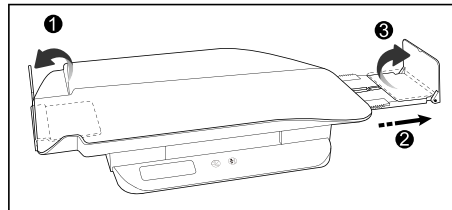


Auto-turn on only works when ruler is fully extended

Turn off

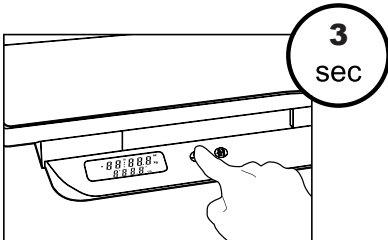


OR

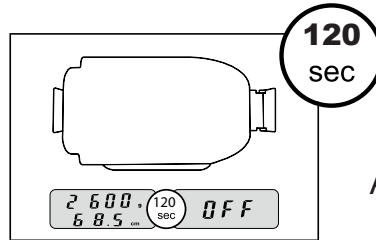


Auto-turn off after 2 minutes


Measure weight + height

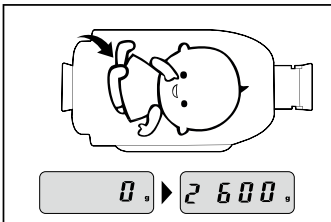

Put baby on scale to measure weight

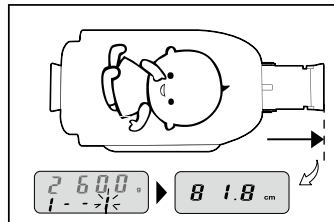


Fully extend the digital ruler

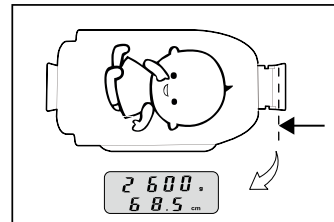


Extend the flap and push in the ruler until it touches your baby's head


Measure height only

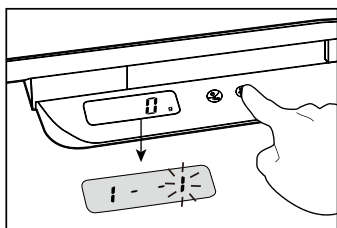

Press once

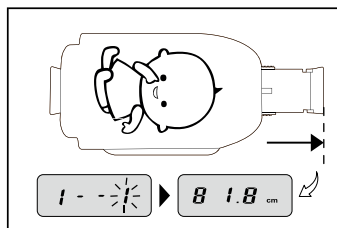


Fully extend the digital ruler


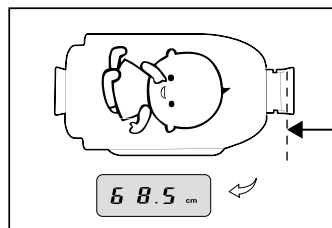


Extend the flap and push in the ruler until it touches your baby's head


Tare (set to 0 even with cloth on the scale)

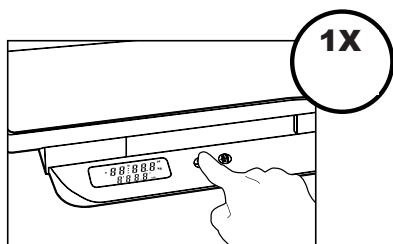


Switch between grams and ounce

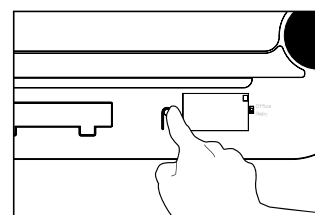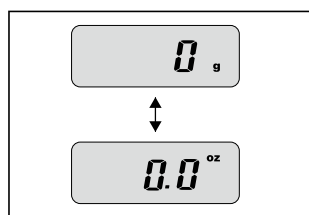

Messages

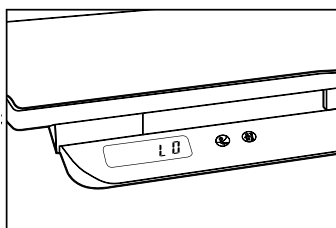

Low battery
Replace the batteries

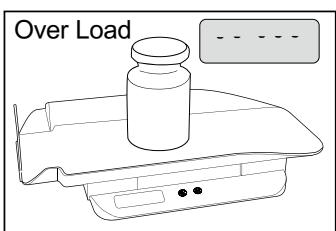

Overload
Maximum 20 kg can be measured

SPECIFICATIONS

Power supply:..... 3x 1,5Volt, size AA
 Battery lifetime:..... > 1 year
 Maximum load:..... 20 Kg
 Gradation:..... 5 grams
 Maximum length: 81,8 cm
 Length gradation: 0,2 cm
 Work temperature:..... 10 - 40°C

Alecto[®]

Service



↓

Help





WWW.ALECTO.NL
SERVICE@ALECTO.NL

 **Hesdo, Australiëlaan 1**
5232 BB, 's-Hertogenbosch
The Netherlands