
medisana

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GB Instruction manual

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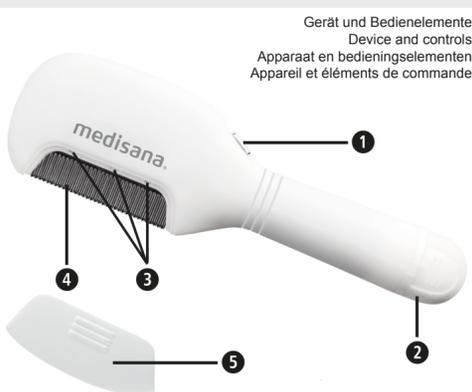GB Safety instructions
Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.
Do not give this device to another person. It is vital that you also pass on these instructions for use.

- Use the unit only as intended in accordance with the directions in the instruction manual.
- Using the unit for any other purpose will void your guarantee.
- This unit is not designed for commercial use.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision and instruction concerning its safe use and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance.
- Cleaning and user maintenance shall not be made by children younger than that specified on words provided.
- The unit is not suited for children under the age of 3.
- If you have a pacemaker, artificial limbs or an electronic implant, please consult your physician before using the unit.
- The unit is not recommended for persons who suffer from epilepsy.
- Use the unit only on dry and combed hair.
- Do not use the unit on dry, cracked or wounded skin.
- When using the electric lice comb, avoid direct contact with your ears, eyes and mouth. Touching the metal teeth can inflict a painful, yet harmless, electrical charge.**

- The comb has rounded pearl-shaped tips on every second tooth. Make sure that no teeth are missing or damaged before you use the unit. Otherwise the lice comb is no longer safe to use.
- In the event of a malfunction, do not attempt to repair the unit on your own; this will void your guarantee. Have all repairs carried out by an authorised service centre.
- Do not use the unit on dry, cracked or wounded skin.
- When using the electric lice comb, avoid direct contact with your ears, eyes and mouth. Touching the metal teeth can inflict a painful, yet harmless, electrical charge.**

- Regular cleaning and good maintenance ensures optimum results and extends the life of your unit.
- Clean the lice comb carefully after each use, particularly if it is used by more than one person.
- Remove the batteries before you clean the unit.
- To clean the housing, use a dry, soft paper or fabric cloth. To clean the comb teeth, use the included cleaning brush.
- Do not use any strong detergents for cleaning, e.g. scouring cream or other cleaning agents, as they may damage the surface of the unit.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Do not expose the unit to direct sunlight and protect it from dirt and moisture.
- Remove the batteries from the unit if you do not plan on using it again anytime in the near future. Otherwise, the batteries may start leaking.
- Store the unit in a clean, dry and cool location.

- Protect the unit from moisture. Should liquid get into the unit, remove the batteries immediately and refrain from using the unit. In such a case, contact your dealer or notify address below.
- Our address and phone number are located on the address page.
- Do not use the unit if the housing or comb teeth are damaged, if the unit does not function properly or if it has fallen down or fallen into water.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Do not expose the unit to direct sunlight and protect it from dirt and moisture.
- Remove the batteries from the unit if you do not plan on using it again anytime in the near future. Otherwise, the batteries may start leaking.
- Store the unit in a clean, dry and cool location.

- This product must not be disposed of together with household waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.
- Please remove the battery before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

Items supplied and packaging
Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the it and contact your dealer or your service centre.

- The following parts are included:

- 1 medisana electric lice comb LC 870 with cleaning brush
- 3 batteries 1.5 V (type AAA)
- 1 instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your dealer without delay.

! WARNING Please ensure that polythene packing is kept away from the reach of children! Risk of suffocation!

What are head lice?
Head lice are small, wingless insects that live off blood on human scalps. Fully adult lice are typically 2 - 3 mm in length. Their colour varies from off-white to brown. The lice eggs are beige-brown in colour and stick to the hairline on your scalp. Head lice hatch after 7 to 10 days; the white shells they leave behind are known as nits. An average infestation is about 10 lice. If left untreated, the number of lice on your scalp can quickly increase to 200. Lice can neither jump nor fly; this is why they can only be transmitted through direct head-to-head contact. However, it is also possible for lice to survive on textiles for a short period of time. Contrary to popular belief, head lice can inhabit all types of hair – not just unwashed or dirty hair.

Special features of the medisana electric lice comb LC 870
The medisana electric lice comb LC 870 effectively eliminates head lice without having to resort to chemical means. It works by creating an electric charge between the teeth that kills off the lice and their eggs. The electric charge is so small that they will not notice it; it is absolutely safe to use. Comb activity is indicated by an audible signal. The unit is powered by conventional AAA batteries. The medisana electric lice comb LC 870 is suitable for adults and children over the age of 3. The comb can also be applied to household pets, such as dogs and cats.

Operating
Inserting / removing the batteries
Remove the battery compartment lid and insert the included batteries (type AAA). Observe the proper polarity (+/-) of the batteries as indicated on the batteries itself and on the inside of the battery compartment. Replace the battery compartment lid. Follow the same procedure when replacing the batteries. We recommend that you do not use rechargeable batteries. Remove the batteries from the unit if you do not plan on using it again anytime in the near future.

! WARNING
BATTERY SAFETY INFORMATION

- Do not disassemble batteries!
- Clean the battery and device contacts if necessary before putting in the batteries!
- Remove discharged batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of this parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Only replace with a battery of the same type!
- Insert the battery correctly, observing the polarity!
- Keep the battery compartment well sealed!
- Remove the battery from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Keep unused batteries in their packaging away from metal objects in order to prevent short circuiting!

Operating
1. The electric lice comb LC 870 may only be used on dry hair. If necessary, detangle your hair and undo any knots or braids prior to using the comb.
2. To switch off the unit, slide the ON/OFF switch forward/upward. A high-pitched humming sound indicates that the unit is in operation. The LED light switches on and illuminates the area to be combed.
3. Position the comb at a 45° - 90° angle and guide it slowly through your hair.
4. While you are combing your hair, the hum will briefly stop whenever the teeth make contact with a louse or egg. If the humming sound does not return, the comb teeth are clogged. Switch the comb off. Use the small cleaning brush that is integrated in the end of the handle, to clean the teeth. Continue combing your hair. Allow enough time to thoroughly clean your hair.
5. Don't forget to switch off the unit after you are done and to clean it thoroughly. When not in use, always put the protective cap back on the unit.
6. Since the medisana electric lice comb LC 870 is particularly effective against adult lice, continue using the comb for about 2 weeks to ensure that you eliminate all of it.
7. To prevent the return of lice, we recommend that you use the electric lice comb once a week as a precautionary measure.

Helpful tips if you are dealing with head lice
You can supplement the application of your electric lice comb by using a chemical agent. For unusual cases, please consult a physician.
• Due to their possible side effects, chemical agents should only be used if you have an actual infestation of head lice and not as a means of prevention. A physician can verify this and prescribe treatment.
Cleaning and maintenance
Regular cleaning and good maintenance ensures optimum results and extends the life of your unit.
Clean the lice comb carefully after each use, particularly if it is used by more than one person.
Remove the batteries before you clean the unit.
To clean the housing, use a dry, soft paper or fabric cloth. To clean the comb teeth, use the included cleaning brush.
Do not use any strong detergents for cleaning, e.g. scouring cream or other cleaning agents, as they may damage the surface of the unit.
When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
Do not expose the unit to direct sunlight and protect it from dirt and moisture.
Remove the batteries from the unit if you do not plan on using it again anytime in the near future. Otherwise, the batteries may start leaking.
Store the unit in a clean, dry and cool location.

Disposal
This product must not be disposed of together with household waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.
Please remove the battery before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

Technical specifications
Name and model: medisana electric lice comb LC 870
Power supply: Dimensions L x W x H: Weight: Storage conditions: EAN number:
4,5Vmax. batteries (type AAA)
193 x 66 x 24 mm
approx. 91 g incl. batteries
dry and clean
41019
40 15588 41019 2

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.
The current version of this instruction manual can be found under www.medisana.com

Warranty and repair terms
Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.
The following warranty terms apply:

- The warranty period for medisana products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear (batteries etc.).
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH
Carl-Schurz-Str. 2, 41460 NEUSS
GERMANY

You will find the service address on the separate enclosure.

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