## VENTURE HEAT 12 MONTH LIMITED WARRANTY

Thank you for your purchase of Venture Heat products. This Limited Warranty applies to physical Venture Heat goods and covers any defects in heated electronic components for one (1) year and a period of six (6) months for lithium batteries from the initial date of purchase from an authorized Venture Heat dealer. This warranty is non-transferrable and proof of purchase is required for warranty service. Damage to the heating system or parts resulting from modification, accident, fire, flood, misuse, or unauthorized service or negligence outlined in the instructions printed herein will void the warranty. **Warranty does not cover wear and tear of the garment; fabric, zippers and other accessories are not covered under this warranty.** Under no circumstances will Venture Heat be liable for the incidental, consequential, or punitive damages resulting from the product, whether or not the product is found to be covered by this warranty.

For more warranty information contact your local Venture Heat distributor:

USA: [service@ventureheat.com](mailto:service@ventureheat.com)  
Canada: [service@ventureheat.ca](mailto:service@ventureheat.ca)  
Australia: [www.ventureheat.com.au](http://www.ventureheat.com.au)  
Europe: [www.varmekle.com](http://www.varmekle.com)

For quickest service, please email the following documents to the Venture Heat Service Department in your region:

- 1) Proof of purchase (receipt, packing list, etc)
- 2) Reason/description for service request.
- 3) Contact information including email address and phone number.

visit [www.ventureheat.com/warranty](http://www.ventureheat.com/warranty)  
to sign up for our complimentary 3 month extended warranty.

## Why choose Venture Heat?

- FDA Registered (USA)
- Medical Device Directive (EU)
  - Targeted area pain relief
  - Safer low-voltage DC power
  - Comfortable infrared warmth
- International travel friendly 100-240v

[www.ventureheat.com](http://www.ventureheat.com)  
[www.ventureheat.ca](http://www.ventureheat.ca)

## EXPERIENCE DRUG-FREE INFRARED PAIN RELIEF

- Increased blood circulation
- Deep penetrating warmth
- FDA registered heat therapy

### Questions?

USA  
[service@ventureheat.com](mailto:service@ventureheat.com)  
877-261-9477

CANADA  
[service@ventureheat.ca](mailto:service@ventureheat.ca)  
778-216-0260

[www.ventureheat.com](http://www.ventureheat.com)

## SOOTHING INFRARED WARMTH

Venture Heat **Far Infrared (F.I.R) Pain Relief Heat Therapy Collection** uses cutting-edge carbon-fiber nanotechnology to provide effective warmth that penetrates deeply into the muscle—far beyond traditional heating pads. Our Heat Therapy Collection stimulates blood circulation, promotes your body’s natural healing processes, provides stress relief, and relaxes your muscles.


## DRUG-FREE PAIN RELIEF

### How do heat therapy wraps provide relief?

For centuries, heat has been used as an effective therapy technique to relax muscles and provide pain relief. Our Heat Therapy Collection utilizes **far infrared ray heat** technology, which puts soothing, sun-like heat **deep under the skin and into the muscles**, easing chronic pain, relaxing muscle spasms, and improving blood circulation. Infrared heat therapy penetrates deeper into the body, aiding in the release of the body’s natural endorphins and making it a great choice for those suffering from pains, injury, or stress.

### Where can I use Far Infrared Therapy?


- Arthritis
- Back pain & lumbago
- Carpal tunnel syndrome
- Muscle stiffness
- Physical therapy
- Post-Op recovery
- Pre-workout warmup
- Raynaud syndrome
- Stress relief
- Sciatica
- Sprains
- Tennis elbow

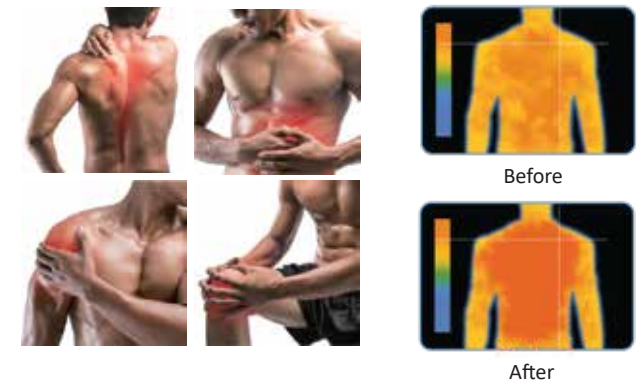
### What are Far Infrared (FIR) Rays?

Far infrared (FIR) rays represent a part of the infrared spectrum that transfers energy in the form of concentrated heat, much like our sun’s rays. FIR can deeply and safely penetrate up to 1.5 inches (4 cm) beneath the skin, stimulating blood circulation and cell repair. Unlike other heat delivery methods, such as traditional heating pads, far infrared heat provides full therapeutic relief at lower temperatures.


Using the newest advancements in carbon fiber nanotechnology, Venture Heat’s Heat Therapy Collection implements this beneficial, far infrared heat into versatile, high quality heating pads and wraps.
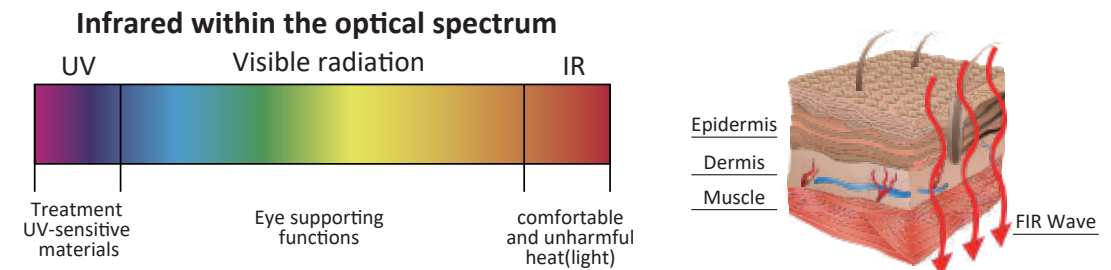
## Operation Instructions

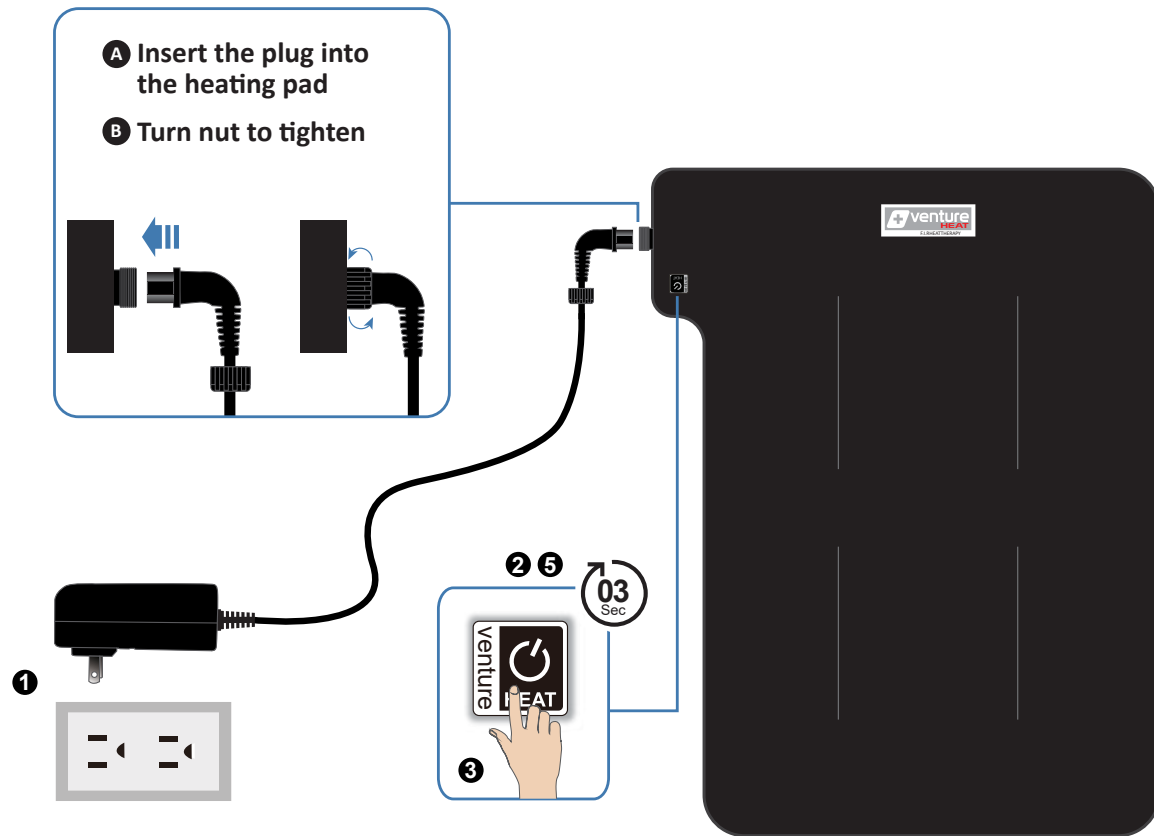1. Connect plug to a power source.
2. Press and hold the power button for 3 seconds to turn on. Pad will start at medium temperature (white).
3. Tap button to select heat level.

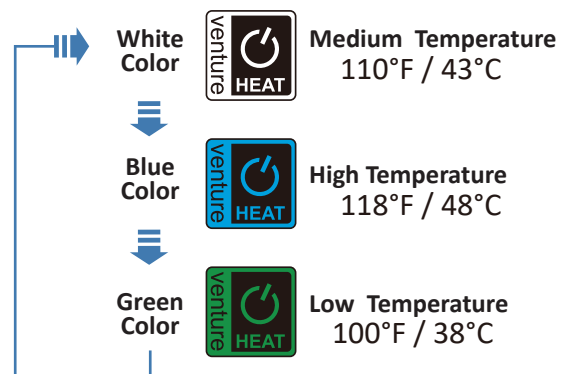4. After 10 seconds of inactivity, the button's Auto Dimming feature will activate. Tap once to wake and change settings.
5. Turning off your heating pad:
  - Press and hold the power button for 3 seconds to turn off
  - Or, the built-in timer will automatically turn off after 1 hour

## Features


### Built-in 1 Hour Timer

Most doctors and physical therapists recommend using heat therapy in 30 minute increments. For your protection, your FIR heating pad features a convenient 1 hour safety shut-off timer.


### Intelligent Heat Therapy

If your FIR heating pad is left on HIGH for more than 30 minutes, the onboard digital controller will automatically reduce your heating pad's temperature to LOW to protect your body from overheating.

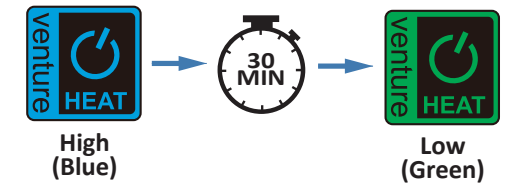
### Automatic Dimming

After 10 seconds of inactivity, the illuminated power button on your FIR heating pad will automatically dim. Tap once to wake. Tap again to change heat settings, or press and hold for 3 seconds to turn off.

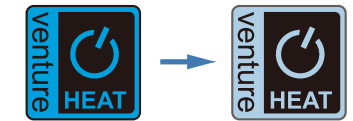
### Designed for the Left and Right Side of the Bed

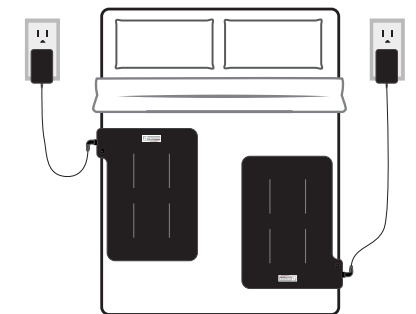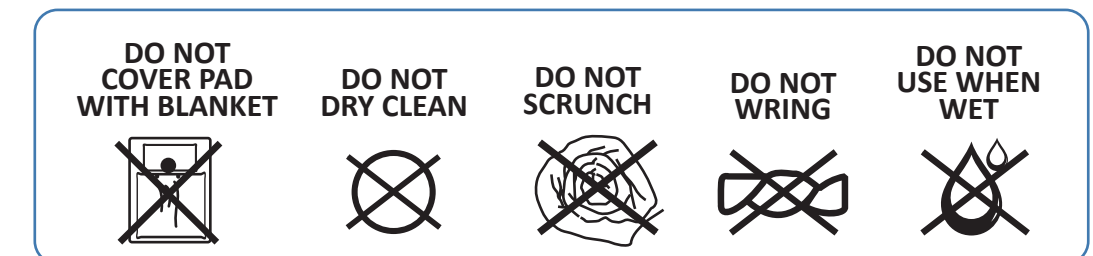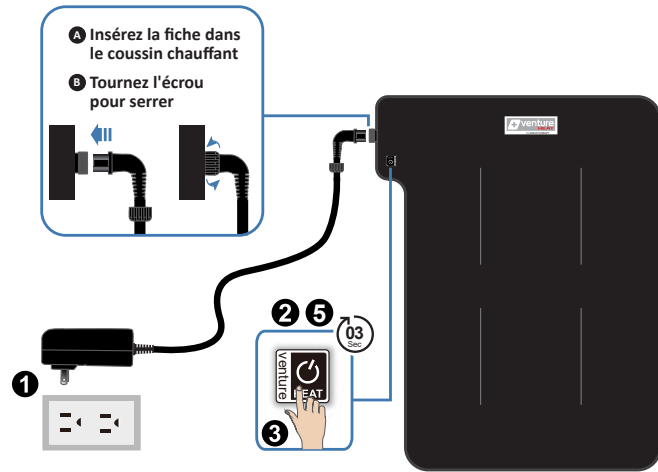
Your FIR heating pad can be placed on either the left or right side of the bed. For optimal results, the logo should face toward the body.

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