

SKINREVERSE<sup>®</sup>  
**ADVANCED DERMALIFT**  

---

**USER'S MANUAL**






## TABLE OF CONTENTS


Get to Know Your Device	4
Technologies Incorporated	5
How to Use	7
Modes of Use	9
Summary Table of Modes	16
Usage Tips	17
Frequency of Use	19
Charging Instructions	20
Cleaning and Maintenance	21
Product Specifications	22
Safety Warnings	23
Warranty	27

# GET TO KNOW YOUR DEVICE


## PRODUCT STRUCTURE

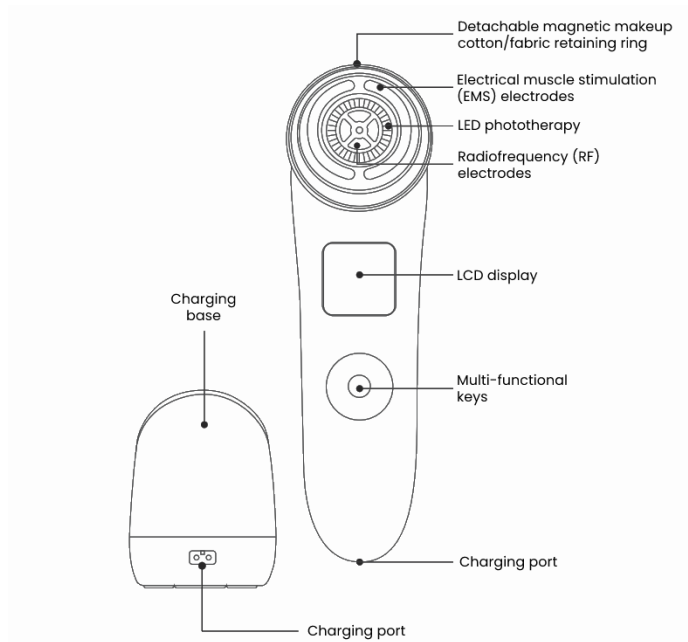


## COMPONENTS & ACCESSORIES

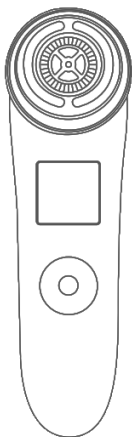


Main device

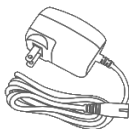

Power adapter

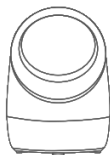

Charging base

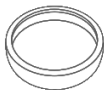

Retaining  
magnetic  
ring



User's manual


## TECHNOLOGIES INCORPORATED



### **MULTIPOLAR RADIOFREQUENCY (RF): 1MHz**

Clinical-Grade **RF 1MHz** to generate collagen stimulation at home.

**Multipolar RF** technology that prevents overheating and tissue damage.

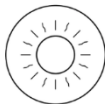

### LED LIGHT THERAPY: RED 623nm & BLUE 465nm

**Red Led** light wavelength of **623nm** to penetrate deep into the skin & increase collagen production & natural hydration, to regulate oil production and improve circulation. **Blue Led** light wavelength of **465nm** to treat acne bacteria & acne scars.



### IONTOPHORESIS (ION): 100Hz

Positive ions are generated to attract negatively charged impurities in pores under the surface layer of the skin.



### SONIC VIBRATION: 10000Hz

In addition to offering gentle yet effective cleansing, sonic vibration stimulates blood circulation and lymphatic drainage, clears away toxins, and enhances the performance of your skincare products.



### ELECTROPORATION ENI: 10Hz

Up to 90% of substances are transferred into the skin due to increased skin permeability. No bruising or reddening of the skin occurs when compared to mesotherapy with needles.



### ELECTRICAL MUSCLE STIMULATION(EMS): 100Hz

It provides deeper muscle stimulation by contracting facial muscles. It's a cost-effective and non-invasive way to firm facial contour and lift sagging skin.



### COOLING: Cooling at 57°F (13°C)

The device reaches 13°C in 30sec. It soothes the skin, calms redness and inflammation, and shrinks pores.

# HOW TO USE

---

## BUTTONS AND DISPLAYS




---

## HOW TO USE THE BUTTONS

- Short press the OK button to turn on/off the device.
- Swipe ▲▼ to change the mode. Once you select a mode press OK to select the mode. While working on a certain mode, press OK to pause mode.
- Swipe ◀▶ to change the level strength.

## DISPLAY MENU

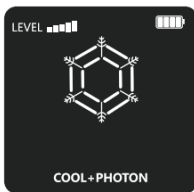

The level indicator on the screen shows the level strength when using a certain mode.



The power icon found on the top-right corner of the screen displays the charge level of the device.



The RF icon in the middle of the screen shows that radiofrequency will be applied. The RF function can be used in a standalone form to activate collagen and reverse the effects of aging, or it can be used in conjunction with other modes, ION, ENI, and EMS.






# MODES OF USE

Skinreverse Advanced Dermalift has 5 modes:

1. RF+PHOTON MODE
2. ION MODE
3. ENI+PHOTON MODE
4. EMS+PHOTON MODE
5. COOL+PHOTON MODE

Select the appropriate mode or combination of modes according to your needs, your skincare routine, or your skin goals. Each mode has a pre-set maximum duration, after which the device will turn to standby mode.

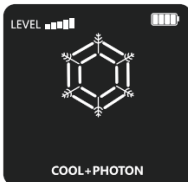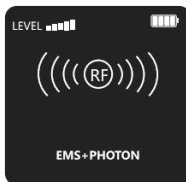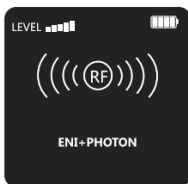

# MODE 1

---

## RADIO FREQUENCY (RF) MODE

RF+RED LED LIGHT



Professional radiofrequency and phototherapy functions are synchronized to enhance anti-aging effects. This mode is mainly aimed at older users with milder signs of aging. RF heats the dermis for 20sec up to 140oF (60°C) after which a constant heat between 104°F -107°F (40-42°C) is maintained. This enables radio waves to penetrate the deeper skin layers where it can stimulate the production of collagen and elastin fibers, promoting skin metabolism thus rejuvenating and firming the skin. Wash your face before use and

use skin care products such as anti-aging gels or water-based products for optimal results.

Pre-set operating time: 5 minutes.

## MODE 2

---

### ION MODE

RF+ION+SONIC MASSAGE



Radiofrequency, ion export, and sonic massage functions are synchronized to deep cleanse the skin. Prior to using the ION mode, wash your face, then you may use a toner to further wet your skin. Remove the magnetic ring on the head of the product, attach a cosmetic cotton, reattach the ring back onto the

device and spray toner on the cosmetic cotton. Clean your face again after operation.

Pre-set operating time: 3 minutes.

## MODE 3

---

### ELECTROPORATION (ENI) MODE RF+ENI+RED LED LIGHT+SONIC MASSAGE



Radiofrequency, electroporation, phototherapy, and sonic massage functions are synchronized to nourish the skin. Prior to using the ENI mode, wash your face and use ION mode to deep clean the skin to prevent the penetration of unknown chemicals into the cells. Apply skincare products to your face (such as essence, lotion, serum, cream, moisturizer, mask,

etc.). Please note that clay mask cannot be absorbed with this product. It is recommended to start with the lowest level for the first time. After gradual adaptation, the appropriate level can be selected according to the human body feeling.


Pre-set operating time: 2 minutes.

## MODE 4

---

### ELECTRICAL MUSCLE STIMULATION (EMS) MODE

RF+EMS+RED/BLUE LED LIGHT+SONIC MASSAGE

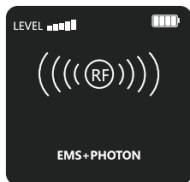

Radiofrequency, electrical muscle stimulation, phototherapy, and sonic massage functions are synchronized to lift and tone facial muscles. Apply a

conductive gel or a water-based product evenly on your face. Depending on your sensitivity level, you may want to adjust the intensity level and choose one that allows you to comfortably use this mode for 8 to 15 minutes (lesser time is recommended for those with milder signs of aging, while those with more prominent and pronounced signs of aging may use the device for longer, but not exceeding 15 minutes per session).

Pre-set operating time: 5 minutes.

## MODE 5

---

### COOL MODE

COOLING+BLUE LED LIGHT



Cooling and phototherapy are synchronized for soothing and pore tightening. After completing






## SKINREVERSE® ADVANCED DERMALIFT

"ENI+PHOTON" or "EMS+PHOTON," clean the cosmetics residue with a cotton cosmetic pad from the head of the device. Choose this mode to calm redness and inflammation, soothe your skin and shrink the appearance of pores. Intensity selection is unavailable for this mode, as this mode has only one pre-set intensity. Since this mode does not use RF, you can hold the device on the same spot for 3-10 seconds. The recommended usage time is 3 to 5 minutes.

Pre-set operating time: 2 minutes.




## SUMMARY TABLE

MODE	TECHNOLOGIES INCORPORATED	SKINCARE GOAL	USE WITH	PRE-SET OPERATION
<b>RF</b> 	RF RED LED LIGHT	ANTI-AGING, STIMULATE COLLAGEN, SKIN REJUVENATION	ANTI-AGING GELS, WATER BASED PRODUCTS, SERUMS, MOISTURIZERS	5 min.
<b>ION</b> 	RF ION SONIC VIBRATION	DEEP CLEANSING	TONER, MAKE UP REMOVER	3 min.
<b>ENI</b> 	RF RED LED LIGHT ENI SONIC VIBRATION	ENHANCE PRODUCT ABSORPTION, DEEP PENETRATION OF ACTIVE SUBSTANCES	ANTI-AGING GELS, WATER BASED PRODUCTS, SERUMS, MOISTURIZERS, MASKS	2 min.
<b>EMS</b> 	RF RED/BLUE LED EMS SONIC VIBRATION	LIFT AND FIRM FACIAL MUSCLES	CONDUCTIVE GELS, ALOE VERA GEL, WATER BASED/OIL FREE PRODUCTS	5 min.
<b>COOL</b> 	COOLING BLUE LED LIGHT	SOOTHE AND CALM, SHRINK PORES	-	2 min.



## USAGE TIPS

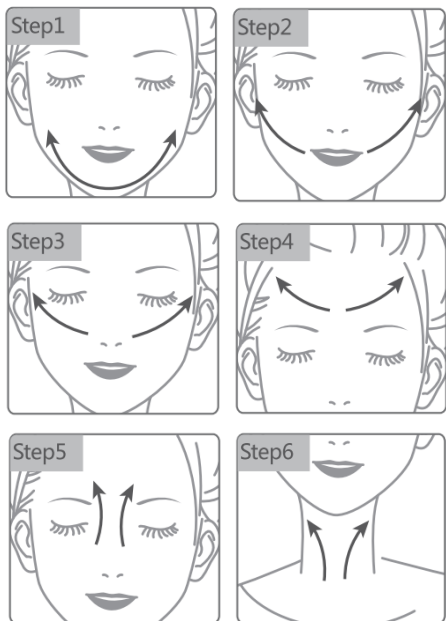

### **IMPORTANT NOTES:**

- **Always keep the device in motion when in use, to avoid overheating, skin damage and discomfort (except cool mode).**

- **Always ensure that the area to be treated has been thoroughly cleaned before use.**
- **If you plan to use all the modes in one session, remember to use first, the ION mode, to deep cleanse the skin, prior to use any other mode.**
- **This device has been designed to work without any need for conducting gels, however a conductive gel helps with conduction-of RF waves and electrical currents-and gives a good glide without dragging the skin.**  
**Tip: During an ENI or EMS treatment, whenever a gel or water-based product starts to dry out, you can spray on, a thermal spring water product, to make it slippery again.**
- **Gently press the device on skin and massage in circular motion or upward direction.**
- **Treat neck skin by moving the device slowly from the lower neck area towards the jawline.**
- **Use lower intensities (up to level 3) around the eyes. The device is not suitable to treat the eyes, the upper or lower eyelid.**
- **Do not use the device over the Thyroid/Adam's apple area.**
- **The device will switch off automatically if not used.**

# FREQUENCY OF USE

MODE SELECTION/ COMBINATION OF MODES	TOTAL SESSION DURATION	FREQUENCY OF USE	
		AGE BELOW 40	AGE OVER 40
ION	3-5 MINUTES	EVERY TWO TO THREE DAYS	EVERY TWO TO THREE DAYS
ENI	3-5 MINUTES	EVERY NIGHT AFTER APPLYING SKINCARE PRODUCTS	EVERY NIGHT AFTER APPLYING SKINCARE PRODUCTS
RF	5-8 MINUTES	ONCE A WEEK OR EVERY THREE TO FOUR DAYS	EVERY TWO TO THREE DAYS
EMS	5-8 MINUTES	ONCE A WEEK OR EVERY THREE TO FOUR DAYS	EVERY TWO TO THREE DAYS
COOL	2-5 MINUTES	NO LIMIT	NO LIMIT
ION+ENI+EMS+COOL	8-15 MINUTES	EVERY THREE TO FOUR DAYS	EVERY TWO TO THREE DAYS
ION+ENI+RF+COOL	8-15 MINUTES	EVERY THREE TO FOUR DAYS	EVERY TWO TO THREE DAYS
ION+ENI+COOL	8-10 MINUTES	EVERY TWO TO THREE DAYS	EVERY TWO TO THREE DAYS
INTENSITY PER SESSION		UP TO LEVEL 4	UP TO LEVEL 5


# HOW TO CHARGE

- This device is portable after charging.
- Please charge the battery for the first time before use.
- Full charge time: about 2 hours

1. Power off the device.
2. Connect the USB cable at the back of the charging base or directly at the bottom of the device.
3. Connect the USB cable to an electrical outlet.
4. The power icon on the screen of the product



lights up to indicate that charging is in progress. After charging, the power icon shows full grid. Remove the adapter from the outlet. If the product and charger are slightly heated at this time, this is a normal phenomenon.

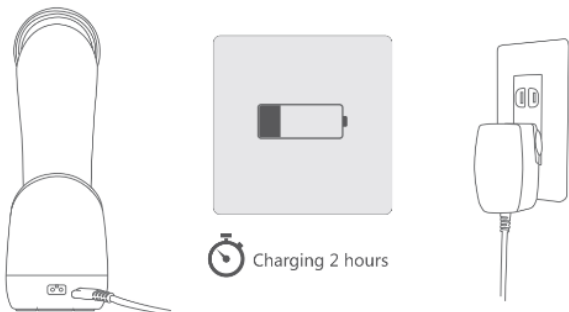

## CLEANING AND MAINTENANCE

- Before use, please clean the head of the apparatus with a wet paper towel or cotton pad.
- After use, wipe the magnetic ring and the head with a cotton pad or soft cloth. Do not flush this product in water.
- Before cleaning, make sure the product is powered off and the power adapter is disconnected from the power supply.
- When cleaning the adapter, make sure to pull it out of the power supply and wipe it only with a dry cloth.

- After each use, please clean the product to avoid the accumulation of grease and other residues to ensure optimal long-term product performance.
- Make sure all components are dry before reconnecting to the adapter.

## PRODUCT SPECIFICATIONS

PRODUCT NAME	:	SKINREVERSE ADVANCED DERMALIFT
DIMENSIONS	:	H 190 × W 55 × D 64.5 (mm)
MATERIALS	:	ABS AND 316 STAINLESS STEEL
BATTERY	:	7.4V, 650mAh
LCD	:	1.5 inches
ADAPTER	:	INPUT 100-240V, 0,35A OUTPUT 9V 1,5A
RATED POWER	:	10W
WEIGHT	:	ABOUT 280gr
WATER-RESISTANT	:	IP65
FULL CHARGE TIME	:	2h

# SAFETY WARNINGS

This product has been designed for facial skincare.

Stop using if you suspect there is a problem with the device.

Contact us at **contact@myskinreverse.com** or by calling us at (+1) 3025202505.

Continued use after unusual performance is not recommended and can lead to secondary failure of the apparatus or damage to the skin. That may lead to the manufacturer's warranty being voided.

Read the safety warnings carefully before using the device. Follow these directions for safety purposes.

The manual categorizes the hazards into three categories:



**DANGER**

Misuse may cause death or severe injury



**WARNING**

Misuse may cause severe injury



**NOTICE**

Misuse may cause injury or property damage which may require a costly repair

 **DANGER**

- Continued use of defective product can lead to user harm.
- When charging this product, please use the adapter in the packing box, and appropriately connect it to the household outlet to avoid malfunction, fire or fire caused by heat.
- Since this product has a built-in rechargeable battery, do not place the product in or near a fire. Avoid using it in high temperatures and superheated environment as this may lead to irreversible damage, malfunction, fire, or other unintended outcomes.
- This product's battery has a certain lifetime. If you continue to use it when it has surpassed the service life of 10,000 hours, chemical reactions may lead to expansion and deformation. Additionally, the device may be deformed at high temperatures or under extreme pressure such as drops from significant heights. If the product is deformed, please do not use it to avoid the risk of electric shock and fire.



**WARNING**

- Improper use that does not comply with the instructions contained in this manual and/or other advice provided by Skinreverse® on its official website and other means, can lead to harm and damage for which the Manufacturer and/or Seller won't be liable.
- Do not use a broken charger or a loose socket attachment because this may trigger electric shock, fire among other bad outcomes.
- Do not use the adapter with high voltage, other than 100V or 240V, to avoid fire, electric shock, etc.
- Do not disassemble or modify the product by yourself, so as this can result to serious harm. If repair is required, please contact Skinreverse®.
- Do not use this device for purposes other than directed as malfunction may occur.
- This product is intended to be used by adults only. Please do not allow the children or those who are physically unable to use this device. Adolescent should use it under the supervision and instruction of adults.

- Keep the device away from children.
- Please stop immediately using the product and disconnect the power cable in case of abnormal heat when in use or during charging.

**NOTICE**

- Do not flush this product in water. After using, remove any cosmetic residue with a tissue or cotton cosmetic pads.
- Do not drop the device. Malfunction may occur, and additional cost for repair may apply even during the warranty period.
- This appliance is not intended for use by persons (including children) with reduced physical or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instructions regarding the use of the device by a person responsible for their safety. Children should be supervised to ensure that they do not play with the device.

## BATTERY INFORMATION

This product has built-in rechargeable lithium batteries. Please pay attention to the following when discarding the product.

- the battery must be removed from the device before it is scrapped.
- the device must be disconnected from the power supply when removing the battery.
- the battery is to be disposed of safely.
- The device contains batteries that are non-replaceable.

To recharge the battery, only use the detachable supply unit provided with this device.


## WARRANTY

This product, by default, is covered by a standard one-year manufacturer warranty. The manufacturer warrants this product to be free from defects in workmanship and materials under normal use and conditions, for a period of 1 (one) year from the original invoice date. Shipping and handling fees are to be paid for by the customer.



---

Innovative, Hi-Tech  
Beauty Devices

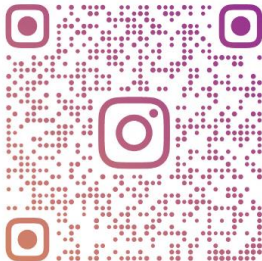

Website: <https://myskinreverse.com>

Email: [contact@myskinreverse.com](mailto:contact@myskinreverse.com)

Tel: (+1) 3025202505

©2022 SKINREVERSE®