

www.thelumia.com

LUMIA

— TECHNOLOGY MEETS BEAUTY —



PENETRADO

Read all directions before use.

Specifications

1	Name	LED Infrared Device
2	Input current	1A
3	Input voltage	5V
4	Input power	5W
5	Type of heat source and number	26pcs 625nm LED,8pcs 850nm far infrared LED
6	Design Principle	Red light+far infrared light+heating therapy
7	Therapeutic temperature range	53°C±5°C
8	Wavelength	Red LED 625 ±5nm, far infrared LED 850 ±5nm
9	Storage condition	Temperature: -10°C to 60°C / 14°F to 140°F Relative humidity: 5%-95%
10	Usage time	5-7 minutes on each treatment area
11	Treatment area	18cm ²
12	Size	51.5*49.8*203.5mm
13	Weight	229g
14	Lithium ion battery	750mah

Step 4: When you are done with your treatment, turn OFF the APPLICATOR by pressing the ON/OFF SWITCH. The device will beep once and the WORKING INDICATOR LIGHT on the APPLICATOR will turn off.

Step 5: Always apply a small amount of a moisturizing cream on the treated skin areas after each session and massage it gently in circular motions.

Step 6: It is imperative to clean the glass probe after each use. Unplug the device from its charging base to clean it. Clean the probe thoroughly by using dry cloth, cotton ball or paper towel. Refrain from using sanitizing substances to avoid paint damage. Only the device's glass probe is suitable for direct contact with creams and/or serums.

●A mild redness - After each use you may notice a slight reddening of the skin (mild erythema) of the treated area. This is normal and should fade within 24-48 hours. Discontinue use if you encounter any discomfort or skin redness that lasts for longer than 48 hours.

Usage Time

First use: In the first month (4 weeks), use at least 3 times a week, 20-28 minutes each time (5-7 minutes on each treatment area).

Problem	Possible cause	Solution
Device suddenly turns off	The rechargeable batteries are empty	Place the device back in the charging station for 3 hours
Yellow ERROR INDICATOR Light is ON	Malfunction of the APPLICATOR	Press the ON/OFF SWITCH twice to switch the APPLICATOR off and on again

Table of Contents

Function.....	3
Protect Your Machine.....	3
Important Safty Information.....	4
Product Description.....	6
Battery Recharging.....	7
Treatment Tip.....	7
Usage Time.....	9
Specifications.....	10

Thank you for choosing us.

Before using this infrared device for the first time, read this User Manual.

Function

The red light heating beauty device uses infrared LED to provide continuous heat therapy. The machine consists of a handle, an electronic control panel, 26pcs 625nm red lights, 8pcs 850nm far-infrared lights, and glass sensor probes. The glass sensor probe can transfer 50.5°C heat to deep skin.

- It's an over the counter LED infrared device that emits energy in the red & infrared region of the spectrum. The device is used for treatment of multiple signs of aging to reduce facial lines, improve the skin's appearance and relieve joint pain.

- It's a normal phenomenon if the red light is not bright when you touch the probe for first use, then please turn it on again. The glass head is activated by touch sensing. Heating and two kinds of lights work at the same time.

Protect Your Machine

Do not store the device when it is powered on. Store the device in a clean dust-free area. Please keep it out of the reach of children, do not let minors use it.

To ensure the service life of the glass head, suggest to use clean soft cloth to clean the glass probe, make sure it is dry and clean.

Important Safety Information

Do not use this device if you have a condition that makes your skin highly sensitive to light, such as systemic lupus erythematosus, porphyria, polymorphic light eruption and solar urticaria.

Step 3: Place the TREATMENT PROBE on the area that you want to treat first. The LED light will automatically activate. Apply light pressure while slowly moving the APPLICATOR in small circular motions as shown in the image.



The Treatment probe surface temperature will slowly increase. Remove the APPLICATOR from the skin if either the skin or the APPLICATOR is too hot. Perform the treatment for 5-7 minutes on each treatment area.

IMPORTANT:

Keep the device in constant motion when in contact with the skin! Holding the device still while treating can cause redness or a mild burn.



Move gently from chin to cheek.



Move on the forehead from left to right.



Gently move the neck from bottom to top.




Move gently from top to bottom.

Battery Recharging

- Make sure the device is turned off.
- Plug the charging cable into the base socket, connect with a charging head to the electrical outlet.
- When charging, the charging indicator flashes, the charging time is about 3 hours.
- Please note your LED infrared device has an inactivity timer that will automatically shut the device off after 2 minute if it is not in contact with the skin.

Treatment Tip

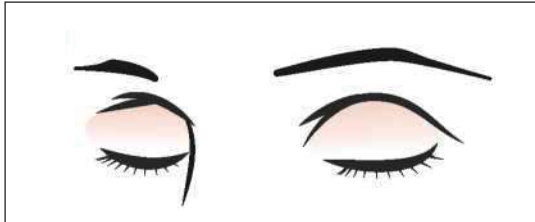

THIS DEVICE IS NOT FOR USE ON EYELIDS OR OPEN EYES!

For your safety, the APPLICATOR only works when there is full skin contact and will stop automatically if it is lifted from the skin.

Step 1: Wash and dry your face. The skin should be clean, dry and free of any powder, gel or cream.

Step 2: Lift the applicator from its cradle in the base unit. Press the On/Off Switch. The device will beep once and the Indicator Light will lit one after another for a short time and then the Green indicator will turn steady on. Your machine is ready for treatment.

- You have a history of light triggered seizures, especially if you have photosensitive epilepsy, as the pattern of bright light could induce seizures.
- You are undergoing radiation or chemotherapy treatments. Discontinue use if you encounter any discomfort or skin redness that lasts for more than two days. You may have sensitivity to light or another condition. Consult with your doctor.



WARNING! TO AVOID ALL SERIOUS INJURY



WARNING! THIS DEVICE IS NOT FOR USE ON EYELIDS AND OPEN EYES



WARNING! THIS DEVICE IS NOT FOR USE BY PERSONS UNDER AGE 18 YEARS



WARNING! DO NOT USE THIS DEVICE ON AREAS OF POTENTIAL SKIN CANCER SUCH AS DARKBROWN OR BLACK SPOTS, LARGE FRECKLES, BIRTH MARKS, MOLES OR WARTS. TO DO SO MAY RESULT IN DELAYING THE BEST MEDICAL TREATMENT

Important Safety Information

WARNING! - TO REDUCE THE RISK OF ELECTRIC SHOCK THAT MAY CAUSE SERIOUS INJURY OR DEATH:

- Do not place or drop the machine into water or any other liquid.
- Do not use the machine if it is damp or wet.

Important Safety Information

WARNING!-TO REDUCE THE RISK OF ELECTRIC SHOCK THAT MAY CAUSE SERIOUS INJURY OR DEATH:

- Do not place or drop the machine into water or any other liquid.
- Do not Use the machine if it is damp or wet.
- Do not attempt to open or repair your machine.
- Do not operate the device if it has a damaged cord or plug.
- Use the device only for its intended purpose and as described in its manual.

PRECAUTIONS-Take the following precautions to avoid minor or moderate injury and damage to your device:

- Do not use the device if you have an active skin condition in the treatment area including sores, psoriasis, eczema, sunburns, herpes simples, open wounds or a rash.
- Do not use the device on skin that is infected, burned, or cut as this could cause discomfort or irritation.
- Never look directly at the light coming from the machine APPLICATOR.
- Never use the device more than 40 minutes each time, to avoid over-heating and damage to the device. Let the device cool down for at least 30 minutes before the next use.

Product Description



Package included:

