# VITORI

## Premium Kristallmatte

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# Commissioning and safety instructions

- Please read the instruction manual carefully. Place the mat on the floor and do not use the product folded.
- First, plug the controller's connector into the mat.
- Then connect the power plug to a power outlet.
- Press the power switch to turn on the mat.
- Press the timer key to set the heating time.
- Turn the temperature control knob to set the desired temperature. When the mat is sufficiently warm, adjust the setting via the control knob to your comfort.
- Turn off the power switch and unplug the power cord from the wall outlet when not using the mat.
- For safety reasons, please monitor the temperature carefully while using the product.
- You can use a cotton sheet while using the product, but not thick padding.
- Do not pierce sharp objects or damage the product.
- The power plug and the control panel must not come into contact with liquid. Make sure that the product is dry before use.
- Heat may cause skin irritation or burns. Please lower the temperature if it becomes too hot.
- Do not use the product folded or crumpled. Do not use the product near fire, heaters or flammable objects. This may cause a fire or electric shock.
- Please do not allow children to play on or near the product.
- Please stop using the product and contact us if there is any problem.
- Please consult your doctor if you have any medical questions.
- Please drink enough water at all times.

Occasionally, a so-called "initial worsening" may occur. In this case, there is initially a subjective intensification of the symptoms before the desired improvements occur. The initial worsening effect, which is also known from other types of therapy, is by definition a "healing reaction" or removal of blockages. This effect, which may occur, is not a criterion for excluding the use of the product. If pain occurs in connection with previous illnesses, please reduce the application intensity first. If, in the course of the initial worsening, persistent severe pain occurs, please discontinue the therapy and consult a doctor immediately to clarify whether and under what circumstances it can be continued.

# VITORI® Premium Crystal Mat

## Forms of therapy :

- Pulsating Magnetic-Field Therapy
- Far-Infrared Heating
- Negative Ions
- Biophoton Light Therapy
- Healing-Stone Effect


## Data:

Name:	VITORI® Premium Crystal Mat
Dimensions:	190cm × 80cm × 3cm
Weight:	12 kilogram
Power consumption:	280 W
Temperature range:	30°C - 70°C
Magnetic field strength:	150 $\mu$ T
Power supply:	AC220V 50Hz
CE EMC No.:	UAK180726216S
CE LVD no:	HUAK180726641E
FDA approval:	3015008734
Protection type:	Class II Type B

## Structure:

1. vegan PU leather
2. Thermoplastic silicone urethane layer
3. healing stones - amplify the long-wave infrared rays  
*Jasper | Jadeite | Agate | Aventurine | Lapis lazuli | Chalcedony | Amethyst*
4. first waterproof layer
5. 36 biophoton dots emit coherent light
6. carbon layer against electromagnetic radiation
7. second waterproof layer
8. ionized mineral tissue generates negative ions
9. hyperfiber layer generates long-wave infrared rays
10. copper silicon mesh for shielding electrosmog
11. CALOR heat layer with thermostat and bimetal
12. aluminum layer to reflect the infrared rays
13. gossypium cotton for thermal insulation
14. magnetic field coils for pulsating magnetic field therapy
15. tile layer for thermal insulation
16. Polyester bottom cover layer

# Control panel and program settings

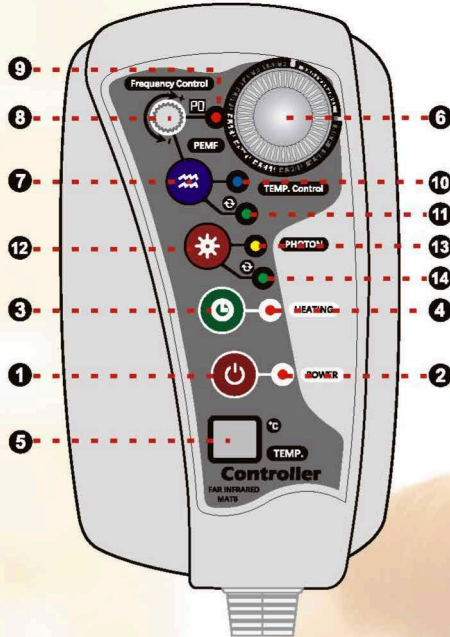**1 Power button:** Turns on/off the power supply for all functions of the mat.

**2 Power LED:** Lights up when the power is on

**3 Infrared timer button:**

0H - no timer (default 12 hours)

1H - 1 hour | 4H - 4 hours | 8H - 8 hours | -H - IFR OFF

**4 IFR LED:** lights up when the mat heats up.

**5 Temperature °C and function display**

**6 Temperature control knob:** Controls the temperature of the mat.

**7 PEMF program key:** 14 PEMF programs

**8 PEMF frequency adjustment dial:** To increase or decrease the PEMF frequency.

**9 P0 PEMF LED:** lights up when you can select the PEMF frequency manually.

**10 PEMF LED:** lights up when PEMF is switched on.

**11 PEMF loop LED:** Lights up, when PEMF loops are running.

**12 Biophoton key:**

-P biophotons off | 1P biophotons on 30 minutes | 2P biophotons in 8 loops-one loop = biophotons on (30min), biophotons off (30min)

**13 Biophoton LED:** Lights up when the biophotons are active.

**14 Biophoton loop LED:** Lights up when the biophoton loop program is on

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# Temperatures and duration of use Infrared-Heating

We recommend an application time of 15 minutes up to 8 hours.

The VITORI® crystal mat can also be used overnight.

In this case, the temperature should not exceed 40°C. Please always pay attention to your well-being and slowly feel your way to the favored warmth.

## **35°C - 40°C - For 120 minutes and longer**

- general support of health
- ideal temperature for headaches and migraines
- immunostrengthening, vitalizing
- mild detoxification (for sensitive people)
- Strengthening of the regeneration phase during sleep, supports cell renewal
- deepens sleep, can help with sleep disorders

## **45°C - 50°C - For application duration of 120 minutes and longer.**

- promotes relaxation as well as function of the nervous system
- Gives off blood pressure as well as blood sugar levels
- can reduce muscle, spine and joint pain
- strengthens the cardiovascular function
- supports breathing (e.g. in asthma & allergies)
- supports regulation of menstrual cycles and hormonal balance

## **55°C - 60°C - 15 minutes to 90 minutes**

- Relief from stress, anxiety and trauma
- possible weight reduction through dissolution of fat cells
- reduces muscle and joint pain
- Regeneration after physical tension
- mood elevating
- stimulates libido and sexual power
- deepest relaxation
- supports tissue regeneration
- supports regeneration after operations

## **65°C - 70°C - Intensive sessions 20 - 60 minutes**

- Defense against acute infections (viral, bacterial and in case of hyperacidity)
- in case of general weakness and blockages
- maximum detoxification of all cells and tissues
- maximum metabolic stimulation
- possible weight reduction through dissolution of fat cells
- session corresponds to sauna application

# Magnetic field (PEMF) programs

Press and hold the PEMF button for 3 seconds to return to P0 while the PEMF programs are running.

Program	Duration	Focus
P-		magnetic field off
P0	7h 20 min	manual program (1 - 30Hz)
P1	40 min	immune system
P2	60 min	detoxification
P3	60 min	concentration
P4	60 min	relaxation
P5	60 min	pain program
P6	60 min	connective tissue
P7	7h 10 min	restful sleep
P8	7h 6 min	regeneration
P9	60 min	General strengthening
PA	60 min	muscle relaxation
PC	60 min	energy
PE	60 min	raising awareness
PF	60 min	Schumann program

The controller has safety software that prevents harmful electromagnetic fields and protects the user from overheating and dehydration during the night. If the temperature is above 50 °C for more than 3.5 hours, the protection function will lower the temperature to 45 °C and the display will show [ ] . It is not recommended to turn off the automatic thermal protection system, but if you need to use the mat at a temperature higher than 50°C for more than 3.5 hours, you can manually deactivate this safety function. Press and hold the timer button for 5 seconds until 5- appears on the display. Press and hold the timer key for 5 seconds until the display shows 5A to reactivate the automatic protection function.

## Intensive programs

Note: The VITORI® crystal mat can be focused, for a particularly intensive application of individual complaints, as follows:

### **Immune System Intensive Program:**

Temperature: 35°C - 45 °C

PEMF: P1

Biophotons: 1P

Duration: 120 minutes

### **Intensive detoxification:**

Temperature: 60°C - 65°C

PEMF: P2

Biophotons: 1P

Duration: 10 - 20 minutes

### **Intensive headache/migraine program:**

Temperature: 35°C - 40°C

PEMF: P5

Biophotons: 1P

Duration: 30 - 60 minutes

### **Intensive program connective tissue**

Temperature: 45°C - 50°C

PEMF: P6

Biophotons: 1P

Duration: 60 minutes

### **Intensive program for muscles, bones and joints**

Temperature: 55°C - 60 °C

PEMF: PA

Biophotons: 1P

Duration: 15 - 90 minutes

### **Night program for regeneration**

Temperature. 35°C - 40 °C

PEMF: P8

Biophotons: -P

Duration: 6 - 8 hours

## Far-Infrared Heating

Infrared radiation energy is a form of energy that heats objects directly without heating the surrounding air. For optimal effect, we recommend using the crystal mat unclothed and covered. For hygienic reasons, you can place a white cotton sheet underneath.

For many millions of years, nature has been using the heat radiation emitted by the sun in the range around 9  $\mu\text{m}$  for evolution, without which no organic life would exist.

The VITORI® crystal mat generates far-infrared heat rays in exactly this wavelength range (4-16  $\mu\text{m}$ ).

Infrared deep heat can have a variety of positive effects:

- kill bacteria and viruses
- strengthen the immune system
- increase the release of endorphins
- permanently increase the body temperature
- act against electromagnetic fields
- eliminate toxic heavy metals
- consume calories
- quench pain

Already the physician Parmenides, who lived 2,000 years ago, said:

**"Give me the ability to create a fever and I will cure any disease."**

You can find more information here:  
[vitori.de/pages/infrarot-tiefenwarme](http://vitori.de/pages/infrarot-tiefenwarme)

# Pulsating magnetic field therapy

Every cell in our body has its own magnetic field between the cell membranes. Just like biophotons, the electrical impulses serve as information mediators between the cells. However, external influences negatively affect the voltage between the membranes.

The cells can no longer work properly and the person tires, there is a drop in performance or recurring pain.

Pulsating magnetic currents can trigger healing processes after just a few minutes of body contact.

In naturopathy, this form of therapy has been used for centuries against physical ailments.

Hippocrates already described the use of magnetic stones. Also the ancient Romans used the positive effects of magnets and the ancient Egyptians wore magnetic jewelry to strengthen your health.

Magnetic field therapy can do the following:

- acceleration of healing processes
- blood circulation promotion
- unclumping of erythrocytes agglutinated by electrosmog
- increased endorphin production
- increased oxygen transport
- release of tension
- pain relief for migraine
- pain relief for rheumatic diseases
- faster regeneration
- metabolism optimization
- improved oxygen supply
- improved blood circulation
- improved nutrient uptake
- improved sleep quality
- improved cell metabolism
- vitalization
- more resistant cells

You can find more information here:

[vitori.de/pages/pulsierende-magnetfeldtherapie](http://vitori.de/pages/pulsierende-magnetfeldtherapie)

# Biophotons

In 1975, the German physicist Prof. Dr. Fritz-Albert Popp provided scientific proof that our cells contain light and communicate with each other via it.

The 36 biophoton points of the crystal mat emit monochromatic, coherent light with a wavelength of 660nm. This can stimulate diseased cells by the coherence of the emitted light to restore the communication between the cells and activate the self-healing powers.

In addition to this main effect, the application may result in some other health benefits:

- activation of the metabolism
- acceleration of healing processes of the body
- increase body resistance to pathogens
- harmonization of electrosmog
- harmonization of hormonal cycles
- Reduction of allergy symptoms
- order of intestinal bacteria and improvement of the microbiome
- stress reduction
- strengthening the immune system
- mood enhancement
- stress Reduction
- improvement of detoxification functions
- improving the quality of sleep

You can find more information here:  
[vitori.de/pages/biophotonen](http://vitori.de/pages/biophotonen)

## Negative ions

If sufficient negative ions flow into the body, channels open in the cell walls. Toxins and heavy metals can be eliminated. This allows new healthy cells to emerge.

The VITORI® crystal mat automatically emits 1500 negative ions/cm<sup>3</sup> when in operation. To stay healthy we need at least 1000 per cm<sup>3</sup>.

Negative ions can:

- improve sleep
- support the weight reduction
- increase the concentration
- lighten the mood

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## Healing stone effect

The crystal mat owes its name to the processed healing stones adapted to the respective chakra. The individual vibrations of the chakras can thus be stimulated and the effect of the long-wave infrared rays can be strengthened.

- red jasper - root chakra
- golden jadeit - sacral chakra
- yellow agate - solar plexus
- green aventurine - heart chakra
- light blue lapis lazuli - throat chakra
- blue chalcedony - third eye chakra
- purple amethyst - crown chakra

You can find more information here:

**[vitori.de/pages/heilsteineundchakren](http://vitori.de/pages/heilsteineundchakren)**

## FAQ

### **How do I clean the VITORI® crystal mat?**

It is best to use the crystal mat with a cotton sheet.

Should it's still get dirty, clean it with a cotton cloth, water and a mild detergent.

### **Do harmful electromagnetic fields occur when using the crystal mat?**

A special layer of carbon acts as an EMF catcher and significantly reduces harmful electromagnetic radiation. The VITORI® EMF catcher functions similarly to a lightning rod by safely dissipating electromagnetic energy into the earth.

### **How can I test if the magnetic field is active?**

With the delivery of the Premium Crystal Mat you get a handy, transparent cylinder including magnet, which moves more or less depending on the set frequency when the magnetic field is active.

### **Can the VITORI® crystal mat influence the intake and effect of medication?**

The effects of drugs or dietary supplements are always dependent on the absorption and utilization in the body. Due to the improved blood circulation, but also due to the better functionality of the organs (gastrointestinal tract), medications are not only better and more utilized, but thus also have a stronger effect. With a permanent use of the VITORI® crystal mat, the medication dosage should therefore be adjusted, food supplements can often be reduced with the same effect.

Please always consult with your primary care physician.

### **How much water should I drink when using the VITORI® Crystal Mat?**

When using the crystal mat, please ensure that you drink enough fluids on a regular basis. Please drink at least 30ml per kilogram of body weight daily. We recommend drinking a full glass of water before and after each use.

### **When should I not use the magnetic field therapy?**

The pulsating magnetic field therapy of the VITORI® crystal mat is not suitable for persons with pacemakers. All other functions can be used and applied to the full extent.

# Warranty VITORI® Premium Crystal Mat

VITORI® provides warranty for the above device according to the rules of the Consumer Protection Act.

Date of purchase:
Customer name:
Address:

If the customer is not in possession of the warranty or has not filled it in, he is entitled to free or chargeable service according to the rules of the Consumer Protection Act regarding the crediting of the quality warranty period.  
If the product fails, contact customer service for assistance.

## Customer Service Service Standards:

1. The warranty on this product is 2 years.
2. The product will be repaired free of charge if it fails within the warranty period (two years) for normal reasons during use by the customer in accordance with the user manual.  
However, under the following circumstances, the customer will be charged for services even within the warranty period:
  - failure or damage due to improper handling of the product
  - failure or damage due to natural disasters
  - Failure due to non-compliance with the operating and safety instructions in the operating manual
  - failure due to unauthorized repair or modification by the customer
  - In case of wear of wearing parts
  - failure in an environment outside the range specified in the product warranty
  - the product fails as a result of its combination with products of a third-party manufacturer
3. The warranty is only valid domestically. Since it is not reissued, you must keep it carefully

This warranty promises you free service under the conditions and for the period described above. The warranty does not limit the customer's statutory rights.

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