

## IMPORTANT FOR OPTIMUM SAFETY: READ ALL INSTRUCTIONS BEFORE USE.

1. Consult your dermatologist prior to using DERMAFLASH® ONE device (ONE or "Device") if you have medical concerns.
2. DO NOT use ONE on unintended areas of the body including on surface of lips, nose and eyelids. ONE is designed ONLY for use on the face, chin and décolleté.
3. DO NOT use ONE on active acne, psoriasis, rosacea, open sores, skin with a rash, sunburned or distressed skin or on skin surfaces where cuts, abrasions or any kind of healing is present. DO NOT use Device on moles, cancerous areas or lesions or on skin which has recently been treated with a laser or recently chemically treated.
4. Some users report a temporary tingling feeling after using the Device this may be due to increased microcirculation. DO NOT use ONE if you experience a worsening of the skin (skin feels tense, uncomfortable, irritated, red, burning, etc.) STOP using ONE and consult your doctor.
5. Children should not use or play with the ONE. The Device has not been tested on children.
6. DO NOT reuse an Edge. For hygiene, safety, and best results ALWAYS use a NEW Edge. Discard Edge directly into wastebasket after use and place ONE in charging base. Ensure Edge tray is on flat surface when loading device with new Edge.
7. DO NOT use ONE device, charging base or Edges that are damaged or appear damaged and never attempt to open or repair the device. Doing this could result in injury or possible fatal electric shock. This product contains no serviceable parts.
8. DO NOT use, place or store ONE device or charging base where it can fall or be pulled into a tub or sink. Neither the charging base nor the device is intended for immersion or use in the shower. NEVER reach for the ONE charging base or plug that has fallen into water. Unplug immediately from the wall socket. Evidence of liquid infiltration into any component of the device will void the warranty.
9. Your ONE is rechargeable with its included USB cable.
10. Use this product only for its intended use as described in the owner's manual. Call DERMAFLASH at 888-652-9579 with any questions or concerns.

CAUTION: Before charging, make sure the charging port and charging base are free from water. DO NOT use your ONE device while charging. Never use the device or charger if they are damaged, and only use the charger supplied by ONE with the device.

DISCLAIMER: Users of ONE do so at their own risk. Neither DDK DRTV LLC nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this product. DERMAFLASH reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes. For complete, up-to-date terms, conditions and warnings please consult the full user manual at [dermaflash.com](http://dermaflash.com) on the Terms & Conditions page.

NOTE: This device complies with Part 15 Subpart B of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## WARRANTY TERMS AND CONDITIONS

Our limited, 1-year warranty covers defects in materials and workmanship of your device as long as you purchased from an authorized retail partner or directly at [dermaflash.com](http://dermaflash.com). The limited warranty begins on the date of original purchase and expires one year later. The warranty period is not extended by any period even if replaced by the Manufacturer, during the warranty period. The ONE device, inclusive of charging base ("Base") and handheld unit and exclusive of Edges, is warranted to be free from manufacturing defects for a period of 1-year. The warranty DOES NOT cover cosmetic deterioration caused by fair wear and tear, or damage caused by accident, misuse or neglect. Any attempt to open or take apart the device (or its accessories) will void the warranty. If you discover a defect and notify DERMAFLASH during the warranty period, DERMAFLASH will, at its discretion, replace the device free of charge. Claims under warranty must be supported by reasonable evidence that the date of the claim is within the warranty period. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period. To claim your warranty, you must call Customer Service at 888-652-9579. Shipping costs are nonrefundable. This undertaking is in addition to your statutory rights as a consumer and does not affect those rights in any way. Please call Customer Service at 888-652-9579 for questions.

This product contains a rechargeable lithium battery. It is installed inside the device. The battery is not a self-serviceable part. Do not attempt to repair or replace the battery sealed inside this product. Products containing an internal battery supply should be recycled or disposed of per state and local guidelines. ©DDK DRTV LLC 2020

## Cautions/Contraindications

1. ONE is designed for use on specific areas of the face, neck and décolleté. It is not intended for use on other areas of the body. Avoid using ONE on skin with active cold sores or warts, wounded, sunburned, excessively sensitive skin, healing problems, dermatitis, inflammatory (couperose) rosacea, eczema, psoriasis or where any kind of healing is present.
2. Do not use ONE on skin which has recently been treated with a laser or recently chemically treated.
3. Discontinue use of prescription strength exfoliants such as Tri-Retinol Complex, Retin-A, Renova®, Differin® (Adapalene 0.1%) 48 to 72 hours before using ONE, you may resume use 48 to 72 after.

## Getting to Know Your DERMAFLASH Device



## DERMAFLASH 3-Step Process: PREFLASH, FLASH, POSTFLASH

The DERMAFLASH 3-Step process instantly unveils a smoother, radiant complexion, and improved tone and texture as it prepares, exfoliates, and soothes the skin.

1. Your ONE includes:
  - ONE Device
  - Charging Base
  - USB Cable
  - Tray(s) of single-use Edges
  - PREFLASH cleanser (volume fill varies)
  - POSTFLASH moisturizer (volume fill varies)
2. The ONE facial exfoliator uses subtle sonic vibration and a specially designed single-use exfoliating Edge to help remove dead skin cells, debris and vellus hair.
3. For best results, ONE should be performed once a week.

**CAUTION:** We do not recommend using ONE more than once a week. The whole treatment takes approximately 5 to 10 minutes. For safety and ultimate efficacy, a new Edge must be used for each facial treatment.

## Charging

1. Before using, always make sure you use the correct plug adapter which corresponds to the country you are in. Insert large end of USB charging cord into a power source (like a power adapter or computer) and small end into charging base.
  2. Place silver side down, without Edge, into charging base. Clear side of device to face up. A full charge allows up to 30 minutes of use.
- IMPORTANT!** Any used Edge must be removed before placing in the charging base.
3. While your ONE is charging, the LED charging indicator will slowly pulse. When fully charged the light will be solid.
  4. To alert you the battery is low, the LED charging indicator at the top of your ONE device will begin to pulse during use and may shut off. If this occurs, return your ONE to the charging base.
  5. If your ONE stops while in use, follow the detailed instructions described in this owner's manual, remove the Edge, then return the ONE to the charging base.

## Universal Charging Base For World Wide Use

Your ONE is rechargeable with its included Charging base and USB cord. The USB cord can be plugged into a power adapter or computer. Power adapter not included.

### READ ALL INSTRUCTIONS BEFORE USE

#### Before Using for the First Time

1. Your ONE is rechargeable with its included Charging base and USB cord. Place your ONE in its charging base overnight to fully charge.

**NOTE:** A full charge allows up to 30 minutes use.

#### Loading A New Edge

**WARNING!** Edges are single-use only. NEVER attempt to load or exfoliate with a previously used Edge. Discard and recycle all Edges after each use.

**WARNING!** Edge is sharp! To avoid personal injury, NEVER attempt to load the Edge into the ONE by hand!

**WARNING!** When loading an Edge in device, be sure to place Edge tray on flat surface before gliding device onto Edge.

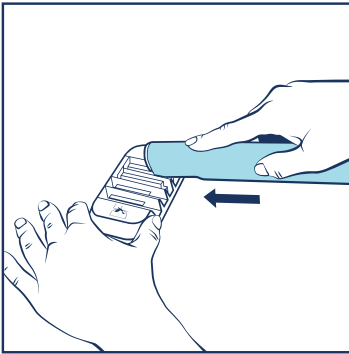


Figure 1

1. Place the Edge tray on a clean, flat, stable surface. To load, slide ONE into any unused Edge in the Edge tray (See Figure 1).

An audible click can be heard as the Edge locks securely into place.

2. Gently lift the device straight up from the Edge tray.

**WARNING!** NEVER slide the ONE BACKWARD in the Edge tray. Doing so may damage the Edge and make it unusable.

3. Once the Edge is securely loaded into the ONE device, press the Power button to turn your ONE ON. The light will illuminate. You will begin to feel a subtle, sonic vibration. This is normal.

## DERMAFLASH Treatment Step 1: PREFLASH

1. Always remove all makeup before cleansing the skin.

2. Wet your face and hands. Using a dime-sized amount of PREFLASH cleanser, create a lather in your hands and massage onto face.

3. Rinse thoroughly with water. Pat skin dry.

**CAUTION:** PREFLASH is for external use only. Avoid contact with eyes.

## DERMAFLASH Treatment Step 2: FLASH

1. The DERMAFLASH treatment takes about 5 to 10 minutes. For safety and ultimate efficacy, a new Edge must be used for each facial treatment.

**WARNING!** To avoid personal injury, NEVER hold your DERMAFLASH at a 90 degree angle or perpendicular to your skin.

2. Before using your ONE, make sure your hair is pulled away from your face, secured with clips or an elastic hairband.

**NOTE:** ONE is intuitively designed to lie on the skin at a 45 degree angle (See Figure 2.)

3. Use the ONE in patterns suggested in the illustration (See Figure 3.)

4. You can DERMAFLASH your cheeks, your entire chin area (including just below chin), above and below your lips and your forehead. DERMAFLASH each area while supporting your skin with a finger from your other hand.

5. Use short feathery strokes. Begin at hairline just in front of the ear and move inward following the contours of your face.

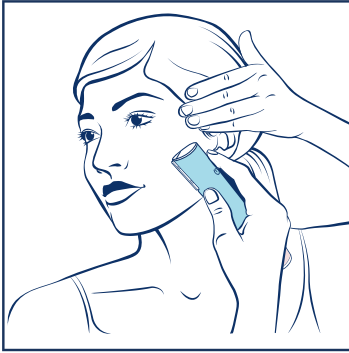


Figure 2

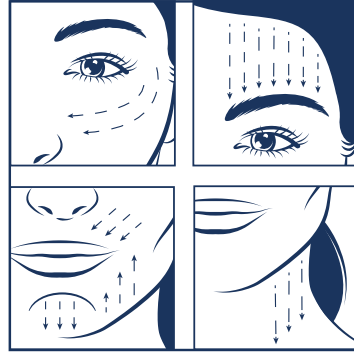


Figure 3

### Forehead Area

6. NEVER USE ONE DIRECTLY ON YOUR NOSE.
7. To avoid unintentional hair loss, NEVER use ONE directly on the brow or close to your eyebrows, sideburns, or anywhere inside hairline.
8. Make short, feathery, gentle downward strokes to the area above your eyebrows.

### Eye and Cheek Area

9. NEVER USE ONE INSIDE THE EYE AREA, OR CLOSE TO EYELASHES OR EYELIDS.
10. Begin the treatment in front of your ear at the top of the cheekbone. With your free hand holding your skin taut and using short, feathery strokes, glide ONE inward towards the nose, over the contours of your cheek and entire side of your face. Take care to avoid the eye socket area and nose.

### Above and Below the Lip and Chin

11. NEVER USE ONE DIRECTLY ON YOUR LIPS.
12. Above the top lip: Begin at the outer corner of the nose area, make short downward strokes following the contours of the cheek down to top lip.
13. Below the bottom lip: Working in short upward strokes, begin at the outside corner of the mouth and work toward the center. (See Figure 3.)
14. Chin: Begin at the center of the chin. Using short downward strokes work your way over the jaw line.
15. When you are finished exfoliating your face, press the ON/OFF button to turn your ONE OFF. The LED light will turn OFF.

### Used Exfoliating Edge Removal

1. To remove your Edge, direct your ONE toward a trash receptacle and press the Eject button.
2. Place your ONE in the charging stand/cradle with the used Edge removed, silver side down. Recharge before your next treatment.

**REMEMBER!** For safety and hygiene purposes, Edges are designed for single-use only.

**WARNING!** KEEP EDGE TRAY AND USED EDGES AWAY FROM CHILDREN.

**WARNING!** The Edge is sharp! Never handle the Edge portion of this part with bare hands.

**NOTE:** All packaging materials and certain components of ONE are recyclable

### DERMAFLASH Treatment Step 3: POSTFLASH

1. Gently massage a small amount of POSTFLASH onto your face.
2. After using ONE you may experience a temporary redness or tingling after treatment. This should be considered normal and will subside.

**CAUTION:** POSTFLASH is for external use only. Avoid contact with eyes.

### Storage

1. Make sure device is unplugged.
2. Store your ONE in the charging base with the used Edge removed.
3. Store device in its original box or in a clean, dry area.
4. Never wrap cord tightly around charging base.
5. Do not put any stress on cord where it enters charging base, as this could cause cord to fray and break.

### Environmental Protection/Disposal

1. All packaging materials and certain components of your DERMAFLASH are recyclable. Please dispose of recyclable materials responsibly.
2. Information on recycling can be obtained from your local authority.