

GEEMARC

# CLEARSOUND®

HEAR THE DIFFERENCE

## CL100

English



!!!User manual to read carefully before start to use!!!

!!!Manuel d'utilisation à lire impérativement avant la mise en service!!!

!!! Bevor Ihr neues Telefongerät zu benutzen, bitte lesen Sie dies Benutzerhandbuch sehr sorgfältig!!!


**CAUTION:** this telephone is specifically designed for the people with hearing difficulties and is provided with a handset earpiece volume control for individual requirements. Due care must be taken by all users that the handset earpiece volume control is set to the lowest level acceptable by each respective user. Care should therefore be exercised to ensure that any naive users are familiar with this requirement.


**ATTENTION:** Ce téléphone peut être destiné aux personnes ayant des problèmes auditifs. il possède un écouteur a volume variable. A la première utilisation, assurez vous que le volume soit réglé au minimum. Les personnes malentendantes pourront par la suite augmenter le volume pour l'adapter à leur besoin.

**WARNUNG:** Dieses Telefon wurde speziell für schwerhörige Personen konzipiert und besitzt eine auf die individuellen Bedürfnisse einstellbare Hörerlautstärke. In Folge dessen sollten alle Nutzer darauf achten, dass der Lautstärkeregler des Hörers, in Rücksicht auf andere Nutzer, auf das niedrigste Niveau gestellt wird. Es ist darauf zu achten, dass unwissende Benutzer mit der Handhabung dieser Apparatur vertraut gemacht werden.


<b>Description</b>	<b>4</b>
<b>Setting up</b>	<b>7</b>
- Connecting the telephone	7
- Ringer volume setting	7
- Tone/Pulse switch	8
- Time break recall setting	8
- Wall mounting	8
<b>Volume and tone adjustment</b>	<b>10</b>
- Receiving volume and tone adjustment	10
- Sending volume adjustment	12
<b>Using the phone</b>	<b>13</b>
- Answer a call	13
- Making a call	13
- Last number redial	14
- Recall, * & # buttons	14
- Mute button	14
- Storing telephone numbers	15
- Dialling stored number	16
<b>General informations</b>	<b>17</b>

# Description

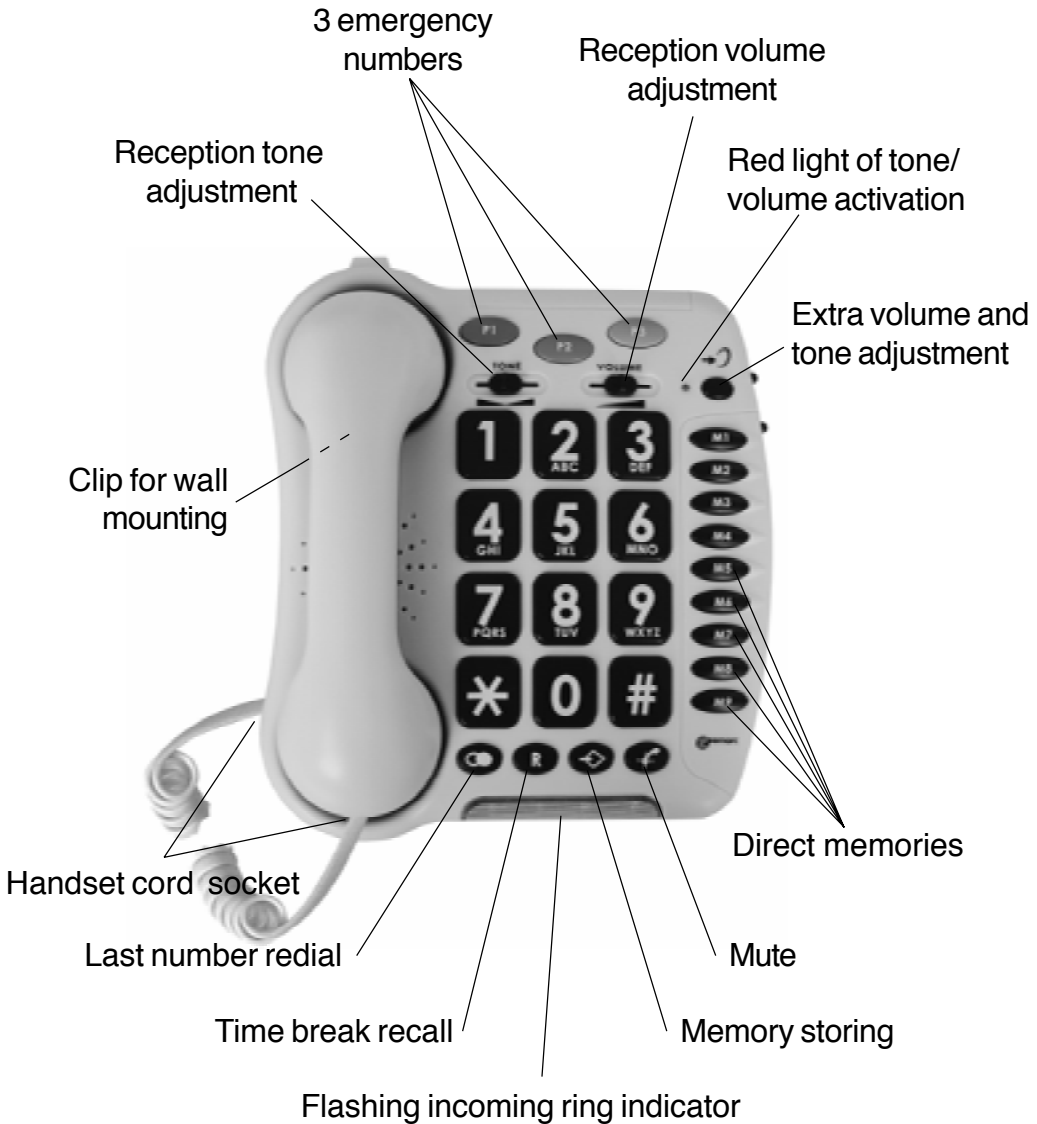

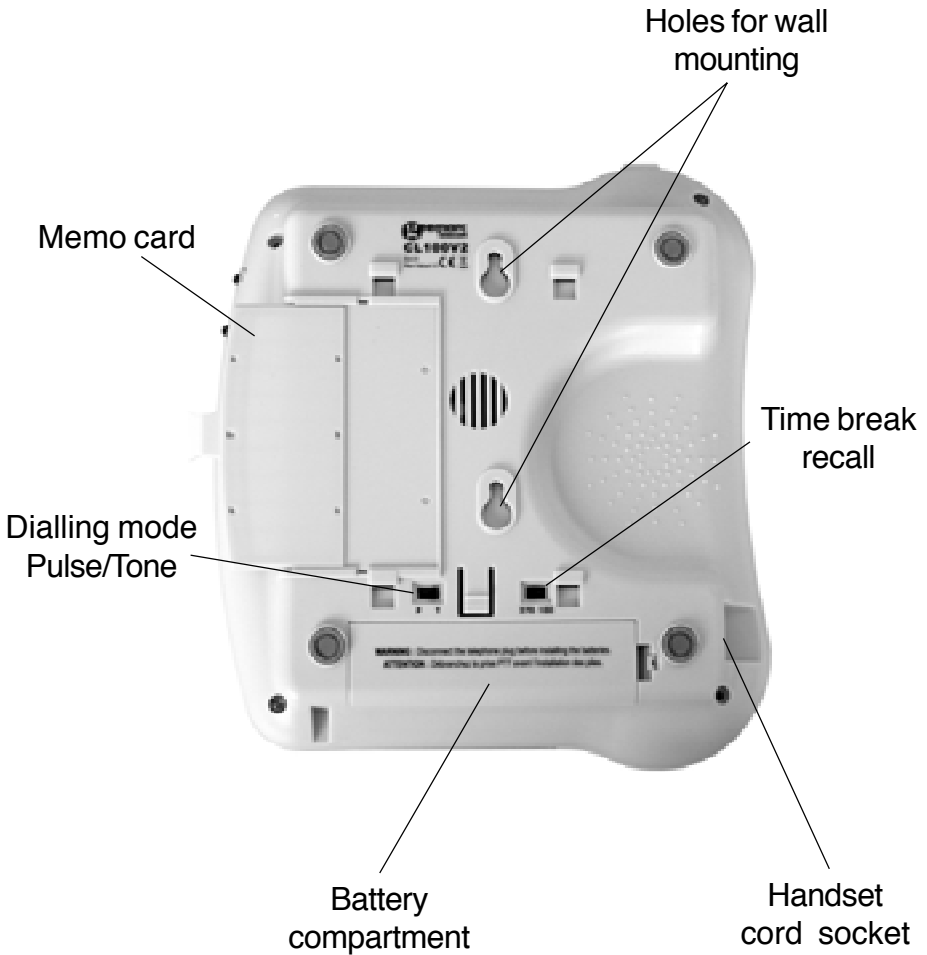


# Description

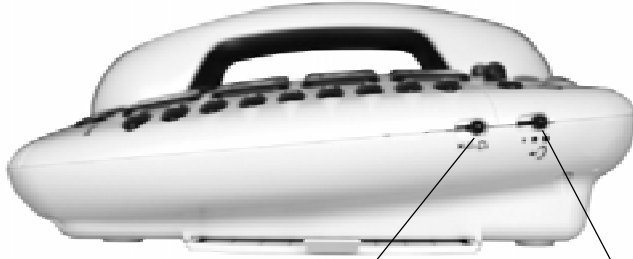

Ring volume  
setting

Sending volume  
setting

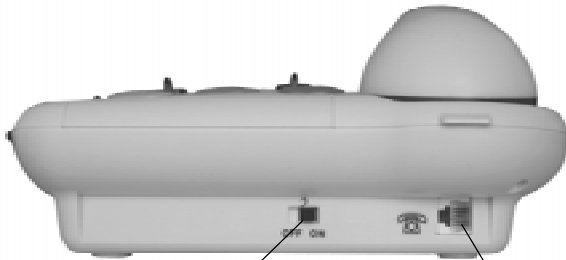

Volume and tone  
adjustment ON/OFF switch


Line socket

## 1 - Connecting the telephone


- Connect the handset cord.
- Open the battery compartment.
- Insert 4 alkaline batteries (type AAA 1.5V, not supplied) in the compartment (\*). Batteries are required if you wish to increase the ringer level up to 10 dB.
- Snap the battery door back in place.
- Connect the line cord plug into the line socket on the rear of the unit (\*).
- Connect the other end of the line cord into the wall PTT socket (\*).

## 2 - Ringer volume setting

A switch located on the right side allows the ringer to be switched to Hi, Low or Off:

 : Off,

 : Low,

 : Hi

**Note:** *The ring level will be increased of 10 dB if you insert batteries.*

## **3 - Dialling mode (Tone/Pulse) settings**

The Tone/Pulse switch is factory pre-set to Tone position. In the UK, all telephone exchanges now use Tone (T) dialing. If your phone does not dial out, it is probably being used from an older private switch board (PBX). In this case, the Tone/Pulse switch must be changed to Pulse (P).

## **4 - Time Break Recall setting (270/100)**

Your phone works with a time break recall for United Kingdom (R=100ms). It can be switched to 270ms for other locations or PBX's.


## **5 - Wall mounting**

To wall-mount your telephone, turn the "clip" (located in the base unit, just below the ear piece).

This will ensure that the handset stays firmly in position.




Drill 2 holes in the wall (80mm apart).  
Insert the wall plugs and the screws leaving about  
6 mm of the screws projecting from the wall.  
Slide the telephone downwards to secure.

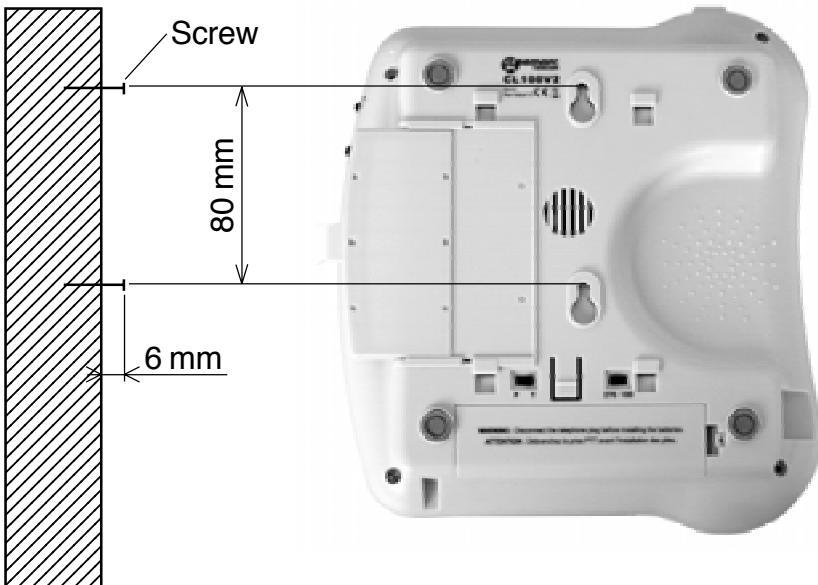


# Volume and tone adjustment


## 1 - Receiving volume and tone adjustment


You can increase the volume and adjust the receiving tone according to your own hearing requirements.


The button  allows extra receiving and tone functions to be activated or not during the conversation.


When the function is activated, the red light of extra receiving and tone function located near the button  is lit.


Adjust the **TONE** and **VOLUME** buttons to fit the levels to your own hearing requirements.

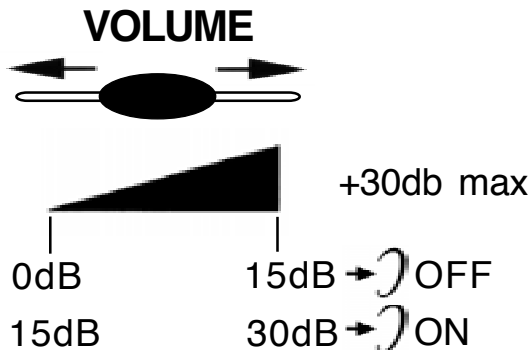

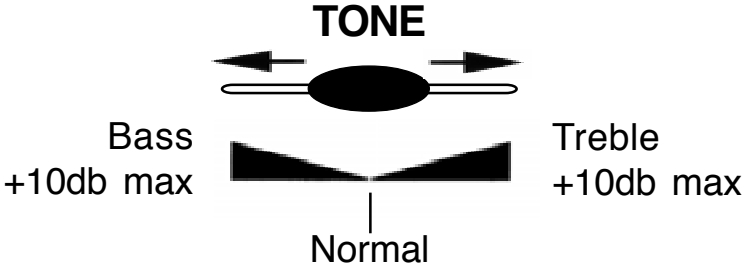


**Note:** When the tone (treble), earpiece volume and sending levels are set to the maximum, whistling may occur. Then the sending level settings should be decreased with the switch .


When the switch  located behind the phone is **ON**, the extra-receiving and tone functions are automatically activated each time you use the phone.

When this switch is **OFF**, you will need to activate the function by pressing the button .

Tone and volume (15dB to 30dB) adjustment is only activated when the button  is set to **ON**. When this button is set to **OFF** only the volume can be adjusted from 0dB to 15dB.

## 2 - Sending volume adjustment

If your voice is quiet or called parties have difficulty hearing you, you can increase the sending level volume with the switch located on the right side. You can also decrease the sending volume if required.

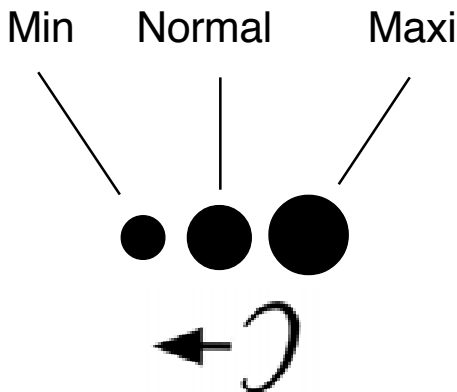

***If you are using this telephone with a hearing aid, set your hearing aid to the “T” setting.***


## **1 - Answer a call**

- When an incoming call is received the phone will ring and flash, to answer the call lift the handset and speak.
- On completion of the call, carefully replace the handset in the cradle to release the line.

## **2 - Making a call**

- Lift the handset.
- Await the dial tone and dial required number.
- On completion of the call, carefully replace the handset in the cradle.



## **3 - Last number redial**

- Lift the handset.
- Await the dial tone and press  button.
- The previously dialled number will be automatically redialled (not in use for memory numbers).

## **4 - Recall (R), \* & # buttons**

These are used with the new services provided on digital exchanges. For details please contact your network operator.

## **5 - Mute button**

If you wish to speak privately to someone else in the room without your caller overhearing, press and hold the  button. You will still be able to hear the callers but they cannot hear you during the mute operation. To resume normal conversation, release the  button.

## 6 - Storing telephone numbers



Your telephone can store 12 direct memories including 3 priority numbers (24 digits maximum for priority memory and 32 digits for M1, M2...M9). For example, the 3 priority **P1**, **P2**, **P3** memories allow you to store the 3 numbers below:

**P1**=Fire

**P2**=Police

**P3**=Doctor

To store a number please follow instructions below:

- Lift the handset.
- Press  button.
- Dial the number.
- Press  button.
- Press **M1**, **M2**, ..., or **M9** or one of 3 priority memories to store your number.

**Note:** *Entering a new number in memory will automatically erase the previous one.*

## **7 - Dialing stored numbers**

- Lift the handset and await the dial tone.
- Press appropriate button **M1**, **M2**, ..., or **M9** or one of 3 priority memories, the stored number will be dialled automatically.

To register your stored phones numbers, use the index located under the base.



# General informations

## Recycling directives:

The WEEE (Waste Electrical and Electronic Equipment) has been put in place for the products in end of life are recycled in the best way.



When this product is out of order, please do not put it in your domestic waste bin.

Please use one of the following disposal options :

- Remove the batteries and deposit them in an appropriate WEEE skip.

Deposit the product in a appropriate WEEE skip.

- Or hand the old product to the retailer. If you purchase a new one, they should accept it.

Thus if you respect these instructions you ensure human health and environmental protection.

## Guarantee:

From the moment your Geemarc product is purchased, Geemarc guarantee it for the period of one year. During this time, all repairs or replacements (at our option) are free of charge. Should you experience a problem then contact our helpline or visit our website at [www.geemarc.com](http://www.geemarc.com). The guarantee does not cover accidents, negligence or breakages to any parts. The product must not be tampered with or taken apart by anyone who is not an authorised Geemarc representative. The Geemarc guarantee in no way limits your legal rights.

**Important:** YOUR RECEIPT IS PART OF YOUR GUARANTEE AND MUST BE RETAINED AND PRODUCED IN THE EVENT OF A WARRANTY CLAIM.

**Telephone connection:** Voltages present on the telecommunication network are classified TNV-3 (Telecommunication Network Voltage) according to the EN60950 standard.

The telephone does not operate if the line current is lower than 18 mA.

**Declaration:** Hereby Geemarc Telecom SA declares that this product is in compliance with the essential requirements and other relevant provisions of the Radio and Telecommunications Terminal Equipment Directive 1999/5/EEC and in particular article 3 section 1a, 1b and section 2.

*The declaration of conformity may be consulted at [www.geemarc.com](http://www.geemarc.com)*





UGCL100V2\_Ver:03

**geemarc™** United Kingdom  
Telecom SA

5 Swallow Court  
Swallowfields  
Welwyn Garden City  
Hertfordshire, AL7 1SB  
For product support:  
01707 384438  
[www.geemarc.com](http://www.geemarc.com)

**geemarc™** FRANCE  
Telecom SA

Parc de l'Etoile  
2, rue Galilée BP20061  
59791 Grande Synthe CEDEX  
Tél. service après vente :  
03 28 58 75 99  
[www.geemarc.com](http://www.geemarc.com)

**geemarc™** Deutschland  
Telecom SA

Siemensstrasse 10,  
53121 Bonn  
Germany  
Tel: +49(0)228/ 74 87 09 0  
Fax: +49(0)228/ 74 87 09 20  
[www.geemarc.com](http://www.geemarc.com)