


owner's manual
medvibe mini massage ball

medmassager.com
189 Continuum Drive, Suite A, Fletcher, NC 28732

charging method

1. Connect the power adapter charging plug to the product charging jack on the side of the massage ball.
2. Plug the power adapter into a 100-240V, 50/60Hz power socket/outlet.
3. The red LED light will flash during charging.
 - *If the massage ball is continuously charged for three hours, and the battery indicator lights are still not fully illuminated (indicating the MedVibe Mini is fully charged), please unplug the power adapter and contact MedMassager™ for further assistance.
 - *The massage ball cannot be on and in use while charging.
4. After the charging is completed, the LED light will continue to illuminate to indicate that the battery is fully charged.
 - *It will take about 1 1/2 hours for the battery to fully charge.
5. When charging is completed, disconnect the AC adapter from the wall and charging plug from the massage ball.

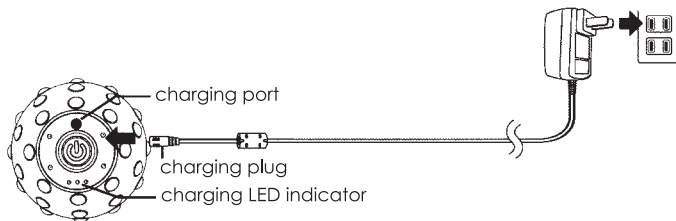


vibration modes

The MedVibe Mini Massage Ball offers three vibration settings and one cycle mode. You can make a selection by pressing the power button.

- Press the power button to turn on the massager. Once turned on, the automatic vibration starts at the lowest setting. One intensity indicator light will be displayed.
- To increase the vibration speed, press the power button again. Two intensity indicator lights should be displayed.
- Press the power button a third time to enter into cycle mode. The vibration speed will automatically rotate through the three intensity levels. Three intensity indicator lights will be displayed.
- Pressing the power button for a fourth time will change the vibration speed to the highest level and four intensity indicator lights will be displayed.

shut down

- Pressing the power button for a fifth time will turn off the massager.
- To shut down from any mode, press and hold the power button for two seconds.

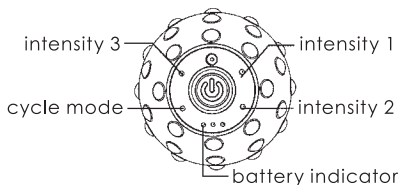

The vibration function will automatically shut-off after 10 minutes. Using this massager for longer periods of time may result in injury.

maintaining your device

properly cleaning your device

- Before cleaning the product, make sure the power switch is turned off and unplug the power cord.
- Wipe the device with a clean, soft cotton cloth. **DO NOT** use corrosive cleaning agents such as benzene or thinners to clean the product.
- **DO NOT** wipe the switch control surface with a damp cloth to avoid machine malfunction. Wipe it with a dry cloth.
- **DO NOT** dip the product into water or other liquids. If water is spilled on the product, it may cause electric shock.

storage method

- To keep your MedVibe Mini Massage Ball clean, keep the device away from high humidity.
- Ideally, you should store the device in its original packaging and store in a clean, dry place.
- **DO NOT** place your device in direct sunlight or store in high temperatures.
- If not used for a long time, cover with a dust cover and remove the power plug from the power outlet.

correct handling



This symbol indicates that this product cannot be disposed of with household waste in any country. In order to avoid possible damage to the "environment" or human health, the device should be recycled to increase the reusability of the remaining parts. Recycle used equipment, follow the recycling process, or contact MedMassager.™ They can recycle products for environmental safety.

troubleshooting

problem	troubleshooting
The massager does not vibrate even if the power button is pressed.	<ul style="list-style-type: none"> • The massager and power adapter are still connected. • The battery needs to be recharged. Note: The vibration function cannot be used during charging.
Vibration intensity has weakened.	<ul style="list-style-type: none"> • There is insufficient power and the battery needs to be recharged.
Vibration only lasts a few minutes, even after a full charge.	<ul style="list-style-type: none"> • Extremely low or high ambient temperatures can cause the product to malfunction. Please use it within the recommended ambient temperature range of 50° F to 104° F. • The battery needs to be replaced. Consult MedMassager.™ Note: DO NOT replace the battery yourself. Doing so may result in damage to the product.
The massager suddenly stopped working.	<ul style="list-style-type: none"> • The product stops automatically after 10 minutes.
There is an abnormal noise during use of the massager.	<ul style="list-style-type: none"> • The product will emit an action sound when vibrating.

If you have tried all troubleshooting options and nothing is working, please contact MedMassager™ for further assistance.

specifications

product name	medvibe mini massage ball
product number	JF.005
battery specification	DC 7.4V 1200mAh
charging time	approx. 1 1/2 hours to fully charge
quiescent current	less than 80uA
full discharge time	50 minutes or 5 cycles of 10 minutes
vibration revolutions per minute	1st gear motor speed: 2000RPM 2nd gear motor speed: 2500RPM 3rd gear motor speed: 2100-3600RPM 4th gear motor speed: 3600RM
auto shut-off time	10 minutes
ambient temperature	50° F to 104° F
product weight	.3kg or 10.6 ou.
material	silica gel
product accessories	power adapter, charging plug and manual
manufacturing origin	China

specifications

power specifications

Input	AC 100-240V, 50/60 Hz
	DC9V to 1000mA

Executive Standard: GB4706.1-2005/GB4706.10-2008

MedMassager™ reserves the right to modify the product and the manual.
Notice that the appearance of the product is subject to the actual product.



medmassager.com

877.501.6334

1.828.475.6022 (International)

Sales@MedMassager.com

189 Continuum Drive, Suite A, Fletcher, NC 28732