KAHUNA  MASK

KAHUNA  **PREMIUM MASK V2**

- ENGLISH
- KOREAN
- CHINESE
-

CERTIFICATION

KAHUNA PREMIUM MASK V2

High standard of quality in manufacture

Premium NIR LED MASK company


100% Korean LED Brand

100% MADE IN KOREA

LED NEAR INFRARED RAY

Miracle of skin regeneration discovered by NASA!

Why Near Infrared Ray is so special?

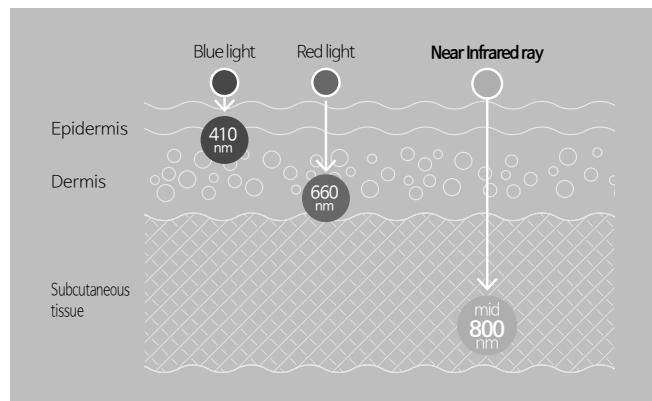Technology used by NASA

LED skin care treatment was researched and developed by NASA to treat astronauts' skin damage caused by space flight.

According to NASA research, 'LED near-infrared light is effective for the recovery and regeneration of skin such as burns, trauma and skin damage caused by strong sunlight.'

Penetration depth

NIR is the only wavelength that can reach up to 4.5mm of dermal layer and subcutaneous tissue.

KAHUNA PREMIUM MASK V2

KEY Points!

SKIN & SCALP

Near-infrared total care solution for skin and scalp!

Skin aging begins with top of head, and facial elasticity must be preceded by scalp care. Now, treat your face and scalp at once with effective Near Infrared Ray Therapy-

Ergonomic design

The ergonomic curved mask is designed to cover the entire face from the chin to the top of the head.

Moreover, the black translucent part of the mask is the material with the highest LED transmittance and is designed to satisfy both design and effect at the same time.

Patented 3 Staged Care System

Since each person's skin condition and thickness are different, staged light intensity control care system enables customized skin care from sensitive to normal skin.

Scalp improvement

Improve skin brightness and tone

Improved skin elasticity and density

Skin texture improvement

Skin moisturizing

Optimized light exposure

Applied optimized near-infrared light dosage for the growth of new collagen fibers and synthesis, and to remove aging collagen fibers.

The uniqueness of the near infrared mask

Unlike red and blue light, near infrared ray has no dazzling effect and has been reported as a harmless light to the eyes.

Premium custom-made near infrared chip optimized for skin improvement guarantees certain effect.


YOUTH ZONE THERAPY

Independent near-infrared therapy on each skin zones with different characteristics

CLINICAL TEST

Effective verification by professional institutions

- Test subjects : Women aged 30–55 who meet the selection criteria / EAST SKIN LED Shine Mist / Scalp Mist
- Test Agency: P & K Skin Research Center Co., Ltd.
- Trial period: 4 weeks

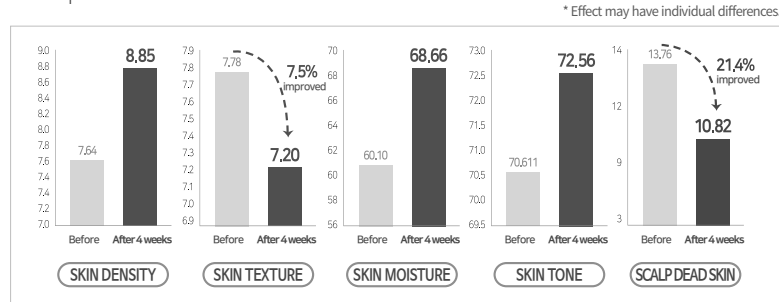KAHUNA PREMIUM MASK V2 Competitiveness:


NIR YOUTH THERAPY

Providing independent near-infrared ray therapy optimized for thin, thick and sensitive skin with different characteristics.


STAGED LED SKIN CARE SYSTEM

- KAHUNA PREMIUM MASK V2 is the one and only LED mask in the market with staged LED skin care system.
- You can control the intensity of the LED near-infrared ray power according to your skin type as well as your skin condition.
- Stage mode1 is recommended for LED mask initial users, thin skin and those who are sensitive. Allow some time for skin to adapt to near-infrared care to experience greater improvement.
- We recommend that you use MODE 2, MODE 3 sequentially, after your skin is adapted to MODE1.

Component Description

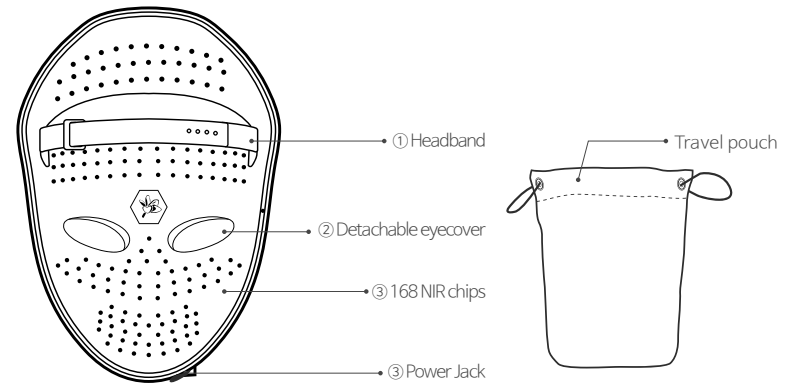Package Components

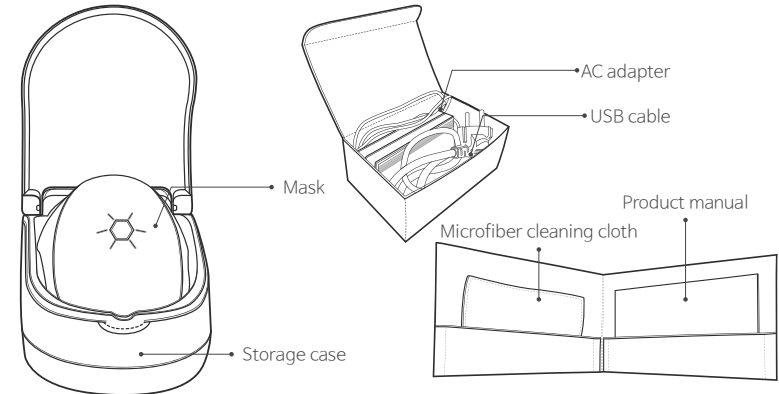Specifications


Number of LEDs	168	Input Voltage	5V/2A
Dimensions	30cmx20cmx15cm	Weight	650g

Power connection

- Connect the usb cable or plug the power adapter to start the mask.
- There will be a "beep" sound on proper connection.

Instruction

- POWER ON & MODE1 : To begin MODE1, touch and hold the sensor on the forehead for 2 seconds.
- Mode change : Touch and hold the sensor for 2 seconds to select the desired mode.
- POWER ON & MODE1 ► MODE2 ► MODE3 ► OFF

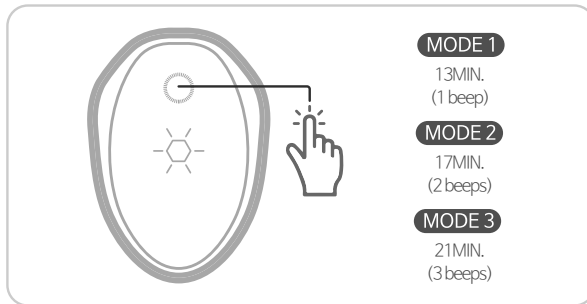Steps to Use

Wash your face.

Apply skin care product.

Select mode and wear the MASK.

Stay comfortable until power off.


KAHUNA PREMIUM MASK V2 uses only near-infrared single wavelength. We recommend to use on a daily basis right before bed: choose the mode, lie on your back and close your eyes while its use

How to use mask / Mode change

STEP 1

Connect USB cable and adjust headband 2 or 3 notches loosened than head size.

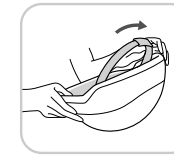STEP 2

Raise the headband all the way up.


STEP 3

Wear headband facing down.

STEP 4

Put the mask on your forehead, hold down the headband and pull down the mask.

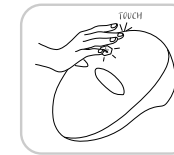STEP 5



Touch and hold the sensor for 2 seconds.
* Select the mode according to your skin type and condition.

STEP 6

After choosing the mode, you may sit or relax on your bed.
*The eye cover can be detachable on your preference.

How to Detach Eye Cover

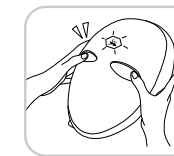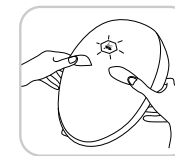To detach the eye cover, force push the back portion of the cover. Check the eye cover's inner side labelled with 'left(L)' and 'right(R)'. Eye cover is detachable so you may use it to your preference. Please keep the eye cover in the safe place to prevent loss.

Step 1. Installation Guide

Search "HiddenTag" and install the app from below app markets.

Download "HiddenTag" for free ▶

Step 2. User Guide

1. Let's verify the authenticity quickly and easily in 3 seconds!

① Check HiddenTag label.

② HiddenTag scan.

③ Check the authenticity.

2. User registration.

① User registration.

② User information registration.

③ Genuine product is successfully registered.

Warnings

Time and frequency of Usage

- People who had facial plastic surgery or cosmetic treatment (e.g. filler or Botox) should consult with doctor first about the usage of near-infrared device, and recommended to wait for the skin condition to fully recover.
- Patients with heart disease or hypertension should consider health status and adjust their usage time and frequency accordingly.
- Using multiple times per use is NOT recommended. To minimize the stress on your skin, maintain single use each day on a regular basis.

Important Safety Warnings

- DO NOT drop the product on floor.
- DO NOT disassemble, self-repair or alter the product.
- DO NOT attempt to clean the product with water, alcohol or benzene.
- DO NOT place the product or power adapter near flammables.
- DO NOT place the product in damp or dust.
- KEEP AWAY from vibrations and magnets.
- DO NOT keep in direct sunlight or expose to heat. (e.g. heater)
- DO NOT place heavy objects on the product.
- DO NOT damage or twist the power plug.
- This product emits near-infrared, which induce slight heat sensation from the device, and the sensation level may differ according to user's body temperature or surrounding environment, including time of season.
- The heat may intensify with continuous operation, so use intermittently.
- Defective power adapter may lead to damage to the product or even explosion.
- Proper use of power adapter will provide longer service life.
- Power adapter's service life may decrease when left unused for extended period or in extreme environment (e.g. below minus 10 degrees Celsius)
- DO NOT use in a flammable environment.
- Using electronic devices in atypical situation may result in fire or explosion.
- If there is smoke or burning smell from the device, shut off and unplug the device and contact the retailer or the head office.
- Automatic temperature control device: It has the function to control the internal temperature according to the external environment, and there may be a slight time difference in use.

Customer Support

- For technical support and repair inquiries, please contact the retailer or use the corporate website.
- The warranty will be void in any of the following events:
 - Damage results from usage other than its normal and customer manner.
 - Damage results from misuse or negligence caused to product.
 - Damage from spillage of toxic or corrosive substances.
 - Damages from accident, unauthorized disassembly, repair, alteration or natural disaster.
 - Warranty period expires.
- If you have dermatitis, sensitive skin or skin allergies, please consult your doctor before use, and stop use and ask a doctor when PLEASE CONSULT YOUR DOCTOR BEFORE USE, and in case of redness on face, itching, or suppurative inflammation, STOP USE AND CONSULT YOUR DOCTOR.
- Please refer to the product manual for all other information.

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DISCLAIMER

In order to achieve the optimal performance and to prevent any damage or misuse of the KAHUNA PREMIUM MASK V2, please read all the contents of the user manual carefully before using the KAHUNA PREMIUM MASK V2. Any amendments of this manual due to typographical errors or inaccurate descriptions shall be made only by the KAHUNA manufacturer. The KAHUNA manufacturer does not assume any liability for the user's personal and property losses caused by the operations not in accordance with this manual. This user manual does not include any form of warranty, subjective expressions or other implied information. The KAHUNA manufacturer and its staff do not assume any liability for any direct or indirect personal and property losses caused by descriptions included in this manual.

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KAHUNA MASK

Address: AJX INC 11125 Knotts Ave. Cypress, CA 90630

Tel: +714-388-2625 Fax: +714-379-6164 Email: info@kahunachair.com