


## Using TePe Interspace™

TePe Interspace™ is a tapered single pointed tufted brush with removable brush tips to connect to either side of the handle, allowing better access. This brush is great for furcations, behind the last tooth, around implants, fixed braces, or any hard-to-reach area. It is important to use and care for your TePe Interspace™ properly—please read below for instructions and helpful tips.

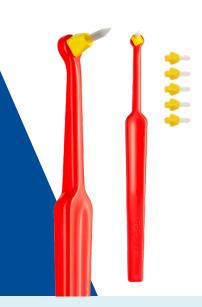

## 5 Simple Steps to Using Your TePe Interspace™:

- 1. Locate the area you want to clean.
- 2. Look in the mirror and decide what angle would work best for you to access the hard-to-reach area.
- 3. Decide which side of the handle you would like the brush tip.
- 4. Screw the brush tip onto the desired side.
- 5. Gently move the bristles back and forth for plaque disruption.

Scan the QR code to watch a video about TePe Interspace™



## Tips for using your TePe Interspace™

- This brush comes with 12 replaceable brush tips.
- The brush tip can be attached to either side of the handle, offering different angles for better accessibility.
- Seristles textures are either soft (yellow) or medium (blue).
- Angle the neck in hot water to adjust to your liking.
- Replace brush tip when bristles start to fray or become worn.







