




## The Rejuvenator

All-In-One Beauty System with Red LED Light Beauty Wand & Activating Lifting Collagen Boosting Serum

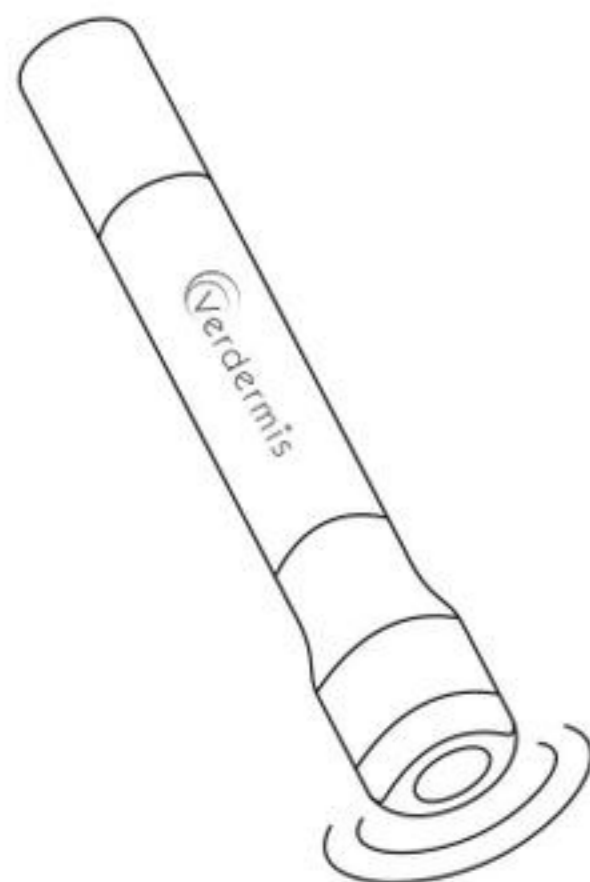

### USER MANUAL

Please read this manual carefully before use and keep it properly for future reference


#### How to Use

**Step 1**—Remove the cap at the top end of the wand. Over clean skin, apply the serum over the face, neck, around the eyes, and other areas of concern. Let the serum absorb. Replace the Cap.  
**Step 2**—apply a 5-10 minutes of Red LED Light Therapy and Vibration over the same areas of concern by gliding the other end of the wand. Note- this beauty tool operates on "touch" technology and will work only when you hold it (like a pen) where your fingers touch the Finger Conduction area (Silver Area) and your skin is in touch with the Skin Conduction area - Silver Ring at the end of the wand where the red light is. (see image below). Enjoy a youthful, revitalized appearance. Use as needed.

#### Notes:

- \*For all skin types.
  - \*Recommended application 5-10 minutes at a time.
  - \*Battery operated - no charge needed.
  - \*Replace batteries when they deplete.
  - \*Serum container can be refilled at any time. if serum is out you can still use the tool for a light therapy.
- TOOL OPERATES ON TOUCH TECHNOLOGY TO PRESERVE ENERGY. TO OPERATE YOU HAVE TO HAVE YOUR FINGERS TOUCH THE FINGER CONDUCTION AREA AND YOUR SKIN HAS TO TOUCH THE SKIN CONDUCTION AREA (END SILVER RING OF THE WAND WHERE THE RED LIGHT IS) - SEE IMAGE BELOW.

#### Product Structure

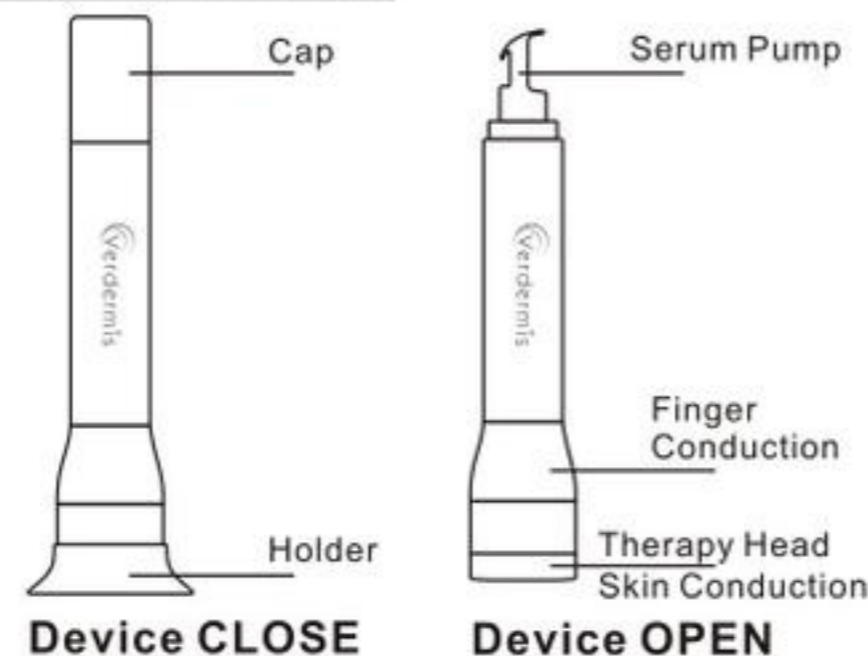

#### Product Specifications

Product	The Rejuvenator	Device Size(CM)	D5*17.5
Model No.	DG7	Rated Current	3V 1A
Weight	68g	Vibration RPM	15000

#### Charging Method & Battery Replacement

This product is operated by 2X3V lithium batteries (included). Once the batteries are depleted, exchange them by turning the bottom end on the wand (where the light is) counterclockwise and replace with new batteries. See pictures below.



#### Maintenance

1. If applying the serum before the light treatment, spread the serum evenly and let it absorb before using the tool.
2. Please wipe the product clean with a dry towel after each application.
3. This product is not waterproof or water resistant, do not rinse or wash with water or other liquid. Keep this product clean and dry.
4. Please don't use corrosive to clean the device.
5. Please store it in a cool dry place.
6. Please don't use the device if it is damaged.

Frequent Questions	Possible Causes	Solutions
The product doesn't turn on or does not operate	<ul style="list-style-type: none"> <li>*Depleted batteries</li> <li>*Batteries don't set up well</li> <li>*Circuit doesn't close because the user doesn't hold the product correctly</li> </ul>	<ul style="list-style-type: none"> <li>*Replace batteries with new ones</li> <li>*Reset batteries (change sides)</li> <li>*This tool operates on "touch" technology and will work only when your fingers touch the Finger Conduction area and your skin touches the Skin Conduction area (see image).</li> </ul>
The vibration is not strong or inconsistent	<ul style="list-style-type: none"> <li>*Batteries getting depleted</li> </ul>	<ul style="list-style-type: none"> <li>*Replace batteries with new ones. Set up correctly</li> </ul>
The Serum doesn't dispense	<ul style="list-style-type: none"> <li>*Serum is out / Serum container is empty</li> </ul>	<ul style="list-style-type: none"> <li>*Replace with a new cartridge of serum OR use the tool without serum</li> </ul>

#### Caution

1. If you have skin disease, skin issues, or a very sensitive skin please consult a dermatologist before use. If you are pregnant, have epilepsy, tumors, cancer, patients with cardiac pacemakers, implanted with cardiac defibrillators or other implanted devices - please consult your physician before use.
  2. Avoid direct contact with the eyes.
  3. For first time users or people with sensitive skin: Apply for a short duration to ensure there are no sensitivities.
  4. For hygiene reasons, please do not share this product with others.
  5. Place the device out of the reach of children to prevent accidents.
  6. Do not use corrosive chemicals to wipe the product.
  7. Please do not disassemble this product (except for changing batteries or replacing a new serum cartridge).
  8. For external use only.
- IN THE RARE CASE YOU HAVE ANY DISCOMFORT WHILE OR AFTER USING THIS PRODUCT, STOP USING IMMEDIATELY AND CONTACT THE MANUFACTURER.

#### Warranty

To ensure your satisfaction, Verdermis offers a comprehensive warranty on its Beauty Tools (check website for details)



To Maximize Your Results Scan Here



WWW.VERDERMIS.COM

Dist by: Y&E Ent. LLC  
 Atlanta, GA 30338 USA  
 Ste 88757 Hauptstraße 22 Neugersdorf 02727 (Saxe) Allemagne Germany  
 71-75 Shelton Street London UK