currentbody.com

CURRENTBODY|skin™ LED Light Therapy Mask is manufactured by
CurrentBody.com Ltd, D6, Stanley Green Business Park, Commercial Avenue, Cheadle Hulme, SK8 6QH, United Kingdom.
EU Authorised Representative: Certification Company, Veluwezoom 42, 1327 AH, ALMERE, The Netherlands

CURRENTBODY|skin™

LED LIGHT THERAPY MASK

USER GUIDE AND MANUAL

Includes important safety information.

Read all instructions prior to use.

PLEASE KEEP THESE INSTRUCTIONS SAFE

CONTENTS

| | |
|--|----|
| CURRENTBODY Skin™ LED Light Therapy Mask | 02 |
| What is in the box | 03 |
| Important safety information | 04 |
| Warnings | 05 |
| Cautions | 06 |
| Contraindications | 07 |
| Device overview | 07 |
| How to use | 09 |
| Warranty | 12 |
| Specification | 13 |
| Symbols explained | 13 |
| Safety and emc | 15 |
| Manufacturer information | 18 |
| Français | 19 |
| Italiano | 24 |
| Deutsch | 29 |
| Japanese | 34 |
| Chinese | 39 |
| Korean | 43 |

Welcome to CurrentBody Skin™ LED Light Therapy Mask

Established in 2009, CurrentBody is the world's only online retailer dedicated to home-use beauty devices. As the first company to recognise devices as one of the fastest growing trends in the world of beauty, we have unrivalled expertise in selecting at-home devices that are science-backed and clinically proven, offering our customers the ultimate choice of devices that treat a wide range of concerns, including anti-ageing, cleansing and hair removal.

Underpinning our range of innovative beauty technology is expert education to help our customers choose the right device for them and excellent customer service.

Our highly trained beauty device experts are on hand to offer impartial advice and answer your questions in person, as well as how-to videos and top tips available throughout our website and social media channels.

What is in the box

System Includes:

CurrentBody Skin™
LED Light Therapy
Mask

Detachable,
adjustable
head strap

Storage bag

USB
charging cord

Comfort goggles

Important Features:

USB port connector

ON/OFF button;

Power indicator light

Type-C Charging port

IMPORTANT SAFETY INFORMATION

- Do not operate the device where it could come in contact with water.
- Do not use or charge the device with a power adapter that is damaged or faulty. Do not harshly bend or abuse the USB charging cord.
- Do not attempt to take the device apart. The device is not serviceable by the user. This will void the warranty.
- Stop using the device immediately if you think it is defective.

Customer Care

Please email sales@currentbody.com if you have any questions or you need any advice.

WARNINGS

- This device is designed for use on adults over 18 years of age and not recommended for use on infants or children.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant.
- Please consult your doctor prior to use if you have, or suspect you have, cancer.
- Do not use the device if it is wet. If the device has fallen into water, you must unplug it first before you touch it.
- Do not look directly into the lights as this may cause temporary discomfort. Ensure that you place the device on the face so that you can see through the eye holes and that the device is not covering the eyes.
- Do not use if your skin is red or swollen from a previous treatment.
- Do not use the device to treat skin that is painful to touch, raised (higher than surrounding skin), crusting, or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- This device contains a Lithium Ion Battery. Before disposal of this device, please review your local laws and requirements surrounding Lithium Ion Battery disposal. The preferred method of disposal is through recycling.

CAUTIONS

- Discontinue use if the device overheats or becomes uncomfortably hot.
- Always unplug the device from the outlet after charging is complete.
- Do not store the device in direct sunlight or in a location with temperatures over 100 Fahrenheit as this may cause damage to the plastic coating.
- Use only the supplied Universal USB charging cord. Use of an incorrect power cord could damage the device.
- Do not attempt to open or repair the device or power supply.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in skin injury.
- If you have any concerns about the skin you are treating, you should cease use and check with your doctor.
- Contact your doctor if you have any medical questions regarding use of this device.

CONTRAINDICATIONS

Ask your doctor prior to using if you:

1. Are taking prescription or over-the-counter drugs known to cause light sensitivity. You can find more information and a list of drugs that can cause light sensitivity at:

<https://www.currentbody.com/pages/led-contraindications>

2. Have any suspicious or cancerous lesions.
3. Have any photosensitive disorders or any genetic conditions of the eye.

Test yourself for light sensitivity prior to using on your face:

1. Power on your device.
2. Place the device on your arm and conduct a full 10-minute treatment.
3. After the treatment is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Return the device to CurrentBody within 14 days for a full refund.

DEVICE OVERVIEW


CURRENTBODY Skin™ LED Light Therapy Mask is a Light Emitting Diode (LED) device that harnesses energy to help reduce the signs of ageing. The device incorporates natural light energy in the, Red (630nm), and Infrared (880nm) spectrums that can help to reduce fine lines and wrinkles.

How it works

Youthful skin is rich in collagen, which gives the skin lift and support. Collagen works deep within the skin like columns of tissue that support the skin's surface. Over time, skin produces less collagen, and with exposure to sunlight, healthy collagen becomes damaged. This can result in wrinkles, especially around the mouth and eyes (crow's feet and elevens). The lights from this device penetrate deep within the layers of your skin to stimulate collagen production.

CURRENTBODY Skin™ LED Light Therapy Mask is programmed to automatically switch off after the recommended 10 minute treatment so it is very easy to use. With daily treatments you can achieve professional results in the comfort of your own home.

Light Wavelength Diagram

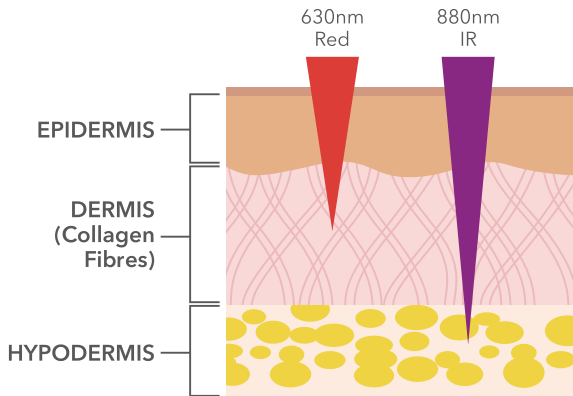HOW TO USE

Set up:

1. Unpack your CURRENTBODY Skin™ LED Light Therapy Mask.
2. Insert the USB power cord into the controller and then into a suitable USB charger (not supplied).
3. Insert the USB charger into a mains outlet and switch on.
4. The power indicator lights will flash to indicate charging. The device will take 4 hours to charge. When the battery is fully charged the green light will stop flashing and turn solid green.
5. Place the device in the protective storage bag when not in use.
6. Ensure you unplug the device from the USB charger before using. For your safety, this device will not work when it is charging.
7. When the charge is low the indicator light will flash. To recharge the device, follow the instructions above.

Instructions for Use:

1. You must conduct a sensitivity test before use (page 7).
2. Wash and fully dry your face prior to use. You can apply a serum or an eye cream (without spf) around your eyes to enhance the LED treatment.
3. Attach the strap and place the device onto your face and adjust the strap to fit snugly around your head.
4. Switch on the device. CURRENTBODY Skin™ LED Light Therapy Mask features an automatic shutoff function after completion of a 10-minute treatment.
5. Remove device from your face and wipe the device clean (see cleaning instructions).
6. Place the device into the protective storage bag.
7. If you find the LED lights too bright, we have included goggles for your comfort. The device can safely be used without goggles - it has been tested to all relevant eye safety standards.

Helpful Tips For best results:

- Use daily.
- Use the device on clean, dry skin only.
- Do not use the device after the application of creams that contain SPF as this could block light penetration.
- Wipe device clean after each use (see cleaning instructions).
- Immediately following treatment, apply a serum or moisturiser of your choice.

Maintenance:

1. This device has no serviceable parts. The battery is not accessible and not serviceable. Users should not attempt to access or service the battery as this will damage the device and void the warranty.
2. Device includes only an ON/OFF button. All settings are set to the factory setting.
3. Charge the device only with the supplied USB power charging cord and a universal USB power adapter (not included).
4. The power indicator light will begin to flash when battery levels are low. This indicates the device should be recharged.
5. Compatible Power Supply: 120V/240V

Cleaning:

1. Turn device off and unplug from any power source.
2. Wipe the tops of the LEDs with alcohol wipes or damp cloth to remove any residue from skin oils after each use.
3. Clean the device body by using damp or microfibre cloth.
4. Allow the device to dry completely prior to use.

Troubleshooting:

System is not turning on:

- a. Is the device fully charged? (power indicator light is on and solid)
- b. Is the power indicator light illuminated?
- c. Is the device disconnected from the mains? The device will not work while charging.

System is not charging:

- a. Is USB power cord firmly inserted into the device?(does power indicator light turn on?)
- b. Try a different wall outlet to determine if there is power to the outlet.

System shuts off abruptly, prior to completion of 10-minute treatment:

- a. System is designed to shut off automatically after 10 minutes.
- b. System may shut off after multiple repeated uses if the device exceeds desired temperature. Allow the device to cool for several minutes and try again.

Some lights appear turned off:

The infrared (IR) light spectrum is not visible to the human eye. These lights will look as if they are not on, but they are operating properly. If the red, visible lights are on, then the device is operating as designed.

The device does not switch off after a 10 minute treatment:

The device is faulty. Please return to CurrentBody.

The LEDs flicker during use:

The device is faulty. Please return to CurrentBody.

WARRANTY

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship. This warranty extends only to the original purchaser. If the product should become faulty within the warranty period, contact CurrentBody.

- CurrentBody reserves the right to repair a defective product.
- CurrentBody reserves the right to replace a defective product with the most compatible product currently available.
- This warranty does not cover damage caused by the following:
 - Accident, misuse, abuse or alteration
 - Servicing by unauthorised persons
 - Connection to incorrect power sources
 - Use with unauthorised accessories

CurrentBody shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you. This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed. If the product is defective within the warranty period please contact CurrentBody.

SPECIFICATION

| Specification | |
|-----------------------|--|
| Produce name | LED Light Therapy Mask |
| Model | MK66R-B |
| Power source | Adapter: Input: 100-240V~50/60Hz, 0.5A Output: 5.0V=2.0A |
| Weight | 350g±10g |
| Dimensions | 429 x 197 x 5 mm (W x L x H) |
| Wavelengths | Red630nm+Infrared880nm |
| Energy Source | LEDs |
| Power Density | 30mw/cm2 |
| Treatment time | 10 minutes, 4 times per week |
| Location for Use | Face |
| Environment of Use | OTC |
| Safety and EMC | IEC 60601-1, IEC 60601-1-2, IEC 60601-1-6, IE C60601-1-11, IEC 60601-2-57 |
| Biocompatibility | ISO 10993-1, ISO 10993-5, ISO 10993-10 |
| Operating Environment | Temperature: -10~35 C Humidity: 30-80%RH |
| Storage Environment | Temperature: -20~60 C Humidity: 30-90%RH |

THE SYMBOLS AND DESCRIPTION

| No. | Symbols | Description |
|-----|---|---|
| 1 |  | Manufacturer. |
| 2 |  | "WEEE (Waste Electrical and Electronic Equipment)". This means that the battery should not be disposed of in the household waste. |
| 3 |  | Consult instructions for use. |

THE SYMBOLS AND DESCRIPTION

| No. | Symbols | Description |
|-----|--|---|
| 4 |  | The date of manufacture. |
| 5 | IP22 | Protection against ingress of solids, it means the volume of the solids should less than 12mm. Protection against ingress of water, it means the device could protected against vertical falling water drop 15°. |
| 6 |  | Type BF Applied Part. |
| 7 |  | This symbol means: Conforms to certain EC Directives on health, safety, and environmental protection standards. CE stands for 'Conformité Européenne'. |
| 8 |  | This symbols implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed. |
| 9 |  | This symbol means: Class B, Civil or household products using light sources |
| 10 |  | Batch code. |
| 11 |  | Symbol for "SERIAL NUMBER". This symbol shall be accompanied by the manufacturer's serial number. |
| 12 |  | Temperature limit. |
| 13 |  | Humidity limit. |
| 14 |  | This symbol means keep dry. |
| 15 |  | This symbol means keep the device away from direct sunlight. |

SAFETY AND EMC

1) It complies with Medical Electrical Safety Standards (IEC 60601-1, IEC 60601-1-11, IEC 60601-2-57).

2) This device also complies with Medical EMC Standard (IEC 60601-1-2).

| Guidance and manufacturer's declaration – electromagnetic emissions | | |
|---|---|---|
| LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should ensure that it is used in such an environment. | | |
| Emissions | Compliance | Electromagnetic environment-- guidance |
| RF emissions CISPR 11 | Group 1 | LED Light Therapy Mask uses LED energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. |
| RF emissions CISPR 11 | Class B | LED Light Therapy Mask is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. |
| Harmonic emissions IEC 61000-3-2 | Not applicable | |
| Voltage fluctuations/ flicker emissions IEC 61000-3-3 | $P_{st} \leq 1.0$; $P_{Tl} \leq 0.65$; $d(t) \text{ sec} \leq 0.5$; $d_{max} \leq 4\% \text{ dc} \leq 3.3\%$ | |

| Guidance and manufacturer's declaration – electromagnetic immunity | | | |
|---|---|---|--|
| LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should ensure that it is used in such an environment. | | | |
| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment --guidance |
| Electrostatic discharge (ESD) IEC 61000-4-2 | $\pm 8 \text{ kV}$ contact $\pm 2 \text{ kV}$, $\pm 4 \text{ kV}$, $\pm 8 \text{ kV}$, $\pm 15 \text{ kV}$ air | $\pm 8 \text{ kV}$ contact $\pm 2 \text{ kV}$, $\pm 4 \text{ kV}$, $\pm 8 \text{ kV}$, $\pm 15 \text{ kV}$ air | Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %. |
| Electrical fast transient/burst IEC 61000-4-4 | $\pm 2 \text{ kV}$ 100 kHz repetition frequency | $\pm 2 \text{ kV}$ 100 kHz repetition frequency | Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT |
| Surge IEC 61000-4-5 | $\pm 0.5 \text{ kV}$, $\pm 1 \text{ kV}$ Line-to-line | $\pm 0.5 \text{ kV}$, $\pm 1 \text{ kV}$ Line-to-line | Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT |

| | | | |
|--|---|---|--|
| <p>Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11</p> | <p>0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°)</p> <p>0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°)</p> <p>0 % UT; 250/300 cycle</p> | <p>0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°)</p> <p>0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°)</p> <p>0 % UT; 250/300 cycle</p> | <p>Mains power quality should be a HOME HEALTHCARE ENVIRONMENT. If the user of LED Light Therapy Mask requires continued operation during power mains interruptions, it is recommended that LED Light Therapy Mask be powered from an uninterruptible power supply or a battery.</p> |
| <p>Power frequency (50/60 Hz) magnetic field IEC 61000-4-8</p> | <p>30 A/m</p> | <p>30 A/m</p> | <p>Power frequency magnetic fields should be at levels characteristic of a HOME HEALTHCARE ENVIRONMENT</p> |
| <p>NOTE UT is the a.c. mains voltage prior to application of the test level</p> | | | |

| Guidance and manufacturer's declaration – electromagnetic immunity | | | |
|--|--|--|--|
| <p>LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should assure that it is used in such an environment.</p> | | | |
| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment –guidance |
| <p>Conducted RF IEC 61000-4-6</p> | <p>3 Vrms 0,15 MHz and 80 MHz</p> <p>6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz</p> | <p>3 Vrms 0,15 MHz and 80 MHz</p> <p>6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz</p> | <p>Portable and mobile RF communications equipment should be used no closer to any part of LED Light Therapy Mask, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance</p> <p>$d = 1,2 \cdot \sqrt{P}$ 80 MHz to 800 MHz $d = 1,2 \cdot \sqrt{P}$ 800 MHz to 2.7 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation Distance in meters (m).</p> |

| | | | |
|---|---|---|---|
| <p>Radiated RF IEC 61000-4-3</p> | <p>10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz</p> | <p>10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz</p> | <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. b Interference may occur in the vicinity of Equipment marked with the following symbol:</p>  |
| <p>NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p> | | | |
| <p>a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which LED Light Therapy Mask is used exceeds the applicable RF compliance level above, LED Light Therapy Mask should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating LED Light Therapy Mask. b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.</p> | | | |

| Recommended separation distances between portable and mobile RF communications equipment and the LED Light Therapy Mask | | | |
|--|---|---|--|
| <p>LED Light Therapy Mask is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of LED Light Therapy Mask can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and LED Light Therapy Mask as recommended below, according to the maximum output power of the communications equipment.</p> | | | |
| Rated maximum output power of transmitter (W) | Separation distance according to frequency of transmitter (m) | | |
| | 150 kHz to 80 MHz $d = 1,2 \cdot \sqrt{P}$ | 80 MHz to 800 MHz $d = 1,2 \cdot \sqrt{P}$ | 800 MHz to 2.7 GHz $d = 2,3 \cdot \sqrt{P}$ |
| 0.01 | 0.12 | 0.12 | 0.23 |
| 0.1 | 0.38 | 0.38 | 0.73 |
| 1 | 1.2 | 1.2 | 2.3 |
| 10 | 3.8 | 3.8 | 7.3 |
| 100 | 12 | 12 | 23 |
| <p>For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p> | | | |

MANUFACTURER INFORMATION

CurrentBody.com Ltd
D6 Stanley Green Business Park Commercial Avenue,
Cheadle Hulme SK8 6QH

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