[currentbody.com](https://currentbody.com)

CURRENTBODY|skin LED Hand Perfector is manufactured by  
CurrentBody.com Ltd, D6, Stanley Green Business Park, Commercial Avenue, Cheadle Hulme, SK8 6QH, United Kingdom.  
CurrentBody BV, Naritaweg 2, 1437EL Rozenburg, Netherlands

Version 2 March 23

# CURRENTBODY|skin™

## LED HAND PERFECTOR

### USER GUIDE AND MANUAL

Includes important safety information.

Read all instructions prior to use.

**PLEASE KEEP THESE INSTRUCTIONS SAFE**

## CONTENTS

CURRENTBODY Skin™ LED Hand Perfector	02
What is in the box	03
Important safety information	04
Warnings	05
Cautions	06
Contraindications	07
Device overview	07
How to use	09
Warranty	12
Specification	13
Symbols explained	13
Safety and emc	15
Manufacturer information	18
Français	19
Deutsch	24
Italiano	29
Español	34
日本語	40
简体中文	45
繁體中文	49
한국인	53
Русский	58
हिन्दी	63

# Welcome to CurrentBody Skin™ LED Hand Perfector

Established in 2009, CurrentBody is the world's only online retailer dedicated to home-use beauty devices. As the first company to recognise devices as one of the fastest growing trends in the world of beauty, we have unrivalled expertise in selecting at-home devices that are science-backed and clinically proven, offering our customers the ultimate choice of devices that treat a wide range of concerns, including anti -ageing, cleansing and hair removal.

Underpinning our range of innovative beauty technology is expert education to help our customers choose the right device for them and excellent customer service.

Our highly trained beauty device experts are on hand to offer impartial advice and answer your questions in person, as well as how-to videos and top tips available throughout our website and social media channels.

## What's in box

### System Includes:

CurrentBody Skin™  
LED Hand  
Perfector

Storage bag

Power supply and  
country specific adaptors

## IMPORTANT SAFETY INFORMATION

- Do not operate the device where it could come in contact with water.
- Do not use the device with a power adapter that is damaged or faulty. Do not harshly bend or abuse the charging cord.
- Do not attempt to take the device apart. The device is not serviceable by the user. This will void the warranty.
- Stop using the device immediately if you think it is defective.

### Customer Care

Please contact us via the website if you have any questions or you need any advice.

## WARNINGS

- This device is designed for use on adults over 18 years of age and not recommended for use on infants or children.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant or are breastfeeding.
- Please consult your doctor prior to use if you have, or suspect you have, cancer.
- Do not use the device if it has fallen into water. Do not use the device if it is wet.
- Do not look directly into the lights as this may cause temporary discomfort.
- Do not use if your skin is red or swollen from a previous treatment.
- Do not use the device to treat skin that is painful to touch, raised (higher than surrounding skin), crusting, or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- Do not use the device if the top is not secured to the base.

## CAUTIONS

- Discontinue use if the device overheats or becomes uncomfortably hot.
- Always unplug the device after use.
- Do not store the device in direct sunlight or in a location with temperatures over 100 Fahrenheit as this may cause damage to the plastic coating.
- Do not attempt to open or repair the device or plug.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in skin injury.
- If you have any concerns about the skin you are treating, you should cease use and check with your doctor.
- Contact your doctor if you have any medical questions regarding use of this device.

## CONTRAINDICATIONS

### Ask your doctor prior to using if you:

1. Are taking prescription or over-the-counter drugs known to cause light sensitivity. You can find more information and a list of drugs that can cause light sensitivity at:  
<https://www.currentbody.com/pages/led-contraindications>
2. Have any suspicious or cancerous lesions.
3. Have any photosensitive disorders or any genetic conditions of the eye.

### Test yourself for light sensitivity prior to use:

1. Power on your device.
2. Wash and dry your hands, and then conduct a full 10-minute treatment.
3. After the treatment is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Return the device to CurrentBody within 14 days for a full refund.


## DEVICE OVERVIEW

CURRENTBODY Skin™ LED Hand Perfector is a Light Emitting Diode (LED) device that harnesses energy to help reduce the signs of ageing. The device incorporates natural light energy in the, Red (633nm), and Near Infrared (830nm) spectrums that can soothe and soften, to help achieve youthful, radiant hands.

## How it works

CURRENTBODY Skin™ LED Hand Perfector is programmed to automatically switch off after the recommended 10 minutes treatment so it is very easy to use. The LED wavelengths from the device penetrate the skin. With daily treatments you can achieve professional results in the comfort of your own home.

Light Wavelength Diagram
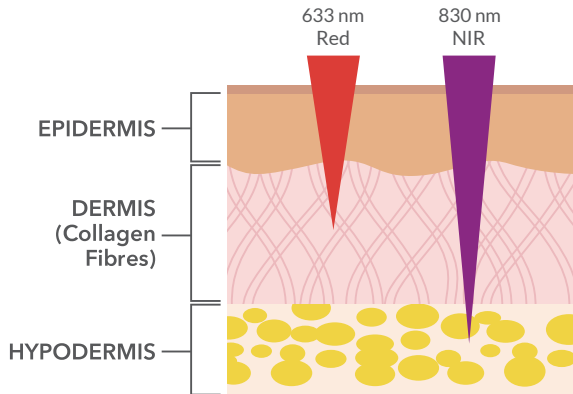
## HOW TO USE

### Set up:

1. Unpack your CURRENTBODY Skin™ LED Hand Perfector.
2. Select the appropriate plug for your territory and plug into a mains outlet.

### Instructions for Use:

1. You must conduct a sensitivity test before use (page 7).
2. Wash and fully dry your hands prior to use.
3. Switch on the device and place your hands in the device.
4. CURRENTBODY Skin™ LED Hand Perfector features an automatic shutoff function after completion of a 10-minute treatment.
5. Wipe the device clean (see cleaning instructions).
6. Place the device into the protective storage bag.

## Helpful Tips For best results:

- Use daily.
- Use the device on clean, dry skin only.
- Do not use the device after the application of creams that contain SPF as this could block light penetration.
- Immediately following treatment, apply a serum or moisturiser of your choice.

## Maintenance:

1. This device has no serviceable parts. Users should not attempt to service the device as this will damage it and void the warranty.
2. Device includes only an ON/OFF button. All settings are set to the factory setting.

## Cleaning:

1. Turn device off and unplug from any power source.
2. Wipe the tops of the LEDs with alcohol wipes or damp cloth to remove any residue from skin oils after each use.
3. Clean the device body by using damp or microfibre cloth.
4. Allow the device to dry completely prior to use.

## Troubleshooting:

### System is not turning on:

- a. Please ensure the device is firmly plugged in to a mains outlet.
- b. If the device still does not work, please get in touch with CurrentBody via the website.

### System shuts off abruptly, prior to completion of 10-minute treatment:

- a. System is designed to shut off automatically after 10 minutes.
- b. System may shut off after multiple repeated uses if the device exceeds desired temperature. Allow the device to cool for several minutes and try again.

### Some lights appear turned off:

The infrared (IR) light spectrum is not visible to the human eye. These lights will look as if they are not on, but they are operating properly. If the red, visible lights are on, then the device is operating as designed.

### The device does not switch off after a 10 minutes treatment:

The device may be faulty. Please get in touch with CurrentBody via the website.

### The LEDs flicker during use:

The device may be faulty. Please get in touch with CurrentBody via the website.

## WARRANTY

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship. This warranty extends only to the original purchaser. If the product should become faulty within the warranty period, contact CurrentBody.

- CurrentBody reserves the right to repair a defective product.
- CurrentBody reserves the right to replace a defective product with the most compatible product currently available.
- This warranty does not cover damage caused by the following:
  - Accident, misuse, abuse or alteration
  - Servicing by unauthorised persons
  - Connection to incorrect power sources
  - Use with unauthorised accessories

CurrentBody shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you. This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed. If the product is defective within the warranty period please contact CurrentBody.

## SPECIFICATION

Specification	
Produce name	LED Hand Perfector
Model	HS-70
Power source	Adapter: Input: 100-240Vac/50-60Hz Output: 5.0V $\bar{=}$ 1.0A
Weight	1050g $\pm$ 10g
Dimensions	403 x 275 x 82 mm (W x L x H)
Wavelengths	Red633nm+Near Infrared830nm
Energy Source	LEDs
Number of LEDs	70 Dual LEDs: 70 x 830nm, 70 x 633nm
Power Density	30mw/cm <sup>2</sup>
Treatment time	10 minutes, For best results use daily.
Location for Use	Hands
Environment of Use	OTC
Safety and EMC	IEC 60601-1, IEC 60601-1-2, IEC 60601-1-6, IEC 60601-1-11, IEC 60601-2-57
Biocompatibility	ISO 10993-1, ISO 10993-5, ISO 10993-10
Operating Environment	Temperature: -10~35 C Humidity: 30-80%RH
Storage Environment	Temperature: -20~60 C Humidity: 30-90%RH

## THE SYMBOLS AND DESCRIPTION

No.	Symbols	Description
1		Manufacturer.
2		"WEEE (Waste Electrical and Electronic Equipment)" This means that the product including its battery must not be disposed of in household waste.
3		Consult instructions for use.

## THE SYMBOLS AND DESCRIPTION

No.	Symbols	Description
4		The date of manufacture.
5	IP22	Protection against ingress of solids, it means the volume of the solids should less than 12mm. Protection against ingress of water, it means the device could protected against vertical falling water drop 15°.
6		Type BF Applied Part.
7		This symbol means: Conforms to certain EC Directives on health, safety, and environmental protection standards. CE stands for 'Conformité Européenne'.
8		This symbols implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
9		Batch code.
10		Symbol for "SERIAL NUMBER". This symbol shall be accompanied by the manufacturer's serial number.
11		Temperature limit.
12		Humidity limit.
13		This symbol means keep dry.
14		This symbol means keep the device away from direct sunlight.

## SAFETY AND EMC

1) It complies with Medical Electrical Safety Standards (IEC 60601-1, IEC 60601-1-11, IEC 60601-2-57).

2) This device also complies with Medical EMC Standard (IEC 60601-1-2).

Guidance and manufacturer's declaration – electromagnetic emissions		
LED Hand Perfector is intended for use in the electromagnetic environment specified below. The customer or the user of LED Hand Perfector should assure that it is used in such an environment.		
Emissions	Compliance	Electromagnetic environment-- guidance
RF emissions CISPR 11	Group 1	LED Hand Perfector uses LED energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	LED Hand Perfector is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/ flicker emissions IEC 61000-3-3	$P_{st} \leq 1.0$ ; $P_{It} < 0.65$ ; $d(t) \text{ sec} \leq 0.5$ ; $d_{max} \leq 4\% \text{ dc} \leq 3.3\%$	

Guidance and manufacturer's declaration – electromagnetic immunity			
LED Hand Perfector is intended for use in the electromagnetic environment specified below. The customer or the user of LED Hand Perfector should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment --guidance
Electrostatic discharge (ESD) IEC 61000-4-2	$\pm 8 \text{ kV}$ contact $\pm 2 \text{ kV}, \pm 4 \text{ kV},$ $\pm 8 \text{ kV}, \pm 15 \text{ kV}$ air	$\pm 8 \text{ kV}$ contact $\pm 2 \text{ kV}, \pm 4 \text{ kV},$ $\pm 8 \text{ kV}, \pm 15 \text{ kV}$ air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	$\pm 2 \text{ kV}$ 100 kHz repetition frequency	$\pm 2 \text{ kV}$ 100 kHz repetition frequency	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT
Surge IEC 61000-4-5	$\pm 0,5 \text{ kV}, \pm 1 \text{ kV}$ Line-to-line	$\pm 0,5 \text{ kV}, \pm 1 \text{ kV}$ Line-to-line	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT

Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°)  0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°)  0 % UT; 250/300 cycle	0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°)  0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°)  0 % UT; 250/300 cycle	Mains power quality should be a HOME HEALTHCARE ENVIRONMENT. If the user of LED Hand Perfector requires continued operation during power mains interruptions, it is recommended that LED Hand Perfector be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a HOME HEALTHCARE ENVIRONMENT
NOTE UT is the a.c. mains voltage prior to application of the test level			

Guidance and manufacturer's declaration – electromagnetic immunity			
LED Hand Perfector is intended for use in the electromagnetic environment specified below. The customer or the user of LED Hand Perfector should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment –guidance
Conducted RF IEC 61000-4-6	3 Vrms 0,15 MHz and 80 MHz  6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz	3 Vrms 0,15 MHz and 80 MHz  6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of LED Hand Perfector , including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance  $d = 1,2 \cdot \sqrt{P}$ 80 MHz to 800 MHz $d = 1,2 \cdot \sqrt{P}$ 800 MHz to 2.7 GHz  Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation Distance in meters (m).

Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz	10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz	Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. b Interference may occur in the vicinity of Equipment marked with the following symbol:  
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which LED Hand Perfector is used exceeds the applicable RF compliance level above, LED Hand Perfector should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating LED Hand Perfector . b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the LED Hand Perfector			
LED Hand Perfector is intended for use in an electromagnetic environment in which radiated RF disturbances are control RGB. The customer or the user of LED Hand Perfector can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and LED Hand Perfector as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2 \cdot \sqrt{P}$	80 MHz to 800 MHz $d = 1,2 \cdot \sqrt{P}$	800 MHz to 2.7 GHz $d = 2,3 \cdot \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

## MANUFACTURER INFORMATION

CurrentBody.com Ltd  
D6 Stanley Green Business Park Commercial Avenue,  
Cheadle Hulme SK8 6QH  
CurrentBody BV, Naritaweg 2, 1437EL Rozenburg,  
Netherlands

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