

MEGAWISE

HANDHELD DEEP PERCUSSION MASSAGER

Model: RD-A10



USER MANUAL

ACTIVATE YOUR 12-MONTH WARRANTY & GET AN EXCLUSIVE GIFT

Register within 7 days after receiving your new product.



Use your phone camera to scan me
& start chatting

or visit: m.me/MegaWiseUSA



Or Visit:

megawise.com/warranty



A Big Thank You

Dear Customer,

We're thrilled to welcome you to the MegaWise family, and we'd like to thank you for purchasing our product. We're proud to serve you as our customer and we hope that you get the most out of your purchase.

SUBSCRIBE

Access exclusive offers, news, and more.



**LIVE CHAT
SUPPORT**

-  m.me/MegaWiseUSA
-  support@megawise.com
-  megawise.com

CONTENTS

PART NAMES	02
PRODUCT FEATURES	03
PRODUCT FUNCTIONS	04
CLEANING & STORAGE	04
METHODS OF USE	05
ATTENTION	05
USAGE	08
TECHNICAL PARAMETERS	08

PART NAMES

5 Massage Modes

Increase Massage Speed

Decrease Massage Speed

ON/OFF button



Peeling Thermal
Therapy Massage
Head



Seven-Point
Massage Head



Scratching
Massage Head



Local Peeling
Massage Head



Point
Massage Head



SEVEN-POINT MASSAGE HEAD:

Covers a larger massage area; suitable for the back, waist, legs, etc.



PEELING THERMAL THERAPY MASSAGE HEAD:

Applies even pressure on the body.



SCRATCH MASSAGE HEAD:

This unique massage head features "powerful claws" suitable for the back, waist, legs, etc.



LOCAL PEELING MASSAGE HEAD:

The most commonly-used massage head.



POINT MASSAGE HEAD:

Suitable for specific acupuncture points; relieves muscle pain, just like a human hand massage.

PRODUCT FEATURES

1. The 5 different massage speeds, allow you to personalize your massage. To achieve the best results, begin at the 3rd speed setting and gradually increase to the most comfortable massaging speed. Press the "-" button of the massager to decrease the speed.

2. Press "Mode" to cycle through the 5 different massage modes.
 3. The long, non-slip handle allows you to massage any hard-to-reach area.
 4. This massager can Relieve swelling, increase muscle flexibility and ease muscle tension.
 5. The product features an Ergonomic and stylish design.
-

PRODUCT FUNCTIONS

1. Promotes blood circulation.
 2. Relieves muscle tension, reduces pressure and eliminates fatigue.
 3. Massages different levels of muscle tissue to relieve pain.
-

CLEANING & STORAGE

1. Always unplug the massager before cleaning.
2. Simply wipe the product's surface with a dry, soft cloth.
3. Do not immerse the massager in water or use any harsh cleaners or solvents.
4. When storing, ensure that the cord is not twisted or kinked.
5. In order to prevent product damage, refrain from dropping, striking, or administering heavy pressure to the appliance.
6. Please use a neutral cleaning agent to clean the appliance. Do not use abrasive or corrosive cleaning agents. Note: The massager must be unplugged while cleaning.

METHODS OF USE

1. Insert the plug into the power outlet.
2. Press the "power switch" to turn the massager on. Press the "+" and "-" buttons to adjust the different massage intensities.
3. Press the "Mode" button to experience the different massage rhythms.
4. All programs will automatically shut down after 15 minutes.
5. After use, turn off the massager and unplug the appliance.

ATTENTION

1. Please use this product according to the instructions.
2. Before use, check that the voltage rating of the appliance matches that of the household voltage.
3. Before cleaning the product, ensure that the power switch is turned off and the power cord has been removed.
4. Do not use this product around children.
5. Do not use this product while bathing.
6. Do not immerse the appliance in water or use it in damp places.
7. Do not use the product close to radiation or other electrical appliances.
8. If the massager accidentally falls into water, please turn off the power supply and unplug the power cord. Do not immediately remove the product from the water.
9. Do not use the massager on someone who is sleeping.
10. Do not place the appliance on a non-ventilated surface, such as a quilt or pillow. This can damage the device and cause it to overheat.
11. Do not use this product near inflammable and explosive gases.

12. Do not use the massager around the eyes or injured parts of the body.
13. Do not use on the head, chest, knees, and joints for long periods of time.
14. Do not allow children and disabled persons to use the appliance.
15. Do not use this product if you are feeling unwell.
16. If the power cord is damaged or there is a problem with the massager, stop using it immediately. Do not attempt to disassemble or repair the product by yourself.
17. Do not insert foreign objects into the gaps around the massage head.
18. It is recommended that each use of the massager should not exceed 15 minutes.
19. Remove the power plug with your hand, and do not pull the power cord. Do not use within 1 hour of eating a meal, and before using the massager, it is advisable to drink 150 ml of water.
20. If you experience physical pain or discomfort as a result of the massager, please discontinue use and consult a doctor.
21. Do not use with any voltage over than 110V.
22. Do not repeatedly insert and remove the plug from the power socket.
23. When cleaning the massager, the power cord must be unplugged.
24. If the power cord is damaged, it must be replaced by the manufacturer or a qualified appliance service center.
25. If the event that the appliance malfunctions, please contact customer support for assistance.

PLEASE CONSULT YOUR DOCTOR BEFORE USING IF ANY OF THE FOLLOWING ARE RELEVANT TO YOU

1. You have an implanted pacemaker or similar implanted medical devices.
2. Are suffering from heart disease.
3. Are suffering from back strain.
4. Have damaged or fragile bones, especially back bones.
5. Have had recent surgery.
6. In pregnancy or recently gave birth.
7. Are suffering from malignant tumors.
8. Are suffering from osteoporosis.
9. Are suffering from blood circulation problems or thrombosis.
10. Are suffering from skin infections.
11. Are being treated for an illness.
12. Have limb abnormalities.
13. Have spinal injuries caused by an accident or disease.

USAGE



TECHNICAL PARAMETERS

Product Name: Handheld Deep Percussion Massager

Model: RD-A10


Rated Voltage: AC 110V

Rated Frequency: 60Hz

Rated Power: 25W

CUSTOMER SERVICE

- 🗨 Please scan the QR code to start chatting and receive an instant reply.

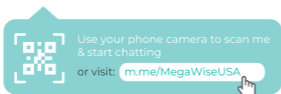

- 🗨 m.me/MegaWiseUSA
- ✉ support@megawise.com
- 🌐 megawise.com



Made in China