

 **LUSTRE®**
CLEARSKIN

REVIVE
MASK

Instructions for Use

REVIVE

MASK

Welcome to LUSTRE®ClearSkin REVIVE 10 minutes LED Beauty Therapy Mask. Discover the revolutionary power of light to transform the appearance of your skin.

Our revolutionary innovations are developed by scientists for beautiful results. We aim to challenge conventions and reinvent skincare through innovative, safe and proven technologies to deliver amazing skin benefits.

LUSTRE®ClearSkin REVIVE 10 minutes Light Beauty Therapy Mask has been designed to work fast and fight against dull and tired looking skin. It helps reveal a fresher looking skin and revived complexion with a visible youthful glow.

The LUSTRE®ClearSkin REVIVE mask has two settings to allow you to create a personalised skincare session. Select the BLUE setting to fight against dull looking skin and reveal a fresher and revived complexion. Select the RED setting to fight against tired looking skin and promote a visibly youthful glow. Alternate between BLUE and RED sessions for even better results.

Get ready for a relaxing experience, let your skin benefit from a proven technology from the comfort of your own home. Simply choose your program, put on the mask and enjoy! It only lasts 10 minutes.

LUSTRE®ClearSkin REVIVE 10 minutes Light Beauty Therapy mask is easy to use, lightweight and flexible for ultra-comfort and a perfect fit. Enjoy your session safely thanks to its unique integrated eye protection.

Benefits of LUSTRE®ClearSkin REVIVE 10 minutes Light Beauty Therapy mask



BLUE
LED LIGHTS



ENERGISES



REFRESHED
AND REVIVED
COMPLEXION



RED
LED LIGHTS



VISIBLE YOUTHFUL
GLOW

Important Information

Important Safety Information

- **DO NOT** operate the device where it could come in

contact with water.

- **DO NOT** use or charge the device with power adapter that is damaged and/or faulty.

- **DO NOT** harshly bend or abuse wire or the USB charging cord.

- **DO NOT** fold/bend or abuse your LUSTRE®ClearSkin REVIVE 10 minutes Light Beauty Therapy mask which damage the LED lights.

- **DO NOT** attempt to take the device apart. This will void the warranty.

- Stop using the mask immediately if you think it is defective. Please contact our customer care team at info@lustreskin.com

Disclaimer

The product is a beauty appliance and is not intended to diagnose, cure or prevent any disease. LUSTRE®ClearSkin makes no claims, representations, or warranties regarding the ability of this product to cure any physical, skin, or mental condition. A qualified healthcare professional should be consulted with regards to any condition requiring medical attention.

Product Representation

LUSTRE®ClearSkin is constantly improving products with the latest technology. As a result, the actual product may differ slightly from the products shown in the images on the packaging and the user guide.


Cautions:

- **DO NOT** look directly into the lights. Doing so may cause temporary vision discomfort. Use the goggles to cover the eyes.
- This device is designed for use on adolescents and adults 13 years or older. Adolescents between 13 and 18 years old shall use the product under supervision of an adult aged 18 years or older. This product shall not be used by children under the age of 13.
- Please consult your general practitioner prior to use if you are or suspect you are, pregnant.
- Please consult your general practitioner prior to use if you have or suspect you have, cancer.
- **DO NOT** place the device near water or any wet surface.
- We highly recommend users perform a light sensitivity test before their first use of the device.

Table of Contents


1. What's in the Box	8
2. Before You Start:	9
- Charging	9
- Light Sensitivity Test	9
- Get Your Mask Ready	10
3. How LUSTRE®ClearSkin REVIVE Facewear mask works	11
4. How to use the device	12
- General use	12
- Cleaning	13
- Storing	14
5. Warnings	14
6. Troubleshooting	20
7. Specifications	21
8. Disposal	22
9. Warranty	22
10. Symbol Glossary	24
11. Warranty Limitations & Disclaimers	24
12. Contact info	26

What's in the Box

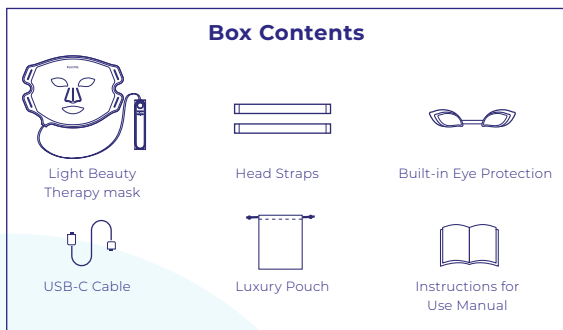

Please contact customer support if any of the mentioned components are missing or damaged.

Product Overview

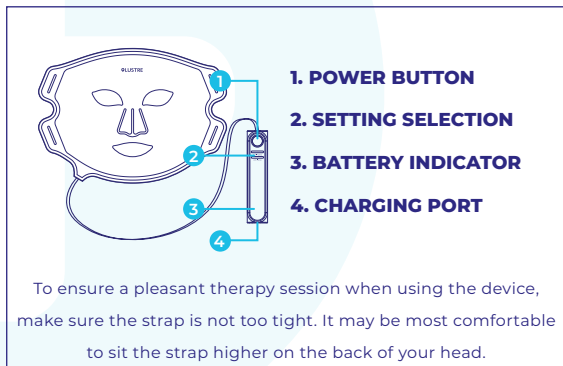

Before You Start

1. Charging

To charge the mask insert the USB-C cable into the mask. Insert the USB-A side of the cable into a USB power supply. The power indicator on the controller will start to blink indicating it is charging.


Note: Always unplug the device before use, the device will not work while being charged.

2. Light Sensitivity Test

1. Place the device on the inner forearm. Switch on the device for a full 10 minute cycle.
2. Don't look directly into the blue light, since it might be harmful to your eyes. Cautiousness is advised.
3. If your skin feels hot, uncomfortable, or itchy during the session stop immediately! You might be sensitive to light.
4. Wait till the session is done. The session will take 10 minutes to complete.
5. The device will automatically turn off after the session is done. Remove the device and wait for at least 6 hours.


Check your skin regularly over the six hours, if your skin appears red and/or itchy the device is unsuitable for you. In the case of light sensitivity please contact your medical advisor or doctor.

3. Get Your Mask Ready

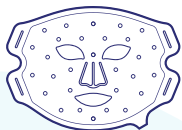


Step 1

Make sure that the inside of the mask is facing you.

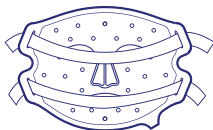


Step 2

Slide both straps through the slots of the mask. Make sure the velcro sides of the strap are facing towards you.

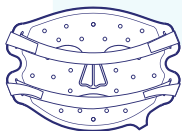

Step 3

Stick the velcro tips onto the strap. Make sure the straps are not too tight.

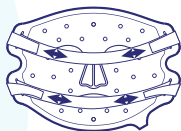

Step 4

Place the mask onto your face and tighten or loosen the straps by adjusting the velcro tips.

How LUSTRE® ClearSkin REVIVE mask works

Discover the revolutionary power of light to transform the appearance of your skin

Light Therapy

LED stands for light-emitting diodes. We use a selection of two different lights, that are the perfect combination to achieve a fresher looking skin and reveal a youthful glow. LUSTRE®ClearSkin REVIVE mask delivers personalised sessions thanks to its two settings:

Choose the RED setting to fight tired looking skin for a visible youthful glow.

Choose the BLUE setting to fight against skin dullness and reveal a refreshed and revived complexion.

Overview

Session time	10 minutes
BLUE LED setting Select this setting to fight against skin dullness and reveal a refreshed and revived complexion	3 times a week for 8 weeks
RED LED setting Select this setting to fight against tired looking skin for a visible youthful glow	3 times a week for 8 weeks
ALTERNATE RED LED and BLUE LED SESSIONS For even better and faster results	3 times a week for 8 weeks

How To Use the Device

General Use

Step 1

Start by cleaning your face with a gentle cleanser. Fit the mask comfortably and switch it on by pressing the power button on the controller for a few seconds until the power indicator lights up. The mask is now in

standby mode.

Step 2

You can select the desired setting by pressing the power button several times until the desired setting is selected.

1. RED Light Setting

Select this setting to fight against tired looking skin for a visible youthful glow

2. BLUE Light setting

Select this setting to fight against skin dullness and reveal a refreshed and revived complexion.

Step 3

After 10 minutes, the power indicator will blink twice before turning off automatically.

Once the session is completed remove the mask and apply your usual serums and creams.

Cleaning

It is important to care of your device. Clean the device with a damp cloth etc. Stubborn stains can be re-moved by gently wiping the stain with a mixture of 70% alcohol.

- **DO NOT** clean the device while it is turned on.

- **DO NOT** clean the device while charging.
- **DO NOT** use any other cleaning products as they might damage the device.

Storing

Store the device in the supplied pouch or the packaging. Make sure the device is not bent, stretched or crushed while being stored or when travelled with.

Warnings

- Fitzpatrick Skin Types V and VI can use the mask safely, however we strongly recommend to assess your skin carefully after each session, there is an increased risk of pigmentation on those skin types.

The Fitzpatrick Scale



- **DO NOT** shine the light from the LEDs directly into your eyes or someone else's eyes. This may cause temporary loss of sight.

- **DO NOT** use the mask under the age of 13. Adolescents over the age of 13 years old must use the mask under the supervision of a parent/guardian aged 18+.

- **DO NOT** let children under the age of 18 clean, charge or use the mask without adult supervision.

- **DO NOT** let children play with the mask. Keep out of reach of children.

- **DO NOT** perform maintenance or take apart the product. In case of malfunctioning or damaged please contact customer support.

- **DO NOT** use the mask outdoors to prevent damage to the product.

- **DO NOT** bend or twist the mask. Handle the mask with care.

- **DO NOT** sit, lay down or step on the mask. This will break the mask.

- **DO NOT** use solvents or abrasive cleaners on the mask unless explicitly specified in this manual or by customer service.

- **DO NOT** use the mask if it is not working properly, or if it has been dropped or damaged. For warranty and

repairs contact customer support.

- **DO NOT** use the mask if you are taking any medication that can cause photosensitivity. If you are not sure about your medication contact your health professional.

- **DO NOT** use the mask if you are pregnant, planning on becoming pregnant, or breastfeeding.

- **DO NOT** use the mask in any other way than specified in this manual.

- **DO NOT** use the mask if you suffer from Lupus erythematosus, photosensitive eczema, Albinism, or other skin conditions. If you use the mask to treat Lupus erythematosus, photosensitive eczema, or Albinism you may cause a severe skin reaction.

- **DO NOT** use the mask if you suffer from any photosensitive disorder (sensitization to light).

- **DO NOT** use the mask over a bacterial infection as this may result in the worsening of symptoms.

- **DO NOT** place the mask directly over an open wound.

- **DO NOT** use the mask if you have epilepsy.

- **DO NOT** use the mask over areas with phlebitis,

thrombophlebitis, and varicose veins.

- **DO NOT** use the mask if you suffer from light-induced headaches.

- **DO NOT** use the mask in any other way than specified in this manual.

- **DO NOT** continue to use if it feels hot or your skin feels uncomfortably warm.

- **DO NOT** use or store the mask near hot surfaces.

- **DO NOT** charge the mask when it is close to or might get in contact with water.

- **DO NOT** store the mask in extreme heat or cold.

- **DO NOT** place the mask in direct sunlight.

- **DO NOT** use the mask if you have a known allergic reaction to silicone. The mask contains silicone and plastics that come into contact with the skin during use. When in doubt of allergies first contact your healthcare provider before using the mask.

- **DO NOT** use the head strap if you have a known allergic reaction to polyester, nylon, and lycra®. The head strap contains a mixture of polyester, nylon, and lycra® that

come into contact with the skin during use. When in doubt of allergies first contact your healthcare provider before using the mask.

- **DO NOT** attempt to modify the mask in any way. Modifying the mask will void your warranty. The mask contains no serviceable parts.

- **DO NOT** use the mask if damaged or when the charging cable is damaged.

- **DO NOT** use the mask with non-manufacturer approved accessories.

- The product is a beauty appliance and is not intended to diagnose, cure or prevent any disease.

Note : The mask should only be supplied with 5V DC low voltage. The voltage is also marked on the mask. If any part of the your main adapter or cord is damaged, discontinue use. Contact Lustre customer support to arrange a replacement.

Note: Consult your healthcare provider before using the mask if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anesthetic, or desensitized areas.

Note: We highly suggest taking caution when taking

any medications. If in doubt, consult your healthcare provider.


Note: Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the mask and contact your healthcare provider.

Note: We highly suggest taking caution when taking any skin related or prescribed medications. If in doubt, consult your healthcare provider.

Troubleshooting


Problem	Cause	Solution
The mask doesn't switch on.	Charging cable is not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	Charging cable is not properly inserted into the device.	Insert the charging cable properly into the device.
	Power outlet is not working.	Try a different power outlet.
	On/off button has not been pushed one full second.	Push the on/off button for 1 second straight.
The device switches off unexpectedly.	Battery level is too low.	Charge the device.

Specifications

LED wavelengths	465nm (Blue), 640nm (Red)
Max. output optical power density	30MW/CM2
LEDs	66 dual-core Red+Blue light LEDs
Operating and storing temperature	0°C to 40°C 32°F to 104°F
Operating and storing humidity	30% – 95%
Power supply	5V  2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	1 meter/3.28 feet
Charging adapter power output	2A
Battery	2600 mAh
Battery type	Lithium-ion
Charging time	3-5 hours
Session time	10 minutes

Dimensions of the device	195 X 290 x 5 mm
Weight of the device	330 ± 30 gram

Disposal

The Waste of Electrical and Electronic Equipment (WEEE ) Directive (2012/19/EU) has been put in place to recycle products using the best available recovery and recycling techniques to minimize the impact on the environment, treat any hazardous substances, and avoid the increasing landfill. This product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of correctly following local authority requirements.

DO NOT dispose the device with household waste.

Warranty

The mask is warranted to be free from defects in material and workmanship for 12 months from its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the mask is used in conjunction with authorized LUSTRE® ClearSkin accessories. If the product should

become defective within the warranty period, contact the retailer where you bought the mask. This warranty does not cover products damaged by:








- Accident, misuse, abuse, or alteration.

- Use with unauthorized accessories or use other than as instructed.

- Connecting to incorrect current and voltage.

To register your product and be kept up-to-date with new and exciting developments from the LUSTRE® team, please take a moment to sign up on our website www.lustreclearskin.com. To complete your registration online you will need to enter your LUSTRE®ClearSkin batch number.

Symbol Glossary

Symbol	Description
! WARNING	WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.
! CAUTION	CAUTION indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.
	Refer to operating Instructions
	UK Conformity Assessed
	European Conformity
	Batch/Lot/ID Number This is the unique identifying number for each batch
	Catalogue number
	Manufacturer
	WEEE Compliance

Warranty Limitations & Disclaimers

This warranty gives you specific legal rights, and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, Lustre Skin Ltd disclaims all implied warranties. If applicable law does not permit the disclaimer of implied

warranties, then the duration of implied warranties are hereby expressly limited to the duration of the express limited warranty set forth above.

In no event shall Lustre Skin Ltd be liable for any special, incidental, consequential, punitive or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability) or otherwise, even if Lustre Skin Ltd is advised or aware that such damages may arise. Such damages include, but are not limited to: loss of Lustre Skin Ltd cost of any substitute equipment, the claims of third parties or damage to property.

If the applicable law does not allow limits on warranties, or on remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between Lustre Skin Ltd and you, the purchaser. This allocation is recognized by both parties and is reflected in the price. By using LUSTRE®ClearSkin REVIVE you acknowledge that you have read these terms, conditions, and limitations, understand them and are bound by them.

Contact Information

Lustre Skin Ltd

Alba Innovation Centre, Alba Campus, Livingston
Scotland, UK, EH54 7GA

T: +(44)1506 592239

E: info@lustreskin.com

W: www.lustreclearskin.com

Lustre Skin Ltd is an ISO 13485:2016 certified company



International Associates Auditing & Certification Limited
The Black Church, St Mary's Place,
Dublin 7, D07 P4AX Ireland
EUAR@ie.ia-net.com

 **LUSTRE**
CLEAR**SKIN**

CE **UK**
CA

SEPT.2023-900015.3