OMNILUX™  
— CLEAR —

English    Deutsch    Español    Français    Nederlands

User Guide

# English

Users of this device can obtain benefit from using it only if they read and comprehend all information in the user guide.

You must follow all the instructions in the User Guide.

Please feel free to contact us if you have any comments about the Omnilux Clear™ device or this User Guide.

## GlobalMed Technologies Co.

163 Camino Dorado Suite B  
Napa, CA 94558  
Tel: +1-707-224-0200  
www.omniluxled.com

The information in this document is proprietary to I-SMART Development Ltd, and supplied on the understanding that it shall be kept confidential to recipients and that it shall not be altered, reproduced and/or used for commercial purposes without the written consent of I-SMART Developments Ltd.

**EC REP** Med-Di-Dia Limited Galway Technology Centre Mervue  
Business Park, Galway, Ireland

**AU Sponsor** Emergo Australia, Level 20 Tower II, Darling Park  
201 Sussex Street, Sydney, NSW 2000, Australia

## Contents

### Section Header

	Section number
Introduction	1.0
Who should use Omnilux Clear?	2.0
Device description	3.0
Intended Purpose	3.1
When should Omnilux Clear NOT be used?	
Contraindications, Warning and Precautions	4.0
Known Side Effects and Benefits	5.0
First time Setup	6.0
Unpacking and Charging your Omnilux Clear	6.1
Connecting the Omnilux Clear to the Controller	6.2
Sensitization Test	6.3
Using Omnilux Clear	7.0
Treatment Recommendations	7.1

Wearing the Omnilux Clear	7.2
Carrying out a Treatment	7.3
Care and Maintenance	8.0
Cleaning the Omnilux Clear Device	8.1
Disposal	8.2
Troubleshooting	8.3
Labelling	9.0
Recycling	10.0
Technical specifications	11.0
Explanation of symbols	12.0
Warranty	13.0

## 1.0. Introduction

Thank you for choosing Omnilux Clear.

The Omnilux Clear LED mask is intended to emit light in the blue and red region of the light spectrum and is indicated for the treatment of mild to moderate acne vulgaris of the face.

(See “Who should use Omnilux Clear” section 2 on page 4, for more information about how to tell if you have mild to moderate acne and if the Omnilux Clear is right for you).

Users of this device can obtain benefit from using it only if they read and comprehend all information in the User Guide.

You must follow all the instructions in the User Guide.

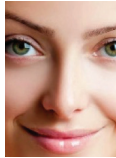
Please feel free to contact us if you have any commented about Omnilux Clear or this User Guide.

Before using the Omnilux Clear you must carry out a sensitization test and wait 6 hours. (See section 6.3).

## 2.0. Who should use Omnilux Clear?

The typical appearance of acne is a mixture of the following: a greasy skin, red pimples or spots, yellow pus-filled spots, and scars. Use this chart to determine whether Omnilux Clear™ is suitable for you.

Acne type	Description		Suitable for Omnilux Clear
Clear Skin	No pimples or spots or yellow pus-filled spots. No whiteheads or blackheads		<b>NO</b> Omnilux Clear is NOT indicated for Clear skin. Omnilux Clear should NOT be used to prevent pimples or spots forming.
Mild acne	A few pimples or spots, a few pus-filled spots, some whiteheads, or blackheads		<b>YES</b> Mild acne is suitable for treatment with Omnilux Clear
Moderate acne	Some pimples or spots, some pus-filled spots, only a few small nodules Many whiteheads or blackheads		<b>YES</b> Moderate acne is suitable for treatment with Omnilux Clear
Severe acne	Many pimples or spots, many pus-filled spots. Some nodules. Many whiteheads or blackheads		<b>NO</b> Omnilux Clear is NOT indicated for severe acne. Please consult your doctor for a more appropriate treatment
Whiteheads or Blackheads – Comedonal acne	Many whiteheads or blackheads. No pimples or spots. No pus-filled spots		<b>NO</b> Omnilux Clear is NOT indicated for whiteheads or blackheads. Please consult your doctor for a more appropriate treatment.

To find out if the Omnilux Clear is right for you and to find a full list of who can and who can't use the Omnilux Clear scan the QR

If Omnilux Clear is not suitable for you, return the device to the retailer.

## 3.0. Device Description

1. Flexible silicone CLEAR device with integrated eye-wear
2. Controller (containing a lithium battery) & power supply.
3. Country Specific adapters
4. 2x adjustable Straps
5. Carry Bag

Omnilux Clear is a home-use wearable LED light therapy device. It produces a cool, narrow band of light which reduces mild to moderate acne vulgaris on the face. It consists of a flexible silicone device that contains light emitting diodes (LED's) and a controller. The LEDs generate the light. The device is worn on the face and is held in place by adjustable Velcro straps. The controller turns the LEDs ON/OFF and controls power to the mask. The controller contains a rechargeable lithium-ion polymer battery and uses a visible display comprising of 3 micro-LEDs to show the user the battery charge status of the device. The power supply is used to charge the lithium battery and is connected to a suitable mains outlet via a 2 or 3 pin input socket and wall plug. The power cable is connected to the controller by a standard micro-USB A-C connector.

### Omnilux Clear cannot be operated whilst charging.

The device is not used to make measurements of any sort, or to draw any conclusions regarding the indication to treat. The device does not require checks on the light output as the LEDs DO NOT dim with age to any practical extent.

### 3.1. Intended Use

Omnilux Clear is an over-the-counter LED device that emits energy in the red and blue region of the light spectrum, and is specifically indicated to treat mild to moderate acne vulgaris on the face through a non-thermal mechanism called endogenous Photodynamic Therapy.

The device is intended to be used in the home by people  $\geq 14$  years of age.

## 4.0. When should Omnilux Clear NOT be used?

### 4.1. Contraindications, Warnings and Precautions

#### Contraindications

- **DO NOT** share the Omnilux Clear Device with any other person. The Omnilux Clear Device is for single person use only.
- **DO NOT** use the Omnilux Clear Device to treat any other conditions apart from those listed in the intended use (mild to moderate acne vulgaris). The Omnilux Clear Device has not been tested for any other conditions than those listed and the risk is unknown.
- **DO NOT** use Omnilux Clear Device if you are pregnant, trying to conceive, or breastfeeding. The Omnilux Clear Device has NOT been tested on pregnant or breast-feeding women and therefore the risk to the foetus, neonate or pregnant woman is unknown
- **DO NOT** use the Omnilux Clear if you are a child under 14 years of age. The safety of the device for children younger than 14 years of ages has NOT been tested and the risk is unknown.
- **DO NOT** use the Omnilux Clear Device if you suffer from light induced headaches.
- **DO NOT** use the Omnilux Clear Device if you suffer from any genetic conditions of the eyes. The Omnilux Clear device has NOT been tested on subjects who suffer with genetic conditions of the eyes and therefore the risks are unknown. If you are unsure about any related eye condition and using the Omnilux Clear contact a health care professional.
- **DO NOT** use Omnilux Clear Device if you suffer from any photosensitive disorder (sensitization to light) of the eye. The Omnilux Clear device has NOT been tested on subjects who suffer with photosensitive disorders of the eye and therefore the risks are unknown. If you are unsure about any related eye condition and using the Omnilux Clear contact a health care professional.
- **DO NOT** use Omnilux Clear Device if you suffer from lupus erythematosus, photosensitive eczema or albinism. If you use the Omnilux Clear Device to treat lupus erythematosus, photosensitive eczema or albinism you may cause a severe skin reaction.
- **DO NOT** use Omnilux Clear Device if you suffer from any photosensitive disorder (sensitization to light) of the skin. If you use the Omnilux Clear Device and you suffer from a photosensitive disorder you may cause a severe reaction.
- **DO NOT** use Omnilux Clear Device if you are taking any medication that can cause photosensitivity. If you use the Omnilux Clear Device and you are taking any medication that can cause photosensitivity you may cause a severe skin reaction.

- For a full list of drugs and other substances that can cause photosensitivity, please visit: [www.omniluxled.com](http://www.omniluxled.com)

**Photosensitivity** is a common side effect of various medications. These can include certain antibiotics, chemotherapy drugs, and diuretics. If you are unsure about any medication, you may be taking consult your healthcare provider.

**Other substances not listed above can also cause photosensitivity.**

Common examples of these substances are: St John's Wort, coal tar, deodorants, antibacterial soaps, artificial sweeteners, naphthalene (mothballs), petroleum products, brightening agents found in laundry detergent, and cadmium sulphide (a chemical injected into the skin during tattooing).

#### Warnings

- Before using the Omnilux Clear you must carry out a sensitization test (see section 6.3).
- **DO NOT** use the Omnilux Clear Device in any way that is not exactly as specified in these instructions for use. Use of the Omnilux Clear other than exactly as specified in these instructions for use has not been tested and the risk is unknown.
- **DO NOT** use the Omnilux Clear if your acne pustules are oozing or bleeding. Stop treatment and resume when your pustules have formed a plug or scab.
- **DO NOT** use the Omnilux Clear in an environment that is NOT specified by the manufacturer, use of the Omnilux Clear in an environment NOT specified by the manufacturer may cause damage to the device or harm to the user.
- **TURN OFF** and **REMOVE** the Omnilux Clear immediately if it feels hot or your skin feels uncomfortably warm. STOP using the device and contact your supplier.
- **STOP** using the Omnilux Clear during a treatment if you experience severe headache, changes in vision or pain in your eyes. Should this occur, contact a healthcare professional.
- **DO NOT** use the Omnilux Clear without first removing sunscreens, topical creams, lotions, or make-up BEFORE treatment. Some ingredients in these products may cause light sensitization or may reduce the effect of the treatment. If you use the Omnilux Clear without removing these products, you may cause a skin reaction, or may reduce the effect of the treatment.
- **NEVER** use the Omnilux Clear near water or where water can be spilt onto the equipment.
- **DO NOT** use the Omnilux Clear if you think it has come into contact with water.
- Keep the Omnilux Clear away from hot surfaces or direct sunlight.
- The Omnilux Clear contains cables that could have the potential to cause strangulation – Keep out of the reach of infants.
- **DO NOT** attempt to modify the Omnilux Clear in any way. Modifying the Omnilux Clear will void your warranty. Modifying the Omnilux Clear could cause an electric shock when the system is used again. The system contains no serviceable parts.
- **DO NOT** use the Omnilux Clear if the mask, controller, charging plug or cord are damaged in any way.

## General Safety

The Omnilux Clear Device is an electrical device, which under certain circumstances could present an electrical shock hazard to the user. Please read this User Guide thoroughly and follow directions stated in the manual to assure maximum safety during operation.

Omnilux Clear may only be charged with the power supply that is supplied with the device. Use of other power supplies may not be safe and could lead to electric shock.

## Precautions

- **DO NOT** pick or rub acne pustules that have scabbed or blocked. Picking acne pustules may lead to acne scarring or worsening.
- The Omnilux Clear straps contain a mixture of polyester, nylon, and lycra. If you are allergic to any of these materials, you may not be able to use the head strap provided.
- The Facemask is made from Medical Grade Silicone, if you have an Allergy to Silicone, you will not be able to use the Omnilux Clear device.
- **The Omnilux Clear Device cannot be used in charging mode.**
- Before each use, with the device switched off, clean the face mask. All surfaces of the Omnilux Clear device may be cleaned by following the cleaning instructions in Section 8.1.
- **DO NOT** use any other cleaning products other than those stated in the cleaning instructions. Use of other cleaning products may damage the Omnilux Clear
- Store the Omnilux Clear at room temperature in dry conditions, away from direct sunlight.

## 5.0. Side Effects and Clinical Benefits

### 5.1. Potential Side Effects

- Immediately after treatment skin may appear red. Skin should return to normal within 1-2 hours.
- Some users have reported pigmentation, skin flaking or peeling, associated pruritis and dryness. This is transient and Clears within a few days.
- Temporary headaches immediately following the light treatment. Headaches can be treated with an appropriate over the counter analgesic and rehydration.

If you experience a side effect that is greater than that described above or a serious side effect not listed, STOP using the device and seek advice from a healthcare provider and report such an incident firstly to the manufacturer Globalmed Technologies, info@omniluxled.com, who will contact the Competent Authorities if deemed applicable. Country specific contact numbers for Globalmed Technologies can be found at [www.omniluxled.com](http://www.omniluxled.com)

## 5.2. Clinical Benefit

People who have undergone a course of 4 weekly treatments for 6 weeks using the Omnilux Clear report a reduction of mild to moderate acne vulgaris.

## 6.0. First Time Setup

### 6.1. Unpacking and Charging your Omnilux Clear

- The Omnilux Clear device uses a controller which contains a lithium-ion polymer battery. **Omnilux Clear cannot be used in charging mode.**
- The Lifetime of the lithium battery can be maximised by allowing it to fully charge and then fully discharge during regular use.
- Unpack the Omnilux Clear and check off the parts described in Section 3.
- Visually inspect the parts and connections. Return the device to your supplier if you think the device or connecting parts or components are damaged, missing, or have been previously used.
- Before use, fully charge the battery.
- To charge the Omnilux Clear battery insert the USB-C cable into the controller. Insert the power supply into a suitable electrical outlet and press the ON/OFF button on the controller.
- The 3 LEDs in the controller will start to flash to indicate charging. The Omnilux Clear battery is fully charged when the three indicator LEDs on the controller stay illuminated.

### 6.2. Connecting the Omnilux Clear to the Controller

- Ensure the controller is no longer charging by unplugging the controller from the power supply.
- Connect the Omnilux Clear mask by inserting the USB A connector of the mask into the USB port of the controller.

### 6.3. Sensitization test

**Before carrying out your first treatment you must carry out a sensitization test. Follow the step-by-step guide below.**

1. Connect the Omnilux Clear mask to the controller (see section 6.2)
2. Cleanse your inner arm to remove any skincare such as sunscreens.
3. Place the Omnilux Clear over your inner arm.
4. Press the ON/OFF button on the controller until lights are illuminated and allow the device to run for its full 10 minute treatment program.
5. At any time if your skin feels HOT or uncomfortable, immediately turn OFF your device and remove from skin.
6. After the treatment wait at least 6 hours. If your skin on your inner arm is red or itchy then the Omnilux Clear is unsuitable for you. Return the device to the retailer.

## 7.0. Using the Omnilux Clear Device

Regarding electromagnetic compatibility (EMC) and electromagnetic fields (EMF), in and around the home environment the Omnilux Clear should be kept away from radio frequency (RF) equipment to avoid interference.

### 7.1. Treatment Recommendations:

- People who have undergone a course of 4 weekly treatments for 6 weeks using the Omnilux Clear device report a reduction of mild to moderate acne vulgaris.

**DO NOT over treat using the Omnilux Clear. There is no additional benefit of exceeding the treatment time or treatment frequency.**

### 7.2. Wearing the Omnilux Clear Device

- Omnilux Clear comes with two adjustable straps that allow you to comfortably wear the mask.
- To fit the mask, thread the straps through the inlet holes on either side of the mask.
- Loop the straps around your head and secure firmly.
- Adjust if necessary.

### 7.3. Carrying out a treatment

1. Cleanse and exfoliate the skin to remove any traces of sunscreen, skincare, lotions, or topical creams.
2. Place Omnilux Clear on the face, secure with the adjustable straps.
3. Connect the Omnilux Clear Device to the controller and press the ON/OFF button on the controller. The device switches off after its 10-Minute treatment. You may stop the treatment at any time by pressing the ON/OFF button on the controller.
4. Remove the mask, apply your normal skincare.
5. Clean the device using a water-based wipe or damp, clean cloth. Ensure that the device is NOT being charged when cleaning.
6. For best results treat for 4 days a week for 6 weeks.

## 8.0. Care and Maintenance

- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the Omnilux Clear device flat and at room temperature in dry conditions. Use the carry bag provided to store your device.
- Store the Omnilux Clear in conditions set out in the labelling provided (Storage and Transport Temperature: -5°C to 55°C).
- Avoid storing the Omnilux Clear device in dusty conditions, in direct sunlight or in damp and humid conditions.
- The Omnilux Clear device contains no serviceable parts.
- Before each use, with the device turned off, clean the device following the cleaning instructions below.

- Before each use visually inspect the Omnilux Clear, controller, and cables for damage. **DO NOT** use the device if you notice any signs of visible damage.

If you follow the care and maintenance for your Omnilux Clear, then the Omnilux Clear facemask has a minimum lifetime of 2 years. The lithium battery within your controller has a minimum lifetime of 1 year.

### 8.1. Cleaning the Omnilux Clear Device

Before cleaning ensure that the Omnilux Clear device controller is unplugged from the Mains outlet.

**DO NOT use any other cleaning products as these may damage the plastics.**

- Clean the front and back of the device with a water-based wipe or damp, clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.

### 8.2. Disposal

This symbol indicates: **DO NOT** dispose of any part of the Omnilux Clear with normal household waste (2006/66/EC)

This symbol indicates: The Omnilux Clear controller contains a built in rechargeable Lithium-ion polymer battery. We advise you to take your product to an official collection point to have a professional remove the rechargeable battery. Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent damage to the environment and human health. Please contact customer service at [info@omniluxled.com](mailto:info@omniluxled.com) or visit our Contact page at [www.omniluxled.com](http://www.omniluxled.com) if you would like further information on proper disposal of this device.

### 8.3. Troubleshooting

Problem	Possible Cause	Action
No blue or red light from Omnilux Clear	Not charged	Plug in the controller to the charger making sure that the connections are attached and leave charging until the three indicator LEDs on the controller remain lit and <b>DO NOT</b> flash
No blue or red light from Omnilux Clear	Device connected to an electrical outlet	Unplug USB connector from controller
LEDs <b>DO NOT</b> turn OFF after 10 minutes	Device is faulty	Return to your supplier
LEDs flicker continuously during treatment	Omnilux Clear is not fully inserted into controller	Ensure the controller and Omnilux Clear device are connected correctly

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to your provider.

**DO NOT** attempt to investigate any further.

**DO NOT** attempt to modify the device in any way. The Omnilux Clear Device has NO parts you can repair yourself.

### 9.0. Labelling

The labelling supplied with the device contains the following information.

- Date of manufacture
- Catalogue Number and Batch Code.
- Operating environment
- Manufacturers Information

### 10.0. Recycling

The outer packaging is recyclable, the device components (mask, strap, controller, power supply and connecting cables) are not and should not be disposed of in household waste, see section 8.2 on Disposal for more information.

### 11.0. Technical specification

- Model Number [REF]: TN2112
- Program time: 10 minutes +/- 11 Seconds
- Wavelengths: 630nm and 415nm +/- 10nm
- Maximum Output intensity 44 mW/cm<sub>2</sub>
- Operating Temperature: 0°C to 30°C
- Storage and Transport Temperature: -5°C to 55°C
- Operating Humidity: 30% to 95% without condensation
- Storage and Transport Humidity: 20% to 85% without condensation
- Electrical Input from power supply: 100-240Vac, 50/60Hz. 0.25A
- Electrical Output; 5V, 1A
- IP22

## 12.0. Explanation of symbols

There are several safety symbols on the Omnilux Clear and associated packaging and in the User Guide. Use this guide to interpret them.

Not following these symbols may result in one or more of the following hazards: electrical shock and/or fire.

Caution/Warning

It is Mandatory to read the IFU Prior to using the Omnilux Clear device.

Please Read the IFU

Applied part type BF

Single patient multiple use

Conforms to Medical Device Regulation (MDR)- 2017/745. CE stands for 'Conformité Européenne'

Legal Manufacturer

Manufacturing Date

Keep Dry

Serial number

Model number for the Omnilux Clear Device - **TN2112**

Batch Code

Unique Device Identifier

Medical Device

Made in China

Keep away from Direct Sunlight

Temperature limit

Humidity limitation

The cardboard packaging that the Omnilux Clear came in can be recycled.

The device is protected against solid objects over 12.5mm in diameter (such as fingers) and vertical water droplets when the device is tilted 150C.

Product shall not be disposed of with normal household waste (2012/19/EU).

*This symbol indicates: Input for the power supply. This symbol means that the device may only be charged with the adapter that is supplied with the device. Usage of other adapters may not be safe and could lead to electrical shock. The adapter has the following electrical specification: **100-240Vac, 50/60Hz.***

**0.25A. Electrical Output: 5V.**

European Authorized Rep (EU only)

Australian Sponsor (AU only)

Importer (EU only)

Distributed by (EU only)

### 13.0. Warranty

Please register your warranty at <https://omniluxled.com/pages/product-registration>

#### 13.1. Standard Warranty:

Under normal use the device is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship.

#### 13.2. Warranty exceptions:

The lithium battery is guaranteed for 1 year after its original purchase date against manufacturing defects in materials or workmanship.

This warranty extends only to the original purchaser. This warranty does not cover damage caused by the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorized persons
- Connection to incorrect power sources
- Use with unauthorized accessories

#### 13.3. What to do if your product is faulty?

GlobalMed Technologies Co. or its authorized distributor shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to one year from date of original purchase.

This warrant gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries **DO NOT** allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you.

This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model, or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed.

Please register your warranty at [www.omniluxled.com](http://www.omniluxled.com).  
If the product is defective within the warranty period, please contact your provider.

**I-Smart Developments**  
**129 Green Lanes**  
**Sutton Coldfield**  
**West Midlands**  
**B73 5LT**  
**United Kingdom**

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