

STEP 1


Before using your Mini, charge the device using the supplied magnetic USB cable. A red light will appear to indicate that the device is now charging.

NOTE: To ensure a proper charging connection, first plug in the USB end of the charger into a compatible power source. Next, carefully align the two magnetic prongs on the other end of the cable with the corresponding ports on the LED device.




STEP 2

The indicator light color will change from red to green when the device is fully charged and ready for use. Detach the device from the USB charging cable to begin treatment.

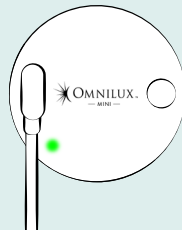

STEP 3

Peel back the Hydrocolloid Patch adhesive side 1 (marked with a device icon) and adhere to the Mini device. Peel back adhesive side 2 (marked with a face icon) and place device directly to skin where desired.

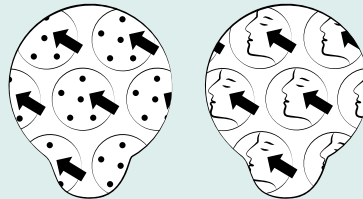


Side 1

Side 2

STEP 4

Press the power button on your Mini device to conduct your 10 minute treatment. Each session is 10 minutes and your device will automatically shut off when the treatment is complete.

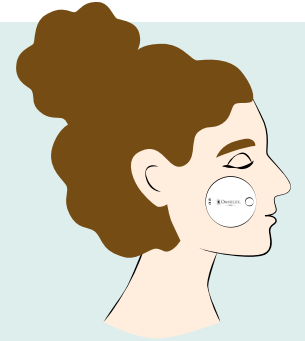

STEP 5

Store your device in a cool, dry environment and clean the treatment surface as needed with a nonirritating water-based wipe. For deeper cleaning, wipe gently with 70% alcohol.



Any questions? Scan the QR code on your smartphone or visit omniluxled.com/apps/help-center for more detailed instructions on how to get the most from your Omnilux Contour Mini.



omniluxled.com