



EVENSKYN  
*ECLIPSE*



---

## USER MANUAL

The Everyday Skin Maintenance, Toning and Lifting Handset



# TABLE OF CONTENTS

About The Eclipse Handset	4
EvenSkyn® Eclipse At A Glance	5
What's In The Box	6
Technologies	7
The Eclipse's (Backlit) Display Modes	8
Operational Overview	9
Mode Explanations & Sequence	10
Selecting Mode Intensity	11
Using The Eclipse	12-14
Usage Tips	15-17
Duration And Frequency Of Use	18
Treatable Areas	19
Results Expectation	20
Product Specifications	21
Technology Deployment By Mode	22
General Recommendations	22
Appendix I	23
Target Audience(s) & Interoperability	23
Warnings & Guidance	24-25
Warranty	26
Cleaning The Eclipse	26
Troubleshooting	26

---

## ABOUT THE ECLIPSE HANDSET

---

The EvenSkyn® Eclipse handset has been specifically designed as a daily-use facial toning, cleansing and skin maintenance device that can readily be incorporated into a user's daily morning or evening skin care routine. So that it may deliver noticeable skin rejuvenating effects, it leverages a smart array of proven and effective technologies, among which are firming microcurrent pulses, stimulating ultrasound waves, collagen-activating light therapy, and massaging sonic vibrations.

### Common concerns addressed by the Eclipse:

- **Skin laxity and sagging skin:** Helps lift and firm (more lifting & firming when compared to the EvenSkyn® Phoenix, but milder lifting & results when compared to the EvenSkyn® Lumo).
- **Dull or tired skin appearance:** Helps stimulate skin to rejuvenate it using a number of proven technologies.
- **Inadequate absorption / limited effects of skin care products:** Helps with active uptake of skin care products and enhances their deep penetration into skin.
- **Pore appearance:** Reduction in prominence of pores.
- **Excess oils, dirt, blackheads and whiteheads:** Aids in removal of excess build-up.


### The Eclipse device is ideal for:

- Individuals between the ages of 28 to 45.
- Persons with moderate signs of aging.
- Those looking for help with facial contouring and toning.
- Individuals wanting an all-in-one lifting, rejuvenating, massaging and cleansing solution.


# EVENSKYN® ECLIPSE AT A GLANCE




## WHAT'S IN THE BOX

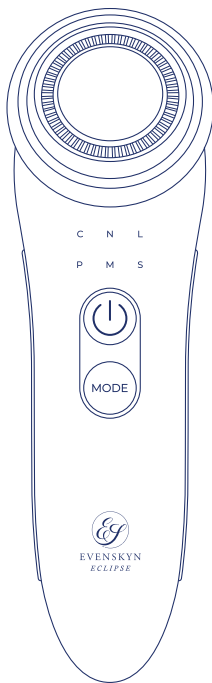


1 x Eclipse Handset

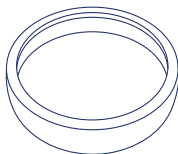


1 x Retaining Ring

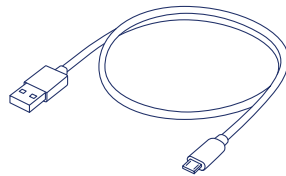

1 x USB Charging Cable

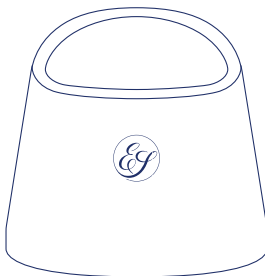

1 x Device Base

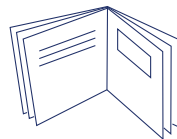

1 x User Manual

# TECHNOLOGIES

The Eclipse features a plethora of technologies that help a user achieve optimal cleansing, lifting and improvements to their skin texture. Ultrasound and microcurrent are the two main core technologies employed by the device. Depending on the mode, red and blue LED light therapy, as well as the sonic vibrations, may also be applied in an assistive capacity to bring about the desired results.

## Quick Eclipse Facts




### Medical Grade Stainless Steel:

Features a medical grade stainless steel (SAE 316) head that is a non-porous surface, meaning, bacteria, viruses and other microbes are unable to adhere or attach to it.



### USB charging:

Quick USB charging makes Eclipse a highly-portable device that one can use on the go.

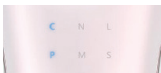

### Level 1



### Level 2



### Level 3

### Three intensity levels:


Three intensity levels for each of the modes allows for leveling up, as and when a user's experience, desire and needs require it.



### Smart contact:

Smart contact technology, adapted from the EvenSkyn® Venus, helps the device conserve energy, while allowing the device to register active usage, and thus allow for full power application / deployment.

## THE ECLIPSE'S (BACKLIT) DISPLAY MODES

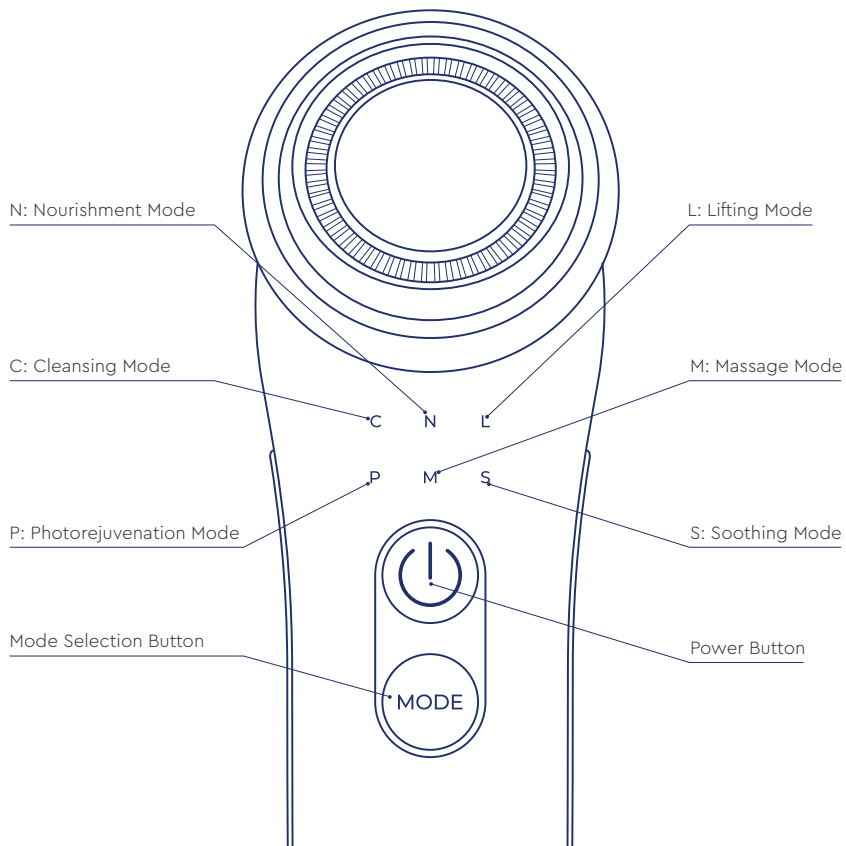



# OPERATIONAL OVERVIEW

## 1 Six Modes:

The Eclipse features the following six modes, whereby each mode is represented by a corresponding letter on the device:

- C: Cleansing Mode
- L: Lifting Mode
- M: Massage Mode
- N: Nourishment Mode
- P: Photorejuvenation Mode
- S: Soothing Mode

## 2 Mode Preview Step:

Once the device is powered on (by holding down the power button for two to three seconds), a backlit light will continuously cycle through all six modes, indicating that it is ready for mode selection.

## 3 Mode Selection:

To select a mode, continue to press the mode button until the letter corresponding to the desired mode is lit up.

## 4 Engaging Mode:

Once the desired mode is reached, as indicated by the backlit indicator light lighting up the letter corresponding to the letter, usage may commence at intensity level one (details on how to use the device follow in the next section).


## 5 Mode Intensity:

Once the usage mode has been selected, select the mode intensity by pressing on the power button. Each mode can be set at one of three intensity levels. A single vertical bar (formed by the lighting up of the letter C and P), represents the lowest intensity level, Level 1. When the first two vertical bars, or left six lights are turned on (letters C, N, P and M) this indicates that the medium intensity level, Level 2, is active. All three vertical bars (whereby all six letters are simultaneously lit up) indicates that the highest intensity level, Level 3, has been reached. Once the desired level is reached, the light will go back to the mode letter to indicate the mode is on. See next page for more detailed illustrations.

# MODE EXPLANATIONS & SEQUENCE

## 6 modes | 2 Methods:

All the modes of the Eclipse can either be used in sequence\*, 1 through 6, or a user may only use select modes based on individual concerns & preferences.

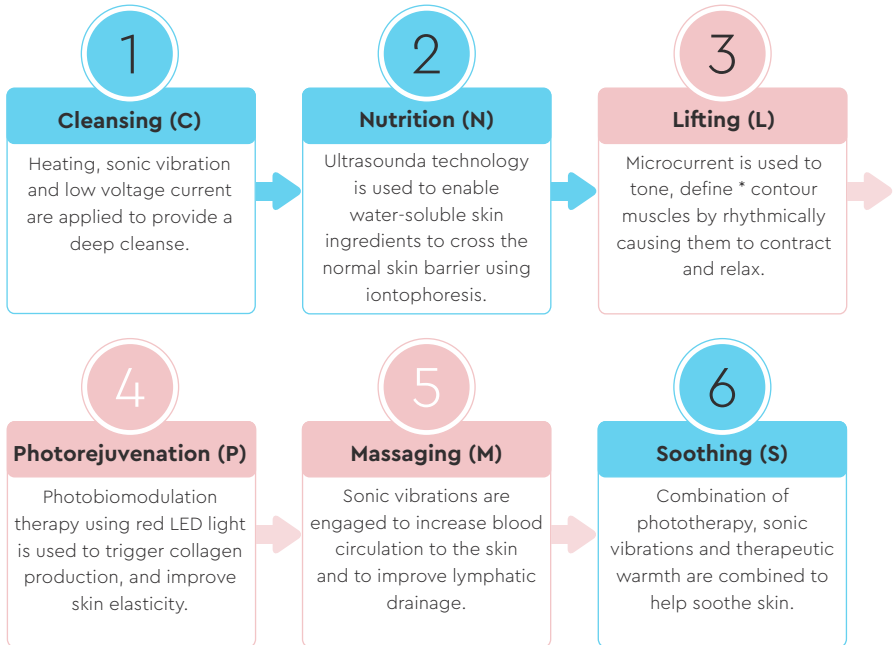


● Anti-Aging focused mode    ● General skin maintenance mode

\*All 6 modes, max. 4-times per week. See "DURATION AND FREQUENCY OF USE" section for details.


# SELECTING MODE INTENSITY

Mode intensity can be chosen by pressing the power button after a mode has been selected. Pressing the power button once to select Level 1, twice to select Level 2, and three times to select Level 3.

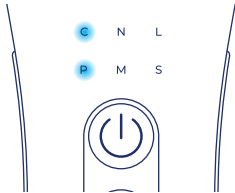


## Low Intensity (Level 1):

Low Intensity (Level 1): The base / lowest intensity is the default setting when a given mode is selected. This will be indicated when the letters and C & P are simultaneously lit up after a given mode is selected.

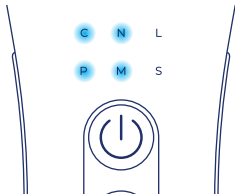

## Medium Mode Intensity (Level 2):

The medium intensity setting deploys more power than the low mode for any of the selectable modes and can be engaged when stronger stimulation is desired. This will be indicated when the letters and C, N, P & M are simultaneously lit up after a given mode is selected.

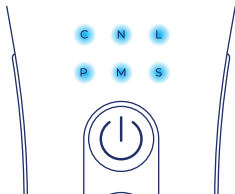


## High Mode Intensity (Level 3):


Level 3 is the highest intensity level that can be chosen for each mode of the Eclipse. When all six letters are lit up, it indicates that Level 3 has been selected and is active.

# USING THE ECLIPSE


## Mastering Use of the Eclipse in 5 Easy Steps:

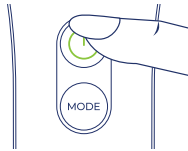


**1** Hold the power button for 2 to 3 seconds or until the backlit light starts to cycle through the different mode letters on the panel above the power button.

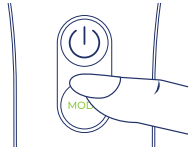

**2** Press the mode button to select the mode you would like to use the device in. The light being stationary over a given mode letter means that mode is ready to be engaged.

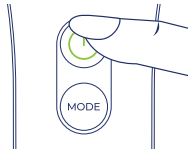

**3** Press the power button once, twice, or three times, to select low, medium or high intensity, respectively. You may now opt to go through all six modes in sequential order, C through S, or you may choose to only use specific modes of the device. If using the cleansing mode, you may attach a cotton pad to the head of the device using the plastic retaining ring around the head of the device, which is removable, and can hold the cotton pad in place. For best cleansing results, use a water-based exfoliating or toning solution. Prior to engaging the nourishment mode, apply the desired skin product to be deeply absorbed into the skin, and similarly, to ensure best results, make sure that it is also water-based. For the lifting, phototherapy and massaging modes, a dedicated conduction medium such as the EvenSkyn® Conduction Gel, or aloe vera gel, may be used for optimal results. An oil-containing skin care product, such as an anti-aging day or night cream may be applied before engaging the last (soothing) mode of device.

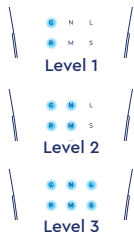


## USING THE ECLIPSE

# 4

Once the mode and intensity selection have been made, and the appropriate skin product applied, a user may begin using the device. Please note that the device will only begin working once the head of the device comes into contact with skin. Circular motions are preferable in the Cleansing, Nourishment, Massage and Soothing modes, whereas linear (outward and upward) motions are recommended for the Photorejuvenation and Lifting modes. Care should also be given to ensure all areas of the face, neck and décolletage are evenly treated. The device should stay in slow but constant motion over the treatment area, such that it is not held in one spot for a duration exceeding more than 10 seconds. This will allow for the avoidance of overstimulation and ensure that heat energy, where applied (as all modes do not engage thermal activation), is evenly distributed over the target tissue. Second and third passes over the same area are perfectly acceptable in all modes, as long as a person does not exceed the total allowable daily/weekly limits\*.

### 1. Lift the Device up to your face

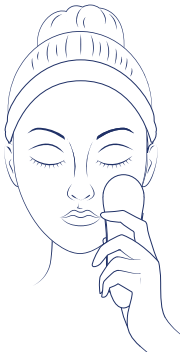

### 2. Circular Motions



### 3 Upwards and diagonally outwards



\*Refer to the 'DURATION AND FREQUENCY OF USE' section for daily/weekly allowable usage limits per mode.

## USING THE ECLIPSE



- 5** When holding the device during active use, ensure that you are either partially or fully touching and making contact with the linear touch-sensitive metal sensor strips on either side of the device. When the device is being used, the metal head of the device will be in direct contact with face and/or neck skin, and if held normally, your hands and fingers will wrap around the device to hold it, thus completing a smart conduction loop and enabling safe use of the device. This intelligent mechanism ensured the avoidance of unintentional triggering of a mode, given that when the device is not in use, since both terminals, i.e., the head and the side strips, will not concurrently make contact with skin, no mode of the device will engage fully. If not in use, this detection system also serves to conserve power and turn off the device after a prolonged period of inactivity.


## USAGE TIPS

The following usage tips are simple easy to follow steps, primarily for the lifting (L) and photorejuvenation (P) modes, in which linear upwards and outwards (from the centre symmetry of the face) motions are recommendable.




### Step 1

Lift upward from the jaw towards the ears gliding the device on your skin slowly in a circular motion.

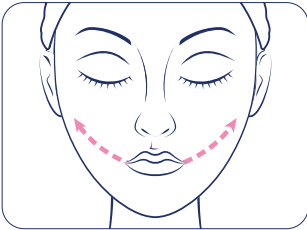

### Step 2

Lift upward from the side of your mouth towards your ears while gliding the device on your skin slowly in a circular motion.

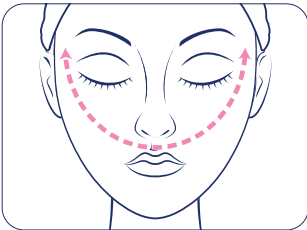

### Step 3

Lift upward from the philtrum to the place above the ears in circular motion.

---

## USAGE TIPS



### Step 4

Massage from the wings of nose to the lower eyelids in circular motion, and lift upward in circular motion to the temple.



### Step 5

Lift upward from the brows in frontal part to the hairline in circular motion.




### Step 6

Lift up slowly back above both sides of the neck but pay attention to avoid the Adam's apple.




## USAGE TIPS

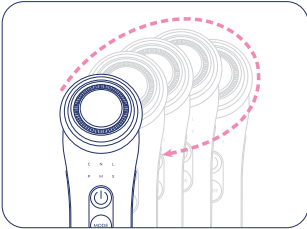


### Constant Motion Required

The Eclipse needs to be in slow but constant motion when being used. This is to ensure uniform dermal stimulation, and to avoid overheating (which can lead to skin damage and discomfort).

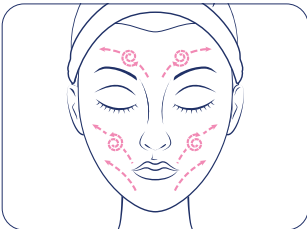

### Circular & Linear Motions

While more benefit will be derived from using the device in linear motions in the photorejuvenation (P) and lifting (L) modes, and while similarly, a user is likely to have a more optimal experience using the device in circular motions in the other remaining modes, a mix of the two motions is permissible in any mode. This is to allow for a range of scenarios of usage where, for example, a user may want to focus on one area, or where they find that a particular motion or usage regimen agrees better with their skin routine and/or skin type.

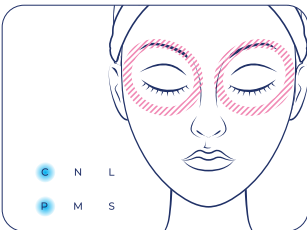

### Low Intensity Around Eyes

Use the lowest intensity (level 1) when treating the areas immediately around the eyes, such as when trying to tighten skin in the Crow's Feet region.

## DURATION AND FREQUENCY OF USE

### Weekly Allowable Use by Mode:


The following table summarizes the number of times each mode can be used per week.

**Table 1. Weekly Use Frequency by Mode:**

MODES	USES PER WEEK	MINUTES PER SESSION	MON	TUES	WEDS	THURS	FRI	SAT	SUN
Cleansing	3	4 Minutes	✓	✓	✓	✓	✓	✓	✓
Nourishment	7	3 Minutes	✓	✓	✓	✓	✓	✓	✓
Lifting	3	4 Minutes	✓		✓		✓		
Photorejuvenation	4	4 Minutes	✓		✓		✓		✓
Massage	7	3 Minutes	✓	✓	✓	✓	✓	✓	✓
Soothing	7	3 Minutes	✓	✓	✓	✓	✓	✓	✓

**Important:** Please refer to Table 2 in Appendix I to confirm that your use of the Eclipse is taking place on days when another radiofrequency or ultrasound-based device, such the EvenSkyn Lumo, is not being used. This is to avoid overstimulation using similar technologies, and thus ensuring that your use of multiple skin rejuvenation-focused EvenSkyn devices in tandem keeps such use well within safe use guidelines.

## TREATABLE AREAS

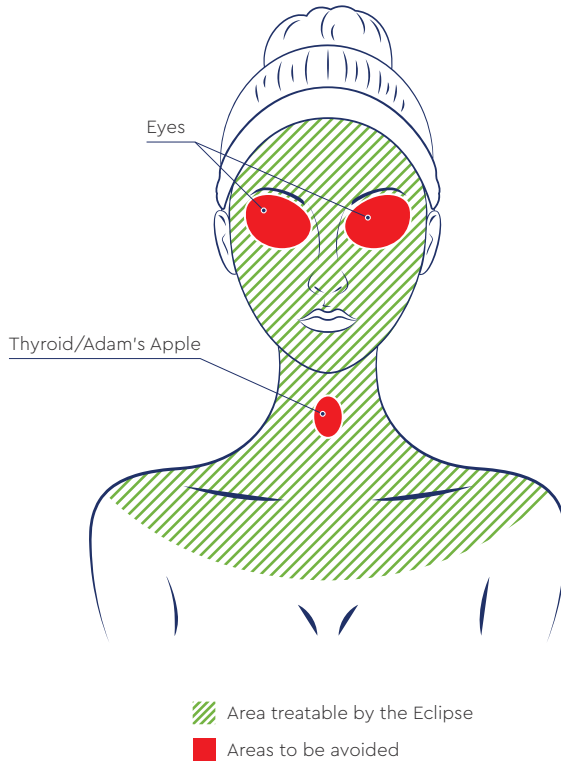

## RESULTS EXPECTATION

- **Cleansing Mode:** Use in this mode should result in the user's skin feeling cleaner, tighter and less oily, as the mode is designed to actively remove dirt, dead skin, excess oils and build-up from the skin.
- **Nourishment Mode:** Use of this mode, which causes active uptake of skin products into deeper layers of the skin, should allow for the realization of more improved results from the skin care products.
- **Lifting Mode:** The application of toning microcurrents in this mode, ideally when used with the help of a conductive medium, should cause the post-use / 'after' appearance to be noticeably lifted.
- **Photorejuvenation Mode:** Overall improved complexion and appearance of skin should result after successive and regular uses in this mode.
- **Massaging Mode:** Increased blood flow that results to the target tissue when this mode is engaged should help with improved recovery of the dermis (post-stimulation) and should contribute to a healthier glow to the face.
- **Soothing Mode:** When used, often typically as the final step, should help lock in moisture and skin products into the skin, while also calming and soothing skin, resulting in the appearance of skin looking less inflamed and reducing pore size.

# PRODUCT SPECIFICATIONS

**PRODUCT NAME:** EvenSkyn® Eclipse

**PRODUCT MODEL:** EV-1009

**BATTERY SPECIFICATION:** 3.7V, 800MAH

**ADAPTER:** Input 100~240V, 0.35A Output 5V 1A

**POWER RATING:** 3W (MAX)

**STANDBY CURRENT:** 20 $\mu$ A

**BARE WEIGHT:** ~145g

**CHARGING DURATION:** ~2 hours

## Global Voltage Compatibility

The EvenSkyn® Eclipse handset, has been designed to operate on all preexisting global voltage ranges (110-127V or 220-240V). The device also ships with the requisite charging apparatus to enable effective use of the device (adapter provision may vary by country).

## TECHNOLOGY DEPLOYMENT BY MODE

Table 2. Technologies Engaged By Mode

MODE	TECHNOLOGY					
	THERMAL ACTIVATION (42°C / 107°F)	SONIC VIBRATIONS	MICROCURRENT	IONIC (+VE / -VE)	LIGHT THERAPY	ULTRASOUND
Cleansing	✓	✓		✓		
Nourishment	✓	✓	✓	✓		✓
Lifting	✓	✓	✓			✓
Photorejuvenation		✓			✓	
Massage	✓	✓			✓	
Soothing	✓	✓			✓	

## GENERAL RECOMMENDATIONS

- The device is best used with water-based products, but other non-water based cleansing solutions will also work with the Eclipse, albeit, likely less effectively.
- For the first five modes, the device is best used with water-based products, such as the EvenSkyn® Conduction Gel. The device can also be used in conjunction with other water-based products, such as water-based toners for the Cleansing mode, water-based serum for the Nourishment mode, water-based conduction mediums such as Aloe Vera or RF / Ultrasound gels in the Lifting, Photorejuvenation and Massaging modes. The soothing mode is the only mode for which a water-based product is not required and a day cream, for example, can be applied prior to use of the Soothing mode of the Eclipse.

## APPENDIX I

Table 3. Weekly Use Frequency by Device:

DEVICES	MON	TUES	WEDS	THURS	FRI	SAT	SUN
EvenSkyn Lumo	✓				✓		
EvenSkyn Venus	✓	✓	✓	✓	✓		
EvenSkyn Mirage		✓	✓	✓		✓	✓
EvenSkyn Phoenix	✓	✓	✓	✓	✓	✓	✓
EvenSkyn Eclipse		✓	✓	✓		✓	✓
Other RF, HIFU or Ultrasound device			✓				

## TARGET AUDIENCE(S) & INTEROPERABILITY

The Eclipse having been designed as a mid-tier skin rejuvenation and maintenance handset, on the continuum of power and results, is situated squarely in the middle between the significantly lesser powered EvenSkyn Phoenix solar-powered microcurrent bar, and on the other hand, the much more powerful flagship at-home skin tightening handset device, the EvenSkyn Lumo.

It is, therefore, ideal for someone with a need for better product absorption (of skin care), skin concerns such as acne requiring deep proactive cleansing, moderate aging-related concerns such as some fine lines and minor sagging, or someone with generally good skin but one who could benefit from an 'all-rounder device' that does everything from massaging and soothing the skin to lifting and cleansing it.

The device may safely be used over the face, neck and décolletage, and can be used in conjunction with the following EvenSkyn devices: Mirage, Venus, Phoenix and Lumo. See the 'Weekly multi-device concurrent usage allowability parameters chart' in Appendix I for more information.

## WARNINGS & GUIDANCE

- Experiencing tolerable heat, as well as the sensation of stimulating microcurrents is consistent with the normal operation of the device. However, if excess discomfort results, discontinue use immediately, consult your physician, and contact our support for more guidance.
- The device has been designed for use on the face, neck and décolletage areas. Use of the device on other parts of the body is not recommended, as safety testing and efficacy have not been established for other body parts.
- Do not use the Eclipse device over areas with an active infection, areas lacking sensation, with conductive metal implants or those implants prone to deforming under moderate heat, inside the orbits of the eyes, immediately above the thyroid in the midline of the neck, over areas with uneven surfaces and those with moles or abrasions.
- Use of the device is not recommended if a user has a known existing medical conditions of the face, neck or décolletage, or other known systemic ailments such as cancer, epilepsy, seizures, rosacea, or hypersensitivity to heat.
- Use of the device is not recommended if a person has a medical device implanted in their body, especially a pacemaker or other ECG device, unless our team and your attending physician, both, have unequivocally cleared you for use of the device.
- As safety research has not been conducted using the device on pregnant individuals, those under the age of 18, and those that have had body and cosmetically-altering procedures, users who fall in any of these categories should not use the device.
- Use of the device by individuals who have had serious cosmetic or other medical procedures should only be made if your attending physician has allowed you to do so. Dermal and subdermal fillers such as Botulinum Toxin fillers, or injectable Hyaluronic Acid temporary fillers, however, do not interfere with the safe use of the device.
- A user should not attempt to open, modify or repair the device, as this may cause the device to malfunction.
- The Eclipse device should not be used while plugged in. A user should ideally charge the device fully before commencing use.



## WARNINGS & GUIDANCE

- Oil in the target treatment areas will limit the conduction, and therefore, effectiveness of the device. A user should engage in cleansing the face to rid the skin of as much oil as possible prior to using the device. The products used to enhance conduction, such as water-based Aloe Vera gel, or our dedicated Conduction gel, are void of oils and should be used with the Eclipse. Similarly, the product absorption assistance provided by the Eclipse will work where a water-based/oil-free serum is used.
- To avoid electrical shocks, and as a precaution against other possible dangerous conditions developing, please avoid submerging the product in water which may create a hazardous situation whereby electrical current may inadvertently transfer to unintended surfaces or persons.
- Battery Notice: This product has built-in rechargeable Lithium Ion batteries. If travelling with the device, it should be kept within reach in-cabin. If discarding the device, please ensure proper disposal of the battery to ensure environmentally safe disposal.
- Discontinuation due to abnormal heating: The body of the device may get warm while charging. The device's body, however, should not get excessively hot when using the device, and should this happen, a user needs to discontinue use of the device and contact EvenSkyn's support immediately. The warming of the device's head is part of the normal operation of the device, but this is limited to 42°C / 107.6°F, ensuring usability. If it exceeds this temperature, causing unbearable heating, use of the device should similarly cease and EvenSkyn's support be contacted.
- Prior to engaging in use of the Eclipse device, a user who is unclear on the safe and effective use of this device, should get in touch with our support (support@evenskyn.com or toll-free at 1-866-386-7596) to obtain the requisite clarifications.

---

## WARRANTY

- This product is covered by EvenSkyn's comprehensive 1 year (365 day) Manufacturer Warranty. The warranty period begins the day the device is purchased.
- The manufacturer, EvenSkyn Corporation, guarantees that this product is free from defects in workmanship or materials. Shipping and handling charges associated with successful warranty claims may be assessed as payable by the customer.

---

## CLEANING THE ECLIPSE

Always clean the Eclipse after each use. The head of the device is waterproof and can be wiped clean with water (and soap if needed), however, the rest of the body, especially the bottom USB port region of the device should be kept clear of water. Alcohol wipes or solution can also be used to clean the device. Once cleaned, pat dry with a microfiber or non-abrasive cloth.

---

## TROUBLESHOOTING

If the power button is held down and the Eclipse handset does not power on, the device may simply be out of charge. Plug in the device for two hours to charge using a working charging cable and then try powering on the device. Trying a different USB cable can sometimes fix the issue where an issue is limited to the USB cable and/or wall charger.


If the device is unresponsive, it could be that the internal microprocessor of the handset is temporarily malfunctioning. Pressing down the power button until the device restarts should reset the device and return it to normal operating status.





EVENSKYN  
*ECLIPSE*

@evenskyn

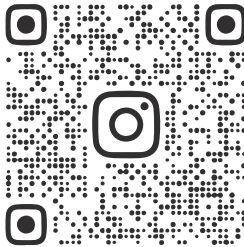

[www.evenskyn.com](http://www.evenskyn.com)

email us at [support@evenskyn.com](mailto:support@evenskyn.com) or call us at

1-866-evn-skyn (386-7596)

© 2023 EvenSkyn®