

# AVIAIR



SINCE 1968

**AVIAIR TECHNOLOGY CO., LTD.**

Designed by AVIAIR in California Assembled in China  
Copyright 2022 AVIAIR



Global



Official

# AVIAIR

## WIRELESS KNEADING MASSAGE CUSHION

Helps to relieve stress,  
pain, and soreness.

User Manual

Read & Keep These  
Important Instructions

AMC-220



## ⚠ Caution

---

Please read the following precautions before use:

People with the following conditions, do not use this product :

- |                                  |                                   |
|----------------------------------|-----------------------------------|
| 1. People with malignant tumors  | 4. People with internal implants  |
| 2. People with bleeding tendency | 5. People with acute inflammation |
| 3. Pregnant women                | 6. People with rhythm regulators  |

Consult a doctor before use to avoid accidents or physical discomfort :

1. Patients who regularly see doctors for treatments.
2. Patients with cardiac and cranial nerve abnormalities.
3. People who have fevers.
4. People with infectious diseases.
5. Patients with dermosensory disorders or skin abnormalities.
6. Patients who are taking medication.
7. People with abnormal blood pressure.
8. Patients who cannot express themselves freely.
9. Patients with perceptual disorders and diabetes.

## ⚠ Caution

---

1. Minors and people who cannot clearly communicate, do not use.
2. If product is not functioning properly, please disconnect power and contact your dealer's customer service.
3. Do not plug or unplug with wet hands. When unplugging, do not pull on the cable.
4. It is forbidden to use this product in flammable and explosive places such as gas stations, high-dust workshops, and gaseous environments.
5. Do not disassemble, modify or repair the product in any way, shape, or form.
6. This product and all accessories are not toys, please do not let children and infants play with product.
7. Do not use this product other than what the product is design for and what is written in this instruction manual.
8. Do not use this product near the heart, head, mouth, genitals, or feet.
9. Do not use in bed by lying on top of product, as excess weight and pressure may damage the product.
10. Do not use or place the product in environments with high temperature, strong electromagnetic field, thunderstorms, and direct sunlight.
11. Do not use this product in places with water, such as bathrooms or other locations with high humidity, to avoid danger and product damage.

## Caution

---

12. Do not use it while driving a car to avoid risk of physical discomfort due to excessive stimulation.
13. Do not pierce this product with sharp objects.
14. Do not use multiple sets of this product at the same time. Do not press against internal organs during use.


## Reminder

---

- > The product will automatically turn off after 15 minutes of usage. it's not recommended to use on the same part for more than 30min.
- > If you feel uncomfortable during the massage, please stop using it immediately and consult a doctor.
- > When not in use, please store in a ventilated and dry place.
- > When the product operates for a long time, the surface may feel heated, this is considered normal.

## ❁ Multi-Purpose

As people deal with high pressure and work fatigue, or extended time being in one posture using computers and mobile phones, it can easily result in muscle stiffness and pains. AVIAIR Wireless Kneading Massage Cushion helps to relieve stress, pain, and soreness.

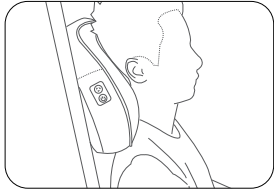


Neck




Back




Waist

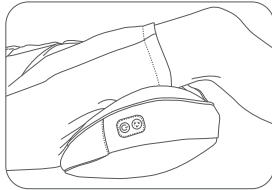

Thighs

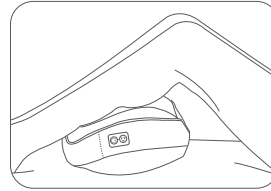


Calf's

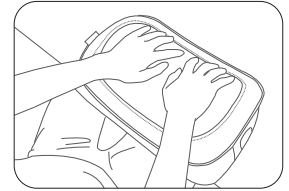

Palm




Upper Arm

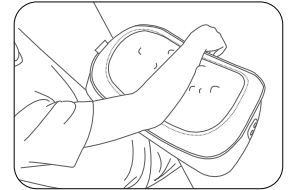

Forearm

## How To Use

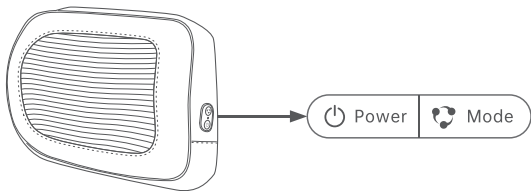

### Power Button

Press and hold for 3 seconds to enter the automatic mode: with warmth light on, kneading massage turns on for 2 minutes, then reverse for another 2 minutes. This cycle repeats for 15 minutes.

Press one time, kneading massage goes in reverse and warmth light kept on.

Press a second time, keep massage direction and warmth light is turned off.

Press the third time to change massage direction and warmth light is back on.

This cycle of functions repeats if user continues to press.

Press and hold for 3 seconds to turn off.


### Mode Button

After turned on, press to choose different modes : Soft, Soothing, or Strength.

## Charging

Insert the cable round plug into cushion's lead wire charging port. Plug in the USB end into any 5V phone charger. Charge up in about 4 hours and can be used for about 100minutes\*. Do not use quick chargers. Do not turn on while charging.

-  During charging red light stays on
-  After charging completes green light stay on

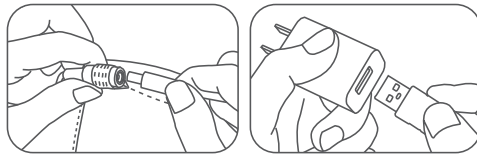


Only use 5V charger

\* Different mode strengths and pressure applied for massage will affect the length of total usage time per charge.


## Cleaning & Maintenance

Before cleaning this product, be sure to disconnect the power supply and do not touch the plug with wet hands.


1. Only the outer cover can be disassembled and cleaned. All other parts cannot get wet.
2. Before cleaning, the cover must be disassembled according to the following steps.
3. The outer cover is not machine washable, please hand-washed with gentle rubbing. water temperature cannot exceed 40°C/ 104°F.
4. The outer cover must be completely dry before it can be fitted on to the main unit.

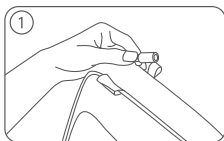


Remove the charging port from the positioning ring

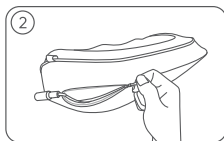

Unzip

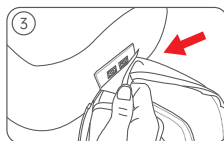

Gently tear off the Velcro

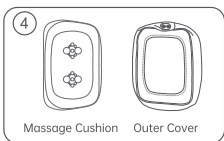


Massage Cushion Outer Cover


Remove the outer cover




Hand wash gently

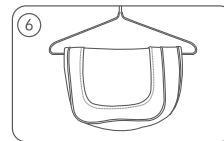


Hang to let dry

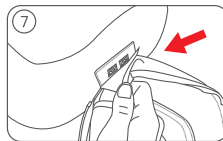


Fit the Velcro correctly to the control panel

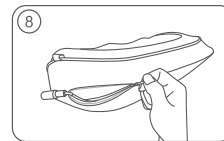

Zip up

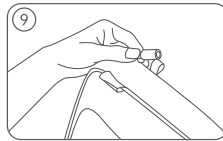

Fit the charging port into the positioning ring

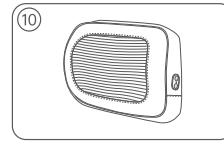

Ready to restart!

## Limited Warranty

---


This product is warranted against product design and manufacturing defects for one year from the date of your purchase. Fill out the form and keep proof of purchase to reserve your rights. This limited warranty excludes USB cable and any misuse or use that does not conform to the instructions in the manual. For any questions, please contact your dealer.

## Warranty

Customer Name	
Phone Number	
Address	
Product Name	
Model	
Purchase Date	
Dealer	
Customer Signiture	

## Specifications

---

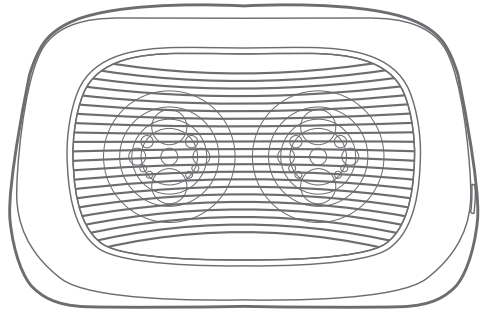

---

Model	Charging Voltage	Power	Battery
-------	------------------	-------	---------

AMC-220	5V	8W	1800mAh
---------	----	----	---------

---

Core Material	Product Dimensions	Weight
---------------	--------------------	--------

Flame Retardant Memory Foam	385 x 135 x 250 mm	1.48kg
-----------------------------	--------------------	--------