


7 Colour Neck and Face Mask

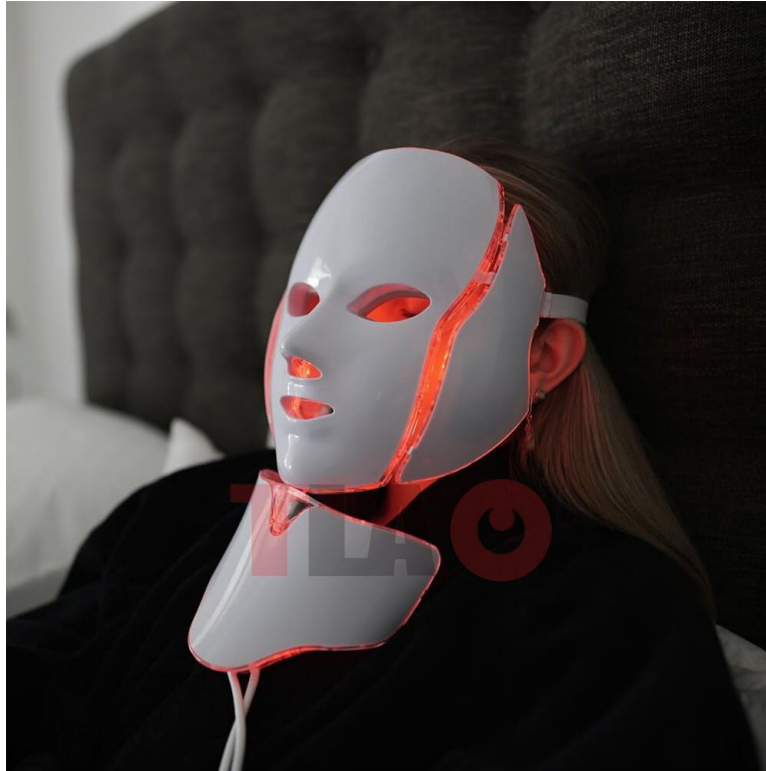

User Guide

Assembly

- Attach the microelectronics nubs before plugging in. *Note these are optional and you do not have to install them
- Install head strap and neck piece as desired.

Follow these steps

1. Plug your mask and neck piece into the remote.
2. Place the mask on in a comfortable position.

3. Turn power on and select colour using your remote
 4. If using microelectronics, chose a comfortable setting
- *Note: set all settings to comfort.
5. Test each wavelength/ colour on your skin for a short time, to make sure no adverse reaction occurs. To change colours, power output, or time settings, simply use the remote
 6. For best results: cleanse skin, before using this product
 7. Use your remote to turn your product off, and unplug between uses
 8. Don't bend the mask in any way. Be very gentle with it and store it in a safe place.
 9. Do not use wet, damp, or use in conjunction with products that cause light sensitivity

Time Guide

You can also use your mask up to 4 times a week for 20 minutes.

WARNING

1. This is an electrical device; to avoid electric shock:
 - a. Do not stand on device.
 - b. Do not touch the plug or device with wet hands or skin.
 - c. Do not use the unit if the band is wet.
 - d. Do not use if cord or plug appear to be damaged.
 - e. Do not use if the unit is damaged, torn, or has come in contact with water. f.

Do not lift or carry unit by power cord.

 - g. Do not use the cord as a handle or pull cord against sharp edges.
 - h. Do not use near hot surfaces or flammable items.
2. Do not allow children to use, this is not a toy.
3. Do not use on open wounds, heat may increase circulation and cause bleeding, and warmth may increase infection in some climates.
4. Do not use if you are allergic to neoprene.
5. Do not use with liniments, salves or ointments that contain heat producing ingredients as skin burns or irritation may occur.
6. Do not use this device if you are sensitive to light, or you are currently taking medications that increase light sensitivity. Doing so may result in a photosensitive reaction. Discontinue use immediately if you experience any discomfort or sensitivity.
7. Contact your doctor before purchasing or using this device if you have any concerns about using this device in general or while taking any current medications.

Eye Protection

EDP advises the use of eye protection with all **EDP** products that may come into your direct vision. **DO NOT STARE INTO THE LED'S FOR ANY REASON** as they are incredibly bright.

How to Test for Light Sensitivity

1. Power on your device
2. Place within 1cm of the skin on your forearm for 3 minutes.
3. Check the area, if the area has turned red, and redness lasts for more than 2 hours, you have symptoms of light sensitivity and should discontinue use of this product.

Care

To preserve the life and functionality of the device, always store in a cool, dry place when possible. For proper use and safety, please be sure to read this manual and follow all instructions and do not do any of the following with the device:

DO NOT let children use the device without adult supervision.

DO NOT use this device on infants.

DO NOT let sleeping, disabled (unable to move away from or operate the light), or unconscious persons use the device.

DO NOT exceed 20 minutes of use in one therapy session.

DO NOT use on open, fresh wounds.

DO NOT use in combination with liniments, salves, ointments, or balms that contain heat producing ingredients. A skin burn could be the result of doing so.

DO NOT modify or attempt to repair this equipment. There are no parts that are serviceable by the user. The **EDP** light is an electrical device. To avoid electric shock and other electricity-related dangers, adhere to the following instructions:

DO NOT wash electrical parts with water or other liquids.

DO NOT drop your device in water or other liquids or place it where it may fall into water or other liquids. Do not expose to running water. If your device comes into contact or becomes submerged into water or other liquid, unplug it immediately from the electrical power outlet. If device is to be used in the vicinity of water, we recommend using an electrical outlet protected by Ground Fault Circuit Interrupter.

DO NOT directly touch the device or touch the plugs or switches with wet hands.

DO NOT apply or touch the device to, or with, wet or damp skin.

DO NOT pull, carry, or lift the device by its cord. If the cord is damaged, don't use the device.

DO NOT use on or near heated surfaces.

DO NOT operate this device in areas where it could be exposed to flammable or combustible products or vapors. Explosion of fire may occur.

DO NOT store your device in hot temperatures exceeding 140 degrees Fahrenheit.

DO NOT use your device if it is damaged. Continuous use of a damaged device may result in electric shock or injury.

DO NOT use extension cords unless they are designed to carry the total wattage of the device being powered.

DO NOT unplug by pulling directly on the cord. Grasp the plug itself, not the cord, when unplugging from an outlet.

Clean plastic cover with a slightly damp cloth, It is best to carefully remove the LED pad when cleaning.

Caution

EDP recommends consulting with your healthcare provider prior to use if you believe that you are sensitive to light. Some cold and allergy medications, pain medications, and medications used to treat infection may cause light sensitivity. We also recommend consulting with your healthcare provider prior to using the device if you are pregnant, have any suspicious or cancerous lesions, or have recently had a steroid injection or use topical steroids. If you experience discomfort or have a concern about the device, stop use immediately and contact your healthcare provider. Unplug from outlet when not in use to avoid tripping hazards.

Disposal directions; found on the device or within its packaging, indicate that this device may not be disposed of with other household waste. This device requires disposal via a designated collection point for the recycling of electrical waste. Please recycle your equipment at the time of disposal to help conserve natural resources and protect human health and the environment.

Disclaimer

This product is not intended to diagnose, cure, or prevent disease. **EDP** makes no claims or representations or otherwise warranties regarding the ability of this product to cure any physical, skin or mental conditions from using this product. A qualified health professional should always be consulted prior to using this product with regards to any condition that requires medical attention