

IPL HAIR REMOVAL SYSTEM

-For Home Use

User Manual



Thank you for choosing our IPL system, for the most effective use of the product, please read the user manual carefully before using .

1.Safety Tips

To protect the skin, it's not suitable for the following customers


- You have a known skin disorder.
- You have a skin disease such as allergic dermatitis , eczema etc.
- Who have photosensitization.
- If you are experiencing Menstruation ,Pregnant or Lactation .
- If you have sunburn.
- Black or sunburn skin (Please refer to the following specific color contrast card).
- Black or fuscous nevus.
- Broken skin, a scab or healed wound, cancer or hemangioma.
- Your eyes and any artificial aesthetic, prosthetic or plastic area.

2.Product Introduction

2.1 Packing list:

IPL machine, Power adapter, Warranty card, Instruction manual.

2.2 Structure Description

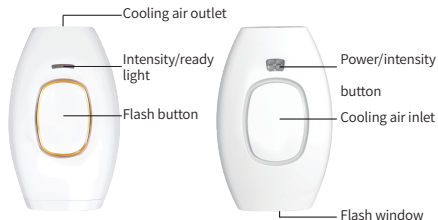

3. Technical Parameters

Lamp tube	Xenon quartz tube
Light source	Intense Pulsed Light
Energy Frequency	1.5-3.9J/cm ²
Focus size	33*10mm ²
Wavelength	470nm~1100nm
Lamp tube life	300000 times
Adaptor	AC100~240V @ 12V DC 2.5~3A
Weight	186g
Size	124*83*48.5mm ³
Operating temperature range	5~30°C
Operating humidity range	Relative humidity 30%~60%

4.Functions Introduction

- Skin rejuvenation
- Hair removal

This product can be used for: excess hair such as lip hair , armpit hair , limb body hair, forehead hair etc.; Painless hair removal is suitable for black hair and dark-hair people, it is not suitable for white, gray or blonde (Refer to the following color contrast card). Skin and Hair color contrast card

<div>Skin color</div> <div>Hair color</div>	White	Yellow	medium brown	light brown	Brown	Dark brown /black
White/grey 	X	X	X	X	X	X
Red 	X	X	X	X	X	X
Light gold 	X	X	X	X	X	X
Light brow 	√	√	√	√	X	X
Brow 	√	√	√	√	X	X
Dark brow 	√	√	√	√	X	X
Black 	√	√	√	√	X	X

5.Instructions For Use

5.1 Use a medical alcohol wipe to clean the glass window.

5.2 Press the "power button" for 1 second, it will start and enter into the working state (The green indicator light is on, and you can hear the system is running) .

Note:






1. You will see the red indicator light is on after starting up.

2. Then you should move the flash window to your skin, you will see the red indicator light flickering, because there is a human body induction sensor built-in the flash window and the sensor will sense skin when you get close to it .
3. The red light indicator must be flickering, otherwise the device will not work. If you find that the red indicator is not flickering, bring it closer to your skin until it does flicker.



Power button

- 5.3 Click the "power button" to select the suitable intensity level, there are five intensity levels indicated with the blue indicator lights. The diagram below shows intensity level from low to high power, level 1 is the minimum power setting, level 5 is maximum power intensity. Any adjustment can be made in the sequence shown with a click of the power button.

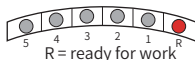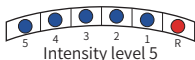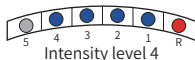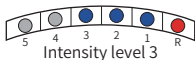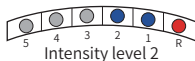

5.4 Once the intensity level is set, move the flash window to the skin which can cover the flash window completely, you will see the red indicator light flickering and then click the flash button, the device will flash one time; after every flash operation the red indicator light will be turned off to charge the power, after 1~3 seconds the light will return to indicating that it is ready for use again.

Note:

- For perfect results, it is suggested the process of hair removal is repeated 2-3 times and ensuring all the body parts should not be repeated and omitted .
- The device is not working correctly if the red indicator light remains on without flickering, in this case, please follow the steps below.
- Ensure the ‘flash window’ is close enough to the skin.
- If the red indicator light is still on ,please move to your skin for 3-5 seconds and check if it has return to a normal function.
- Please restart it if this operation failed.
- If the above operation cannot solve it , please contact our customer service.

5.5 If left idle, after 5 minutes the device will go into a stand-by mode, simply click the power button to wake up it. Note on restarting, the intensity will return back to level 1.

5.6 Press the power button for 2 seconds to shut down it.

5.7 When finished using, please clean the flash window using the suggested medical alcohol wet-wipe.

Note:

After frequent use of the device there will be dew white spot inside the lens, this will not affect the normal use and service life of the instrument.

5.8 Automatic mode/manual mode switching

The default boot is manual mode, need to short press flash button to work; If you need to switch automatic mode, should be long press flash button for 2 seconds, and then the blue led by illuminated steady into flashing, that be switched into the automatic mode, in the automatic mode the flash window began to automatically flash work when detected the skin, and the flash stop after leaving the skin, if again close to the skin will continue to flash; To exit automatic mode, you need to press the flash key for 2 seconds, and then to enter manual mode.

5.9 Overheat protection

If continuous use for a long time, the cooling environment is not ideal, the device may appear overheating protection, at this time all the blue lights flashing together, the red light is often out, the device flash function off, to the internal ambient temperature back to normal, then the equipment began to work again, but at this time the intensity level will switch to the level 1, if need the high-end intensity level, can be reset it by yourself , the internal cooling time is about 2~5 minutes after the device overheating protection, the specific time is depend on the use environment!

6.Hair Management Process

6. 1 In test before use

It is strongly advised to check for adverse reaction by conducting a skin test 2 hours in advance of any hair removal procedure.

Test method:

On a small area of skin on the inside of the upper arm (the medial).From a low energy to high energy apply 1~2 times at the same place. After 2 hours if youfeel redness, burns, stabbing pain or other symptoms consult a doctor before

using again. On slight discomfort, use the product again on a lower intensity setting, however if discomfort continues, consult a doctor before further use.

Test method:

6.2 Operating steps

6.2.1 Clean the flash window, then shave the hair on your body and clean the residual hair.

6.2.2 Start the device , choose the best suitable intensity level to hair removal.

6.2.3 After finish hair removal, you should cool the skin using Ice water or a cold damp towel.

6.2.4 You can use a toner to help your skin retain moisture after hair removal.

6.2.5 Suggestions about intensity level operations for all types of skin and different parts of the body, please refer to the below comparison table.

6.2.6 It is strongly suggested you do not use multiple times on the same part of skin. You should flash less than two times in succession only at one location. The operation frequency times should be reduced when using high energy intensity level.

Intensity Level and skin color contrast table

Skin color Parts	White color	Yellow color	Medium brown	Light brown	Brown	Dark brown/black
Armpit	3-5 level 3-4 Times	3-5 level 3-4 Times	3-4 level 2-3 Times	2-3 level 2-3 Times	2-3 level 2-3 Times	inapplicable
Arm	3-5 level 3-4 Times	2-4 level 3-4 Times	2-4 level 2-3 Times	1-3 level 2-3 Times	1-3 level 2-3 Times	inapplicable
Leg	3-5 level 3-4 Times	2-4 level 3-4 Times	2-4 level 2-3 Times	1-3 level 2-3 Times	1-3 level 2-3 Times	inapplicable
Around Mouth	3-4 level 3-4 Times	2-3 level 3-4 Times	1-2 level 2-3 Times	1-2 level 2-3 Times	1-2 level 2-3 Times	inapplicable
Face	3-4 level 3-4 Times	2-3 level 3-4 Times	2-3 level 2-3 Times	1-2 level 2-3 Times	1-2 level 2-3 Times	inapplicable

Note:

- Using for 6-10 times as a hair removal cycle.
- Use one time every two weeks for up to six times.
- Use one time every four weeks when using 6-10 times.
- If you find new thin hair grow after use it, operate one time every 2-3 months at this place , and continue gaining satisfying results.
- To gain the best results it is suggested that on areas of thick hair you use many times over a period of time.
- Use will vary from person to person, use less frequently when hair area grows by a small amount, about 30%.

7. Device Operation Notes

- 7.1 Do not flash near or directly at the eyes when in use.
- 7.2 After operation, some skin redness is a normal you can cool the area with ice water or cold damp towel.
- 7.3 Do not clean the flash window directly with spray alcohol.
- 7.4 Do not use in conjunction with other cosmetics such as water gel products.
- 7.5 Try not to eat light sensitive food during operation (such as spinach, cole, mustard, etc).
- 7.6 For best results use the higher energy setting. (Warning - only after following the previous instructions).
- 7.7 The darker your skin any skin reaction may not be so apparent and the reaction may be delayed, in this case ensure the 2 hour waiting period is followed before increasing the intensity setting.
- 7.8 Continue to look after your skin moisturizing and using Sun block in a normal way after using this product.
- 7.9 Avoid UV light during use and always protect your skin when out in sun. Also try to Avoid long exposure to computer radiation when working, protecting yourself with sun block or sunscreen.

- 7.10 No not block the Air inlet and the cooling outlet so as not to affect the cooling of the device.
- 7.11 Hair stuck to the device may cause some blank spots but it does not affect the performance of the device and the normal use;
- 7.12 In the late period of life, the inner wall of the Light bulb will appear different degrees of black phenomenon, this is a result of the normal aging process of the tubes.
Don' t worry to continue to use as normal within the specified parameters.
- 7.13 This product is not waterproof, avoid contact with water.
- 7.14 Any faults, immediately unplug the device, If you find the product cover or glass of flash window is broken, return to factory for maintenance.
- 7.15 Do not store the device near Strong magnetic fields, moisture , dust or at high temperatures.
- 7.16 Avoid damage and keep out of the reach and sight of children.
- 7.17 Prohibit using it with other application which is not allowed
- 7.18 Prohibit other people opening this device except our authorized persons otherwise it may be danger electric shock.

