

STERILIZE +

NON TOUCH THERMOMETER

INTRODUCTION


The STERILIZE+™ Non Touch Thermometer is a versatile and must-have device in the new normal. This easy-install thermometer can be used at home, in commercial areas, or in the office, which comes in a simple and classic design. It reads the temperature in 1 second and can be configured in C or F. A quick, reliable tool that keeps you and your environment safe.



FEATURES

- Safe and easy to use
- Yields results within 1 second
- Built with Infrared thermometer and ranging sensor
- Options for Celsius and Fahrenheit for measurement
- Initiates measures to detect symptoms
- Provides alert with red light and audible sound
- Industry leading 1 year warranty
- Includes a rechargeable battery and charging cord

PRODUCT DIAGRAM

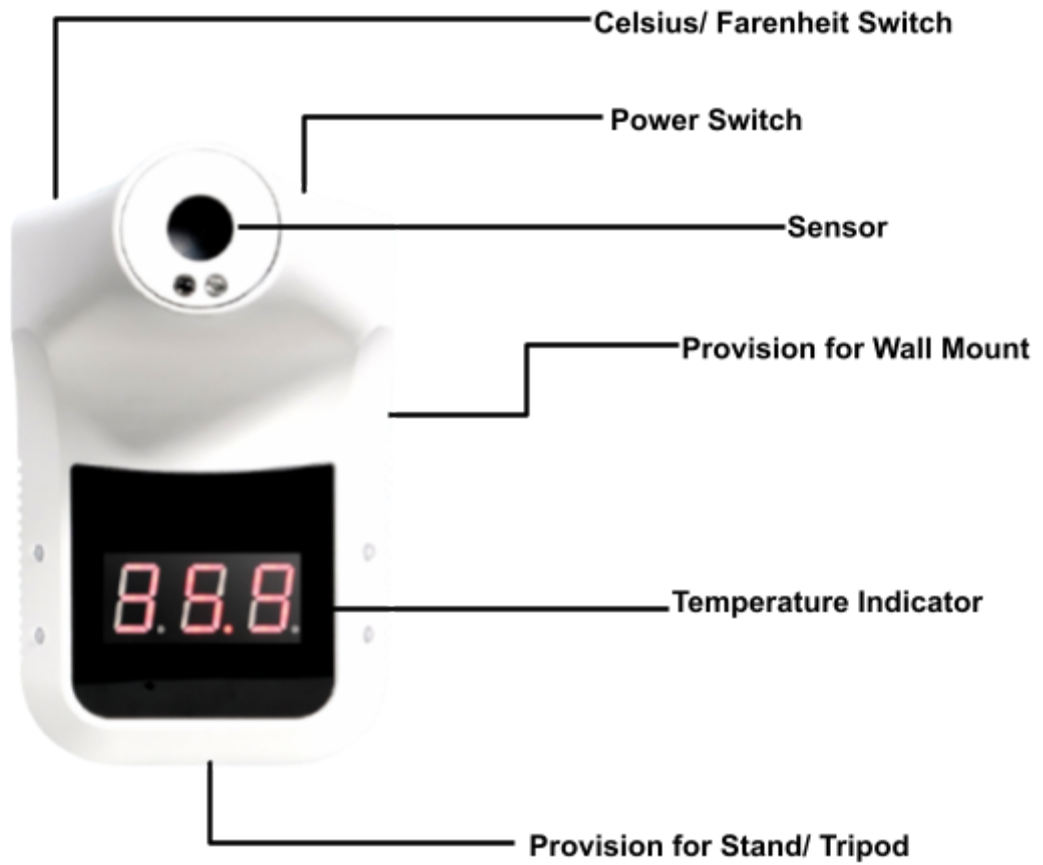

OPERATING THE STERILIZE+™ NON TOUCH THERMOMETER



[Click here to watch the Instructional Video](#)

Charging the Device

Take the USB cable and insert it into the USB connecting port. Connect to a USB power source (such as a powerbank). Charge the device for 1-2 hours depending on the power source. (Note: the battery has 300-500 charge cycles) Fully-charged devices can have a max standby time of 1 week,

Alternatively, the device can be used while plugged in.

Turning On and Off

To turn the device on, turn the power switch on the upper right portion of the device. (Refer to diagram)

THE STERILIZE+™ NON TOUCH THERMOMETER CLEANING GUIDE

- Store in a cool, dry place
- Do not leave outdoors
- Wipe with a dry towel to clean daily

SPECIFICATIONS

Size: 6.88 in x 4.52 in x 5.51 in

Weight: 1.15 lbs

Power Supply: Rechargeable 18650 battery/ USB

WARRANTY

STERILIZE+™ branded items purchased from sterilizeplusuv.com include a 1-year limited replacement warranty. If you have any questions not covered in this guide, please email support@sterilizeplusuv.com. Customer care specialists are available Monday through Friday, 9:00 a.m. to 6:00 p.m. ET.