

The **Aquafacial™** device uses a combination of vacuum suction and specialized cleansing serums to remove dead and dull skin, increase skin absorption, and revitalize your glow. The Clean Face Company™ serum enables a dual approach to cleaning. It removes oils, blackheads, and grease while replenishing the essential nutrients needed for optimal skin function. Used consistently, it delivers both instantly visible results and long-term benefits.

**Table of Contents**

Page **3** ..... Important Safeguards  
Page **4** ..... What's Included  
Page **5** ..... Assembly and Set Up  
Page **6** ..... Cleaning and Maintenance  
Page **7** ..... Using the Aquafacial Device  
Page **9** ..... About Us

# IMPORTANT SAFEGUARDS

## For Household Use Only.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### READ ALL INSTRUCTIONS BEFORE USING

#### **DANGER— To reduce the risk of electrocution:**


1. Always unplug the product immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.

#### **WARNING — To reduce the risk of burns, electrocution, fire or injury to persons:**

1. Product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug. If it is not working properly, if it has been dropped or damaged or dropped into water. Return the product to a service center for examination and repair. **Contact customer support at 1-877-CLN-FACE (1-877-256-3223)**
5. Keep the cord away from heated surfaces.
6. Never block the air openings of the product or place it on a soft surface such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any opening or hose.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Connect the power adaptor to the device and then connect it to an AC power outlet.
11. Unplug this product before filling. Fill the reservoir with water and Aquafacial Cleansing Serums only. Do not overfill or add anything other than water or Aquafacial Cleansing Serums as it may cause damage to the device.

## SAVE THESE INSTRUCTIONS

## WHAT'S INCLUDED

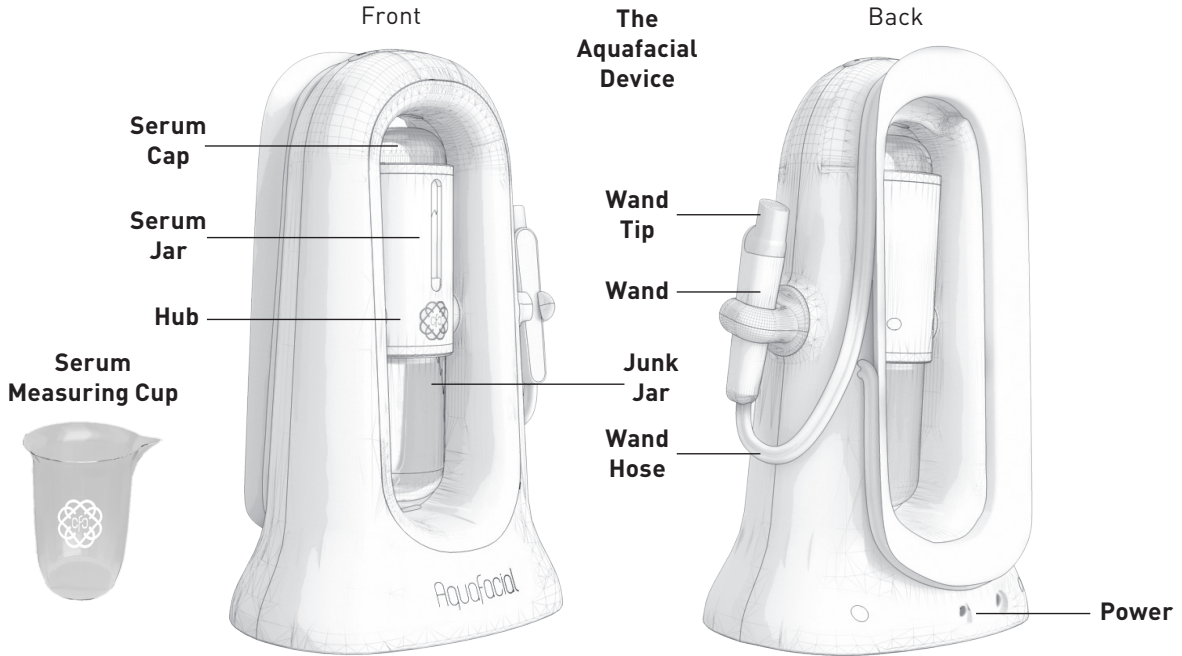

### Tips



**Gentle Cleansing Tip** - Use this tip to clean large areas of the face such as forehead, cheeks, and neck.





**Pro Cleansing Tip** - For a deeper clean to target skin impurities or blackheads. Suggested use as needed or 2-3 times per week.



**Fine Tip** - For those smaller or hard to reach areas of the face such as around the eyes, nose, or mouth.

Power Adaptor: GA-1201000, Huizhou Guoatong Technology Co Ltd, primary: 100-240Vac, 0.6A,50/60Hz, Secondary: 12VDC, 1000mA

## ASSEMBLY AND SET UP

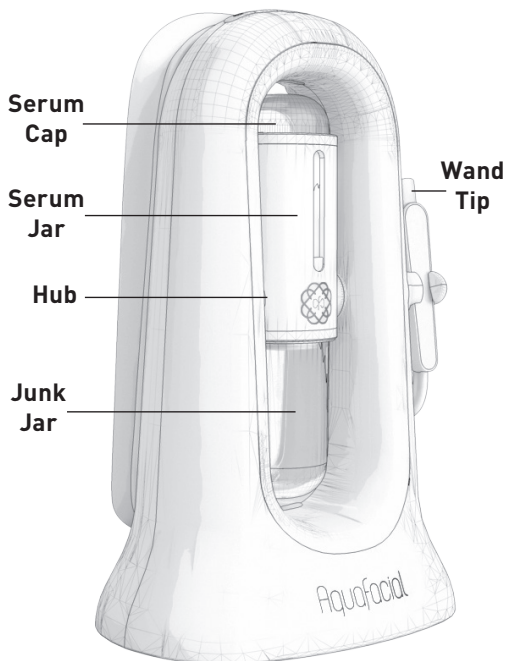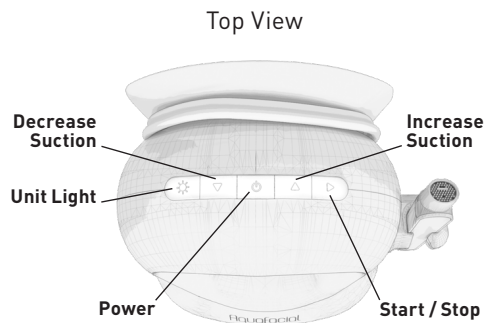

### ASSEMBLY AND SET UP

1. To add serum, rotate the top of the Hub forward and pull the Serum Cap off.
2. Fill Serum Measuring Cup with desired amount of serum and pour into the Serum Jar. Do not exceed 50 ml per use.

**Note:** Only use Clean Face Company Cleansing Serums. Adding any other serums may cause damage to the device.

3. Replace the Serum Cap by firmly pressing it back into the Serum Jar.
4. Rotate the bottom of the Hub forward, and securely screw the Junk Jar into the Hub.

**Note:** An airtight seal is required between the Junk Jar and Hub for proper device suction.

5. Rotate the Hub back to the vertical position.
6. Select your desired Wand Tip and firmly press it onto the top of the Wand.
7. Plug the Aquafacial device in and press the Power button to turn it on.
8. Select your desired level of suction by using the ▲ and ▼ buttons on the top of the device.

**Note:** When initially using the Aquafacial device, start on the lowest suction level and gradually increase over multiple uses to find your optimal setting.



9. To begin the flow of serum to the Wand, press the Start/Stop ► button and place the Wand Tip directly to the surface of your skin.

**Note:** An airtight seal between your skin and the Wand Tip is required for serum to flow. Upon initial use and if necessary, place your thumb over the Wand Tip until you begin to feel the serum flow.

10. Move the Wand across your face as shown on page 7.
11. When you have finished your facial, press the Power button to turn the device off and unplug the Aquafacial from the outlet.
12. Discard any remaining serum in the Serum Jar.
13. Empty the Junk Jar following each use. Do not refill the Serum Jar without first emptying the Junk Jar.
14. It is recommended to clean the device after each use as shown on page 6.

## CLEANING AND MAINTENANCE

Top View

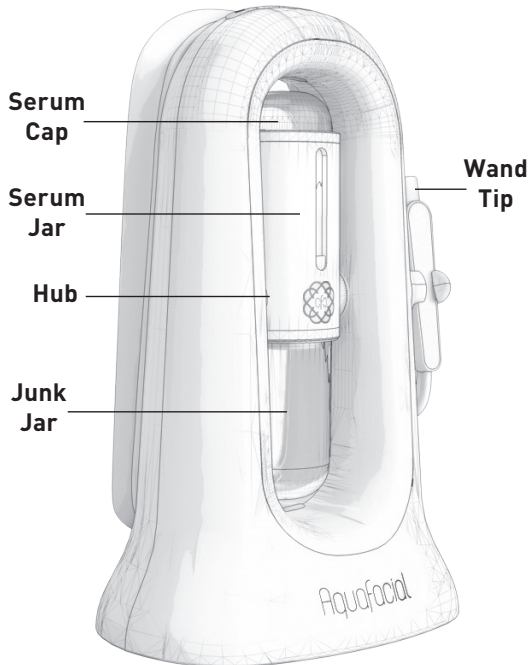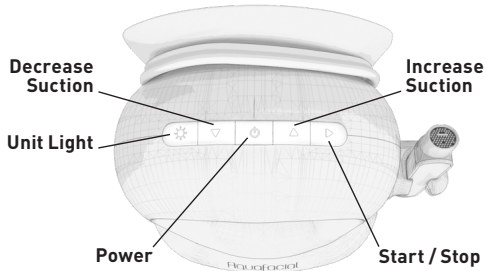

### DISASSEMBLY AND CLEANING


1. Confirm the device is turned off and unplugged from the outlet.
2. Rotate the bottom of the Hub forward and unscrew the Junk Jar. Remove the Serum Cap and the Wand Tip.
3. Wash the Junk Jar, Serum Cap, Wand Tip, and Serum Measuring Cup with warm soapy water.
4. You may use a bottle brush or sponge with warm soapy water to clean the Serum Jar.
5. Thoroughly rinse all washed parts with clean water to remove any remaining soap.
6. If needed gently wipe the Aquafacial device with a damp cloth.
7. Dry all parts and Aquafacial device with a soft, dry cloth.
8. Reassemble the Aquafacial device per Assembly and Usage instructions.

### MAINTENANCE

1. After every week or 3–5 uses, we suggest cleaning the Wand Hose.
2. To clean the Wand Hose, fill the Serum Jar with 25 ml of clean, warm tap water. Press the Power Button to turn the device on and press the Start/Stop ► button.
3. Place your thumb over the Wand Tip for 5–10 seconds, then remove your thumb for 1–2 seconds.
4. Repeat Step 3. 8–10 times to flush the clean water through the Wand Hose.
5. Once you are finished, empty the Junk Jar and any remaining water in the Serum Jar. Then dry all parts and Aquafacial device with a dry, soft cloth.

**WARNING – Do not submerge the device in water.**

## USING THE AQUAFACIAL DEVICE

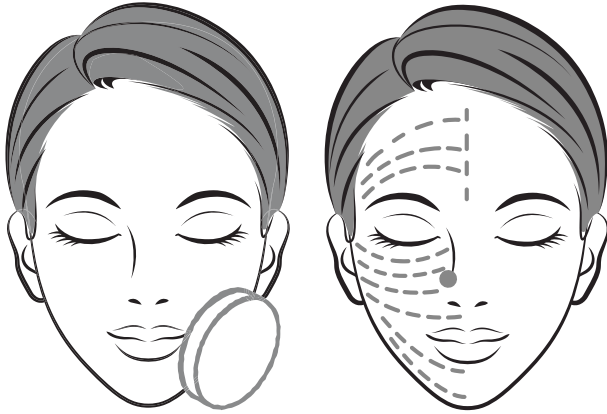

- 1. Begin by washing your face using warm water to remove any make up and pat dry with a towel.**
  - a.** Gently move the tip across the skin using short and even strokes. Always start on the lowest suction level when using the device for the first time.



### 2. Cheeks

- a.** Move the recommended tip a short distance outward from the nose towards the sides of the cheeks.
- b.** Move the tip in a downward motion towards the lips.



### 3. Around the mouth

- a.** Move the tip downward over the upper lip and outward around the lips.
- b.** Cover the chin area in downward strokes and around the lower lip region.

## USING THE AQUAFACIAL DEVICE



### 4. Around the eyes

- a. Lightly move the recommended tip around the eye region with outward strokes only.
- b. Move the tip along the ends of both eyebrows.



### 5. Forehead

- a. Move up the center of the forehead and outwards towards the temples.



### 6. Around the nose

- a. Move down along the ridge of the nose, making sure to maintain suction at all times.
- b. Move the tip sideways across each side of the nose.



## OUR PROMISE TO YOU

At Clean Face Company, our goal is to create a mindful relationship with our customers by providing clean, results-driven skincare products for the face and body. Based in Southern California, we live and love the bond with nature and keep our planet in mind. Our products are sourced and formulated with clinically proven botanical ingredients in order to provide the safest and most effective results.

We want to hear from you!

Please send your insights to:

Clean Face Company

Alpine, California

**[www.cleanfacecompany.com](http://www.cleanfacecompany.com)**

**1 877-CLN-FACE**

**1-877-256-3223**



Model: CFC-DEV-000001-501

©2022 Clean Face Company



# Aquafacial

Welcome to  
the cleanest skin  
of your life

U S E R   G U I D E

