Dermastir Gold 24k Microneedling Stamp-On

Dermastir Gold 24k Microneedling Stamp-On by Alta Care Laboratoires is a quick, effective, hygienic and secure system of delivering active ingredients and anti-age skin rejuvenation solution. In addition to the maximized nourishment process, the Microneedling Stamp-On creates microtraumas or channels, that stimulating a natural response process of the body by repairing these microwounds, triggers skin cells to produce high doses of collagen to construct firmer and vounger-looking skin.



DIRECTIONS FOR USE

Clean the Dermastir Microneedling Stamp-On by first pouring at least 75% alcohol in a clean cup. The Stamp-On should be all covered with the alcohol for about 20 minutes. Take the Dermastir Microneedling out, then immerse in a physiological saline for about 5 minutes, and allow it to air dry on paper towel. Pour the serum inside the sterilized bottle of the Dermastir Gold





24k Microneedling Stamp-On and twist the Stamp-On on top and perform different movements; from left to right. from top to bottom or diagonally.







diagonally

Dermastir Gold 24k Microneedling Roll-On

The Gold 24K Microneedling Roll-On is an effective microneedling device that helps to regenerate the skin and boost the rejuvenating processes. The Dermastir Gold Roll-On is an effective and hygienic microneedling system that delivers active ingredients to the skin and enhances the anti-aging and rejuvenating process. The Gold Roll on creates microchannels that stimulate the cells to start producing higher quantity of collagen. In this way, the skin results firmer and younger-looking.

DIRECTIONS FOR USE

First of all, disinfect the Dermastir Microneedling Roll-On in 75% alcohol solution for at least 20 minutes. Leave the Roll on in a physiological saline for 5 minutes. Attach the Roll-On to the top of one of the Sterile Vials and perform different movements: from left to right, from top to bottom or diagonally.





DERMASTIR®

WWW.DERMASTIR.COM





6, Piazza di Spagna - 00187 Roma, Italia. www.altacare.com



Dermastir Sterile Vials - Lipo

Dermastir Sterile Vials Lipo are one of the few products on the market with clinically proven fat burning ingredients. Dermastir Lipo is unique because it is effective when used topically. In many aesthetic centres, Dermastir Lipo vials are combined with radiofrequency and ultrasound machines, since they are sterile. The Dermastir Lipo vials are presented in sealed glass vials that guarantees the ingredients are present actively in high concentrations to burn down fat.

Dermastir Sterile Vials - White

Dermastir White Vials by Alta Care Laboratoires is an innovative product to treat pigmentation caused by the sun rays. Many whitening and lightening products on the market are acid-based and they cannot be used in summer. Dermastir White vials are acid-free. Alta Care Laboratoires uses peptides to brighten the skin and to treat pigmentation. Dermastir White vials can be used all the year round, even in summer to treat pigmentation.

SERUMS MONODOSE

Medical device - STERILE A



Dermastir Sterile Vials - H53EGF

The exclusive formula of the Dermastir Sterile Vials H53EGF, with 53 active ingredients, stimulates the epidermal growth and elasticity of the skin. Peptides are small enough to penetrate the skin topically. When the peptides are absorbed, the cells receive signals to start epidermal growth.

Innovative technology from Alta Care Laboratoires

A NEW LOOK AT BIOREVITALIZATION

LIQUEFY BODY FAT ACCELERATE LIPOLYSIS


5 % Lecithin

Dissolves fat deposits



1 % Carnitine
Oxidizes fat deposits

2 Biomimetic peptides Stimulate metabolism

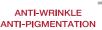





3 Biomimetic Glutathione Peptides Biotech skin lighting







7 Minerals

12 Vitamins

Natural Moisture Factor

24 Amino Acids





