


1

## GET TO KNOW YOUR NEWA

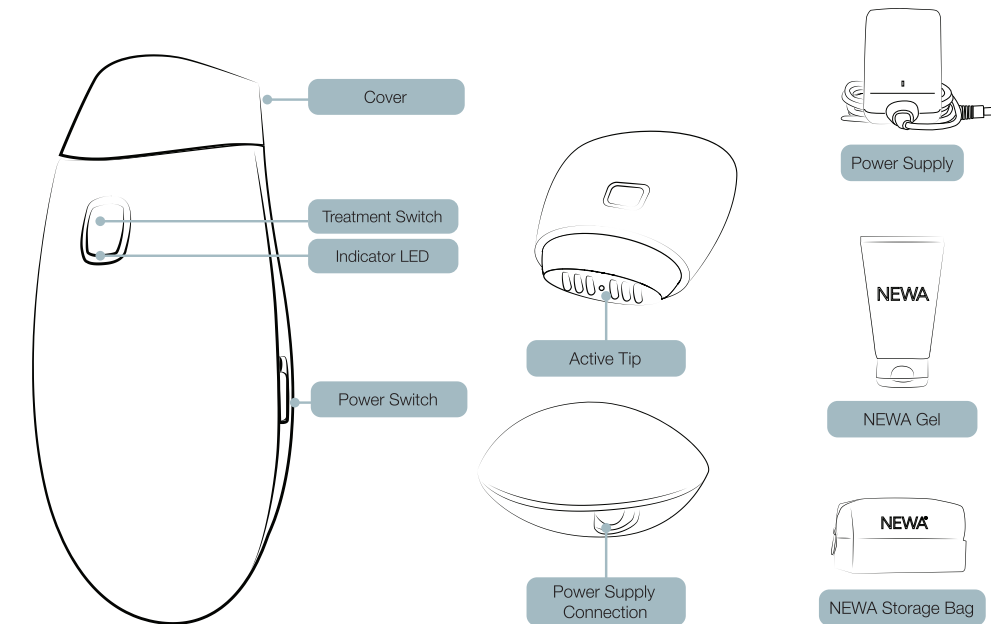

2

## BEFORE YOU START YOUR TREATMENT

## BEFORE YOU START

Take a few minutes to read these important tips to ensure you receive the safest, most effective NEWA treatments. Please note for the full contraindications / warnings we recommend you to take a moment, and review our User Manual.

## 1. Test

Before you begin, follow the instructions below and treat your forearm first to see how your skin reacts to NEWA.

## 2. Choose your treatment areas


Divide your face into 6 "treatment areas": two on each upper cheek, two on each lower cheek/jawline, and two under the chin.

## 3. Prepare yourself for the treatment

- Remove all jewelry from the treatment areas. This includes earrings, necklaces or rings.
- Thoroughly clean and dry your face. Make sure all make-up and creams have been removed from the treatment area.

3

## TREATMENT PREP

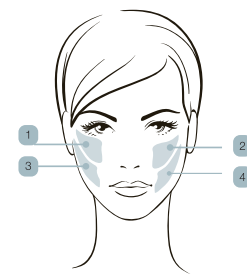

1 **Connect the power supply** to the power supply connector on the device and plug it into an electrical outlet. **Remove the cover.**




## 2 Apply NEWA Gel

Apply two thick "lines" of the NEWA Gel directly onto the active tip.

Gently place the NEWA active tip on the treatment area. With a circular motion, spread the NEWA Gel from the center to the borders of the treatment region.

You should not use any other gel product with the NEWA device.

3 **Turn ON the device** - Move the Power Level Switch to "2". The LED will emit a steady light. We recommend you start with the High heating level (level "2"), but if you find the sensation uncomfortable, you may select the lower power level (level "1").

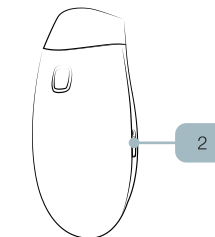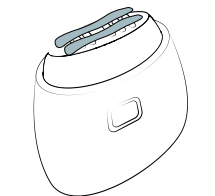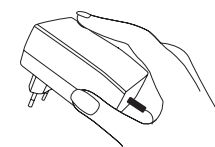

4

## TREAT

1 **Start treatment** - Press and release the treatment switch to start energy delivery. The LED will start blinking.

**Make sure all of the active tip area is in full contact**

with the skin during treatment, while maintaining a reasonable, comfortable circular movement speed. When the optimal treatment skin temperature has been reached, the LED on the NEWA changes from blue to green. If you feel that your skin temperature is too high, simply increase the size of the circular motion during treatment.

Do not stop the circular movement of the device

To stop a treatment before the 4 minute end time, just press the power button.

2 **Move to the next treatment area** - When treatment of the first area is completed (a total of 4 minutes treatment time), the LED will stop blinking and the device will vibrate to indicate that energy delivery has stopped. Apply more NEWA Gel on the active tip area, continue to the next treatment area and repeat the procedure.

3 **Complete your skincare routine** - Gently remove any excess treatment gel from your face. Cleanse your face thoroughly, than apply any skincare cream of your choice.

## Personalize your NEWA® experience

Maximize your skin rejuvenation results with the NEWA® Mobile app - an interactive resource that has everything you need to get the most out of your NEWA® experience.



For further skin care advice or information: [www.mynewa.com](http://www.mynewa.com)

### CONTACT US

**NEWA inc.**  
4400 Route 9 South  
Freehold NJ 07728, USA

Customer Support:  
Tel: 800-314-7508  
Email: [support@mynewa.com](mailto:support@mynewa.com)

Web: [www.mynewa.com](http://www.mynewa.com)

Copyright © 2012 EndyMed Medical Ltd.  
NEWA®, 3DEEP and EndyMed Medical™ are trademarks of EndyMed Medical Ltd.  
**All rights reserved.**  
Patented technology. 15

## TREATMENT FREQUENCY

In order to maximize your results with the NEWA, it is recommended you perform treatments 3-5 times a week for a period of 1 month. After the end of the first month, you may want to continue maintenance treatments twice a week. If you stop using the NEWA, your results will start to lessen. If you restart your treatments after a period of non use, you may want to repeat the full treatment series.

To check out NEWA FAQs and other user information go to [www.mynewa.com](http://www.mynewa.com)

### SAFETY INFORMATION

PRECAUTIONS AND SAFETY INSTRUCTIONS  
Remember, this is an electric device using radiofrequency waves. Please adhere to the following precautions:

- Use this device only for its intended use as described in this user manual.
- Turn the device on only during treatment.
- Do not use the device if it is disassembled, cracked or if the electric cord is damaged.
- Do not use the device in the vicinity of inflammable liquids (alcohol, acetone).
- Always thoroughly dry your hands and face before treatment.
- Keep out of reach of children.
- Do not place, store or use the device near sources of water such as the sink or bathtub. If the device is wet or falls into water don't touch the device. Unplug the power supply from the electrical outlet immediately.
- Immediately after use, or before cleaning the device, always turn the device off and unplug it from the electrical outlet.
- Only use the power supply supplied with the NEWA to use the device. Using different power supply can damage it.

### CONTRAINDICATIONS (WHEN NOT TO USE THE NEWA)

- Have an implanted pacemaker, arrhythmias or any other severe known heart disorder. (Flow of electrical energy through the body may damage the pacemaker and cause dangerous heart arrhythmias).
- Have any implantable metal device or body piercing in the treatment area. (Flow of electrical energy through the implanted metal may heat the object and cause burns).
- Are pregnant or breastfeeding. (Safety and effectiveness of the device was not tested)
- Suffer from autoimmune disorders or diabetes. (Skin healing in these cases may be impaired)
- Have Cancer of any type
- Any heart disorders
- Taking any drugs that affect skin sensitivity (it may affect heat sensation, or increase skin sensitivity to heat and may cause overheating and burns).
- Have blood clotting disorders
- Use blood thinning medications. (The mechanical massage by the device may induces skin hematomas)

### WARNINGS

- Do not use the NEWA® over eyelids, eyebrows, lips, ears, breasts or genital area, – the skin in these areas is delicate and sensitive and might be overheated.
- Do not use the device on red, irritated, inflamed or infected skin – the treatment may worsen these symptoms. Discontinue use if irritation occurs.
- Do not use the NEWA® without the supplied treatment gel. The device may not operate properly if used without the NEWA® treatment gel, or if used with another lotion or gel.
- Do not turn on the Newa device when the device is not in use. This includes instances when the Newa device is resting on a table or any other flat surface.

- During the treatment, while the device is on (LED is flashing) you should continue to move the Newa™ in a circular motion on the skin. Once movement is stopped during treatment, the skin may overheat and burns can occur.

### LABELING AND SPECIFICATIONS

Specifications	NEWA Device
Power Consumption	10W
Input Voltage	100-240VAC
Dimensions	72mm, 30mm, 120mm
Weight	100 gr
Input	100-240VAC, 50-60 Hz 0.8A
Output Voltage	9 VDC
Output Current	2A
Direction of Connector	⊕-●-⊖

Symbol	Meaning
	The user should refer to the User Manual
	CE mark
	Type BF applied part
	Waste Electrical and Electronic Equipment (WEEE) Disposal Waste electrical and electronic equipment must not be disposed of as unsorted municipal waste and must be collected separately.
	The device complies with the regulations established by the Federal Communications Commission relating to the amount of radio frequency emissions allowed such that other devices are not adversely affected.

# NEWA®

## QUICK GUIDE

Everything you need to know to start using your NEWA

To read the full user manual,  
go to [www.mynewa.com](http://www.mynewa.com)