


INSTRUCTION MANUAL

MAGMA
MICROWAVABLE
HEATING WRAP

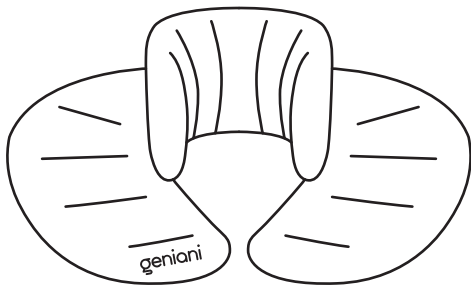

**Inquiries
Help
Issues**

1 (404) 225-3955

Monday to Friday 9 AM to 5 PM EST

info@geniani.com


PLEASE READ BEFORE USING

This step will ensure uniformity of heating and cooling

The Magma Heating Wrap is designed with evenly divided compartments throughout the product. The filling is composed of millets that can shift during transportation and storage.

For best results, we highly recommend shaking the Magma Wrap to distribute the millets evenly before use. The millets can be shifted through the gaps on the outer part of the unit, as marked in the picture.

Since the product contains natural contents, it is recommended to heat the product in the microwave 3 times, with a 20 minute break each time, to remove noxious odor. Put the Magma Wrap into the storage bag and place it in the center of the microwave. Microwave for 1 minute and 30 seconds. Wait 20 minutes. Repeat 2 more times.

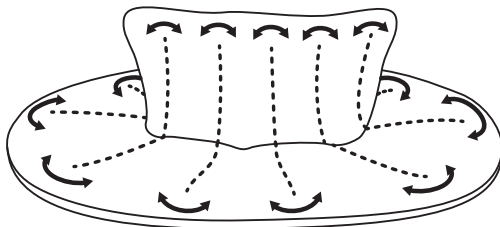



WARNING



Keep out of children's reach

3 MIN



Do not exceed a total heating time of 3 minutes



Hot! Be careful when taking out from a microwave oven



Overheating may cause the cover to ignite



Do not store in non-breathable pockets/plastic pockets



WARNING

Please read carefully and save all of the instructions.

Do not use the Magma Wrap for any other use, except its intended use.

Do not leave your Magma Wrap unattended while heating.

The microwaved item may actually be much hotter inside than it seems by touch - so take it out carefully.

As microwave ovens differ, it may be necessary to adjust temperature and time settings.

The heating of this item should only be carried out by an adult.

Not to be used by children, individuals who are ill or insensitive to heat.

Do not use on sensitive skin or areas insensitive to heat.

Discard if found worn or damaged.

Limit heat applications to 20 minutes at a time

CARE INSTRUCTIONS




Do not get the product wet



Do not wash



Do not iron

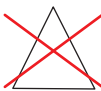

Do not bleach



Spot clean only

CARE INSTRUCTIONS

Magma Microwavable Heating Wrap contains millets and lavender.

Please do not wash or wet the wrap.

You can only clean the visible spots or stains by hand.

Do not subject the Magma Wrap to any dry cleaning, wet cleaning, or hand-washing.


Please make sure that the Magma Wrap and storage bag are clean and dry before heating them in the microwave. Place the Magma Wrap inside the storage bag and then onto a clean, dry plate. Make sure that the storage bag with the Magma Wrap inside is folded so that no part touches the sides of the microwave.

HEATING INSTRUCTIONS

Put the Magma Wrap into the storage bag and place the latter in the center of the microwave.

Please make sure that the microwave is absolutely clean.

Microwave for 1 minute and 30 seconds (1:30). Please flip the unit over after 45 seconds (0:45) and continue heating the remaining time.

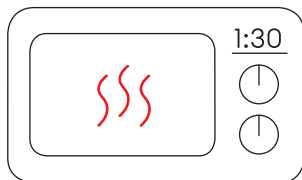

Heating time frames are based on a 1000W microwave. Please adjust according to your microwave wattage.

Before use, please carefully check the temperature of the Magma Wrap. If not warm enough, microwave again for 30 seconds.

Never exceed a total heating time of 3 minutes.

If it is too hot, allow the Magma Wrap to cool.

Shake the product to distribute its contents evenly before carefully placing it around your neck and shoulders.

Your Magma Wrap is now ready to wear! **Expect 20 minutes of warm muscle relaxation.**

Allow your Magma Wrap to cool at room temperature before heating again.

COOLING INSTRUCTIONS

Put the Magma Wrap into the storage bag and place it in the freezer.

Chill the Magma Wrap in the freezer for 2 hours. Leave it in the freezer for a longer period if additional cooling is required.


Take the Magma Wrap out from the freezer. Carefully touch to test. **If it's too cold, allow the product to warm up before using it.**

The Magma Wrap is now ready to wear!

Shake the product to distribute the contents evenly before carefully placing it around your neck and shoulders.

Relax and enjoy it! **Expect 10 minutes of cooling therapy.**

You can put your Wrap back into the freezer any time.

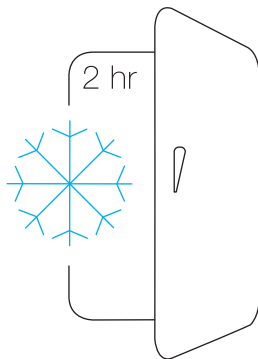

SPECIFICATIONS


Model:	Magma
Size:	Universal
Fabric:	100% polyester
Ingredients:	Milletts, Lavender
Weight:	3.02 lbs \pm 5%

TIPS

Once the Magma Wrap is heated or chilled, it is ready for use. Place the unit over your shoulders and raise the collar around your neck.


Sit back and relax as the Magma Wrap soothes you. The Wrap stays in place while walking, standing, or sitting.

It is best to keep a good posture while using the product.

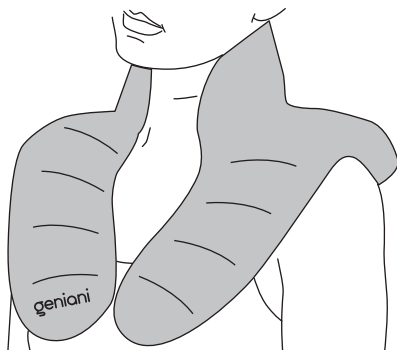



Aromatherapy
Scent for Stress Relief

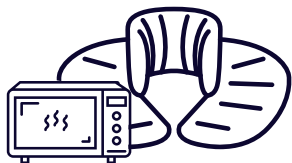

Heats in Microwave
for Hot Therapy



Cools in Freezer
for Cold Therapy



Ultra-soft
Fabric

CUSTOMER SUPPORT



**Inquiries
Help
Issues**

1 (404) 225-3955

Monday to Friday 9 AM to 5 PM EST

info@geniani.com

We are proud of our perfect customer service, and we promise to treat you like family. Feel free to contact us!

To get more information,
please visit our website
www.geniani.com

