



geniani

LAVA

PLEASE READ AND SAVE
THESE INSTRUCTIONS

INSTRUCTION MANUAL

NECK & SHOULDERS HEATING WRAP

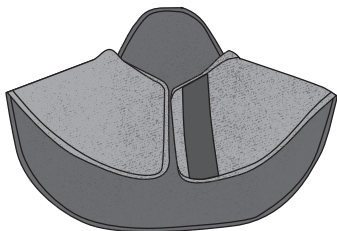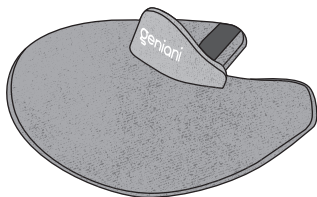

FOR PERSONAL
USE ONLY

Find Us Online

Search for **@GenianiProducts** • **#Geniani** on social media for wellness tips, special deals, giveaways, inspiration, and more.




**QUESTIONS OR
CONCERNS?**

PLEASE CONTACT US

1 (833) 780-5026

Mon-Fri, 9:00 am - 5:00 pm PST/PDT
or at info@geniani.com

CONTROLLER INSTRUCTIONS

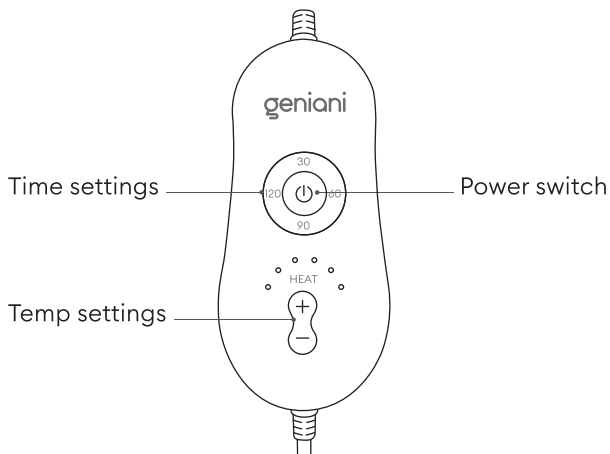

ON /OFF

The digital controller allows to turn the unit ON or OFF, simply by pressing the button 

MEMORY FUNCTION

The default setting for the Lava Heating Wrap is 30 minutes for the switch-off timer and 48°C (118.4°F) for the temperature level. The controller remembers the temperature and time that were set last time. When the heating wrap turns on, it will have the parameters of previous use.

CONTROLLER INSTRUCTIONS

HEAT SETTINGS

The Lava Heating Wrap has six levels of temperature: 40°C (104°F), 44°C (111.2°F), 48°C (118.4°F), 52°C (125.6°F), 56°C (132.8°F), 60°C (140°F). Each has a light indicator from the left to right on the controller. Press “+” button and the indicator lights will flash from the left to right and stop at the chosen temperature. If the temperature need to be decreased, press the button “-”.

AUTOMATIC SWITCH-OFF

The digital controller allows setting the timer automatically to switch off. The Lava Heating Wrap has four switch-off settings: 30 minutes, 60 minutes, 90 minutes, and 120 minutes.

Simply press the suitable interval to set up the timer. Once the timer counts itself down, the heating wrap will automatically switch off. The default setting of the switch-off timer is 30 minutes.

IMPORTANT NOTE

Always unplug this device when not in use.

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY



DANGER! Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire).



Do not sit on the heating wrap



Do not use while sleeping.
Never position a heating wrap underneath the body



Do not bleach



Do not iron




Do not fold the heating wrap

IMPORTANT SAFETY INSTRUCTIONS


DO NOT DESTROY



Not for children under
3 years of age

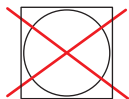


Do not tumble dry

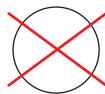

Do not dry-clean



Cold hand wash or gentle cold
machine wash

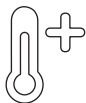

Avoid hot surfaces



Dispose of packaging
in an environmentally
friendly manner

IMPORTANT SAFETY INSTRUCTIONS

DANGER! To reduce the risk of burns, electric shock, and fire, this product must be used in accordance with the following instructions.

1. Do not use while sleeping.
2. Do not use on an infant.
3. This pad is not to be used by or on an invalid, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
4. Do not use on areas of insensitive skin.
5. Burns can occur regardless of control setting, check skin under pad frequently.
6. Never use pad without cover in place (all household pads if provided with supplementary cover).
7. Do not use pins or other metallic means to fasten this pad in place.

IMPORTANT SAFETY INSTRUCTIONS

8. Do not sit on or crush pad - avoid sharp folds.
9. Never pull this pad by the supply cord.
10. Do not use the cord as a handle.
11. Carefully examine inner cover before each use. discard the pad if inner covering shows any signs of deterioration (all household pads if provided with supplementary cover).
12. Loop cord loosely when storing, tight wrapping may damage cord and internal parts.
13. Read and follow all instructions on box or packed with pad before using.

Thank you for purchasing Geniani heating wrap, designed to provide soothing heat therapy to areas of stress in the neck, shoulders, and upper back.

**Ensure this pad is used
only on a standard 120V outlet.**

POLARIZED PLUG

This heating wrap has a polarized plug (one blade is wider than the other) as a safety feature (Figure 1). This plug fits into a polarized outlet only in one way.

If the plug does not fully fit into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.

HOW TO USE YOUR HEATING WRAP

The wrap connector must be connected to the heating wrap connector before plugging the unit into an electrical outlet (Figure 2).


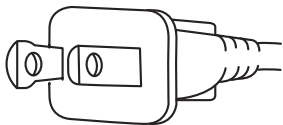


Figure 1

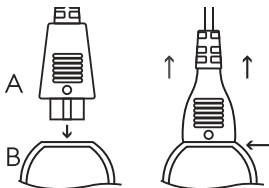


- A. Control connector
- B. Heating wrap connector


Figure 2

OPERATING INSTRUCTIONS

NOTE: If the unit does not turn on when you press the power button, make sure the control connector is securely plugged into the heating wrap connector.

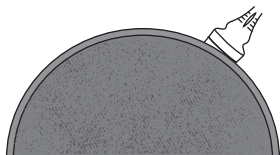

Body side
soft flannel fleece

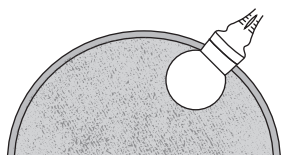

Exterior side
soft flannel fleece

Neck & Shoulders Heating Wrap Lava is designed to stay in place throughout the entire treatment session. The extra-soft fleece fiber surface and the breathable material make it especially comfortable. Use for just 20 minutes daily (or as recommended by a doctor) to help temporarily reduce pain and muscle tension.

CARE AND MAINTENANCE

- Do not dry clean.
- Do not spray this device with any cleaner, alcohol based solvent, or any other chemical.
- Do not bleach, wring, tumble dry, or iron.
- Clean with soft detergent or foam soap.
- Wash by hand or as follows in automatic washer using slow agitation and normal spin speed: Fill washer with warm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure heating pad before washing. Add heating pad and soak five minutes. Agitate two minutes. Spin. Fill for cool rinse, agitate one minute then spin. Dry in automatic dryer on high heat, five minutes only, or line dry. Remove damp heating pad from dryer, block to original size, and drape to finish drying. Do not iron your heating pad.

SPECIFICATIONS

Size:	22" x 24"
Model Number:	H31N1
Voltage:	120 V
Frequency:	60 Hz
Wattage:	75 W
Material:	Flannel fleece top and bottom
Heat Settings:	6
Auto Shut Off Timer:	30, 60, 90, 120 minutes
Cord:	9 ft.

**For more information, visit our website at
www.geniani.com**

WARRANTY

The manufacturer does not bear any responsibilities when personal injury or property loss is caused due to the following reasons:

- Damage to the product due to improper use or misuse of the product.
- Alteration, change, maintenance of the product with other equipment without following the instructions provided by the manufacturer.
- After verification, if the defects are due to improper operation during the transportation of the product.
- After verification, if the problem or dispute is caused by the quality specification or the performance of the parts and components that were produced by other manufacturer.
- In case of operation, repair or maintenance of the unit without following the manual's instructions or related regulations.
- If the damage is caused by natural calamities, improper environmental use or force majeure.

CUSTOMER SUPPORT

If you have any questions or concerns about your new product, please contact our helpful Customer Support Team

*Please have your order invoice and order ID ready before contacting Customer Support.



QUESTIONS OR CONCERNS?

PLEASE CONTACT US

E-mail: info@geniani.com

Phone: 1 (833) 780-5026

Support Hours: Mon–Fri,
9:00 am–5:00 pm PST/PDT



SCAN ME

We are proud of our exquisite customer support service, treating everyone like they're part of our family.

For more information, please visit our website www.geniani.com

Contact us anytime!

GENIANI

800 Progress Center Ct, Ste 100,
Lawrenceville,
GA 30043, United States
Phone: 1 (833) 780-5026

EXTEND YOUR WARRANTY BY **ONE YEAR**

Register your product within **14 days** of purchase at www.geniani/warranty to extend your one-year warranty by additional 1 year.