

geniani

PLEASE READ AND SAVE
THESE INSTRUCTIONS

INSTRUCTION MANUAL



VISTA MICROWAVABLE HEATING PAD

Size 8"x17"

Find Us Online

Search for **@GenianiProducts** • **#Geniani** on social media for wellness tips, special deals, giveaways, inspiration, and more.



**QUESTIONS OR
CONCERNS?**

PLEASE CONTACT US

1 (833) 780-5026

Mon-Fri, 9:00 am - 5:00 pm PST/PDT
or at info@geniani.com

BEFORE FIRST USE


The Vista Heating Pad is designed with evenly divided compartments spread throughout the product. The filling is composed of clay beads, millets, and flaxseed which may shift during transportation and storage.

For best results, we highly recommend shaking the heating pad to distribute the filler evenly before use. This enables the filler to be shifted through the gaps on the outer part of the unit, as shown in the picture.

Since the product contains natural contents, it is recommended to heat the product in the microwave 3 times, with 20-minute breaks in between, to remove noxious odors.

Place the Vista Heating Pad in the center of the microwave. Microwave it for 1 minute and 30 seconds. Wait 20 minutes. Repeat 2 more times.

Heating time frames are based on a 1000W microwave. Please adjust according to your microwave's power.

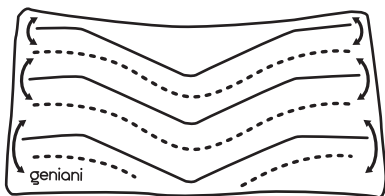



WARNING



Keep out of children's reach



Hot! Be careful when taking out from a microwave oven



Overheating may cause the cover to ignite



Do not store in non-breathable containers/plastic containers




WARNING

1000w



For microwave with
1000w power

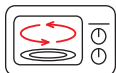

The microwave plate
must be on and spinning
throughout the heating cycle



If the pad comes in contact
with any oil or food, the fabric
will burn



Never microwave without
a plate!



WARNING

Please read carefully and save all of the instructions.

Do not use the Vista Heating Pad for any other use, except its intended use.

Do not leave your Vista Heating Pad unattended while heating.

The microwaved item may actually be much hotter inside than it seems by touch - so use caution when taking it out.

As microwave ovens differ, it may be necessary to adjust temperature and time settings.

The heating of this item should only be carried out by an adult.

Not to be used by children, or individuals with health conditions.

Do not use on sensitive skin, or areas insensitive to heat.

Discard if found worn or damaged.

Limit heat applications to 20 minutes at a time

CARE INSTRUCTIONS



Do not get the product wet



Do not wash



Do not iron



Do not bleach



Spot clean only

CARE INSTRUCTIONS

Vista Microwavable Heating Pad contains clay beads, millets, flaxseed, and lavender.

Please do not wash or wet the pad.

You can only clean the visible spots or stains by hand.


Do not subject the Vista pad to any dry cleaning, wet cleaning, or hand-washing.

Please, make sure that the Vista Pad is clean and dry before heating it in the microwave. Place the Vista pad onto a clean, dry plate. The Vista pad should be folded so that no parts touch the sides of the microwave.

HEATING INSTRUCTIONS

Put the Vista Heating Pad on a clean plate and in the center of the microwave.

Please make sure that the microwave is thoroughly cleaned.

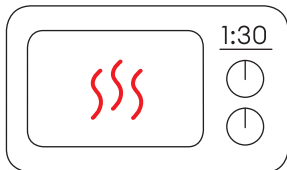

Microwave for 45 seconds (0:45). Turn the pad over and continue heating for another 45 seconds (0:45).

Heating time frames are based on a 1000W microwave. Please adjust according to your microwave's power.

Before use, carefully check the temperature of the Vista pad. If not warm enough, microwave again for 30 seconds.

Never exceed a total heating time of 3 minutes.

If it is too hot, allow time for the Vista heating pad to cool down.

Shake the product to distribute its contents evenly before carefully placing it.

Your Vista pad is now ready to wear! Expect 20 minutes of warm muscle relaxation.

Allow your Vista Pad to cool at room temperature before heating again.

If the heating pad is used 2 or 3 times continuously, a decrease in moisture output may occur. To avoid this, simply allow the pad to update moisture from the environment, and wait 2 to 3 hours before reusing.

COOLING INSTRUCTIONS

Put the Vista Pad in the freezer.

Chill the pad in the freezer for at least one hour. Leave for a longer period if additional cooling is required.


Take the Vista Pad out from the freezer. Carefully touch to test. **If too cold, allow the Vista Pad product to warm a bit before using.**

The pad is now ready to wear!

Shake the product to distribute the contents evenly before carefully placing it.

Now relax and enjoy! **Expect 10-15 minutes of cooling therapy.**

You can put your Pad back into the freezer at any time.

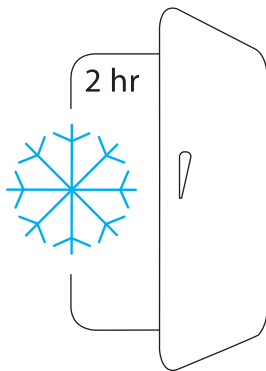

SPECIFICATIONS


Model:	Vista
Size:	8"x17"
Fabric:	100% Polyester
Ingredients:	Milletts, Flax seeds, Clay beads, Lavender
Weight:	3.3 lbs \pm 5%


TIPS

Once the Vista Pad is heated or chilled, it is ready for use. Place the Pad. It can be applied to your neck and shoulders, back and stomach, arms and leg, etc.


Sit back and relax as the Vista Pad soothes you.

It is best to keep a good posture while using the product.

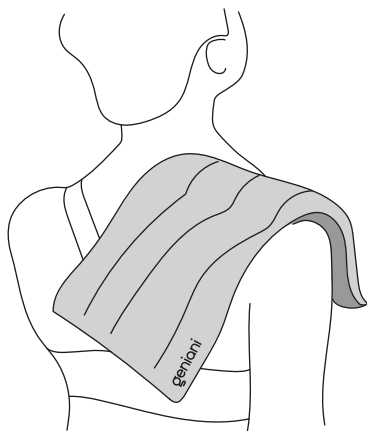


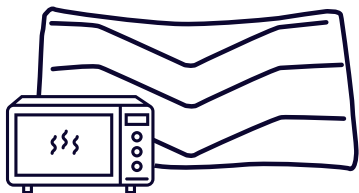


Heats in Microwave
for Hot Therapy

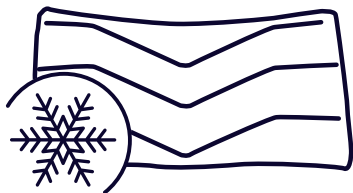

Cools in Freezer
for Cold Therapy

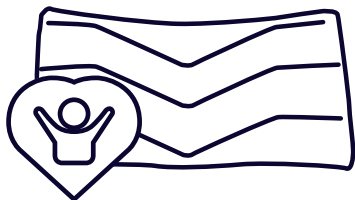

Aromatherapy
Scent for Stress Relief

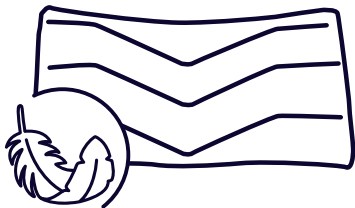

Ultra-soft
Fabric

CUSTOMER SUPPORT

If you have any questions or concerns about your new product, please contact our helpful Customer Support Team

*Please have your order invoice and order ID ready before contacting Customer Support.



QUESTIONS OR CONCERNS?

PLEASE CONTACT US

E-mail: info@geniani.com

Phone: 1 (833) 780-5026

Support Hours: Mon–Fri,
9:00 am–5:00 pm PST/PDT



SCAN ME

We are proud of our exquisite customer support service, treating everyone like they're part of our family.

For more information, please visit our website **www.geniani.com**

Contact us anytime!

GENIANI

800 Progress Center Ct, Ste 100,
Lawrenceville,
GA 30043, United States
Phone: 1 (833) 780-5026