




ANJOU MOLE REMOVAL PEN User Guide

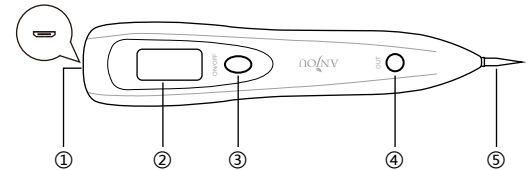

Product Diagram


- ① Micro USB Charging Port
- ② LED Display
- ③ Power Button
- ④ Work Button
- ⑤ Pen Head

Package Contents

- 1 x Anjou Mole Removal Pen
- 2 x Coarse Needles
- 40 x Fine Needles
- 1 x Pinhead Cover
- 1 x Charging Cable
- 1 x User Guide

Charging the Pen

1. Charge the pen before use.
2. Charge via USB cable.
3. The LED display will light up when charging and turn off when fully charged.

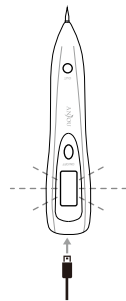

ANJOU
www.anjou.com

NORTH AMERICA

E-mail: support@anjou.com
Tel: 1-888-456-8468 (Monday-Friday: 9:00 – 17:00 PST)
Address: 46724 Lakeview Blvd, Fremont, CA 94538

EUROPE

E-mail: support.uk@anjou.com(UK)
support.de@anjou.com(DE)
support.fr@anjou.com(FR)
support.es@anjou.com(ES)
support.it@anjou.com(IT)
EU Importer: ZBT International Trading GmbH,
Halstenbeker Weg 98C, 25462 Rellingen, Deutschland

ASIA PACIFIC


E-mail: support.jp@anjou.com(JP)
Tel: 03-5542-0907 (10-18時/土日祝除く)
JP Importer: 株式会社SUNVALLEY JAPAN
Address: 〒104-0032 東京都中央区八丁堀3丁目
18-6 PMO京橋東9F

MANUFACTURER

Shenzhen NearbyExpress
Technology Development Co., Ltd.
Address: 333 Bulong Road,
Shenzhen, China, 518129

AJ-PCA025

Replacing / Installing the Needle

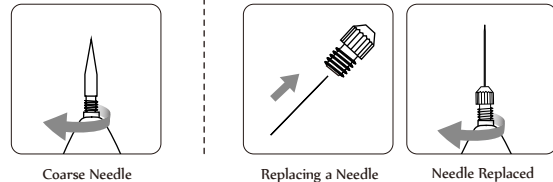

Choose an appropriate needle as per your skin conditions.

1. Installing the Coarse Needle: Insert the coarse needle into the pen head and turn clockwise to tighten.
2. Installing the Fine Needle: Thread the fine needle through the pinhead cover and then insert the pinhead cover into the head of the pen; and then tighten clockwise.
3. Sterilize the needles by boiling them in water or wipe them with alcohol pad before use.

Removing Moles / Freckles / Tattoos

1. Select an appropriate needle to the size of your target spot and install it on the pen.
2. Wipe the target spot with alcohol pad.
3. Apply anesthetic the size of a coin onto the target spot. And clean it after half an hour.
4. Press and hold the Power Button to turn on / off the pen.
5. Press the Power Button to cycle through 9 levels of rotation speed.
6. Press and hold the Work Button to start removing moles, release to stop working.
7. Tighten your skin to relieve the pain.
8. Do not press the Work Button for more than 1 second each time to avoid burning skin.
9. Apply repair cream onto the treated area, avoiding touching with water before scabbing.

Note:

1. Don't squeeze with force, move gently with circular motion.
2. It is recommended to use coarse needle (rotation speed: 1 to 2 level) to treat freckles. While the fine needles are used to remove moles, verruca, fat granules, plane warts, and acne. Refer to the following chart for more information:

Rotation Speed	1-2	2-3	3-4		5-9	
Needle	Fine Needle	Fine Needle	Fine Needle	Coarse Needle	Coarse Needle	Coarse Needle
Applicable Skin Conditions	Fat granules	Small moles / freckles	Skin conditions the size of a mung bean	Skin conditions the size of a mung bean	Skin conditions the size of a soybean	Skin conditions the size of a soybean

Caution

- Test on a piece of pigskin first before using the pen for treatment.
- Keep the pen away from children as it produces high temperature in an instant.
- Do not scratch the treated area to avoid infections, it takes 20 to 30 days to heal.
- Avoid dark pigmented food, such as soy sauce or stimulating food such as ginger or pepper during the 3 month recovery period.
- Fully charge the product every 3 to 4 months to maintain the battery lifespan.
- Do not dismantle the product.
- Keep away from water or other liquids.

尺寸：90*65mm

折页

材质：128克铜版纸

Model:
AJ-PCA025