



GLO SCIENCE

we all speak smile

User Manual

BURSTING WITH BRILLIANCE

Meet Your New
brilliant white smile
Teeth Whitening Device Set

Teeth Whitening GLO Vials®
(5 pack, Up to 20 Applications)

GLO® Vegan
Lip Care



time to take charge

5

It takes 5 hours to fully charge your GLO® device. Charge after each use so you'll always be ready to GLO®!



To charge, plug the small end of the included power cord into the bottom of your GLO® control. Connect the USB end into a power adapter or USB wall plug. On the go? You can also charge your GLO® device using a computer USB port.



Your device will blink a brilliant blue light when charging. When charging is complete, all eight bands will light up continuously. Your GLO® control will only turn on when charged, connected to the mouthpiece, and the GLO® button is pressed.

Scan the QR code to learn about your Brilliant White Smile Set.



heat up, chill out



There's no need to slow down while you work on a brighter smile. While using your Brilliant device, you have your hands free to multitask!

Complete 3 back-to-back 8-minute sessions totaling just 24 minutes. For best results, complete a 24-minute GLO® session for 5 consecutive days for your most brilliant smile. Touch up once a month or as needed.

Depending on your starting shade, desired results, and enamel sensitivity, you can continue whitening for up to 10 consecutive days.

let's get started



1

CHARGE: Your device will blink a brilliant blue light when charging. When charging is complete, all eight blue bands will light up and no longer blink.

1x ON



2

TURN ON: Connect your mouthpiece to the GLO® control. Push the GLO® button once to turn on your device. Your mouthpiece will not illuminate but will warm up while you apply lip care and whitening gel.



3

PREP YOUR MOUTH: Apply lip care to inside and outside of lips. Remove cap from whitening gel vial and squeeze a small amount of gel out of the vial onto the soft bristles. Apply gel in a circular motion to front of teeth that are visible when you smile. Each GLO® Vial will have enough gel for 4 sessions, plus extra for touch up whitening.



4

START GLO SESSION: Push the GLO® button once more to begin treatment. Your GLO® mouthpiece will illuminate bright blue. The 8 blue bands on the control will count down your 8-minute treatment.

6



5

FINISH SESSION: When your session time is almost complete, your GLO® control will begin to blink. At the end of your GLO® session, your control will beep, and your mouthpiece will shut off. This concludes your 8-minute GLO® session.

Each GLO® treatment consists of 3 consecutive 8-minute sessions. Spit out any excess GLO® whitening gel and reapply for the next session. You may also reapply GLO® lip care if desired. The orange heat indicator light will remain on, keeping your mouthpiece warmed and ready in between your GLO® sessions. **Repeat steps 2 and 3 for a total of three applications.**

2x OFF



6

COMPLETE TREATMENT: Double-tap the GLO® button to shut off your device. Your device will automatically turn off after 5 minutes of inactivity.

Disconnect your mouthpiece from your GLO® device and rinse the mouthpiece thoroughly to remove any excess whitening gel. Dry with a clean cloth.



7

smile without compromise

Transparency is at the core of our values.

Our unique whitening formulations are made from the highest quality ingredients. Our team of scientists meticulously test each formula to deliver visible results.



VEGAN



CRUELTY
FREE



GLUTEN
FREE



SENSITIVITY
FREE



HYDROGEN
PEROXIDE

powered by science

Subscribe & Save: Keep glowing all year long.
Save up to 15% with a subscription of enamel-safe,
sensitivity-free whitening vials.



america's smile expert™


Founded by world-renowned dentist, professor, and oral-health expert Dr. Jonathan B. Levine, GLO® Science was created where science meets soul. We believe a smile is worth a thousand words, so we create products using only the highest-quality ingredients to give you exceptional results.

And with each product sold, we donate a portion of proceeds to provide oral healthcare services to those in need through the GLO® Good Foundation. That's something to smile about.

To learn more visit GLOScience.com

10

often imitated, never duplicated

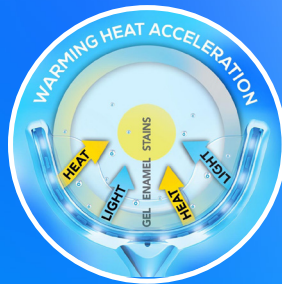

This device is as unique as your smile, because it contains G.L.O. Technology. Guided Light Optics or G.L.O. is a revolutionary smile care technology developed by Dr. Jonathan B. Levine.

This innovative technology combines gentle heat and light to accelerate whitening results for a brighter smile in a fraction of the time, without sensitivity.

To learn more about our technology visit GLOScience.com

11

look at you GLO®

Join the GLO Science® rewards program for insider-only benefits just for shopping with us. Get exclusive discounts, first access to new products, and more.

Sign-up at GLOScience.com/rewards to start rewarding yourself today.



we GLO® together

Save your receipt as proof of purchase. This will activate your one year limited warranty.

Register Your Device @
GLOScience.com/Registration

Drop us a line:
Leave a review,
we'd love to
hear from you!



smiles are best when shared

WE LOVE TO SEE YOU SMILE.

Share your grin using #ASmileThatGLOs for a
chance to be featured. @GLOScience



ingredients

GLO® Whitening Gels (GLO® Vials)

Net wt 5 x 2.9 ml (0.1 oz)

Ingredients: Distilled Water, Polyethylene Glycol, Hydrogen Peroxide, Potassium Nitrate, Glycerin, Edetate Disodium, Sodium Saccharin, Sodium Fluoride, Sodium Hydroxide, Flavor.

Warnings:

- If irritation (such as redness, swelling, soreness) of the gums or the mouth occurs, discontinue use and consult a dentist.
- Products containing peroxides are not recommended for use by children under 12 years of age.
- Use for periods of longer than 14 days is to be only under the supervision of a dentist.
- Avoid contact of the product with the eye.
- Avoid direct contact of the active surface of the tooth whitening product with the gums and/or salivary flow.

GLO® Vegan Lip Care

Net wt 0.1 oz, 3.0 g

Ingredients: Ricinus Communis (Castor) Seed Oil, Myrica Cerifera (Bayberry) Wax, Helianthus Annuus (Sunflower) Seed Wax, Mentha Piperita (Peppermint) Oil, Butyrospermum Parkii (Shea) Nut Butter *, Menthol, Medium Chain Triglyceride, Tocopherol, Triolein, Copernicia Cerifera (Carnauba) Wax *, Glyceryl Dioleate, Ceramide Np, Monoammonium, Glycyrrhizinate, Natural Flavor, Ethyl Alcohol *organic

Cautions: Consult with a dental professional before use if you have untreated dental cavities or any form of gum disease. If gum irritation occurs, stop using the product and contact your dental professional. Avoid getting the gel on fabric or clothing as it may lead to discoloration.



1 year limited warranty

GLO Science®, Inc. commits to the original retail purchaser that we will replace your GLO Brilliant® Personal Teeth Whitening Device and/or components within one year from the date of purchase, for any manufacturing defect in the GLO Brilliant® Personal Teeth Whitening Device or components. (Components include: the GLO® Mouthpiece and Silicone Case, GLO® Control, Lanyard and USB Cord, GLO Brilliant® Whitening Gels and GLO® Vegan Lip Care). If the product you have shipped back is a model no longer being made and is unavailable at that time, we will happily replace the product with the most recent model. To process a return, please call 855-GLO-5976 (855-456-5976) or email info@gloscience.com. We will provide you with a valid Return Authorization Number, which must be accompanied with your return together with a valid itemized receipt for the returned goods. To ship back your Device, for your convenience, we will provide you with a UPS prepaid postage label mailed directly to you for the return of your product.

If you choose to return your Device or components and not use the UPS label we provide, we will not be responsible for the product being lost or damaged.

Our GLO® Limited Warranty does not cover any of the following:

- Defects, damage or other conditions caused, in whole or in part, by misuse, abuse, neglect, negligence, alteration, accident, freight damage, improper repair, improper maintenance, improper, abnormal or excessive use and/or failure to obtain a replacement in a timely manner, as prescribed by GLO Science®, Inc.
- Defects or damages that are a result of an act of nature (such as a flood or earthquake) or acts of war or terrorism.
- Damage resulting from improper use of parts or accessories that are not compatible with GLO Brilliant® Personal Teeth Whitening System and/or components, including but not limited to the breakage or fracture of the GLO Brilliant® Personal Teeth Whitening System and/or components by any reason other than by manufacturing defect.

- Ordinary wear and tear on GLO Brilliant® Personal Teeth Whitening Device and/or components.
- Cosmetic damage, including scratches or dents.
- Accessories or parts that are not those of GLO Science®.
- Products purchased, new or used, from unauthorized distributors or through unauthorized distribution channels.
- GLO Science's only obligation under this warranty is the replacement of defective material that comprises the GLO Brilliant® Mouthpiece and the GLO Brilliant® White Smile Personal Teeth Whitening Device. GLO Science® shall not be liable and hereby disclaims any direct, special, indirect, incidental, exemplary or consequential damages or delays or losses under any other legal theory, including but not limited to, loss of use, loss of revenue, or loss of profits.

Warranty Disclaimer

This warranty is GLO Science's only warranty and is in lieu of all other warranties, express or implied. GLO Science® makes no implied warranties of any kind, including any implied warranties of merchantability or fitness for a particular purpose.

All of us at GLO Science® are committed to making your GLO® experience as positive as it can be. Simply contact us anytime at 855-GLO-5976 (855-456-5976) or info@gloscience.com and we will be happy to assist you. Visit www.gloscience.com for news and updates on GLO Science®.

brilliant white smile

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

- Consult the dealer or an experienced radio/TV technician for help.

Battery Disposal Instructions:

Batteries should be recycled or disposed of as per state and local guidelines.

GLO® (stylized) and GLO Brilliant® are Registered Trademarks of GLO Science®, Inc. GLOing, G-Vial, G.L.O. and Guided Light Optics are Trademarks/Service marks of GLO Science®, Inc. All right reserved.

Electrical Specification:

Input: 100-240 VAC, 50-60 Hz, 0.2A

Output: 5 VDC, 1A

Light output: 10/0.085W 470 nm Blue LEDs

Made in China/Designed in the USA

Distributed by GLO Science®, Inc
10 W. 37th Street, New York, NY 10018

©2022 GLO Science®, Inc
P/N# 6581042/6581042BB