

CORDLESS BACK CUSHION



Before use, please read the manual carefully
and keep it properly.

BELMINT
HOME SLEEK HOME™



CORDLESS BACK CUSHION USER GUIDE

USER GUIDE
for the
Belmint Cordless
Back Cushion

Model: BEL-CBC

Voltage: 13VDC 2000mA

Power: 26 W

Frequency: 50Hz/ 60Hz

Congratulations on your purchase of a Belmint Cordless Back Cushion

Your Belmint Cordless Back Cushion gives you a spa-quality massage in the palm of your hand. The portable massager lets you relieve tired and sore muscles at home, in the office, or while away on travel.

Genuine shiatsu back massage relieves muscle tension and pain with high-definition, selective area targeting. The Belmint Cordless Back Massage Cushion is perfect for anyone with casual or chronic back pain from overuse, exercise, or health conditions. Massage for up to 1.5 hours at a time for the perfect combination of relaxation and pain relief.

Feel the comfort of a soothing spa-quality shiatsu back massage with the Belmint Cordless Massage Cushion. Choose from three massage zones to focus on your upper, middle and lower back or use the spot massage to target specific areas of pain. The rechargeable battery lasts for 90 minutes of blissful massage on a three to four-hour charge. The auto shut-off allows you to set 15-minute increments for the perfect massage session.

Portable and lightweight, the Belmint Cordless Back Cushion is perfect for use at home, in the office, or while you are traveling. The easy-to-use controller gives you the ability to change the massage zone or time with the push of a button. With the Belmint Cordless Back Cushion, relief is at your fingertips.

To celebrating your comfort,
The Belmint Team

IMPORTANT INFORMATION

IMPORTANT:

Complete the information below and retain for your records. You may need this information when calling Belmint Customer Services at:

845-388-1270

MODEL NUMBER:

BEL-CBC

.....

PURCHASED AT:

.....

DATE PURCHASED:

.....

BEFORE YOU GET STARTED:

WHAT'S IN THE BOX?

Let's do a quick check of what's included in your purchase. You should have the following items:

- 1 – Back Massage Cushion
- 1 – 2000 MAH Rechargeable Battery
- 1 – Controller



PRODUCT DETAILS:

CAUTION:

All servicing of this massager MUST be performed by authorized BELMINT service personnel ONLY.


NOTE:

The Belmint Cordless Back Cushion will automatically shut off after 15 minutes for your safety.

BELMINT
HOME SLEEK HOME™



INTENDED USE:

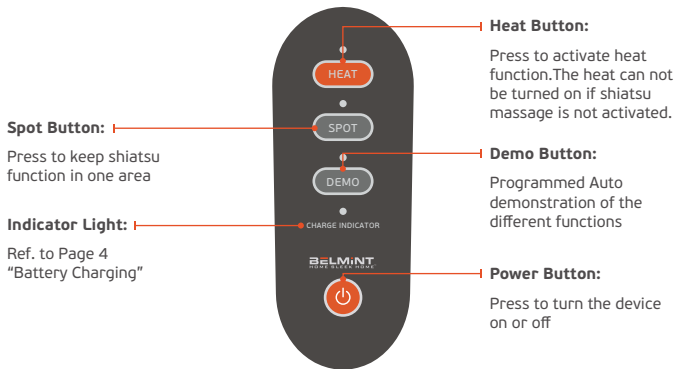


- This device will switch off automatically after 15 minute timer run out.
- Do not set up or use the massager in a bathroom or similiary wet / damp areas.
- Do not use this device while driving.

- For your own safety we equip the massager with a protection against over-heating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk and after a runtime of 15 Minutes. The massager will switch off automatically and will only allow switching on again after an appropriate downtime.
- Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 Minutes.

BATTERY CHARGING:

1. Position the massager on the chair or other support.
2. Connect the adapter cable to the corresponding cable on the cushion.

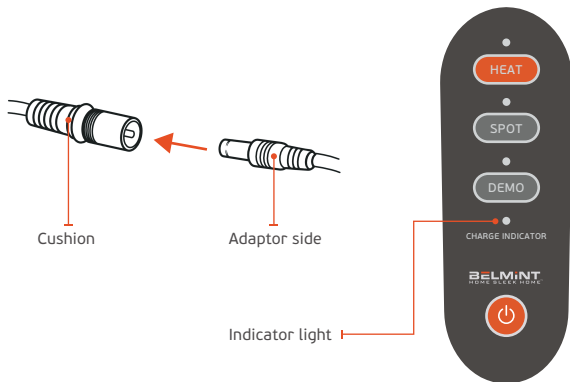

3. Plug the charging adapter into an outlet.
4. Normally it takes 2.5 hours to fully charge the battery.
5. Indicator light :

While charging, the indicator light will turn 'RED' and blink

When fully charged, the indicator light will keep 'RED' and not blink

(If the battery is full, the massager can work continuously for around 80 minutes.)

When the battery is under 20%, the indicator light will turn 'Green' and blink.

When the battery is over 20%, the indicator light will keep 'Green' and not blink.

6. It is safe to use the massager while charging.

NOTE: Attention ! Please only use provided charging adapter.

SAFETY INSTRUCTIONS:



PLEASE READ THE FOLLOWING INSTRUCTIONS CAREFULLY BEFORE PUTTING YOUR MASSAGE DEVICE INTO USE TO ENSURE TROUBLE-FREE OPERATION AND OPTIMUM EFFICIENCY. PLEASE RETAIN THESE OPERATING INSTRUCTIONS FOR FURTHER USE!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- **DO NOT** wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- **DO NOT** use any damaged cables, plugs or loose sockets.
- **NEVER** operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- **DO NOT** use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- **DO NOT** use this device while driving.
- **DO NOT** use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

CARE AND CLEANING INSTRUCTIONS:

1. Do not wash in washing machine or immerse the product or the adapter in water.
2. Always switch off before cleaning and remove the adapter from the socket with clean, dry hands.
3. The massager should be wiped with a slightly moist cloth and dried thoroughly. Do not wipe the product with alcohol or harsh chemicals such as thinner or bleach.
4. Keep away from heat, fire or exposure to direct sunlight.
5. Do not apply excessive pressure or impact onto massager/massage heads. Application of a direct force in excess of 110 lbs. may damage your massager and invalidate your warranty. Refrain from over usage of the device.
6. Store the massager in a cool, dry place out of reach of children.
7. Avoid contact with sharp edges or pointed objects which might cut or puncture the cover of massager.

TROUBLESHOOTING GUIDE:

Fault	Solutions
Massage heads rotating at reduced speed.	The load on the massage heads is excessive. Reduce load and try again.
Device switches off suddenly	Device may have been in use for 15mins. Auto-shutoff function activated. Continue using by pressing the 'power' button after resting the unit for 15mins.
	Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30mins before restarting.
Device fail to start	Check if device has been used for more than 15mins. If so, let device rest for 15mins before continuing. Device may require longer time to cool in warmer environments.
	The battery is low, please charge in time.
	Check if adapter is plugged into socket and adapter output is connected to device power inlet wire .
Jerky motion of massage heads / Vibrations from device	Excessive pressure may have been applied onto the device. Reduce the load and try again.
Infra-red lights are on but no motion from massage heads	Device may be jammed. Remove pressure / load from the device and press the 'Power' button to try again.

TROUBLESHOOT GUIDE:

If you encounter operational issues with your Belmint Cordless Back Cushion, please refer to the Tips and Guidelines listed below.

If the problem persists, please turn the POWER OFF and unplug the unit or call Belmint Customer Service at **845-388-1270**.

DO NOT attempt to make adjustments to the plug or alter the plug, cord, or massage unit. Doing so may nullify the warranty.

CONDITION	REASON	SOLUTIONS
Massager won't start or turn ON.	Unit may be unplugged.	Check that the unit is plugged into the appropriate outlet or power source. Check that the unit is plugged into the adapter. Press the POWER button. You should feel it click. If problem persists, please contact Customer Service at 845-388-1270 .
Massager stops working.	Unit may be unplugged.	Check that the unit is plugged into the appropriate outlet or power source. It may have come unplugged.
	Unit may have automatically shut OFF.	Press the POWER button. The unit may have shut OFF due to the automatic timer.
	You may have accidentally pressed the POWER button and turned the massager off.	Pressing the POWER button will also turn the massager back on if you inadvertently turned the unit off.
	Unit may have overheated.	Unplug unit from wall and let it cool before using again.

