

▼
**PILLOW
MASSAGER**
▲



www.belmint.com



support@belmint.com



1 (844) 423-5646




@belmintproducts



@belmintproducts

PILLOW MASSAGER

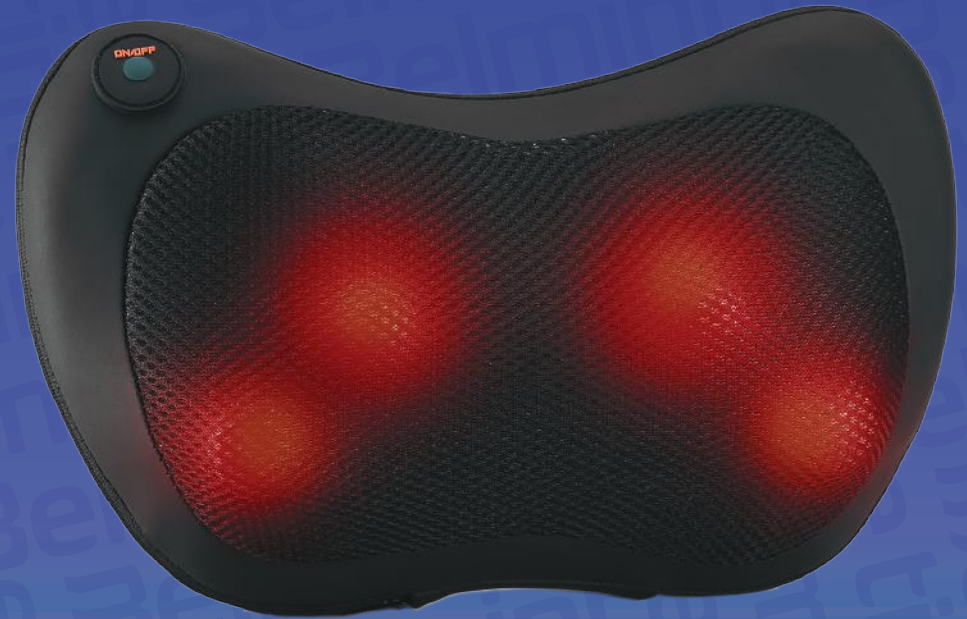

Congratulations on your purchase!

You're on your way to Spa-level experience

Your Belmint Pillow Massager electronic shiatsu massager with heat makes it possible for you to get a spa-level massage anywhere, anytime.

Like all Belmint products, your pillow massager is crafted with high-quality, durable materials to ensure your comfort. Specially designed to ease tension in tired and sore muscles, your massager is equipped with rolling and kneading nodes that enables you to target those hard-to-reach areas of your body. For a truly custom massage experience, your belmint pillow massager is equipped with optional heat, directional control and modular speed control because some days we need more kneading than others. Cleanup and storage are easy. Simply wipe with a damp cloth after use and store it in its included carry case.

You won't want another pillow massager. You've got a Belmint

To celebrating your comfort,

The Belmint Team

Sometimes peace is hard to find

But good service shouldn't be

At Belmint, we strive to provide our customers with ultimate relaxation.

If you encounter any issues with the performance of your device, our authorized service personnel are here to help. Check back cover for information.

Massage Position Guide



Legs



Abdomen



Back



Neck

Parts & Features



1 User Manual



1 AC/DC Plug Adapter



Power & Heat Button

Technical Details

Model	Power Consumption	Auto Shut-Off	Voltage
BEL-CPM	20 W	15 minutes	12V

Massager Functions

ON

Press once to turn unit on.

DIRECTION

Press again to change node direction.

HEAT

Press again to turn heat off.

OFF

Press again to turn unit off.

When Things Go Wrong

ISSUE	REASON
Unit won't start or turn on	<ol style="list-style-type: none"> 1. Check that the unit is plugged into the appropriate outlet, or power source. 2. Check that the unit is plugged into the adapter. 3. Press the POWER button. You should feel it click
Unit Stops Working	<p>You may have inadvertently pressed against the POWER button thus turning it off.</p> <p>15 minute Auto Shut-Off may have engaged</p>

Words to the Wise

How to prevent damage and injury

Please do not disassemble or open the product. Do not touch the inner parts of the product to avoid any trouble or more damage.

Unauthorized repair of the product will VOID THE WARRANTY.

Power Adapter and Power Cable

Please check if the power adapter is matched with the specified voltage with the output.

When plugging out, hold the power plug and NOT the power cable. The power cable should not be pushed down by any other heavy materials or pressure.

Controller Safety Information

1. Don't use your wet hand to operate the controller.
2. Don't drop the controller to the floor or in water.
3. Don't press or drag the cable of the controller.
4. Don't put any heavy goods on the controller.

For your Safety

1. When using the massager, the massager head should be near the two sides of neck. Otherwise, it will result an ache.
2. Please don't use massager when sleeping.
3. Please don't use this product when you're under the influence of alcohol.
4. Please don't put the massager wheel directly near the hand, elbow parts, knee, or breast.
5. Please don't use the massager for more than 30 minutes.

