

GlotechTM Mask Lite

QUICK START GUIDE

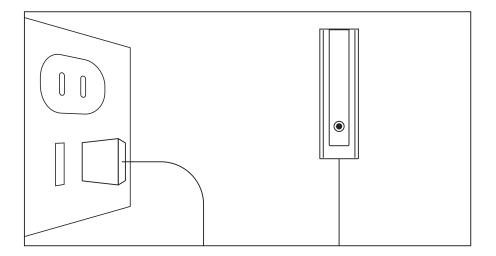


Charge the battery pack using the supplied cable. Insert the small end of the cable (mini USB) into the battery. Insert the large end (USB-A) into a wall plug.

Once fully charged, the indicator light will be illuminated steady blue.

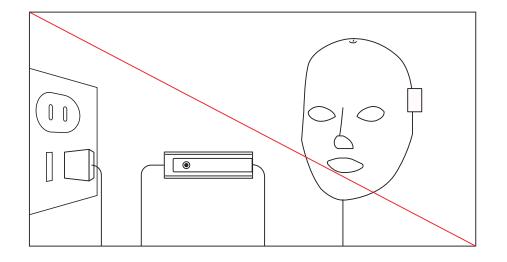


NOTE: The battery must be charged on its own (not attached to the mask) as this may cause damage to the battery.


Detach the battery from the charging cable, before plugging the battery into your mask.

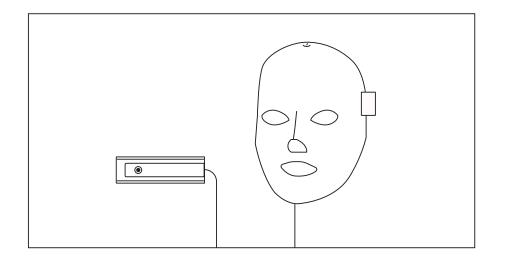


Connect the mask to the battery using the built-in cable on the mask. Your mask is operated by the supplied battery pack, making it completely portable.

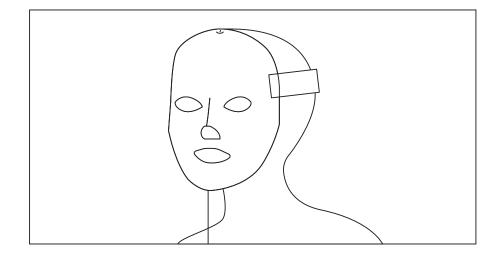


Adjust the strap to fit comfortably but snug on your face. You may find placing the strap higher on the crown of your head makes it more comfortable and secure.

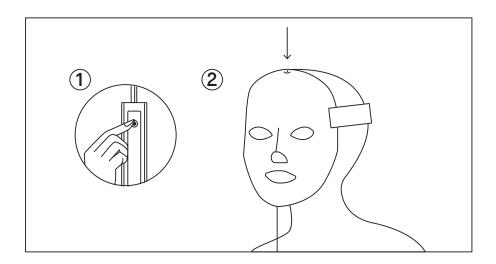

- 1 Activate the battery by pressing the button on the battery so the indicator light illuminates blue.
- 2 Turn on your mask by pressing and holding your thumb to the power button on the mask for 3+ seconds.

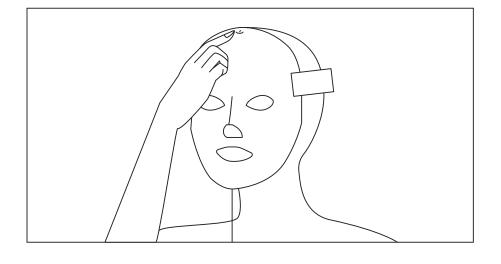

The default mode is red light. To switch to blue light, short-press the button on your mask. To turn off your mask before the auto 10-minute timer, hold the power button on the mask for 3 seconds.

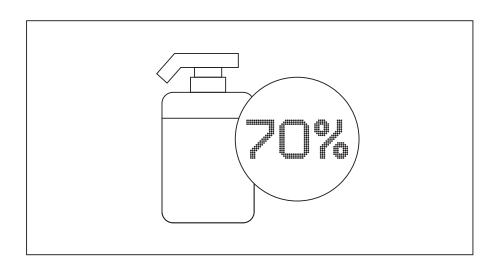

Clean your mask after each use with a damp cloth or 70% diluted alcohol